

VERSION 2

NOV 08, 2022



WORKS FOR ME 1

## 🌐 Implant Surgery: Chronic recoverable Neuropixels in mice V.2

**This protocol is published without a DOI.**

Emily A Aery Jones<sup>1</sup>

<sup>1</sup>Stanford University



Emily A Aery Jones  
[Stanford University](#)

COMMENTS 0

### ABSTRACT

This protocol collection explains how to build a low-cost, lightweight system to implant Neuropixels 1.0 probes into mice, record during freely moving behavior, then recover the probe for future use. This protocol explains how to implant a Neuropixels probe, previously prepared by following the *Assembly* protocol instructions, on a mouse.

Chronic recoverable implant based off of Juavinett et al, 2019 (<https://elifesciences.org/articles/47188>), then modified by Janna Aarse, then modified by Emily Aery Jones.

See also other chronic recoverable designs:

\* for rats: Luo & Bondy et al, 2020 (<https://elifesciences.org/articles/59716>)

\* Neuropixels 2.0: Steinmetz, Aydin, Lebedeva, Okun, & Pachitariu et al, 2021  
(<https://www.science.org/doi/10.1126/science.abf4588>)

\* overview: Neuropixels 2021 course lecture 3.7 by Yoh Isogai ([https://www.youtube.com/watch?v=7cZqYliGvBQ&list=PLfhWmWntvjI64ti\\_a-MzHlwqwEU0ZIALb&index=18](https://www.youtube.com/watch?v=7cZqYliGvBQ&list=PLfhWmWntvjI64ti_a-MzHlwqwEU0ZIALb&index=18))

### PROTOCOL CITATION

Emily A Aery Jones 2022. Implant Surgery: Chronic recoverable Neuropixels in mice. **protocols.io**  
<https://protocols.io/view/implant-surgery-chronic-recoverable-neuropixels-in-cizbuf2n>  
Version created by [Emily A Aery Jones](#)



## KEYWORDS

electrophysiology, silicon probe, entorhinal cortex

## LICENSE

————— This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

## CREATED

Nov 08, 2022

## LAST MODIFIED

Nov 08, 2022

## PROTOCOL INTEGER ID

72451

## GUIDELINES

- Durations are estimates given an experimenter with some stereotactic surgery experience. Once practiced, a skilled surgeon can complete this procedure in 3 hours (excluding set up and clean up). Given the duration and difficulty of this surgery, it is recommended that experimenters with no surgical experience first become comfortable with headbar attachment surgery and craniotomy surgery before attempting.
- This surgery is adaptable for headfixed recordings. First, perform this surgery, stopping right before the *Implant the probe* section and skipping the *Drill the craniotomy* section, performing the *Ending the surgery* section as normal except skipping the dexamethasone. Train mice until performance criteria are reached. For second surgery, repeat this surgery, mounting directly in headbar stereotax adapters. Skip the shaving and scalp removal steps, skip the *Level and mark skull* section, drill craniotomy through the metabond, and skip the *Mount skull components* section except for mounting the well. Continue through *Ending the surgery*.

## MATERIALS TEXT

### Tools:

- Skinny electric razor
- Compressed air
- Dovetail holder, stereotax adapter, 0.9mm hex key (Sensapex uMp-NH & uMp-NPR-200, McMaster-Carr 6958A22)
- Headbar holder, headbar stereotax adapters, cannula stereotax adapter, headbar screws (18-8 Stainless Steel Socket Head Screw 0-80 Thread Size, 1/8" Long, McMaster 92196A052) & hex key for headbars & stereotax
- Forceps
- Bent forceps
- Scissors
- Serrefines
- #3 scalpel handle
- Tiny flathead screwdriver (Wiha, 0.8 x 40)
- Sharpie pen size 005
- Dental well mixer plate (store in freezer)
- Dental cement powder scoop, solvent dropper bottle, & rubber well
- 20uL pipet
- UV glasses (Grainger)

### Consumables:

- All materials from the chronic Neuropixel assembly protocol
- Isoflurane
- Absorbent pads
- Drill bits (0.5mm, Meisinger, HM1-007-FG 0.7mm for Robodril & HM1-005-HP 0.5mm for air drill; 0.1mm carbide burr, Beavers Dental BWFG9903 for air drill)
- Nair
- Iodine (povidone-iodine solution, 10%)
- 70% ethanol
- Small weigh boats
- Cotton swabs (Puritan)
- Precision applicator brushes (Parkell S379)
- Metabond (Parkell S380, clear L-powder)
- Dental cement powder & solvent (Ortho-Jet, Lang Dental, 1303CLR)
- Kim wipes
- KwikCast
- Toothpick
- Refresh Lacrilube Ointment
- Dow Corning 3-4680 Sil Gel Kit, pipetted into 1.5mL Eppendorf tubes
- Light Cure Acrylic (Vivid Flow Syringe 1g w/Tips A1, D33-0110)
- Electrical tape
- Copper tape (McMaster-Carr 76555A724)

#### *Sterile Consumables:*

- Sterile pads (Busse 696)
- Insulin syringes
- Scalpel blades (Henry Schein, 11)
- Sugi Eyespear Cleaner (Kettenbach 30601)
- 0.5mL Eppendorf tube
- 20uL tips
- Isopropyl alcohol
- DiD Cell-Labeling Solution (Invitrogen V22885)

#### *Injections:*

- Saline
- Opioid pain reliever (e.g. Buprenorphine)
- NSAID (e.g. Rimadyl)
- Antibiotic (e.g. Baytril)
- Steroidal anti-inflammatory (e.g. Dexamethasone)

#### *Animal support:*

- Hydrogel (ClearH2O, 70-01-5022)
- DietGel (ClearH2O, 72-07-5022X)
- Sprinkles or forage mix

#### *Equipment:*

- Benchmark Scientific Autoclave
- Germinator Dry Bead Sterilizer
- UV light
- Smart Weigh 500g Scale
- FHC Heating Pad
- VetEquip Isoflurane Machine with tubing to nosecone and chamber
- pressurized O2 and air
- Leica IC90E camera, M80 stereoscope, & KL2500 light
- Neurostar Stereodrive
- David Kopf Stereotax with attachments for Neurostar drill and stereotactic rod holder
- Schott gooseneck lamp
- NSK Pana-max Pax-Tu-M4 air drill
- Insignia Mini-fridge & freezer (for injectables and well mixer plate)
- Vivid Light Cure System (003-0042)
- Stryker T/Pump hot water pad

## BEFORE STARTING

- At least 1 day prior to surgery, singly house mouse in post-operative conditions to help mouse overcome neophobia of supporting materials. Prepare a fresh cage with food on the floor, half packet of gel food, half packet of gel water, and sprinkles/forage mix to encourage retaining weight during recovery. Add a nestlet and additional nesting material to support heat retention in single housing. Add enrichment such as chewing materials. Remove any overhead structures like huts or hoppers to prevent implant getting caught.
- Machine headbar stereotactic adapters from files at [https://github.com/emilyasterjones/chronic\\_NPX\\_mouse/tree/main/headforks](https://github.com/emilyasterjones/chronic_NPX_mouse/tree/main/headforks)
- Plan your probe trajectory using this tool: [https://github.com/petersaj/neuropixels\\_trajectory\\_explorer](https://github.com/petersaj/neuropixels_trajectory_explorer)

## Prepare surgical tools and field

- 1 Sterilize tips of metal instruments, ground screw, and headbar in autoclave ⌚ 00:25:00 or hot bead sterilizer ⌚ 00:00:05 . Disinfect cotton swabs, toothpick, and a weigh boat under UV light ⌚ 00:30:00 .

### Note

Sterilize all components that might contact the brain. Silicone sealants cannot be sterilized. Only contact these with sterile instruments and keep containers closed when not in use.

- 2 Disinfect surgical field with 10% bleach. Lay absorbent padding across field. Lay down a sterile drape and place tips of metal instruments on this field. Load scalpel blade and drill bits and keep tips over sterile field.

### Note

Always return sterile tips to sterile field when not in use. Do not place other items on this field.

- 3 Draw injectables into syringes and place needles on the sterile field or cap them. Place iodine and ethanol each in a weigh boat with 3 cotton swabs. Place remaining tools and consumables and implant assembly on nonsterile pad or within easy reach.

#### Safety information

Only recap needles using the 1-handed method



Sterile drape (top) with sterile tips, needles, drill bits, and implant components. Non-sterile drape with nonsterile components and capped needles.

- 4 Turn on surgical heating pad, lights, scope, O2, and air flow.

30m

## Prepare animal

- 5 Weigh mouse and record weight.

- 6 Set O2 flow rate to 1-1.5L/min and isoflurane to 3%. Place mouse in anesthetic chamber.

#### Safety information

Use vacuum systems connected to isoflurane lines and work on a downdraft table to move uninhaled isoflurane away from the surgeon.

- 7 When mouse is immobile, remove from chamber and begin to shave. Shave back to front from behind the ears to just behind the whiskers, and side to side just above the eyes and from the outer edges of the ears. Return to chamber when moving and repeat until area is completely shaved. For final pass, apply Nair to entire shaved area, wait up to 1 minute, then wipe clean with a cotton swab. This ensures no follicles remain in the incision area.
- 8 When breathing has slowed to 1Hz, move animal to the toothbar. Switch isoflurane from chamber to nosecone. Wait until unresponsive to pedal reflex test, then lower to 1.5% isoflurane.

#### Note

Monitor and record anesthetic depth throughout procedure. Ears and toes should be pink, respiration should be once per 1-2 seconds, and pedal reflex should be lost. If respiration slows or extremities start to lose color, lower isoflurane levels. If respiration speeds up or reflex is present, increase isoflurane levels and halt the procedure until animal is fully anesthetized again.

- 9 Place lacrilube ointment on eyes. Inject dexamethasone and buprenorphine.

#### Note

Giving dexamethasone pre- and post-operatively and maintaining sterility will prevent microglia from encapsulating the probe and reducing cell yield.


- 10 Immobilize skull in earbars. Adjust nosecone AP and DV position as needed to ensure head is flat and centered. Set earbar on dominant hand side (at 5th tickmark), then hold neck scruff with dominant thumb and middle finger with pointer finger on the nose to manipulate the head. Position and secure earbar on nondominant side (at 5th tickmark). Tighten on both sides, then test by pressing down next to each ear. Confirm immobile and level before proceeding.



- 11 Sterilize scalp by passing iodine and ethanol with a cotton swabs alternately three times.



- 12 Cut and remove the scalp from the skull. Lift the middle of the scalp well above the skull and make a small horizontal cut. Then, make small circular cuts around the first incision. The anterior boundary is a couple millimeters behind the eyes. The posterior boundary is where interparietal and occipital bones meet (suture behind lambda suture). Whenever possible, move skin away, rather than cutting. Hold remaining skin open with serrefines. Tighten earbars, which are often slightly loose following incision.
- 13 Clean skull with scalpel. Scrape away any connecting tissue and muscle connections to the back of the skull. Superficial bleeding will occur; clear with cotton swabs.

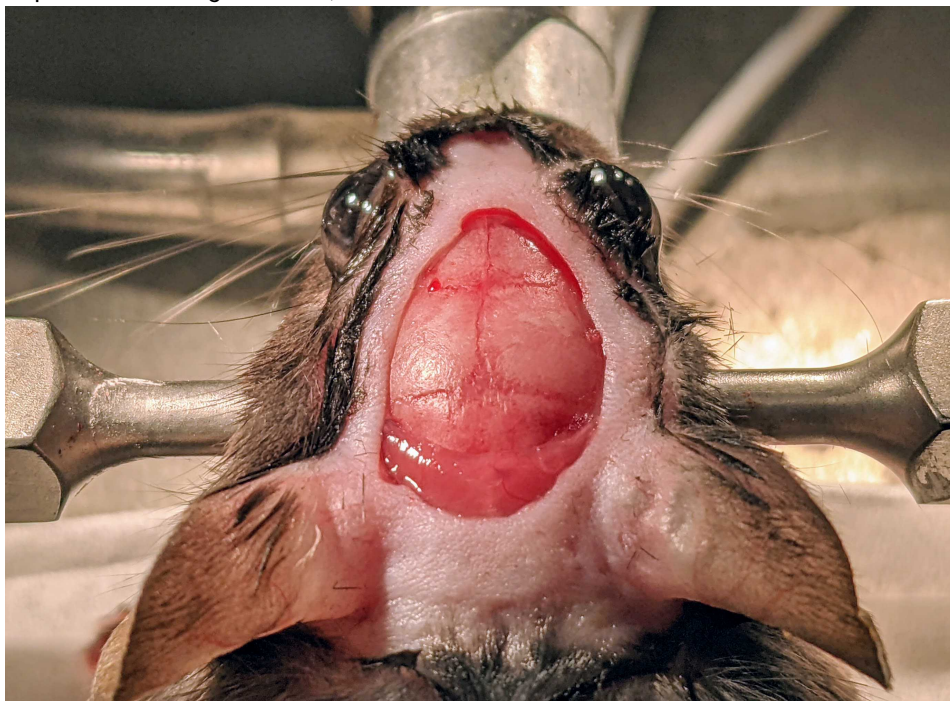

Mouse immobilized in earbars, scalp shaved to eyes and ears, scalp excised to eyes and most posterior suture within shaved boundary, and skull scraped clear.

20m

## Level and mark skull

- 14 Find bregma, lambda, 2mm left and right of midpoint along midline suture with drill. Record changes in DV between these points. Level midline axis by adjusting nosecone DV, then repeat measurements. If using Neurostar system, use Adjust for tilt and scaling feature to adjust coordinates.
- 15 Make two marks to indicate the ML distance for the two MECs. Go to coordinates AP -3.9, ML +/-3.3. Bring the drill down to the surface of the skull and drill a superficial hole. Remove shavings with air duster, then fill holes with Sharpie.

### Note


This step is specific to MEC recordings. Adjust coordinates as needed for your target structure. You can choose to implant into left or right MEC in the next steps.

- 16 Go to behind the sinus, midline (approx. coordinates AP -6.15, ML 0) to drill a hole over cerebellum for the reference screw.

#### Note

Referencing subtracts off the signal at the reference site. If analyzing LFP with an external reference, the reference screw needs to be somewhere electrically quiet, like cerebellum, so as not to introduce a new signal. If using an internal reference, this screw is just for grounding and thus can go anywhere. If not shorting reference and ground pads, add another screw anywhere on the skull and attach the wire from the ground pad (on the left, if the Sensapex holder is facing you during insertion) to that screw.

- 17 Score the skull with the scalpel to increase the surface area for metabond to bond to in later steps.

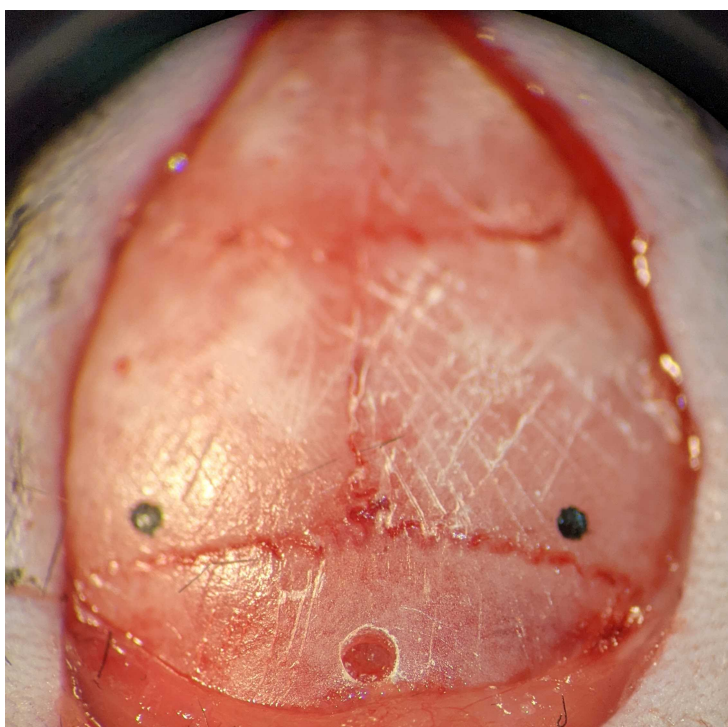

Right and left MEC marked, reference screw hole drilled, and skull scored.

## Drill craniotomy

- 18 Switch to cool light (e.g. from the scope) to better visualize the sinus. Using the air drill with the 0.5mm bit, shave down the skull in a ~1mm diameter circle behind (not including) your selected dot. Smooth down edges. Regularly stop to remove skull shavings with air duster and to clear and cool the area with saline to improve



sinus visualization and to cool the brain to reduce swelling.

#### Note

It is vital to visualize the sinus for good targeting, but not to drill into it, as this will cause severe bleeding. Once the sinus is visualize, shave down directly anterior to it to try to make the craniotomy as close to it as possible without including the sinus.

- 19 Once the skull is very thin above the craniotomy, use the fine-tip drill bit to carve a moat (approx. 0.5mm diameter) around the area that you will insert the probe (posterior to the dot, right in front of the sinus). You can test whether the skull is thin enough for this step by pressing gently on the skull with forceps and seeing how much it gives. If it's rock-hard, keep shaving it down.

#### Note

Making a small craniotomy will reduce the amount of brain exposed post-operatively, but also reduces options for probe insertion.

- 20 Once the moat has broken through the skull all the way around carefully lift the piece of skull off the brain using fine forceps. Some bleeding will probably occur when you lift off the piece of skull. Stop with Sugi spears. Use the fine-tip drill bit to smooth down edges, and use the bent forceps to look for hidden skull ledges to be removed. Try not to remove dura, and don't push the skull fragments towards the sinus.

#### Note

If you rupture the sinus, severe bleeding will occur. Administer fluids and absorb with Sugi spears until a clot forms. Terminate the surgery, as it is you cannot use this craniotomy and it is unlikely you will be able to securely cover the clotted craniotomy to prevent future rupture under the implant.

- 21 When the craniotomy is clear and finished, mark the sinus with a Sharpie and extend the dot anteriorly with a Sharpie, to improve visualization once the well is added in next steps. Apply KwikCast to seal the craniotomy from glue in subsequent steps.

45m

## Mount skull components

- 22 Affix screws in the drilled hole. Grasp with the bent forceps and place it on top of the hole. Use a screwdriver to slowly screw it in. Screw to the point where you can't pull it out with your bent forceps (about 3 rotations; no need to screw it all the way into the skull).

23

\*


OPTIONAL: Attach the probe to the stereotax (see next section) and tax to just above the dot. Tax up, remove the probe holder, and store affixed in a clamp out of the way. When you insert the probe in the next section, tax posterior by 1.4mm to get to the target location (AP -5.1 for MEC).

24

Select a well size that will allow you to access the full craniotomy with an angled probe, visualize the sinus and center Sharpie marks, but not be too large. Apply a thin layer of UV cure glue, position around the craniotomy, and cure with UV light. This allows rapid curing so the well can be held in place and won't move around during gluing. Once secure, apply and cure a second layer around the outside of the well to make sure it can hold liquid.

#### Safety information


Wear UV-filtering glasses when using the UV light.

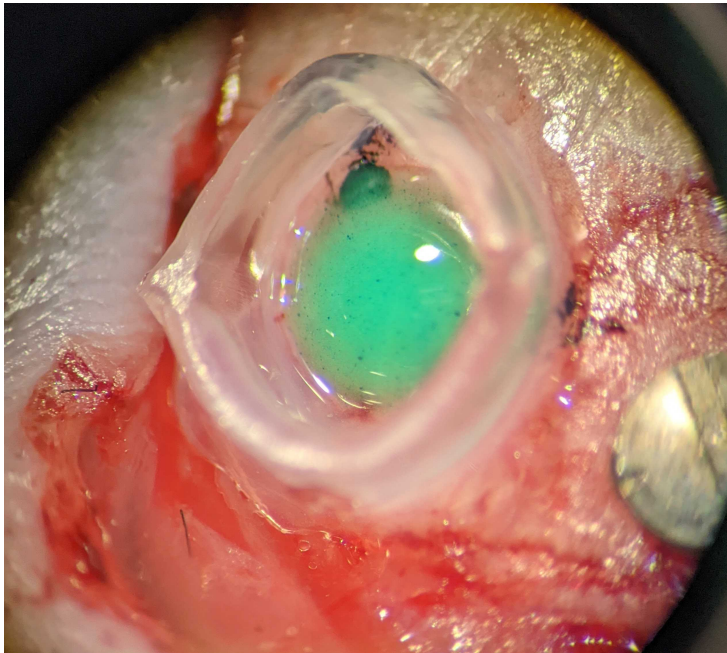

Craniotomy surrounded by a well. Dot anterior to center of MEC (top) and sinus mark (middle left) are visible. Craniotomy sealed with Kwikcast (green) to protect UV glue (clear) from invading craniotomy while well was secured. Ground screw visible in bottom right.

25

Place headbar into holder and screw it in. Lower onto the skull. Position just behind the eyes. Center the marker line over the central suture. Lift up 1/2 turn.

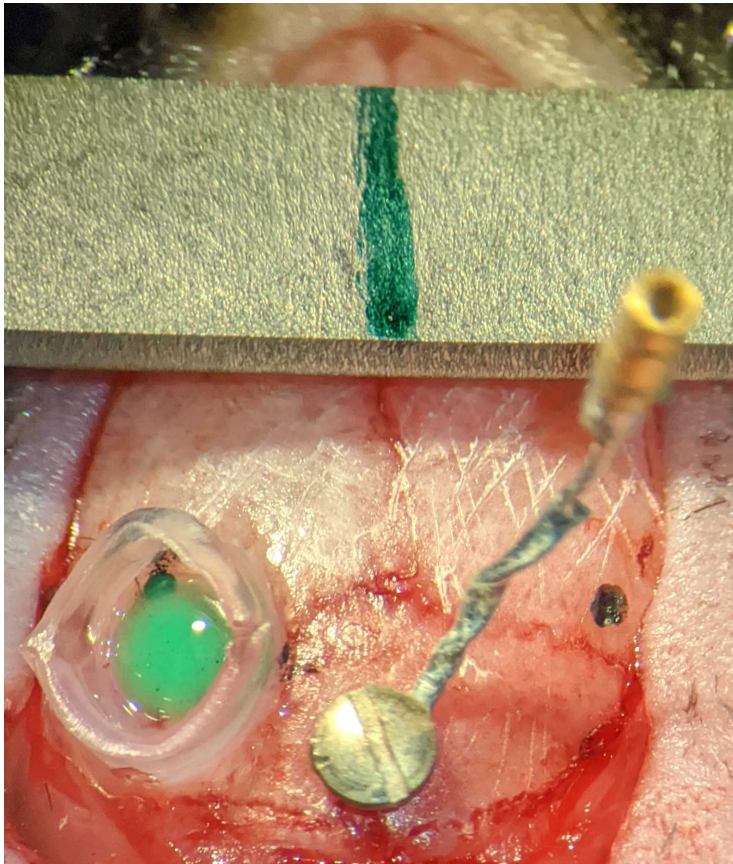


Positioning of skull components prior to metabonding. Headbar (top), skull screw (bottom), and craniotomy with well (bottom left).

- 26 Prepare metabond in a mixing well that has been chilled to -20C. Add 1 scoop of metabond powder, 4 drops of solution, and 1 drop of catalyst. Mix well and apply with precision applicators. Let metabond run under the headbar, then lower it down 1/2 turn to put it back on the skull. Use metabond to cover the ground screw (including the soldering junction) and the entire skull. Make a bridge over the headbar with metabond. Apply a thin layer to allow it to cure completely quickly.

#### Note

Metabond makes a strong connection to secure the headbar for headfixing and bonds well to bone, but is too runny to build up. Dental cement bonds well to metabond, but not bone, and can be used for construction. Covering the skull with metabond allows building atop with dental cement in future steps.

- 27 Allow metabond to cure completely for  00:10:00 .

10m

#### Note

Glues do not cure if wet. During gluing steps, constantly dry the margin (skin and muscle) and pull it away from the glued area with cotton swabs. While glue is curing, check the margins and pull off any excess glue so that none of the implant is attached to the skin, which can disrupt the implant when it is pulled during normal mouse behavior.


- 27.1 While waiting: Mount rod to Sensapex holder, then to stereotactic arm. Angle the stereotax to 8-10 degrees (for MEC) and tax to just above the target location, either by eye or using coordinates from the earlier measurements from the dot.
- 27.2 While waiting: Remove rod from stereotax and affix probe. Slowly and carefully dip probe into IPA 3 times to sterilize, then DiD 10 times to dye. Store on a clamp out of the way.
- 28 Release headbar from headbar holder and mouse from earbars. Mount headbar stereotactic adapters into earbar slots. Move nosecone anteriorly and ventrally until headbar is positioned snugly in headbar stereotactic adapters. Affix with headbar screws. Record nosecone position for explant surgery.

#### Note

Moving to headbar mount means that the insertion angle can be replicated exactly during explant surgery.

25m

## Insert the probe

- 29 Remove Kwikcast and wash craniotomy with saline. Run bent forceps around the rim of the craniotomy to ensure no dried tissue covers the brain; if so, remove it.
- 30  Mount probe in Sensapex holder to stereotax with wings facing rostrally and lower to brain surface. Record DV start coordinate. Make sure entry is unimpeded, advancing 800um quickly (e.g. 0.5mm/s) to clear dura. If the probe bends, stop immediately and try a different insertion point.

#### Note

If your target region is more anterior, you can flip the orientation so the wings are mounted to the back of the skull instead. Position the stereotax at a steep angle so you can visualize the probe through the middle slot of the body piece, watching to confirm that you don't hit the probe mount when positioning the scope. It is not recommended to rotate the probe 90 degrees (sites facing the ears instead of the nose/tail), as

the probe will experience AP shearing forces during animal movement and the probe cannot bend in this direction.

- 31 Drive probe at slowest speed (e.g. 3.3microns/s) to 3400 angled depth (for MEC), then retract 200 microns to ensure the probe will not hit skull when the brain settles. If the probe bends during insertion, retract until bend is relieved and either stop driving the probe or try another insertion. Driving to this depth takes 15 minutes.

#### Note

Move scope to view the probe from the side. If the probe bends during insertion, it will bend along the axis visible from the side.




Probe at final depth inside craniotomy.

1h

## Secure the probe


- 32 Load blue & clear Dowsil components into 1.5mL Eppendorf tubes for easy access. Mix 20uL each of blue & clear Dowsil gels in a 0.5mL Eppendorf tube. Draw up slowly as it is viscous. After mixing well, draw up 20mL, and apply gently to the inside wall of the well. If it spills out, wait to harden then scrape away.
- 33 Plug ground/reference pin into screw pin. Apply UV cure acrylic to the joint and cure. Bend the wires inward so they're flush along the back of the probe.



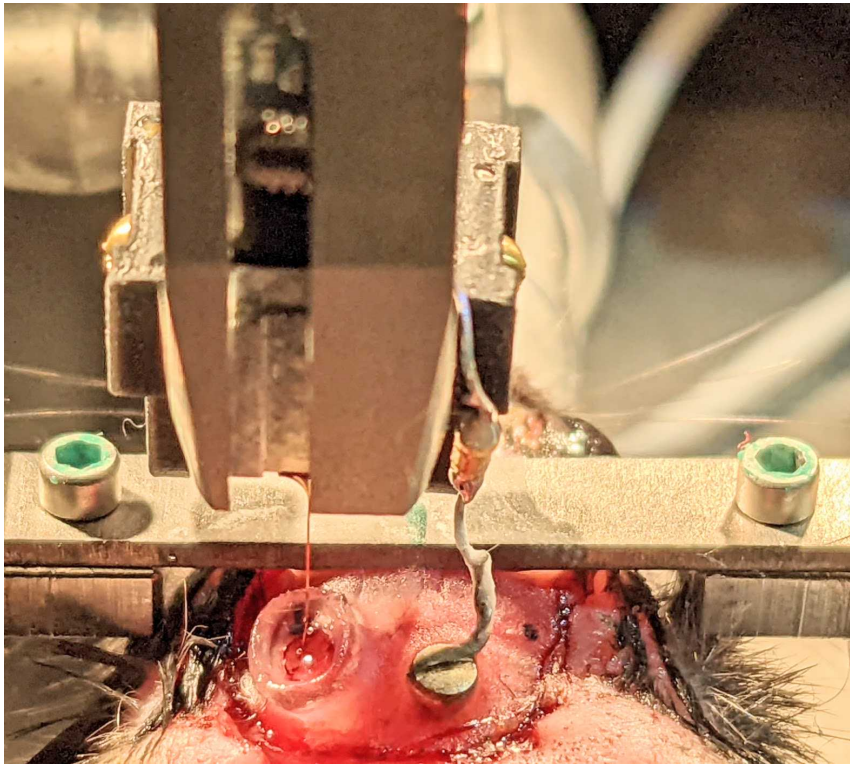


Mouse mounted in headbar stereotax adapter with skull components metabonded and probe at full depth. Craniotomy filled with Dow-Sil and ground/reference wire plugged into screw pin and secured with UV glue.



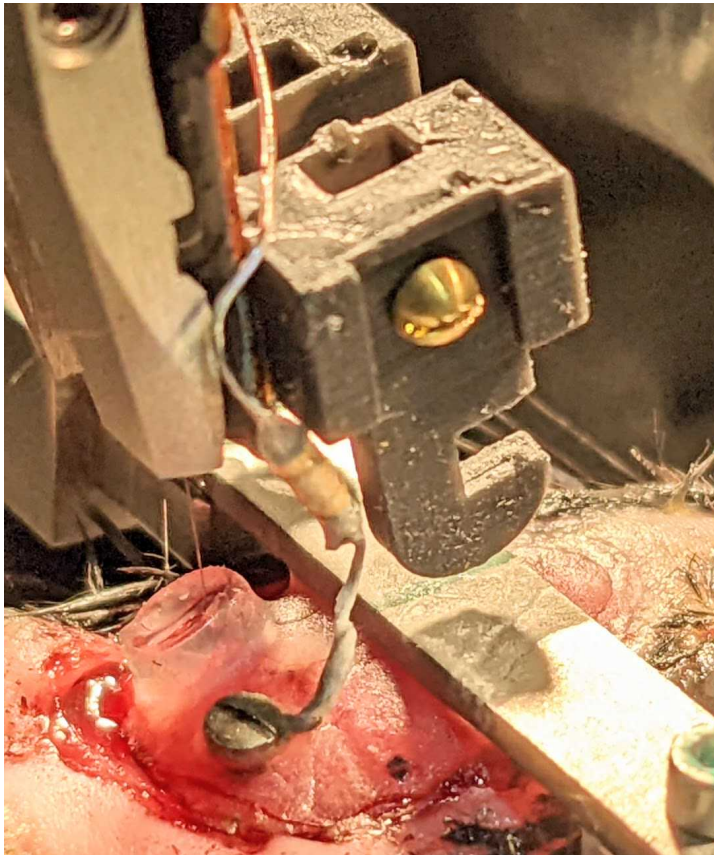


Same as previous figure, from the side.

- 34 Build up dental cement on the skull, covering the screw. Using thicker cement, build up from the headbar to the towards wings, covering up to the notch. Do not get any dental cement on the body piece. Build a cement bridge between the wings and from the wing to the screw. Apply cement to fill the gap between the well and the headbar on the left edge. Your goal is to create a ring of dental cement around the probe so that once the outer tape layer is applied, the bottom ring of dental cement will prevent anything from entering the probe zone from below while the tape secures anything from above.

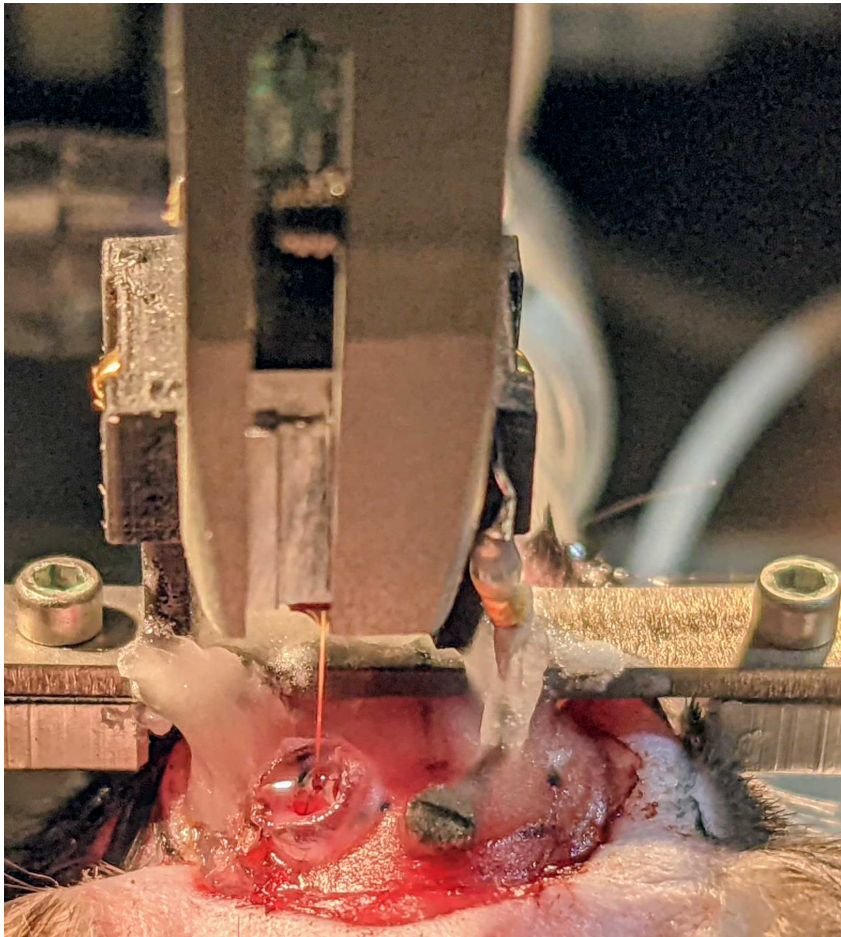


Dental cement connects the wings to the headbar and the skull and creates a ring extending from the lateral edge of the well around to the ground/reference wire. The ring will be completed with the dome piece mounted in the next step.

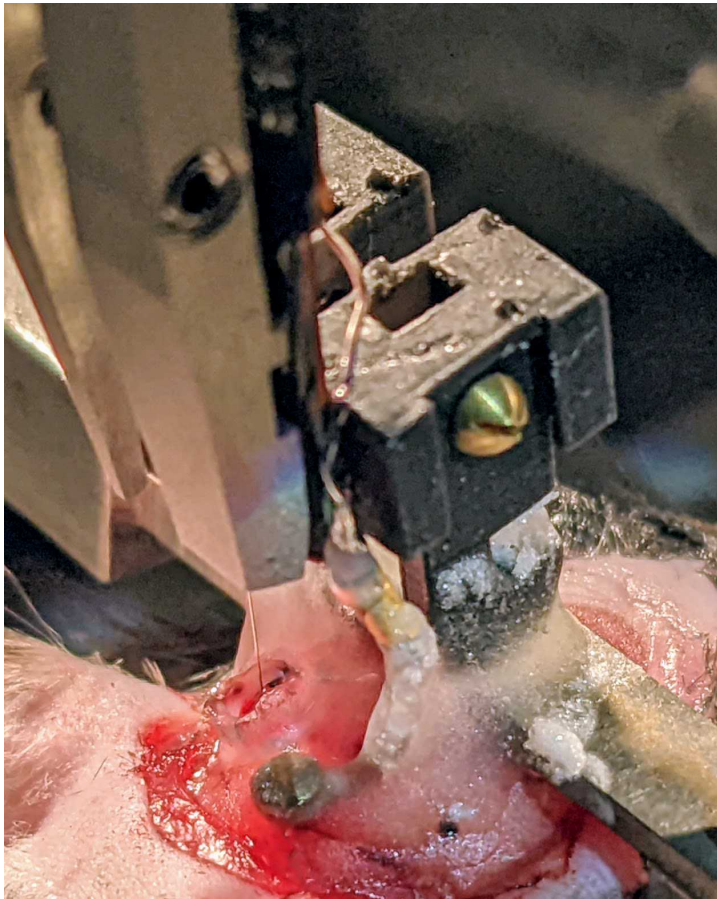

Same as previous figure, from the side.

- 35 While waiting for the dental cement to cure: apply a thick layer of UV cure glue and affix the dome just behind the well and screw. Make sure there is enough room for the Sensapex holder to open in the next steps. Cure the glue, then apply additional layers to secure it, pulling the neck muscle away from the glue so it securely attaches to the skull. Complete the dental cement ring by using dental cement to attach the dome to the left wing and to the screw.
- 36 Release the Sensapex holder from the probe and tax up.



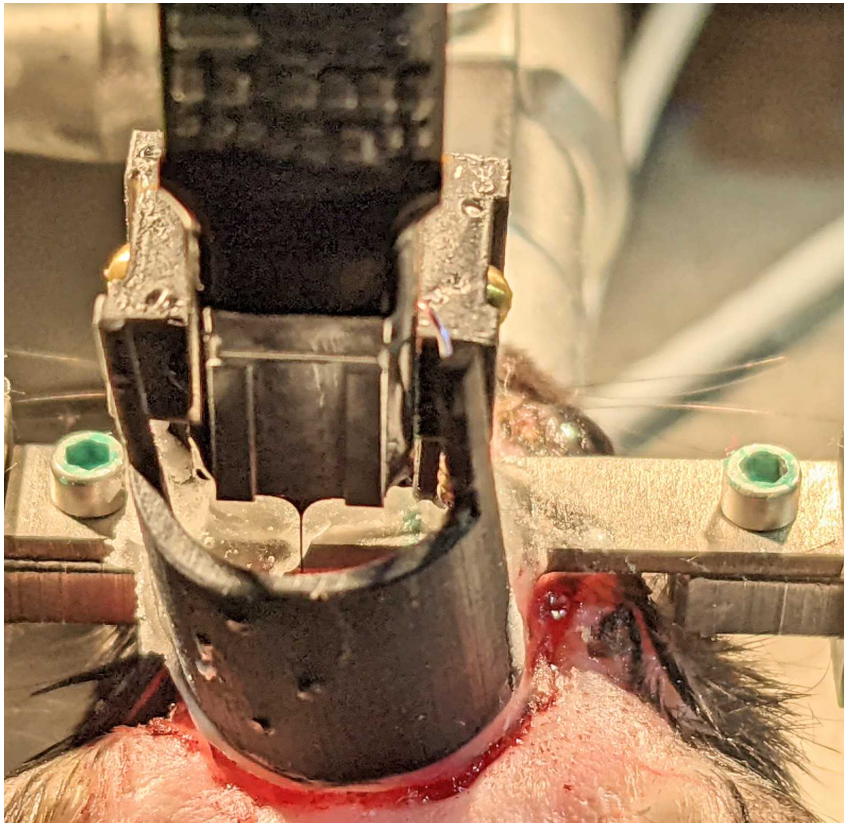


Dome is attached, dental cement ring is complete, and Sensapex probe holder is removed.



Same as previous figure, from the side.

- 37 Using electrical tape to close, snap the flex cable holders on, tab slot facing towards the body piece. Connect the dome antennae to both holder pieces with dental cement.

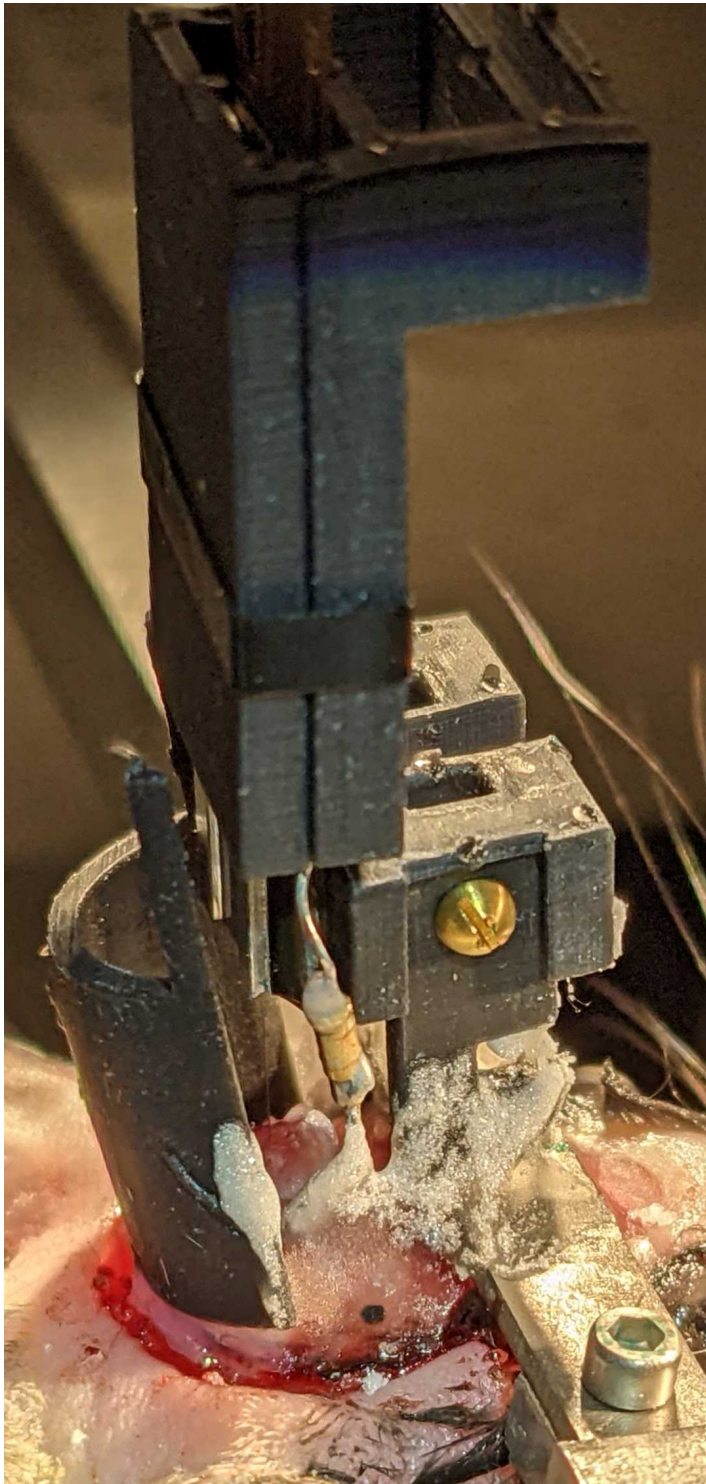

Flex cable holders cover the flex cable, held closed with electrical tape, and dental cemented to the left dome antenna. Design is symmetric, so that right MEC can be targeted (and right dome antenna attached instead).

- 38 Wrap in a single piece of copper tape. Make sure the bottom margin doesn't touch skin, but is low enough to not have any gaps above the dental cement. Press to secure to dome, body piece, and sides of flex cable holder. Make slits in front of and behind the flex cable holder and fold the tape in. Press in the remaining larger piece. Put a thin piece of electrical tape along the top margin to prevent it from releasing. Fold over flex cable



into tab slot and secure with a similarly dimensioned piece of electrical tape.

#### Note

Electrical tape is flexible and secures well, but has no structure. Copper tape is structured and provides a faraday cage around the probe shank, but isn't very sticky.



Headfixed implanted mouse with probe folded into tab slot.

5m

## Ending the surgery

- 39 Inject saline, rimadyl, and baytril. Turn off isoflurane and remove mouse from stereotax. Weigh mouse and record weight. Place mouse on heated pad in a recovery cage. Monitor mouse until walking normally.
- 40 Daily for 3 days, weigh mouse and administer saline, rimadyl, baytril, and dexamethasone. Keep free feeding for at least 3 days post-op to allow weight to stabilize. Monitor for any signs of pain or failure to recover. If needed, follow humane endpoint procedures as detailed in your animal protocol.

**41** Recording can begin 1 day post-op. Number of detectable neurons is highly variable for first 5 days, then remains steady.

Step 41 includes a Step case.

**Freely moving recordings**

**Headfixed recordings**