

Version 4

Apr 16, 2021

# nCoV-2019 sequencing protocol for illumina V.4

Version 1 is forked from [nCoV-2019 sequencing protocol](#)

Kentaro Itokawa<sup>1</sup>, Tsuyoshi Sekizuka<sup>1</sup>, Masanori Hashino<sup>1</sup>, Rina Tanaka<sup>1</sup>, Satsuki Eto<sup>1</sup>,  
Risa Someno<sup>1</sup>, Makoto Kuroda<sup>1</sup>

<sup>1</sup>National Institute of Infectious Diseases, Japan

1

*Works for me*

This protocol is published without a DOI.

Kentaro Itokawa  
National Institute of Infectious Diseases, Japan

## ABSTRACT

This protocol was originally forked from "[ARTIC amplicon sequencing protocol for MinION for nCoV-2019](#)" by Josh Quick to adapt for any illumina sequencers (iSeq, MiSeq, NextSeq, etc.).

Because the PCR products are fragmented and ligated with adapters, this protocol is compatible to many sequencing kits with various read lengths (75PE, 150PE, etc.) for Illumina machines.

While the library preparation uses QIAseq FX by Qiagen and is basically straight forward (as per kit instruction but modified to 1/4 scale), some tweaks for much of simplicity and speed were added.

## Change histories (V2):

1. The amount of cDNA input to multiplex PCR is increased (now almost same amount to the ARTIC Network's original protocol) (Step 10). This change gives better results for samples with extremely low RNA copy.
2. Amount of adapter solution input was corrected (Step 19).

## Change histories (V3):

1. RT step has been changed to use of LunaScript ([Tyson et al., 2020](#)).
2. Corrected typos and wrong descriptions about amount of reagents.
3. Added a link to library quantification protocol.
4. Other minor changes.

## Change histories (V4):

1. Added description about the primer set ver. [N1/N2](#).
2. Added suggestion to double the concentration of nCoV-2019\_74\_LEFT and nCoV-2019\_74\_RIGHT.
3. Default annealing & extension temperature of multiplex PCR was changed to 64 °C.
4. Optional normalization step after PCR clean-up was omitted.

## PROTOCOL CITATION

Kentaro Itokawa, Tsuyoshi Sekizuka, Masanori Hashino, Rina Tanaka, Satsuki Eto, Risa Someno, Makoto Kuroda 2021. nCoV-2019 sequencing protocol for illumina. [protocols.io](#)  
<https://protocols.io/view/ncov-2019-sequencing-protocol-for-illumina-btjqnkmw>  
Version created by Kentaro Itokawa

## FORK NOTE

## FORK FROM

Forked from [nCoV-2019 sequencing protocol](#), Josh Quick

## KEYWORDS

SARS-CoV-2, Genome Sequencing, illumina

## LICENSE

— This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

CREATED

Mar 23, 2021

LAST MODIFIED

Apr 16, 2021

PROTOCOL INTEGER ID

48464

GUIDELINES

The important issue you have to consider first is how many samples you will multiplex in single run. This all depends on the capacity of your sequencer and the sample's viral loads represented by Ct-values in qPCR clinical test.

If a sample contains a relatively high copy number of virus genome (say, Ct < 25), the obtained reads usually distribute evenly across the genome. In such cases, only 10 Mb (330x) per sample is enough to cover the whole genome with sufficient coverage (>30x).

On the other hand, coverage bias increases as a sample's Ct-value increases. Hence, you will need more data to recover relatively weak regions. For samples containing only a low copy number of virus genome (Ct ~ 32), at least 100 Mb (3300x) is desirable to sequence the most part of the genome.

Also, if you are going to multiplex many samples in a low-throughput model (e.g. iSeq100), intensive normalization for each sample index will be needed.

MATERIALS TEXT

STEP MATERIALS

 [QIAseq FX DNA Library](#)

[Kit Qiagen Catalog #180475](#) Step 14

 [QIAseq FX DNA Library](#)


[Kit Qiagen Catalog #180475](#) Step 14

## cDNA preparation

1

This protocol uses 5/8 reagents per sample compared to the original LunaScript protocol.

Mix the following components in an 0.2mL 8-strip tube or 96 well PCR plate;

Component	Volume
LunaScript RT SuperMix	 1.25 µl
Template RNA (purified)	 5.0 µl
<b>Total</b>	 6.25 µl

A mastermix should be made up in the **mastermix cabinet** and aliquoted into PCR strip tubes. Tubes should be wiped down when entering and leaving the mastermix cabinet.

- 2 Gently mix by pipetting and pulse spin the tube to collect liquid at the bottom of the tube.

Incubate the reaction as follows:

⬆ 25 °C for ⌚ 00:02:00

⬆ 55 °C for ⌚ 00:20:00 (\*)

⬆ 95 °C for ⌚ 00:01:00

Hold at ⬆ 4 °C

\*Incubation for 20 min on 55 °C is used instead of 10 min of the kit protocol to make sure cDNA is synthesized efficiently. However, we have not confirmed the benefit of this modification, yet.

#### Primer pool preparation (not everytimes)

- 3 If concentration of your primer stocks are 50 µM, generate primer pool stocks by adding each primers equally or as described in tables below to Eppendorf tubes labelled either "Pool 1 (50 µM)" or "Pool 2 (50 µM)". These are your 50 µM stocks of each primer pool.

If concentration of the primer stocks are 100 µM, dilute this primer pool **1:1** in molecular grade water, to generate 50 µM primer stocks.

It is recommend that multiple aliquots of each primer pool are made to in case of degradation or contamination.

Primers should be diluted and pooled in the **mastermix cabinet** which should be cleaned with decontamination wipes and UV sterilised before and after use.

## 4 About N1/N2 primer set

[N1 \(NIID ver.1\)](#) primer set is a modified version of the [ARTIC Network's V1 primer set](#). The N1 primer set includes replacement of 12 primers which involved in dimer formation with other primers ([Itokawa et al., 2020](#)).

Since March 2021, we have added another small modification on the N1 primer set (ver. [N2](#)). The modification is addition of one primer "nCoV-2019\_72\_RIGHT\_C22A" to pool2 mixture. This primer contains an alternative nucleotide corresponding to G22017T mutation seen in [R.1 \(B.1.1.316.1\)](#) lineage to prevent amplicon 72 dropout.

Also, we recently confirmed that doubling the concentrations of nCoV-2019\_74\_LEFT and nCoV-2019\_74\_RIGHT primers improves the low coverage of amplicon 74.

Here is the recommended formulation of N2 primer set.

Primer	Volume
nCoV-2019_1_LEFT	5 µL
nCoV-2019_1_RIGHTv2	5 µL
nCoV-2019_3_LEFT	5 µL
...	...
nCoV-2019_97_LEFT	5 µL
nCoV-2019_97_RIGHT	5 µL
<b>Total</b>	<b>490 µL</b>













Example of primer mixing for Pool1 of the N2 primer set

Primer	Volume
nCoV-2019_2_LEFT	5 µL
nCoV-2019_2_RIGHT	5 µL
...	...
nCoV-2019_72_LEFT	5 µL
nCoV-2019_72_RIGHT	5 µL
<b>nCoV-2019_72_RIGHT_C22A</b>	<b>5 µL</b>
<b>nCoV-2019_74_LEFT</b>	<b>10 µL</b>
<b>nCoV-2019_74_RIGHT</b>	<b>10 µL</b>
...	...
nCoV-2019_98_LEFT	5 µL
nCoV-2019_98_RIGHT	5 µL
<b>Total</b>	<b>505 µL</b>


Example of primer mixing for Pool2 of the N2 primer set

## Multiplex PCR

- In the mastermix hood set up the multiplex PCR reactions as follows in 0.2 mL 8-strip PCR tubes or 96-well PCR plate:

Component	Pool 1	Pool 2
5X Q5 Reaction Buffer	 5 µl	 5 µl
dNTPs (10 mM each)	 0.5 µl	 0.5 µl
Q5 Hot Start DNA Polymerase	 0.25 µl	 0.25 µl
Primer Pool 1 or 2 (50 µM)	 0.72 µl	 0.72 µl
Nuclease-free water	 16.53 µl	 16.53 µl
<b>Total</b>	 23 µl	 23 µl



A PCR mastermix for each pool should be made up in the **mastermix cabinet** and aliquoted into PCR strip tubes. Tubes should be wiped down when entering and leaving the mastermix cabinet.

- In the **extraction and sample addition cabinet** add  2 µl diluted cDNA to each tube and mix well by pipetting.

The **extraction and sample addition cabinet** should be cleaned with decontamination wipes and UV sterilised before and after use.

- Pulse centrifuge the tubes to collect the contents at the bottom of the tube.

- Set-up the following program on the thermal cycler:

Stage	Temperature	Time
Heat Activation	 98 °C	 00:00:30

30 cycles\* of ...


Denaturation  $\uparrow$  98 °C 00:00:15

Annealing  $\uparrow$  64 °C \*\* 00:05:00

Hold  $\uparrow$  4 °C Indefinite

\*) We recommend 30 PCR cycles regardless of Ct-value because comparing PCR product yield and Ct-value gives a simple safeguard against sample swapping. However, one could increase this number to 35 to achieve more even PCR product yields.

\*\*) Previously, 65 °C was used for the annealing/extension step in multiplex PCR. However, inspecting many results in different thermocyclers and institutes, we now consider setting this temperature to  $\uparrow$  64 °C rather than 65 °C consistently yields better result with N1/N2 primer set. Nevertheless, fine-tuning of this temperature may still recommended to obtain the best result. Generally, amplicon 64 tends to show low coverage if annealing/extension temperature is higher than optimal.

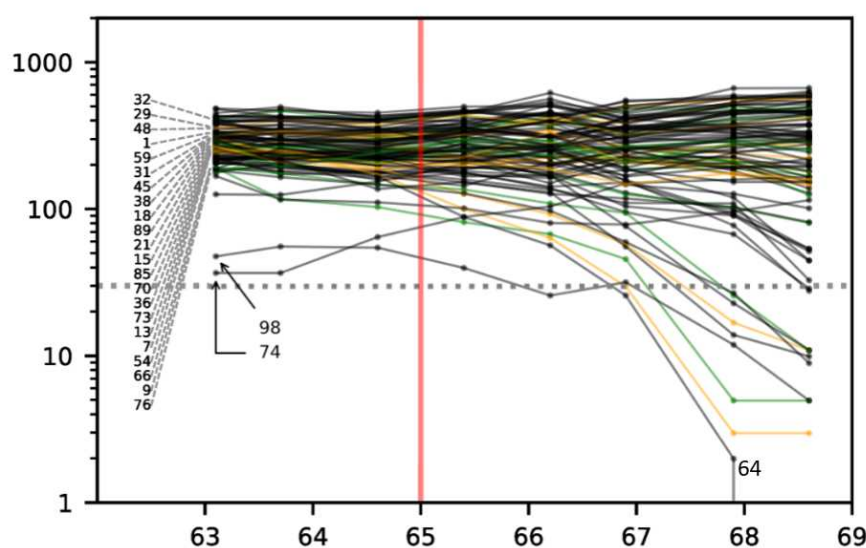

Temperature and coverage of amplicons with N1 primer set and the protocol version 1 in TaKaRa Thermal Cycler Dice Touch. Modified from [Ito et al., 2020](https://doi.org/10.1101/2020.03.10.289111) (CC BY 4.0).

#### PCR clean-up

- 9 Combine each 10  $\mu$ l contents of "Pool 1" and "Pool 2" PCR reactions for each biological sample into a new 8-strip PCR tubes or 96-well PCR plate.

Clean-up the amplicons using the AmpureXP using 1x volume.

10

Elute in **20 µl** of low-TE buffer (10 mM tris-HCl pH8.0, 0.1 mM EDTA).

Transfer eluted DNA to a new 8-strip PCR tubes or 96-well PCR plate.

Amplicon clean-up should be performed in the **post-PCR cabinet** which should be cleaned with decontamination wipes and UV sterilised before and after use.


#### Quantification and normalisation (optional)

11 Quantify the amplicon pools using fluorescent based method such as Qubit dsDNA HS kit (Invitrogen), Qaunt-it High-Sensitivity dsDNA Assay Kit (Invitrogen), or QuantiFluor® dsDNA System (Promega).

12 **!!! We do not recommend normalizing library concentrations at this step but after adapter ligation. It is more immune to sample swapping. !!!**

\* In our experience, samples with DNA concentration less than **2 ng / µL** in 30 cycles at this point do not generate meaningful results.

You can expect the quantity of DNA input will be directly proportional to the amount of data you will obtain if no further adjustment takes place.

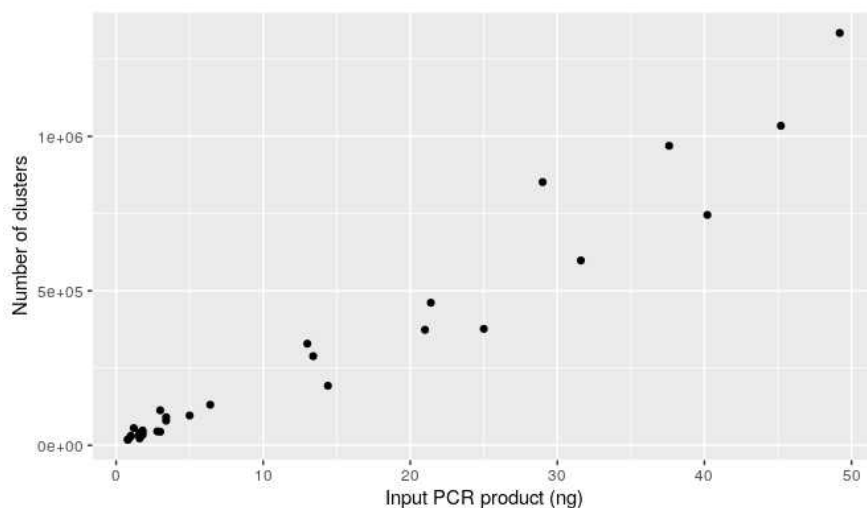

An example of input DNA vs number of read-pairs in single NGS run

## 13

This protocol uses **1/4** reagents per sample compared to the original QIAseq FX DNA library kit (Qiagen: 180473 for 24-plex, 180475 and 180477 for 96-plex, and 180479, 180480, 180481 and 180482 for 384-plex).

Transfer **2 µl** of pooled & purified PCR products to each well of a new 8-strip PCR tubes or 96-well PCR plate. Then, place the tubes or plate at 96 well aluminum block **On ice**.

Set the thermal cycler with a program below and start. Keep the heat-lid at **80 °C**.

**4 °C pose**

**32 °C** **00:06:00**

**65 °C** **00:30:00**

Prepare a master mix per one sample below.

Component	Volume / sample
FX Buffer, 10x	<b>1.25 µl</b>
H <sub>2</sub> O	<b>6.75 µl</b>
FX Enzyme Mix	<b>2.5 µl</b>
<b>Total</b>	<b>10.5 µl</b>

Add **10.5 µl** of the above master mix to each well of the 8-strip PCR tubes or 96-well PCR plate. Mix well by pipetting.

Make sure the plate is always **On ice** during this procedure to avoid fragmentation reaction proceeds at this step.

Place the library plate into the thermal cycler posing at 4°C and immediately skip to the next step ( **32 °C** ).

## 14 Remove the tube or plate from the thermal cycler after finishing the thermal program.

Take the adapter plate out of the kit box and thaw the content.

[QIAseq FX DNA Library](#)

[Kit Qiagen Catalog #180475](#)

Keep both end-prepped DNA mixture and thawed adapter solution **On ice** during this step.

Add **1 µl** adapter solution to each end-prepped DNA mixture.

Prepare a master mix per sample below.

Component	Volume / sample
DNA Ligase Buffer, 5x	<b>5 µl</b>
DNA Ligase	<b>2.5 µl</b>
H <sub>2</sub> O	<b>4 µl</b>
<b>Total</b>	<b>11.5 µl</b>

Add **11.5 µl** of above master mix to each end-prepped DNA mixture mixed with adapter **On ice**.

Set a thermal cycler with the following program with heat lid at **80 °C**.

- 20 °C** **00:15:00**
- 65 °C** **00:20:00** (ligase inactivation)

Start the thermal program, and place the tubes or plate immediately.

Seal and store remaining adapter plate in -20 °C. Extreme care should be taken to avoid cross-contamination between different indexed adapters.

#### Library pooling & purification

**15** Prepare a 1.5 or 2.0 ml low-binding tube.

Take **5 µl** \* ligated mixture from each well and pool them into the 1.5 or 2.0 ml low-binding tube.

\*You can take different volumes by sample to adjust the DNA quantities for even representation of each index See the **(Guidelines & Warnings)**. Use the DNA concentrations measured after the PCR clean-up.

Briefly measure the volume of pooled mixture in the 1.5 ml or 2.0 ml low-binding tube.

Purify by Ampure XP using **x0.8** volume.

Finally, elute DNA in **50 µl** low-TE or Elution Buffer.

Transfer the eluted DNA to a new 1.5 or 2.0 ml low-binding tube.

Purify again by Ampure XP using **x1.2** volume.

Finally, elute DNA in **25 µl** low-TE or Elution Buffer.

Transfer the eluted DNA to a new 1.5 or 2.0 ml low-binding tube.



Now, the library is ready for sequencing after quantification.

You may use any your routine methods for library quantification.

We recommend qPCR based methods in terms of accuracy and sensitivity.

Here is our homemade protocol.



Illumina TruSeq Library quantification with qPCR probe method


by Kentaro Itokawa,

National Institute of Infectious Diseases, Japan

PREVIEW

RUN

If you do not have a good quantification technique, here is a relationship between DNA mass concentration measured by Qubit dsDNA HS kit (Invitrogen) versus molar concentration measured by qPCR for nine libraries we obtained with this protocol.

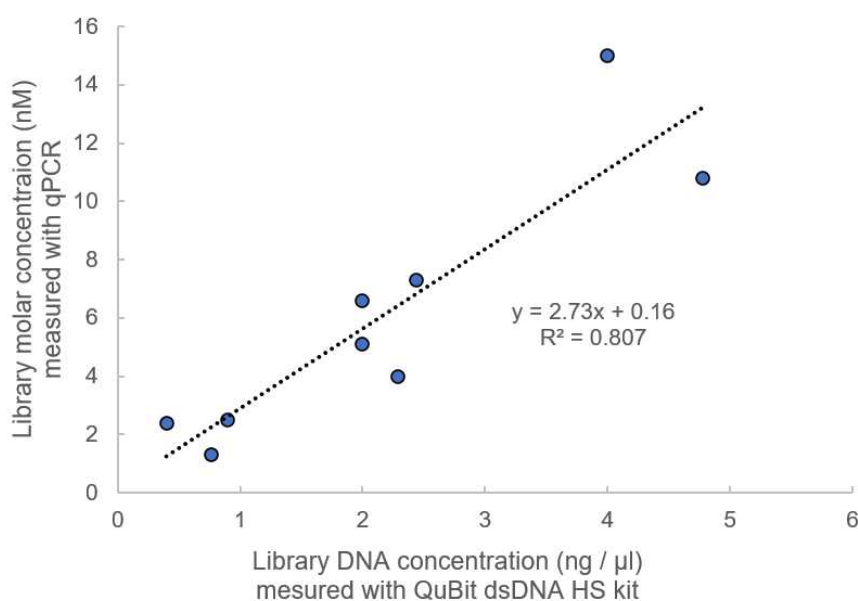

Library DNA mass concentration measured by Qubit versus molar concentration measured by qPCR