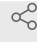


Version 6 ▼

Aug 12, 2022

An X-HTDC method for estimating particulate phosphorus from microalgae V.6

Yingyu Hu¹, Zoe V. Finkel¹¹Dalhousie University1 *Works for me* Sharedx.doi.org/10.17504/protocols.io.kqdg35dq7v25/v6

Marine Microbial Macroecology Lab
Tech. support email: ruby.hu@dal.ca



Yingyu Hu
Dalhousie University

ABSTRACT

Total particulate phosphorus (TPP) is often determined using the High Temperature Dry Combustion (HTDC) method followed by hydrolysis of the ash and then molybdenum colorimetry. Here we show that a higher than traditionally-used combustion temperature, 800 °C vs. 450 - 550 °C, improves phosphorus recovery from several organic phosphorus compounds, marine phytoplankton cultures and particulate samples from the field. The ashing auxiliary MgSO_4 further improves P recovery by improving decomposition, reducing volatilization during combustion, and increasing the efficiency of hydrolysis. A 0.2 M HCl hydrolysis, under 90 °C for 30 minutes yields a higher P recovery compared with hydrolysis at room temperature or 60 °C. In aggregate these improvements to the method double the P recovery from phospholipids to 97%. TPP recovery from laboratory phytoplankton cultures and field samples increased an average of 11%, primarily due to the improvements in P recovery from phospholipids, polyphosphates, and nucleic acids. We refer to this new method as the eXtra high temperature dry combustion ash/hydrolysis method (X-HTDC) and recommend its application for measuring particulate phosphorus from organic compounds in aquatic systems.

The working range of this assay is 1.22 to 500 μM orthophosphate. Minimum sampling biomass is 0.19 μg P/filter.

In order to assess the intracellular phosphorus in microalgae, we recommend an oxalate reagent (Tovar-Sanchez 2003) to wash the microalgae collected on the filter to remove surface adsorbed phosphorus.

P.S. Chen, T.Y. Toribara and Huber Warner. Microdetermination of Phosphorus. Anal. Chem..

<https://doi.org/10.1021/ac60119a033>

Antonio Tovar-Sanchez, Sergio A Sañudo-Wilhelmy, Manuel Garcia-Vargas, Richard S Weaver, Linda C Popels, David A Hutchins. A trace metal clean reagent to remove surface-bound iron from marine phytoplankton. Marine Chemistry.

[https://doi.org/10.1016/S0304-4203\(03\)00054-9](https://doi.org/10.1016/S0304-4203(03)00054-9)

DOI

dx.doi.org/10.17504/protocols.io.kqdg35dq7v25/v6

PROTOCOL CITATION

Yingyu Hu, Zoe V. Finkel 2022. An X-HTDC method for estimating particulate phosphorus from microalgae. **protocols.io**
<https://protocols.io/view/an-x-htdc-method-for-estimating-particulate-phosph-ce8mthu6>
Version created by Yingyu Hu



KEYWORDS

particulate phosphorus, intracellular phosphorus, phosphomolybdenum-ascorbic reduction, orthophosphate, oxalate reagent, adsorbed phosphorus, X-HTDC, High temperature dry combustion

LICENSE

_____ This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

CREATED

Aug 12, 2022

LAST MODIFIED

Aug 12, 2022

PROTOCOL INTEGER ID

68589

SAFETY WARNINGS

Polycarbonate filter can release toxic gas and smoke during combustion. An exhaust system is required for muffle furnace while using the X-HTDC method.

BEFORE STARTING

We have found that crucibles may lose their temperature resistance after acid-washing or long soaks in alkaline detergent. Crucibles tended to shatter in the oven during the initial increase in temperature from room temperature to 500 °C, even when the ramp rate was carefully controlled at 150 °C/h. We recommend not soaking crucibles in acid but instead we suggest the crucibles be filled with 0.2 M HCl and then incubated at 90 °C for 30 minutes as the acid-washing step. It is necessary to inspect the temperature resistance of newly acquired crucibles by combusting them at 500 °C for 6 h (ramp rate: 150 °C/h) after acid-washing. We found that crucibles that pass this inspection do not usually shatter when heated to 800 °C.

Sampling

1 Sampling microalgae for total particulate phosphorus (i.e. intracellular phosphorus and adsorbed phosphorus)

- 1.1 Filter microalgae in liquid media onto polycarbonate filters, using gentle vacuum pressure (130 mmHg).

Filter forceps
blunt end, stainless steel
Millipore XX6200006P

- 1.2 Rinse filter tunnel with filtered artificial seawater (phosphate free) to avoid sample loss.

- 1.3 Place sample filters in 2 mL Cryogenic Vials.

Cryogenic Vials with Closures
Polypropylene, 2 mL
Corning® 66021-974

- 1.4 Filter blank media (without cells) through polycarbonate filter as blank.



- 1.5 Flash freeze filters and store at **-20 °C**.

2 Sampling microalgae for intracellular particulate phosphorus

- 2.1 Filter microalgae in liquid media onto polycarbonate filters, using gentle

vacuum pressure (5 inches Hg).

Filter forceps
blunt end, stainless steel
Millipore XX6200006P

2.2 Add  **5 mL** oxalate reagent onto the filter, and let oxalate reagent sit in the ^{5m} filter funnel for  **00:05:00**




Preparation of oxalate reagent
by Yingyu Hu,
Dalhousie University

PREVIEW

RUN

2.2.1 Add  **50 mL** MilliQ water in a 250 mL beaker.

2.2.2 Weigh  **40 g** NaOH and slowly pour into the beaker.

2.2.3 Use squeeze bottle to rinse the weighing boat and transfer rinse water into the same beaker.

2.2.4 Use glass rod to gently stir and fully dissolve NaOH.

The solution is very hot and corrosive. It can cause skin burns and eye damage.

2.2.5 Carefully transfer NaOH solution into 100 mL volumetric flask by using glass rod.

2.2.6 Rinse beaker with small amount of MilliQ water three times, transfer rinse water into the flask.

2.2.7 Mix the solution by gently shaking the capped volumetric flask and top to 100 mL with MilliQ water.

2.2.8 Transfer the prepared reagent into a 250 mL PP bottle.

2.2.9 Label the bottle with SDS pictogram.



2.10 In a 1000 mL beaker with stir bar, add **600 mL** MilliQ water.

2.11 Add **18.6 g** EDTA, **14.7 g** sodium citrate, **0.74 g** KCl and **5 g** NaCl into the beaker, stir until all ingredients are dissolved. **pH 5.7**

2.12 **10 Molarity (M)** NaOH is added dropwise to bring pH in between 6 to 7 by using a transfer pipet

2.13 Add **12.6 g** oxalic acid to the solution, stir the mixture while heating.

2.14 After oxalic acid is completely dissolved, stop heating and let it cool to room temperature. A water bath filled with tap water can be used to speed up cooling. **pH 3.3**

2.15 Add **10 Molarity (M)** NaOH dropwise to bring pH to **pH 8**

2.16 Top to 1 L in volumetric flask with MilliQ water.

2.17 Filter oxalate reagent by rapid flow to a 1 L PP bottle.

Sterile Disposable Filter Units with PES
Membrane

Thermo Scientific™ Nalgene™ Rapid- Flow™	5964 520
---	-------------

2.18 Label the bottle and keep it at **Room temperature**.

2.3 Drain and then rinse filter funnel with filtered artificial seawater (phosphate free).

2.4 Place sample filters in 2 mL Cryogenic Vials.

Cryogenic Vials with Closures
Polypropylene, 2 mL

Corning®	66021-974
----------	-----------

2.5 Filter blank media (without cells) through polycarbonate filter as blank.

2.6 Flash freeze filters and store at -20°C .

X-HTDC-ing

3 Mark number at the bottom of each crucible with pencil, log the following information:

- (1) The number of crucible
- (2) The code of sample in the crucible

Porcelain crucibles

40 mL

VWR 89037-996

Crucible cover

VWR 71000-146

4 Transfer sample to crucible with clean filter forceps and lay filter at the bottom.

5 0.17 M MgSO_4 reagent:

Dissolve 1.023 g MgSO_4 in 50 mL MilliQ water

[Magnesium sulfate anhydrous Fisher](#)

Scientific Catalog #M65500

6 Add $200\text{ }\mu\text{L}$ 0.17 M MgSO_4 directly onto each sample and blank filter.

Single-use pipet tip to avoid cross-contamination.

- 7 Cover the crucibles and place in the oven at 🔥 **90 °C** until samples are completely dry.

Forced air oven

VWR 89511-410

Remove samples out of the oven as soon as they are dried. If muffle furnace is not available, keep crucibles in the vacuum desiccator.

- 8 Combust dried samples at 🔥 **800 °C** for ⌚ **09:00:00**

9h

Muffle furnace

F30428C

Thermo 10-505-13

Map the location of crucibles in the oven, in case pencil mark disappears under 800°C.

Ramp rate should be controlled at < 🔥 **200 °C** /hour or follow the instruction provided by manufacture, otherwise the crucibles might shatter.

SP.RAT: 150/PAMPU: hour

Or

SP.RAT: 2.5/PAMPU: minute

- 9 Allow samples to gradually cool down in the muffle furnace.

- 10 Pencil mark on crucibles should be still visible, however, it can be easily removed by water. Therefore, when removing samples out of the furnace, label the lid and crucible with sharpie immediately.

Digesting

- 11 **[M]0.2 M** HCl reagent:

In a reagent bottle, dissolve one part of **[M]12 N** HCl in 59 parts of MilliQ water

[M]12 N Hydrochloric acid Contributed by users

Volume of HCl_0.2M_mL = (5_mL) X (#Sample + #Blank)

- 12 Preheat oven to **90 °C**

- 13 Add 5 mL 0.2 M HCl to each crucible.

- 14 Gently swirl the crucible.

- 15 Cover the crucibles and place crucibles in the muffin tin pan for easier-handling.

- 16 Incubate in the oven for **00:30:00**

30m

- 17 Cool samples down to **Room temperature**

- 18 Gently swirl the crucible and then transfer 500 uL solution to 2 mL microtube. Duplicate each sample and blank.

Maxymum Recovery® Snaplock
Microcentrifuge Tube
2.0 mL, Polypropylene, Clear, Nonsterile,
Axygen® MCT-200-L-C

Preparing standard working solutions

2h



19 Standard working solutions and reagents can be prepared during sample digestion.

20 KH_2PO_4 primary standard stock solution (≈ 1 mM)

 Potassium dihydrogen orthophosphate **ACP**

Chemicals Catalog #P-4550

20.1 Transfer about 1 g KH_2PO_4 into a beaker, cover the beaker with foil

20.2 Place the beaker into an oven, dry KH_2PO_4 at  **110 °C** for at least
 **02:00:00**

2h

20.3 Move KH_2PO_4 into a vacuum desiccator, allow KH_2PO_4 to cool to room temperature

20.4 Dissolve around 0.136 g dried KH_2PO_4 in 1 L milliQ water.

- Use 1 L volumetric flask
- Take notes of the actual weight of KH_2PO_4 for final concentration of standard stock solution

20.5 Transfer standard stock solution into a 1 L bottle and store in the fridge.

This stock solution lasts quite a long time, unless there is evidence for growth of algae or other extraneous biotic material.

21 Standard working solution

KH ₂ PO ₄	Primary (ul)	MilliQ (ul)
S1	0	1000
S2	5	995
S3	10	990
S4	20	980
S5	50	950
S6	100	900
S7	150	850
S8	200	800

22 Transfer 500 uL of each standard working solution to 2 mL microtube.

Preparing working reagents

2h

23

All reagents are freshly prepared before colorimetric measurement.

24 [\[M\]6 N](#) (3 M) sulfuric acid reagent:

Carefully add 1 part [\[M\]18 M](#) concentrated sulfuric acid into 5 part MilliQ water

[\[M\]18M sulfuric acid](#) **Contributed by users**

25 [\[M\]2.5 %](#) ammonium molybdate reagent:

Weigh [\[M\]0.25 g](#) ammonium molybdate in a Falcon tube and top to [\[M\]10 g](#) with MilliQ water. Cap and shake until totally dissolved.

[\[M\]Ammonium molybdate](#) **Sigma**

Aldrich Catalog #09878-100G

26 [M]10 % ascorbic acid reagent:

Weigh 1 g ascorbic acid in a Falcon tube and top to 10 g with MilliQ water;
Cap and shake until all dissolved.

⊗ Ascorbic acid Sigma

Aldrich Catalog #A5960-100G

Wrap the tube with foil if the reagent is not used right after prepared.

27 Calculate the volume of molybdate-ascorbic reagent:

Total volume of reagent_mL = (0.5 mL) X (#standard working solution + #samples + #blanks)

28 Mix the reagents into Falcon tube:

A	B
Reagent	Part(s) as in volume
MilliQ	2
6N sulphuric acid	1
2.5% ammonium molybdate	1
10% ascorbic acid	1

Colorimetric measurement

2h

29 Preheat incubator/shaker to 37 °C

SHAKING INCUBATOR

71L

Corning® LSE™ 6753

30 Add 500 µL reagent to each standard, sample and blank, starting from blanks, including blank for standards and blank for samples.

Finntip Stepper Tips

5 mL

Thermo Scientific 9404200

Before dispensing the reagent, wipe or dab the liquid drop on the outside of the tip, avoid wiping the open tip.

31 Vortex each tube.

32 Incubate at **37 °C** for **03:00:00** while shaking at 200 rpm

3h

33 Load microplate with 250 uL reactant from each tube, duplicate.

	<u>1</u>	<u>2</u>	<u>3</u>	<u>4</u>	<u>5</u>	<u>6</u>	<u>7</u>	<u>8</u>	<u>9</u>	<u>10</u>	<u>11</u>	<u>12</u>
<u>A</u>	S1	S1	Samples and sample blanks: 40 with duplicate									
<u>B</u>	S2	S2										
<u>C</u>	S3	S3										
<u>D</u>	S4	S4										
<u>E</u>	S5	S5										
<u>F</u>	S6	S6										
<u>G</u>	S7	S7										
<u>H</u>	S8	S8										

Example of loading the microplate

96-Well Microplates, Polystyrene, Clear,

Greiner Bio-One 655101

34 Read plate in microplate reader

A	B
Shake duration	00:00:05
Shaking type	Continuous
Shaking force	High
Shaking speed [rpm]	600
Wavelength [nm]	820
Use transmittance	No
Pathlength correction	No
Measurement Time [ms]	100

Varioskan LUX Multimode Microplate Reader

Thermo Fisher VL0L00D0

Calculating

- 35 Subtract the average absorbance at 820 nm of the blank standard replicates from the absorbance at 820 nm of all other standard working solutions.
- 36 Subtract the average absorbance at 820 nm of the blank sample (i.e. blank filter) replicates from the absorbance at 820 nm of all other individual samples.
- 37 Prepare a standard curve by plotting the average blank-corrected 820 nm absorbance for each standard working solution versus its concentration in μM .
- 38 Use the standard curve to determine the orthophosphate concentration of each unknown sample by using its blank-corrected 820 nm absorbance.
- 39 $(\text{P per sample})_{\mu\text{g}} = (\text{orthophosphate})_{\mu\text{M}} \times (\text{V_HCl})_{\text{mL}} \times (0.001) \times (30.97)$

