

AUG 08, 2023

OPEN ACCESS



**DOI:**  
[dx.doi.org/10.17504/protocols.io.q26g7py6kgwz/v1](https://dx.doi.org/10.17504/protocols.io.q26g7py6kgwz/v1)

**External link:**  
<https://www.jax.org/jax-mice-and-services/ipsc>,  
<https://www.jax.org/jax-mice-and-services/ipsc/cells-collection>

**Protocol Citation:** michael. ward, Erika Lara Flores, Mark Cookson 2023. iNDI Maintenance protocol of iPSCs Version 1 .  
**protocols.io**  
<https://dx.doi.org/10.17504/protocols.io.q26g7py6kgwz/v1>

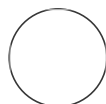
**License:** This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

## iNDI Maintenance protocol of iPSCs Version 1

Erika Lara  
michael.ward<sup>1</sup>, Flores<sup>1</sup>, Mark  
Cookson<sup>1</sup>

<sup>1</sup>National Institutes of Health, NIH Center for Alzheimer's and Related Dementias (CARD)

NIH Center for Alzheimer's and Related Dementias



Erika Lara Flores

### ABSTRACT

IPSC maintenance protocol

- Matrigel procedure for coating plates
- Vitronectin procedure for coating plates
- Thawing iPSC
- Splitting and Passaging iPSC
- Freezing iPSC

**Protocol status:** Working  
We use this protocol and it's working

**Created:** Aug 08, 2023

**Last Modified:** Aug 08, 2023

**PROTOCOL integer ID:**  
86179



**Keywords:** iNDI, iPSC, neurodegeneration, genetic engineering, Jackson Laboratory, freezing, maintenance

## MATERIALS

- ✕ Matrigel hESC-qualified (Corning Cat# 354277) Corning Catalog #354277
- ✕ Vitronectin (VTN-N) Recombinant Human Protein Truncated Thermo Scientific Catalog # A31804
- ✕ Essential 8™ Medium Gibco - Thermo Fischer Catalog #A1517001
- ✕ • Chroman I MedChemExpress Catalog #HY-15392
- ✕ • Phosphate Buffered Saline (1X) without Calcium or Magnesium Lonza Catalog 17-516F
- ✕ 0.5 mM EDTA Gibco - Thermo Fischer Catalog # AM9260G
- ✕ Accutase Gibco - Thermo Fischer Catalog # A1110501
- ✕ Dimethylsulfoxide (DMSO) Corning Catalog #25-950-CQC
- ✕ Knockout™ Serum Replacement Gibco - Thermo Fischer Catalog #10828028
- ✕ KnockOut™ DMEM/F-12 Gibco - Thermo Fischer Catalog #12660012

## Matrigel Coating

1h 30m




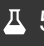
- 1 **Aliquot concentrated Matrigel:**
  - Gradually thaw a 5ml bottle of Matrigel on ice in a Styrofoam container
  - Pre-chill labeled Eppendorf tubes by placing in a cool rack on ice.
  - Before pipetting concentrated Matrigel into pre-chilled tubes, chill a 1 ml pipet tip by pipetting ice-cold [KnockOut™ DMEM/F-12](#) up and down several times, then immediately use the tip to aliquot Matrigel.
  - Prepare aliquots of  500 µL of concentrated Matrigel and freeze down at  -80 °C

### Note

*Matrigel concentrated can polymerize rapidly at room temperature, so it's very important when aliquoting or preparing coating solution do on ice.*




## 2 Coating plates with Matrigel solution:

1h 30m


- Reconstitute a  500  $\mu\text{L}$  aliquot of Matrigel in  50 mL of cold KnockOut™ DMEM/F-12, pipet  500  $\mu\text{L}$  cold media into the aliquot of Matrigel tube and pipet up and down several times, then transfer what has thawed to the tube containing  50 mL of cold KnockOut™ DMEM/F-12, repeat until the frozen aliquot of Matrigel has been completely transferred to the 50 ml of cold media. Mix inverting several times.

### Note

*Matrigel solution can be store at 4°C until it finishes.*

- Add half of the normal culture volume of Matrigel solution to the culture surface (i.e., 1 ml per well of a 6-well plate).
- Gently rock plate to spread the Matrigel solution evenly across the plate.
- Place in  37 °C incubator for  00:30:00 to  01:00:00 before use it.

### Note

*For KOLF2.1 iPSC you can get better results from  Overnight coating*

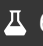

- Aspirate Matrigel and add culture medium (E8).

There is another option for coating plates besides Matrigel. **Vitronectin** is a recombinant human protein that provides a defined surface for feeder-free culture of iPSC. When used with E8 medium, vitronectin has been proven to maintain pluripotency and normal growth characteristics in multiple iPSC lines.





## Vitronectin Coating

30m

## 3 Aliquot concentrated Vitronectin 0.5 mg/ ml 0.5 mg/mL


- Thaw a vial of Vitronectin on ice in a Styrofoam container.
- Prepare aliquots of  60  $\mu\text{L}$  of concentrated Vitronectin and freeze down at  -80 °C

#### 4 Coating plates with Vitronectin solution [M] 5 µg/mL

- Reconstitute a  60 µL aliquot of Vitronectin into a 15 ml conical tube containing  6 mL of PBS, gently resuspend by pipetting up and down, do not vortex.
- Add half of the normal culture volume of Vitronectin solution to the culture surface (i.e. 1 ml per well of a 6-well plate).
- Gently rock plate to spread the Vitronectin solution evenly across the plate.
- Place in  37 °C incubator for  Overnight before use it.
- Aspirate Vitronectin and add culture medium (E8).


### Thawing iPSC



5m

- 5 Remove iPSC stock cryovial from liquid nitrogen and thaw in  37 °C bead bath. Thaw quickly by gently swirling until a small piece of frozen material remains. Spray the vial with 70% ethanol before transferring to a biological safety cabinet.

#### Note

*Since DMSO is toxic to cells at room temperature, perform the following steps in a time-efficient manner to obtain optimal cell viability.*

- 6 Gently add the thawed cell suspension dropwise to a conical tube containing  10 mL culture medium or PBS, rinse cryovial with 1ml of medium and add the rinse to the tube, gently mix cells by swirling.

- 7 Centrifuge tube  00:05:00 at 200 - 300 x g at  Room temperature

5m

- 8 Aspirate the supernatant and gently resuspend cells in culture medium (E8) supplemented with

**1mM** 50 nanomolar (nM) Chroman I and transfer to Matrigel or Vitronectin-coated plates.

#### Note

*After thawing the cells is recommended to maintain high cell density to maximize cell viability. The table below suggests some number of cells and vessel to use.*

Number of cells in cryovial	Vessel
3 x 10 <sup>6</sup> cells	100 mm dish
0.5 x 10 <sup>6</sup> to 1 x 10 <sup>6</sup> cells	a well of 6 well plate
0.2 x 10 <sup>6</sup> to 0.4 x 10 <sup>6</sup>	a well of 12 well plate
0.08 x 10 <sup>6</sup> to 0.15 x 10 <sup>6</sup>	a well of 24 well plate

9 Gently rock plate to evenly distribute cells.

10 Return plate to **37 °C** incubator.

11 Next day replace the media with fresh E8 medium (**2 mL** /well of a 6-well plate).

12 When the well is 80% confluent pass to expand stock.

## 13 EDTA

**13.1** Aspirate culture medium and wash with PBS 1X.

**13.2** Aspirate PBS and add half of culture volume (1ml/well of a 6-well plate) of EDTA [M] 0.5 millimolar (mM) in PBS.

### Note

*I have been using*



Cell Dissociation Buffer, enzyme-free, PBS Thermo Fisher Catalog #13151014

as

well and it works the same as EDTA.



**13.3** Incubate for 00:03:00 at 37 °C or 00:08:00 at Room temperature

11m

### Note

*The time can vary by cell line and density (the optimal density is 70-90%)*

**13.4** Aspirate EDTA solution, the cells colonies should remain attached so be careful not to disturb them.

**13.5** Add  1 mL of culture medium supplemented with  50 nanomolar (nM) Chroman I to cells to dissociate by pipetting two or three times.

**13.6** Typically splitting ratios for 6 well plates are between 1:6 and 1:12, so, add the desire volume of culture medium to the cells and discard any excess of cells or re-plate into a new Matrigel or Vitronectin-coated well.

## 14 Accutase

**14.1** Aspirate culture medium and wash with PBS 1X.


**14.2** Aspirate PBS and add half of culture volume of Accutase.

**14.3** Transfer to  37 °C incubator for  00:08:00

8m

### Note



*The time can vary by cell line and density (the optimal density is 70-90%) and the goal to use accutase is singularize as single cells.*

**14.4** Meanwhile aspirate Matrigel/Vitronectin from plates and add culture medium E8 supplemented with  50 nanomolar (nM) Chroman I.

**14.5** When Incubation is ready, tilt the plate and pipet the accutase solution two to three times up and down the culture surface to break the colonies.

**14.6** Quench the Accutase adding half of the culture volume of PBS.

**14.7** Transfer to a new conical tube and rinse with more PBS the culture surface, combine with the cell solution in the tube.


**14.8** Centrifuge  00:05:00 at 200 - 300 x g at  Room temperature

5m

**14.9** Aspirate supernatant.

#### Note


*Remains of Accutase could interfere with the cell viability after re-plating, make sure you aspirate everything.*

**14.10** Resuspend the cell pellet in culture medium E8 supplemented with  50 nanomolar (nM) Chroman I.

**14.11** Count cells and plate cells the desire amount into a Matrigel-coated plates.

**14.12** Gently rock plate to evenly distribute cells.



**14.13** Return plate to  37 °C incubator.

**14.14** The next day replace the medium with fresh E8 (2ml/well of a 6-well plate)

**14.15** When the well is 80% confluent pass to achieve an assay or to expand stock.

**Note**




*KOLF2.1 iPSC does not behave well when is more than 80-90%, most of the cells will die and it will be very difficult to get them back.*

## Freezing iPSC

**15** Prepare freezing medium as combining Knockout™ Serum Replacement with 10% DMSO.

**16** Leftovers of EDTA or Accutase dissociation procedure could be cryopreserved by centrifuging as in step 14.8 and resuspending the cell pellet with freezing medium.

**17** When freezing cells from a well/plate, prepare the cells as for EDTA split.

- 18 Aspirate EDTA and gently dissociate cells with freezing medium.
- 19 Transfer  1 mL of cell suspension to a 1 mL cryovial and freeze in a CoolCell freezing container at  -80 °C  
 Overnight
- 20 Next day transfer the cryovials to liquid nitrogen for long term storage.