

Jun 20, 2024

## Metabarcoding PCR Protocol

DOI

**[dx.doi.org/10.17504/protocols.io.6qpvr8842lmk/v1](https://dx.doi.org/10.17504/protocols.io.6qpvr8842lmk/v1)**

Mirayana Marcelino Barros<sup>1</sup>, Tiago Pereira<sup>2</sup>, Alejandro De Santiago Perez<sup>1</sup>, Hunter Powell<sup>3</sup>, Jenna Brown<sup>4</sup>

<sup>1</sup>University of Georgia; <sup>2</sup>Department of Marine Sciences, University of Georgia; <sup>3</sup>UGA; <sup>4</sup>University of Georgia- Holly Bik



Jenna Brown

University of Georgia- Holly Bik

OPEN  ACCESS



DOI: **[dx.doi.org/10.17504/protocols.io.6qpvr8842lmk/v1](https://dx.doi.org/10.17504/protocols.io.6qpvr8842lmk/v1)**

**Protocol Citation:** Mirayana Marcelino Barros, Tiago Pereira, Alejandro De Santiago Perez, Hunter Powell, Jenna Brown 2024. Metabarcoding PCR Protocol. **protocols.io** **<https://dx.doi.org/10.17504/protocols.io.6qpvr8842lmk/v1>**

**License:** This is an open access protocol distributed under the terms of the **[Creative Commons Attribution License](#)**, which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

**Protocol status:** Working

**We use this protocol and it's working**

**Created:** May 30, 2024

**Last Modified:** June 20, 2024

**Protocol Integer ID:** 100934


**Funders Acknowledgement:**  
Holly Bik

### Abstract

18S Metabarcoding PCR Protocol using Invitrogen Platinum Hot Start PCR 2X Master Mix (REF 13000014) buffer.


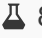


## Preparing the PCR tubes




- 1 Label PCR tubes and place them under a UV light for  00:10:00 .

10m

## Preparing the master mix

- 2 Add  10  $\mu\text{L}$  of **Invitrogen Platinum Hot Start PCR 2X Master Mix (REF 13000014) buffer** and  8  $\mu\text{L}$  of **molecular water** to a 1.5 mL tube PER NUMBER OF SAMPLES *i.e., for 4 samples, add 40  $\mu\text{L}$  of buffer and 32  $\mu\text{L}$  of molecular water to a 1.5 mL tube.*
- 3 Mix well using a vortex, then centrifuge.

## Metabarcoding PCR

- 4 Transfer  18  $\mu\text{L}$  of the **master mix** to each PCR tube.
  - 5 Transfer  5  $\mu\text{L}$  of the **DNA template** to each PCR tube.
  - 6 Transfer  2  $\mu\text{L}$  of the **Metabarcoding Primer** to the assigned PCR tube. Immediately place primers on a PCR Cooler Rack after use.
  - 7 Mix the samples well using a vortex, then briefly centrifuge them using a microcentrifuge or platefuge to ensure all contents are at the bottom.
  - 8 Place samples in a thermo cycler and set to a PCR protocol. Set the thermo cycler to use a heated lid and wait until the lid and block are close to the temperature. See 8.1 for thermo cycler specifications.  
*Tip: In the Bik Lab, select protocol "METAB-HS" on thermo cycler.*
- 8.1 METAB-HS Protocol:
- Lid: 105°C
  - Volume: 25  $\mu\text{L}$
  - 1. 94°C, 3:00
  - 2. 94°C, 0:45
  - 3. 50°C, 1:00



4. 72°C, 1:30

5. GOTO step 2, 34X

6. 72°C, 10:00

7. 4°C, ∞

9 Close the thermo cycler lid well.