




DEC 14, 2023

## POLE TEST PROTOCOL

Michael Lee<sup>1</sup>

<sup>1</sup>University of Minnesota

Team Lee

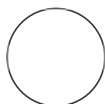

jbalster

### DISCLAIMER

OPEN ACCESS



**DOI:**  
[dx.doi.org/10.17504/protocols.io.n92ldzjkxv5b/v1](https://dx.doi.org/10.17504/protocols.io.n92ldzjkxv5b/v1)

**Protocol Citation:** Michael Lee 2023. POLE TEST PROTOCOL. **protocols.io** <https://dx.doi.org/10.17504/protocols.io.n92ldzjkxv5b/v1>

**License:** This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

**Protocol status:** Working  
 We use this protocol and it's working

**Created:** Jun 04, 2022

**Last Modified:** Dec 14, 2023

DISCLAIMER – FOR INFORMATIONAL PURPOSES ONLY; USE AT YOUR OWN RISK

The protocol content here is for informational purposes only and does not constitute legal, medical, clinical, or safety advice, or otherwise; content added to [protocols.io](#) is not peer reviewed and may not have undergone a formal approval of any kind. Information presented in this protocol should not substitute for independent professional judgment, advice, diagnosis, or treatment. Any action you take or refrain from taking using or relying upon the information presented here is strictly at your own risk. You agree that neither the Company nor any of the authors, contributors, administrators, or anyone else associated with [protocols.io](#), can be held responsible for your use of the information contained in or linked to this protocol or any of our Sites/Apps and Services.

### ABSTRACT

This protocol details pole test.

### ATTACHMENTS

[452-955.docx](#)

### GUIDELINES

**Detailed records should be kept on every experiment, noting date, timing, subject test order, procedural details followed and any relevant observations pertaining to the assay (e.g. specific behavior or procedural errors).**

**PROTOCOL integer ID:**  
63883



MATERIALS

**Keywords:** pole test,  
animals, ASAPCRN

Materials

- Pole
- Timer

**Funders**  
**Acknowledgement:**  
ASAP

Test		22m
1	Place animal on horizontal pole facing the top.	
2	Place pole upright with mouse facing upward.	
3	Record the time to turn (orient facing down), and time to reach the bottom of the pole.	
4	Maximum duration  00:02:00 .	2m
5	All animals tested for 3 trials with a minimum of  00:20:00 between trials.	20m