



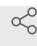
Sep 15, 2022

Low Biomass, high contamination Illumina DNA prep using DNeasy PowerSoil (Pro) Kit

Tania Kurbessoian¹, Jason E Stajich¹, Sonia L. Ghose²¹UC Riverside; ²UC Davis

Sonia L. Ghose: Originally adapted from.;

1 Works for me

 Sharedx.doi.org/10.17504/protocols.io.6qpvr4ey3gmk/v1[stajichlab](#)

Tania Kurbessoian

ABSTRACT

This is an addendum to the already optimized DNeasy PowerSoil Kit. This can also be applied to the PowerSoil Pro Kit.

Use this protocol for Low Biomass or high contamination samples.

DOI

dx.doi.org/10.17504/protocols.io.6qpvr4ey3gmk/v1

PROTOCOL CITATION

Tania Kurbessoian, Jason E Stajich, Sonia L. Ghose 2022. Low Biomass, high contamination Illumina DNA prep using DNeasy PowerSoil (Pro) Kit.

protocols.io<https://protocols.io/view/low-biomass-high-contamination-illumina-dna-prep-u-cgp6tvre>

LICENSE

————— This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

CREATED

Sep 15, 2022

LAST MODIFIED

Sep 15, 2022

PROTOCOL INTEGER ID

70110

GUIDELINES

Using some small tweaks to an already optimized protocol you can get your DNA from contaminated samples!

MATERIALS TEXT

DNeasy Powersoil Kit (Adjusted name to PowerLyzer Kit) -

<https://www.qiagen.com/us/products/discovery-and-translational-research/dna-rna-purification/dna-purification/microbial-dna/dneasy-powerlyzer-powersoil-kit/>

DNeasy Powersoil Pro Kit <https://www.qiagen.com/us/products/discovery-and-translational-research/dna-rna-purification/dna-purification/microbial-dna/dneasy-powersoil-pro-kit/>

100% ethanol

BEFORE STARTING

Make sure you have all the components of the kit ready to go.

1 Add sample into the PowerBead tube.

2 Add 60µL of Solution C1.

2.1 Make sure the sample has not precipitated. If so, heat sample to 60°C until precipitate is dissolved into solution.

3 Vortex to mix and incubate at 65°C for 10 minutes. 10m

4 Bead beat at "homogenize" setting: 90 seconds bead beating, 60 seconds rest, then another^{4m} 90 seconds (3 minutes total of actual beating).

- 5 Centrifuge at 10,000 x g for 30 s.
- 6 Remove the supernatant or upper aqueous layer if you have one, to the new tube.
- 7 Add 100 μ L Solution C2 and 100 μ L of Solution C3. Vortex to mix. Incubate at 4°C or on ice for 5^{5m} minutes.

7.1 *This is for low humic soils (100 μ L of each). If high humic soil, add 150 μ L C2 followed by 150 μ L C3.

- 8 Centrifuge to pellet (1 min at 10,000 x g.) and transfer the supernatant to a new tube. Avoid^{1m} pellet.
- 9 Ideally you will have 650 μ L of lysate, and can add 650 μ L of Solution C4 and 650 μ L of 100% ethanol.

9.1 Make sure you mix C4 solution well. If new kit add ethanol.

9.2 If you have 700 μ L, add 700 μ L of Solution C4 and 600 μ L of 100% ethanol. (Keep it as close to 1:1 as possible)

- 10 Load the lysate to the spin filter 650 μ L at a time and bind in three steps, alternating with^{3m} centrifugation (1 min at 10,000 x g). Discard flow-through each time.
- 11 If the membrane is not stained brown, wash with 650 μ L of 100% ethanol (centrifuge 1 min at^{1m} 10,000 x g and discard flow-through).

11.1 If the membrane is stained brown, per sample, prepare a mix of 300 μ L of

Solution C4 and 370µL of 100% ethanol. Wash the column with this mixture first. Follow this wash with the 100% ethanol wash.

- 12 Wash with 500µL of solution C5 (centrifuge 1 min at 10,000 x g and discard flow-through) ^{1m}
- 13 Dry the spin column for 2 minutes 10,000 x g. Transfer to a clean tube. ^{2m}
- 14 Elute in 60µL of buffer C6 (Heat this buffer to 60°C before eluting). Let the buffer sit on the membrane 5 ^{5m}
minutes before elution.
- 15 Then centrifuge for 1 min 10,000 x g into your storage tube. ^{1m}