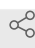




Jul 14, 2022

# Dual Luciferase Assay

vahedilaboratory<sup>1</sup><sup>1</sup>University of Pennsylvania*In Development* Share

This protocol is published without a DOI.

 vahedilaboratory

## ABSTRACT

Dual luciferase assay WIP

## PROTOCOL CITATION

vahedilaboratory 2022. Dual Luciferase Assay . **protocols.io**  
<https://protocols.io/view/dual-luciferase-assay-cdbrs2m6>



## LICENSE

————— This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

## CREATED

Jul 13, 2022

## LAST MODIFIED

Jul 14, 2022

## PROTOCOL INTEGER ID

66641

## MATERIALS TEXT

1. Dual-Luciferase® Reporter Assay System found [here](#).
2. FuGENE transfection reagent found [here](#).

## 1 Cell culture

- If you are reading this, I am assuming that you already have 100,000 cells plated in each well of a 24-well TC plate (or something of equal measure)
- Let them rest overnight because tomorrow you will transfect them

## 2 Transfection (Transient) of Plasmid

- Dual luciferase assays work by expressing and measuring two luciferase plasmids: Firefly Luciferase and Renilla luciferase, which helps normalize transfection efficiencies of the plasmids
- Assemble the following mixture for a 24 well plate:

2uL of FuGENE

1uL plasmid A (200ng)

1uL plasmid B (200ng)

50uL Optimem

- Let sit for max 15 minutes
- Add dropwise to the cells
- Leave for 24 hours

### 3 Cell Lysing and Luciferase Assay

- Remove media from each well and wash with PBS; remove the PBS
- Add 50uL trypsin to each well and coat every cell
- Leave in incubator for 3 minutes
- Pipette all cells off of plate and add cells to 1.5mL tubes
- Spin down cells at 1,500RPM 4C for 5 minutes to pellet cells
- Remove supernatant
- Re suspend cells in 50uL of PLB 1X Passive Lysis Buffer (PLB)
- Incubate lysates on ice for at least five minutes
- Vortex each tube for 10-15 seconds
- Incubate for 5-10 more minutes on ice
- Centrifuge the tubes at 12,000 rcf at 4C for two minutes
- Transfer the supernatant of each tube to a new labeled tube.
- Either store the lysates at -80C or continue to the luciferase assay

### 4 The Luciferase Assay Itself

- Set up a GloMax luminometer to measure for 10 seconds
- For each sample you want to measure, add 20 uL of LARII and 20 uL lysate to a microfuge tube
- Vortex the tube briefly
- Put the tube in the slot in the luminometer, close the lid, and measure
- Usually measure each sample two times