



DEC 08, 2022

WORKS FOR ME

1


🌐 Preparing 1% Eosin Y stock solution

COMMENTS 0

This protocol is published without a DOI.

joshdm¹

¹UC Davis

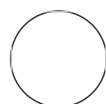

joshdm

ABSTRACT

Preparation of 1% stock solution of Eosin Y for routine H&E staining.

PROTOCOL CITATION

joshdm 2022. Preparing 1% Eosin Y stock solution. **protocols.io**
<https://protocols.io/view/preparing-1-eosin-y-stock-solution-cj9iur4e>



LICENSE

————— This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

CREATED

Dec 08, 2022

LAST MODIFIED

Dec 08, 2022

PROTOCOL INTEGER ID

73738

- 1 Weigh 2g of Eosin Y disodium salt and place it into a 1000mL Earlenmyer flask
- 2 Measure 40mL of dionized water and add it to the flask. Swirl to dissolve Eosin powder.

- 3 Measure 160mL of 95% ethanol and add it to the flask. Swirl to dissolve.
- 4 Using a p1000 Pipette, add 1mL (1000uL) of Glacial Acetic acid to the solution. Swirl until homogeneous.