




MAY 13, 2023

NCEM Drop - Paramyxo confirmation, generation of RNP (TM-014)

sandra.crameri¹

¹CSIRO ACDP

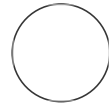

sandra.crameri

ABSTRACT

Generate RNP to confirm Paramyxoviurs

OPEN  ACCESS

DOI:

dx.doi.org/10.17504/protocols.io.x54v9dq8pg3e/v1

Protocol Citation: sandra.crameri 2023. NCEM Drop - Paramyxo confirmation, generation of RNP (TM-014).

protocols.io

<https://dx.doi.org/10.17504/protocols.io.x54v9dq8pg3e/v1>

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working

We use this protocol and it's working

Created: May 13, 2023

Last Modified: May 13, 2023

PROTOCOL integer ID:
81822

HEADER


1 SAN:


SPEC No:

OPERATOR:

Generation of RNP

2 Prepare [M] 0.5 Mass / % volume



 Deoxycholic acid sodium salt Merck MilliporeSigma (Sigma-Aldrich) Catalog #D-6750 in TC Water.

3 Mix 1:1 with the sample for  00:05:00

5m

Conventional

11m

4 Adsorb  10 μ L sample to grid  00:10:00 , inspect to ensure sample does not dry out.

10m

5 Drain excess sample from grid using filter paper, leave wet.

6 Stain  nano-W Contributed by users Catalog ##2018-5ML  00:01:00

1m



7 Drain & dry using filter paper