



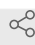
Version 2 ▼

Oct 10, 2022

# MagAttract + Metapolyzyme metagenomic gDNA extraction from skin swabs V.2

Natalie Ring<sup>1</sup><sup>1</sup>Roslin Institute & Royal (Dick) School of Veterinary Studies, University of Edinburgh

1 Works for me

 Share[dx.doi.org/10.17504/protocols.io.q26g7yr19gwz/v2](https://dx.doi.org/10.17504/protocols.io.q26g7yr19gwz/v2) Dogstails Natalie Ring

## ABSTRACT

A protocol for the metagenomic extraction of bacterial DNA from skin swab samples (optimised using canine swabs), for use in a rapid diagnostics pipeline. At the end of the protocol, the DNA is cleaned up and ready for rapid barcoding (SQK-RBK004) library preparation for nanopore sequencing (or whatever other application you want to do).

Unless otherwise stated, all reagents should be included in the listed kits.

DOI

[dx.doi.org/10.17504/protocols.io.q26g7yr19gwz/v2](https://dx.doi.org/10.17504/protocols.io.q26g7yr19gwz/v2)

## PROTOCOL CITATION

Natalie Ring 2022. MagAttract + Metapolyzyme metagenomic gDNA extraction from skin swabs. **protocols.io**

<https://dx.doi.org/10.17504/protocols.io.q26g7yr19gwz/v2>

Version created by Natalie Ring

MANUSCRIPT CITATION please remember to cite the following publication along with this protocol

Coming soon.

## LICENSE

————— This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

CREATED

Oct 10, 2022

LAST MODIFIED

Oct 10, 2022

PROTOCOL INTEGER ID

71097

GUIDELINES

This protocol, an adaptation of Qiagen's MagAttract HMW DNA kit, was developed by Natalie Ring and Alison Low for the Dogstails project, a collaboration between the Roslin Institute and the Royal (Dick) School of Veterinary Studies funded by the Dogs Trust. We are grateful to the dogs (and their owners) who donated samples to the R(D)SVS's Hospital for Small Animals, many of which were used in the development of this protocol.

Please follow on Twitter for latest updates, papers and results:

@NatalieAnneRing

MATERIALS TEXT

### Kits

 [MagAttract HMW DNA](#)

[kit Qiagen Catalog #67563](#)





 [ProNex Size-Selective Purification](#)

[System Promega Catalog #NG2001](#)

 [Qubit® dsDNA HS Assay Kit Thermo Fisher](#)

[Scientific Catalog #Q32854](#) In 2 steps

### Other reagents

- 50 mM Tris, 10 mM EDTA, pH8.0 ("buffer P1")
-  [1X PBS \(Phosphate-buffered saline \)](#)
-  [MetaPolyzyme Sigma](#)
- [Aldrich Catalog #MAC4L-5MG](#)
-  [Nuclease-free Water Contributed by users](#)
-  [Distilled Water Contributed by users](#)

## Equipment

☒ Swabs in tubes (no liquid) **VWR**

**Scientific Catalog # 710-0928**

Any swab tube with no transport/storage medium is fine.

DNA LoBind tubes, 1.5 mL  
Tubes

Eppendorf      022431021      [↗](#)  
1.5 mL

OR

SafeSeal reaction tube, 1.5 ml, PP, PCR  
Performance Tested, Low DNA-binding  
Tubes

Sarstedt      72.706.700      [↗](#)  
1.5 mL

Magnetic Stand  
Magnetic Stand

Thermo Scientific      MR02      [↗](#)

Any magnetic rack that fits your tubes will  
suffice.



### Centrifuge

Benchtop Centrifuge

Eppendorf 5405000441 [Link](#)

Any benchtop centrifuge will suffice



### ThermoMixer

Benchtop Incubator

Eppendorf 5382000023 [Link](#)

Any heat block will suffice



### Mini-centrifuge

Centrifuge

Fisher S67601B [Link](#)

Any standard mini centrifuge with adapters for different tube sizes will suffice



Vortex Mixer

Vortex Mixer

VWR

97043-562



#### BEFORE STARTING

- **"Buffer P1"** is required for the metapolyzyme lysis incubation: 50 mM Tris, 10 mM EDTA, pH 8.0
- **Metapolyzyme** is used here at a concentration of 3.3 mg/ml (resuspend 5 mg lyophilized powder in 1.5 ml PBS pH 7.5)
- We recommend using **low DNA-binding tubes** throughout, but definitely for the elution/storage of DNA

#### Extended pre-lysis spin down

- 1 Bathe swab tip in 3 ml PBS in the swab tube for 10 minutes, with occasional vortexing. 10m  
Remove swab from tube, squeezing the sides as you do.

**3 mL PBS**

**00:10:00**

Bathing the swab tip in PBS overnight yields much more DNA, if you have time.

- 2 Pellet 2x 1.5 ml aliquots of cell-PBS solution in 1.5 ml tubes by centrifuging at maximum speed 20m  
(13,000 RPM) for 20 minutes, then discard supernatant

 **3 mL PBS**

 **16,000 x g, Room temperature, 00:20:00**

We have found that this extended spin at the beginning of the protocol results in much better yield of bacterial gDNA, especially in samples with low bacterial abundance

#### Metapolyzyme & Proteinase K Lysis

- 3** Resuspend cell pellets (which might be invisible) and combine in 160 µl buffer P1 (50 mM Tris, 10 mM EDTA, pH 8.0)

 **160 µL buffer P1**

- 4** Add 20 µl metapolyzyme (3.3 mg/ml, 5 mg resuspended in 1500 µl PBS) and mix by flicking the tube

 **20 µL metapolyzyme (3.3 mg/ml)**

- 5** Incubate on a thermomixer for 60 minutes at 37°C with 900 RPM shaking

1h

 **900 rpm, 37°C, 01:00:00**

- 6** Add 20 µl MagAttract proteinase K and mix by flicking the tube

 **20 µL proteinase K**

- 7** Incubate on a thermomixer for 30 minutes at 56°C with 900 RPM shaking

30m

 **900 rpm, 56°C, 00:30:00**

## MagAttract DNA isolation and washing

- 8 Add 150 µl MagAttract buffer AL and mix by pulse vortexing

 **150 µL buffer AL**

Our standard "pulse vortex" is 10 short (<1 second) pulses per tube

- 9 Add 15 µl MagAttract Suspension G and 280 µl MagAttract buffer MB and mix by pulse vortexing

 **15 µL Suspension G**

 **280 µL Buffer MB**


Make sure the magnetic beads (Suspension G) are really well mixed before adding them! The whole suspension should be black, not separated into a bead layer and a clear layer. We usually resuspended by vortexing for 10 or more seconds.

- 10 Incubate on a thermomixer for 3 minutes at room temperature with 1,400 RPM shaking

 **1400 rpm, Room temperature , 00:03:00**

- 11 Spin down briefly, then pellet beads on magnet and remove supernatant

- 12 Add 700 µl MagAttract buffer MW1 and incubate on a thermomixer for 1 minute at room temperature with 1,400 RPM shaking <sup>1m</sup>

 **700 µL buffer MW1**

 **1400 rpm, Room temperature , 00:01:00**

13 Repeat steps 11 and 12 1m

14 Spin down briefly, then pellet beads on magnet and remove supernatant

15 Add 700 µl MagAttract buffer PE and incubate on a thermomixer for 1 minute at room temperature with 1,400 RPM shaking 1m

 **700 µL buffer PE**

 **1400 rpm, Room temperature , 00:01:00**

16 Repeat steps 14 and 15 1m

17 Spin down briefly, then pellet beads on magnet and remove supernatant



18 Rinse the pelleted beads on the magnetic rack with 700 µl distilled water by pipetting down the opposite wall of the tube, then incubate for 1 minute on magnetic rack

 **700 µL distilled water**

19 Remove distilled water

20 Repeat steps 18 and 19



- 21 Spin down briefly, then remove any remaining supernatant
- 22 Add 50 µl nuclease-free water off the magnet, to resuspend the bead pellet
-  **50 µL nuclease-free water**
- 23 Incubate on a thermomixer for 3 minutes at room temperature with 1,400 RPM shaking 3m
-  **1400 rpm, Room temperature , 00:03:00**
- 24 Spin down briefly, then pellet beads on magnetic rack and **keep supernatant** in a low-DNA binding 1.5 mL tube (e.g. [Eppendorf](#) or [Sarstedt](#))

#### Qubit Pre-clean-up quantification

- 25 Quantify DNA using Qubit dsDNA HS kit. If DNA concentration is an appropriate concentration for your experiment (for us, this means at least 0.2 ng/µl), continue to clean-up steps.

 **Qubit® dsDNA HS Assay Kit Thermo Fisher**

**Scientific Catalog #Q32854**

 **1 µL DNA**

 **199 µL Qubit dsDNA HS working solution**

#### ProNex DNA clean-up

- 26 Add 150 µl room temperature ProNex beads to your entire tube of DNA (49 µl)

 **200 µL ProNex beads**

Like the magnetic beads in Suspension G, make sure the ProNex beads are really well mixed (10+ seconds of vortexing) immediately before you use them.

27 Mix well by slowly pipetting up and down 10 times

28 Incubate at room temperature for 10 minutes (no shaking needed)

10m

🕒 00:10:00

🌡 Room temperature

29 Spin down briefly, then pellet beads on magnet and remove supernatant

30 Rinse the pelleted beads on the magnetic rack by pipetting 200 µl ProNex Wash Buffer down <sup>1m</sup>  
**the opposite wall of the tube**, then incubate at room temperature for 60 seconds (no shaking), then remove Wash Buffer

📄 200 µL Wash Buffer

🌡 Room temperature

🕒 00:01:00

31 Repeat step 26

32 **Air-dry** (lid open) the sample on the magnetic rack for 5 minutes (longer is OK, no more than <sup>5m</sup>  
60 minutes)

🌡 Room temperature

🕒 00:05:00

33 Add 20 µl nuclease-free water off the magnet. Resuspend the pellet by **flicking the tube**, <sup>5m</sup>  
then incubate at room temperature for 5 minutes (no shaking needed)

📄 20 µL nuclease-free water

🌡 Room temperature

- 34 Spin down briefly, then pellet the beads on magnet and **keep supernatant** in a low DNA-binding tube

Qubit post-clean-up quantification

- 35 Quantify DNA using Qubit dsDNA HS kit. If DNA concentration is an appropriate concentration for your experiment (for us, this means at least 0.2 ng/μl), continue to library preparation.

🔗 [Qubit® dsDNA HS Assay Kit Thermo Fisher](#)

Scientific Catalog #Q32854

🧴 1 μL DNA

🧴 199 μL Qubit dsDNA HS working solution