



Nov 20, 2020

Fungal Plate Photography

Craig Bateman¹¹University of Florida

1

Works for me

dx.doi.org/10.17504/protocols.io.bnu6meze

Bark Beetle Mycobiome Research Coordination Network

ABSTRACT

This protocol describes how to photograph fungal plates.

This protocol is part of the Bark Beetle Mycobiome (BBM) Research Coordination Network. For more information on the BBM international network: Hulcr J, Barnes I, De Beer ZW, Duong TA, Gazis R, Johnson AJ, Jusino MA, Kasson MT, Li Y, Lynch S, Mayers C, Musvuugwa T, Roets F, Seltmann KC, Six D, Vanderpool D, & Villari C. 2020. Bark beetle mycobiome: collaboratively defined research priorities on a widespread insect-fungus symbiosis. *Symbiosis* 81: 101–113 <https://doi.org/10.1007/s13199-020-00686-9>.

DOI

dx.doi.org/10.17504/protocols.io.bnu6meze

DOCUMENT CITATION

Craig Bateman 2020. Fungal Plate Photography. **protocols.io**
<https://dx.doi.org/10.17504/protocols.io.bnu6meze>

LICENSE

This is an open access document distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

CREATED

Oct 23, 2020

LAST MODIFIED

Nov 20, 2020

DOCUMENT INTEGER ID


43646

ABSTRACT

This protocol describes how to photograph fungal plates.

This protocol is part of the Bark Beetle Mycobiome (BBM) Research Coordination Network. For more information on the BBM international network: Hulcr J, Barnes I, De Beer ZW, Duong TA, Gazis R, Johnson AJ, Jusino MA, Kasson MT, Li Y, Lynch S, Mayers C, Musvuugwa T, Roets F, Seltmann KC, Six D, Vanderpool D, & Villari C. 2020. Bark beetle mycobiome: collaboratively defined research priorities on a widespread insect-fungus symbiosis. *Symbiosis* 81: 101–113 <https://doi.org/10.1007/s13199-020-00686-9>.

1. Download Canon EOS Utility to your computer: <https://www.usa.canon.com/internet/portal/us/home/support/self-help-center/eos-utility> (Listed as "EOS Digital Solution Disk Software")
2. Ensure photo setup is as pictured (in the flow hood in Adam's office on the 2nd floor):

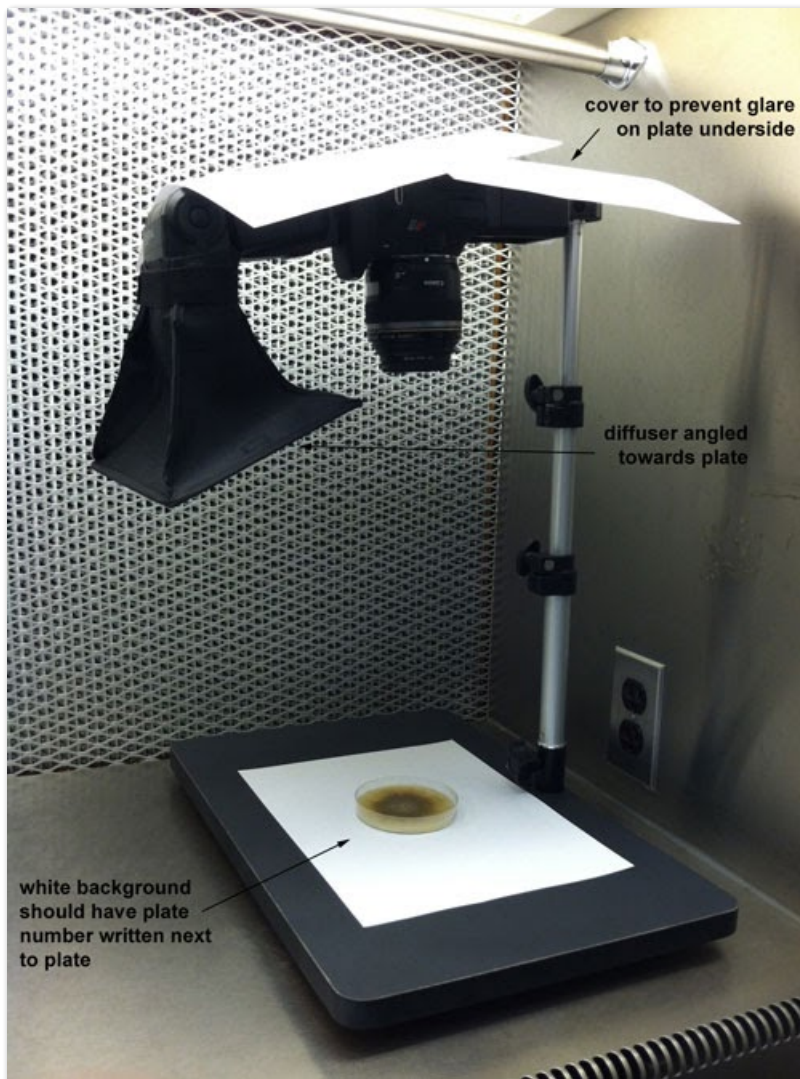

1. Hook your computer up to the camera to enable live shooting and direct photo downloads.
2. Turn on the hood blower and wipe down all surfaces with 70% ethanol.
3. Set the camera dial to M (Manual) and adjust to the following settings:
 - Manual focus
 - Large resolution (use Raw for poster photos)
 - Shutter 1/100
 - F9.0
 - ISO 400
 - Spot metering
 - Flash white balance setting
 - WB Shift: B5 G2
 - As it appears on the EOS Utility:



1. Set a plate open-side-up underneath the camera and manually focus the camera to the top of the agar.
2. Add a piece of paper displaying the underlying plate number in the photo.
3. Proceed taking photographs of the top and bottom of each plate. Settings do not need to be altered between plates/photographs.
4. Save photographs to our lab drive (Forest Entomology:\Plate photos) and name them according to plate number (eg: "plate number XXXX.jpg")
5. Keep the plate at 4-10° C until ready for discard. Plates should be discarded into an autoclave bag.
 **Note on batteries for flash: Only use Panasonic eneloop batteries. All batteries should be equally charged. Disconnect batteries after use. Battery charger is in camera drawer in Hulcr lab.