




Aug 31, 2022

# Proteoform Identification and Quantitation with TopPIC and TDPortal for Human Tissues

James M Fulcher<sup>1</sup>, Yen-Chen Liao<sup>2</sup>, Mowei Zhou<sup>2</sup>, Ljiljana.PasaTolic<sup>2</sup><sup>1</sup>Pacific Northwest National lab; <sup>2</sup>Pacific Northwest National Laboratory

In Development

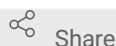

Share

[dx.doi.org/10.17504/protocols.io.3byl4bpj2vo5/v1](https://doi.org/10.17504/protocols.io.3byl4bpj2vo5/v1)

Human BioMolecular Atlas Program (HuBMAP) Method Development Community  
PNNL-TTD



Mowei Zhou  
Pacific Northwest National Laboratory

## ABSTRACT

This protocol describe a workflow for top-down proteomics analysis. Top-down proteomics data are processed with two separate software packages TopPIC and TDPortal. Proteoform identifications were merged from the two software with unified FDR to increase coverage. TopPICR was separately used to cluster TopPIC proteoform to extract abundances for label-free quantitation.

## DOI

[dx.doi.org/10.17504/protocols.io.3byl4bpj2vo5/v1](https://doi.org/10.17504/protocols.io.3byl4bpj2vo5/v1)

## PROTOCOL CITATION

James M Fulcher, Yen-Chen Liao, Mowei Zhou, Ljiljana.PasaTolic 2022.  
Proteoform Identification and Quantitation with TopPIC and TDPortal for Human Tissues. **protocols.io**  
<https://protocols.io/view/proteoform-identification-and-quantitation-with-to-ceqztdx6>



## FUNDERS ACKNOWLEDGEMENT

National Institutes of Health (NIH) Common Fund, Human Biomolecular Atlas Program (HuBMAP)  
Grant ID: UG3CA256959-01

## LICENSE

————— This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

## CREATED

Aug 02, 2022

## LAST MODIFIED

Aug 31, 2022

## PROTOCOL INTEGER ID

68089

## PARENT PROTOCOLS

In steps of

[Overall protocol for MicroPOTS LCMS top down proteomics of kidney tissue sections](#)

Cited in

[Top Down Proteomics Data Collection for Microdissected Kidney Tissue Functional Units](#)

[Top Down Proteomics Data Collection for Microdissected Pancreas Tissue Functional Units](#)

## TopPIC Processing

- 1 Convert Instrument raw data to mzML using MSConvert

**MSConvert** [↗](#)

- 2 Analyze mzML files using the TopPIC Suite (version 1.4.13.1) .

**TopPIC Suite 1.4.13.1** [↗](#)



[source](#) by Xiaowen Liu

## 2.1 TopFD Parameters—————

Spectral data type:	Centroid
Maximum charge:	30
Maximum monoisotopic mass:	50000 Dalton
Peak error tolerance:	0.02 m/z
MS1 signal/noise ratio:	3
MS/MS signal/noise ratio:	1
Thread number:	10
Precursor window size:	2 m/z
Use Env CNN model:	No
Miss MS1 spectra:	No
Generate Html files:	Yes
Do final filtering:	Yes

## 2.2 TopPIC 1.4.13 Parameters-----

\*\*\*\*\* Parameters \*\*\*\*\*

Protein database file:	 <a href="#">ID_008032_8627C6BD.fasta.zip</a>
Spectrum file:	xxxxxxxxxxxxxxxx_ms2.msalgn
Number of combined spectra:	1
Fragmentation method:	FILE
Search type:	TARGET
Fixed modifications:	None
Use TopFD feature file:	True
Maximum number of unexpected modifications:	1
Error tolerance for matching masses:	15 ppm
Error tolerance for identifying PrSM clusters:	0.8 Da
Spectrum-level cutoff type:	EVALUE
Spectrum-level cutoff value:	0.05
Proteoform-level cutoff type:	EVALUE
Proteoform-level cutoff value:	0.05
Allowed N-terminal forms:	NONE,NME,NME_ACETYLATION,M_ACETYLATION
Maximum mass shift of modifications:	275 Da
Minimum mass shift of modifications:	-150 Da
Thread number:	14
E-value computation:	Generating function
Common modification file name:	 <a href="#">Dynamic_mods.txt</a>
MISScore threshold:	0.15
Executable file directory:	
Version:	1.4.13

The protein fasta contains human proteome from UniProt with both SwissProt and TREMBL sequences. Decoy sequences were added as well. Unzip the attachment to use it.

- 3 TopPIC outputs proteoform spectrum matches (PrSMs) as tab-separated files (...toppic\_prsm.tsv) and quantification data within MS1 feature files (...\_ms1.feature). These are both imported into the R environment for post-processing with TopPICR.
- 4 TopPICR is used for post-processing to improve proteoform identification and quantification. All functions are documented within the TopPICR R package.

## TopPICR

[source](#) by Evan Martin

- 4.1 First, result files are read into R using the `read_toppic(file_path = path, file_name = names)` function in TopPICR, where the "path" is the path to the directory containing the TopPIC PrSM files and "names" is a character vector specifying the PrSM files to import. This function can also be utilized to import the MS1 feature files into a separate object.
- 4.2 Next, the data is further processed with the `augment_annotation()` and `rm_false_gene()` functions to account for ambiguity in proteoform identifications
- 4.3 False discovery rate (FDR) filtering is accomplished by finding the appropriate E-value cutoff to filter the results to 1% FDR at the isoform and protein level. This is provided by the `find_evalue_cutoff()` and `apply_evalue_cutoff()` functions.
- 4.4 Proteoform inference is performed with `infer_pf()` function and the proteoform level is determined with `set_pf_level()` function .

Smith LM, Thomas PM, Shortreed MR, Schaffer LV, Fellers RT, LeDuc RD, Tucholski T, Ge Y, Agar JN, Anderson LC, Chamot-Rooke J, Gault J, Loo JA, Paša-Tolić L, Robinson CV, Schlüter H, Tsybin YO, Vilaseca M, Vizcaíno JA, Danis PO, Kelleher NL (2019). A five-level classification system for proteoform identifications.. Nature methods.

<https://doi.org/10.1038/s41592-019-0573-x>

- 4.5 Retention time alignment is processed with the `form_model()` and `align_rt()` functions.
  - 4.6 Mass calibration is accomplished with the `calc_error()` and `recalibrate_mass()` functions
  - 4.7 Clustering and deisotoping error correction is performed with the `cluster()` and `create_pcg()` functions.
  - 4.8 Metadata for each proteoform cluster is generated with the `create_mdata()` function.
  - 4.9 Steps 4.5 and 4.6 are applied to the MS1 feature files as well before features are matched and combined (for MBR) with the `match_features()` and `combine_features()` functions.
- 5 The final table of proteoform identification and quantitation results from TopPIC Suite and TopPICR are exported as comma-separated value (.csv) files.

## TDPortal Processing

- 6 Request TDPortal access and follow their instructions to set up an account.

### TDPortal

by Northwestern University

#### TDPortal search process

##### 6.1 Upload data

##### 6.2 Search on TDPortal

TDportal has an option for label-free quantitation, but it is not used in this workflow.

##### 6.1 Upload data

1. Connect to Northwestern through VPN.

(<https://kb.northwestern.edu/page.php?id=94726>)

2. Copy the files to your user folder. (Eg. [\\resfiles.northwestern.edu\NU-PCEDATA\external\\_users\XXXXX](https://resfiles.northwestern.edu/NU-PCEDATA/external_users/XXXXX))
3. The system will ask you to log in. Please use "ads\your id" with your password to log into your folder.
4. Create a sub-folder under your user folder with each search.
5. Put raw files to the sub-folder accordingly and do not have more folders under the sub-folder. (<https://kb.northwestern.edu/page.php?id=70525>).

## 6.2 Search on TDPortal

([https://portal.nrtdp.northwestern.edu/static/TDPortalSOP\\_043\\_20180301.pdf](https://portal.nrtdp.northwestern.edu/static/TDPortalSOP_043_20180301.pdf))



1. Connect website(<https://portal.nrtdp.northwestern.edu/>)
2. Log in TDPortal with "your email address" and "your password"
3. Your subfolder's name will show as each dataset.
4. Select files into the "Input files" under the selected dataset.
5. Select organism "human".
6. Set parameters as follow:  
User empirical P-score: False  
Filter by FDR: True  
Create SAS input sheet for quant: Select True when we need to.  
Precursor resolution: High resolution  
Fragmentation Type: Auto (or the type we used on MS).  
Code set: Standard 4.0.0  
Include ProSight Error Tolerance Search: False (select "true" when we want to allow one unknown mass shift in the proteoform).

## 7 Exporting TDPortal results

### TDViewer

by Northwestern University

1. Download \*.tdReports file. Note: There can be two separate processes created in the queue. One for ID results in the TDReport. Another is the CSV file for quantitation (if enabled).
2. Click the download icon to download these files.

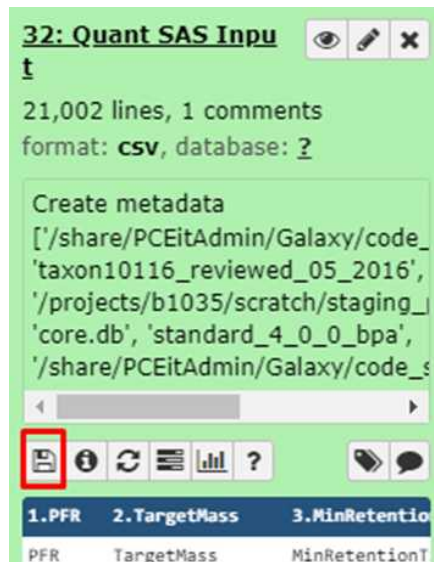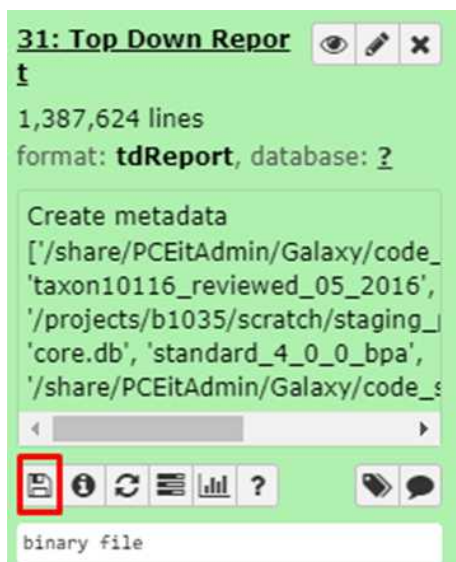

3. Open TDReport by TDViewer 2.0(<http://tdviewer2.northwestern.edu/>)
4. Read and export proteoform ID results from TDViewer with 1% FDR cutoff.

#### Combining Results

- 8 Results (proteoform spectral matches) from TopPIC and TDPortal are then merged using a function written in R that is openly available on GitHub. The input proteoform tables from each software was pre-filtered with FDR cutoff of 1% (adjusted FDR in TopPICR for TopPIC, and the default FDR in TDPortal).

#### TDPortal\_TopPIC\_Join [↗](#)

[source](#) by James M Fulcher

#### Final output

- 9 Results for proteoform spectral matches (merged from TopPIC and TDPortal) and proteoform quantitation (TopPICR) are uploaded to HIVE.