



VERSION 4

APR 19, 2023

OPEN ACCESS

**DOI:**  
[dx.doi.org/10.17504/protocols.io.n2bvj6mdnlk5/v4](https://dx.doi.org/10.17504/protocols.io.n2bvj6mdnlk5/v4)

**Protocol Citation:** Tsu-Chun Hung, Yin-Tse Huang 2023. DNA extraction (BOMB). **protocols.io** <https://dx.doi.org/10.17504/protocols.io.n2bvj6mdnlk5/v4> Version created by Yin-Tse Huang

**License:** This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

**Protocol status:** Working  
 We use this protocol and it's working

**Created:** Apr 06, 2023

**Last Modified:** Apr 19, 2023

**PROTOCOL integer ID:**  
 80088

## DNA extraction (BOMB) V.4

Tsu-Chun Hung<sup>1</sup>, Yin-Tse Huang<sup>1</sup>

<sup>1</sup>KMU



Yin-Tse Huang

### ABSTRACT

DNA extraction (BOMB)

### MATERIALS

1. Lysis master mix (870 uL/sample)

A	B
TE buffer	225 uL
Lysis buffer	375 uL
Ammonium acetate	270 uL

2. TE buffer


A	B
Tris HCl pH8.0	10mM
EDTA	1mM

3. Lysis buffer


A	B
GITC	4M
Tris HCl pH8.0	50mM
SDS	0.5g
EDTA	20mM

## Sample Collection

3m

1 Add  200  $\mu$ L of **1mm beads** to 1.5mL eppendorf tube 30s


2 Add  200  $\mu$ L of **0.5mm beads** to 1.5mL eppendorf tube 30s

3 Add  870  $\mu$ L Lysis master mix to 1.5mL eppendorf tube 30s

#### Note


In 4°C fridge

Lysis master mix: **225  $\mu$ L** of TE buffer + **375  $\mu$ L** of lysis buffer + **270  $\mu$ L** of 10M ammonium acetate

4 Collect  10-20 mg of **sample** to 1.5mL eppendorf tube 1m


## Sample crush

4m

5 Put 1.5mL eppendorf tube in mixmill for sample crush, at this condition: 30 rpm/s, for 4mins 4m  
 00:04:00

## Centrifugation

3m


6 Put 1.5mL eppendorf tube in centrifuge for centrifugation, at this condition: 3m  
 17.0 x g, 25°C, 00:03:00

## DNA purification

37m 30s

## 7 Add 350 $\mu$ L of isopropanol to the 1st well of 96 well plate

30s

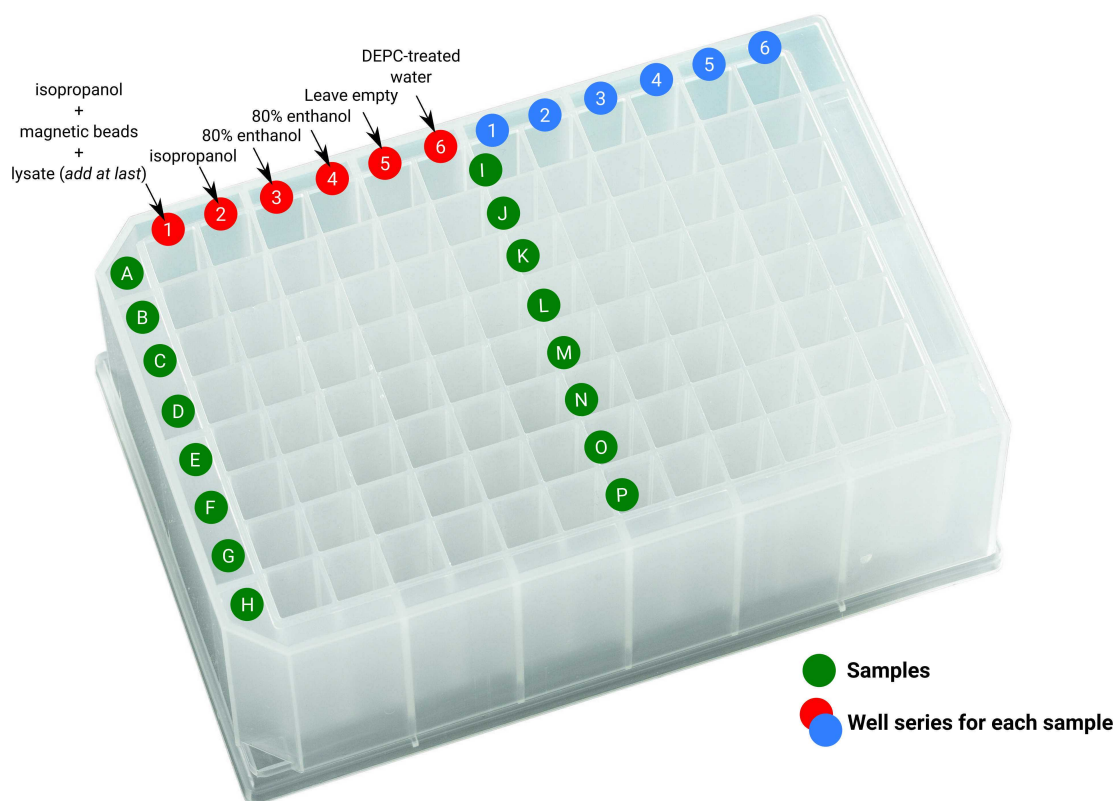

### 7.1 Add 100 $\mu$ L of magnetic beads (10 mg/ml) to the 1st well of 96 deep well plate

30s

#### Note

Vortex the bottle and pipetting before using magnetic beads; re-do vortex after adding to 3 samples to prevent set down of magnetic beads







### 7.2 Add 200-300 $\mu$ L of the sample (lysate) from the 1.5mL centrifuged tube to the 1st well of 96 deep well plate

30s

#### Note

Pipetting **as much lysate as you can**, as long as it's free of any cell debris (no solids in your tip)

**ADD at LAST**

- 8 Add  400 µL of **isopropanol** to the 2nd well of 96 deep well plate 30s
- 9 Add  300 µL of **80% ethanol** to the 3rd well of 96 deep well plate 30s
- 10 Add  300 µL of **80% ethanol** to the 4th well of 96 deep well plate 30s
- 11 Add  300 µL of **DDW** to the 5th well of 96 deep well plate
- 12 Add  100 µL of **DEPC-treated water** to the 6th well of 96 deep well plate 30s
- 13 Put the prepared 96 deep well plate in the automated DNA extraction machine and select the BOMB protocol 34m
- 14 After the extraction is done, collect  50-100 µL of the **eluted sample** as the DNA template for downstream experiments