



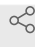
Sep 06, 2022

# 🌐 Maintenance and inactivation of mouse embryonic fibroblasts (MEFs) as feeder cells for human pluripotent stem cell culture

Hanqin Li<sup>1</sup>, Oriol Busquets<sup>2</sup>, Steven Poser<sup>2</sup>, Dirk Hockemeyer<sup>1</sup>, Frank Soldner<sup>2</sup>

<sup>1</sup>University of California, Berkeley; <sup>2</sup>Albert Einstein College of Medicine

1 *Works for me*

 Share

[dx.doi.org/10.17504/protocols.io.b4pbqvin](https://dx.doi.org/10.17504/protocols.io.b4pbqvin)



Devin E Snyder

## ABSTRACT

This collection describes the maintenance and inactivation of mouse embryonic fibroblasts (MEFs) as feeder cells for human pluripotent stem cell (hPSC) culture.

## Collection overview

Thawing of mouse embryonic fibroblasts (MEFs) for hPSC cultures

Expansion of mouse embryonic fibroblasts (MEFs) for hPSC cultures

Freezing of mouse embryonic fibroblasts (MEFs) for hPSC cultures

Harvesting and irradiation of mouse embryonic fibroblasts (MEFs) for hPSC cultures

Mitomycin C inactivation of mouse embryonic fibroblasts (MEFs) for hPSC cultures

Preparation of mouse embryonic fibroblast (MEF) feeder plates for hPSC cultures

A. Starting with frozen irradiated or Mitomycin C inactivated MEFs (optional)

B. Starting with fresh irradiated or Mitomycin C inactivated MEFs

## General notes

1. Throughout these protocols, the term hPSC is used to collectively refer to both hiPSCs and hESC. All described procedures have been tested and work equally well for hiPSCs and hESCs.

2. Either fresh or frozen stocks of irradiated or Mitomycin C inactivated MEFs can be used to prepare hPSC feeder cells.

3. The indicated MEF densities are recommended starting densities and might have to be adjusted for each hPSC line and hPSC media formulation (KSR, serum-free versus serum-containing media).

4. MEFs were obtained as described in *Manipulating the Mouse Embryo: A Laboratory Manual*, Third Edition (ISBN: 0879695919)

Andras Nagy, Marina Gertsenstein, Kristina Vintersten, & Richard Behringer. *Manipulating the Mouse Embryo: A Laboratory Manual*, 3rd ed.. Cold Spring Harbor Laboratory Press.

DOI

[dx.doi.org/10.17504/protocols.io.b4pbqvin](https://dx.doi.org/10.17504/protocols.io.b4pbqvin)

## COLLECTION CITATION

Hanqin Li, Oriol Busquets, Steven Poser, Dirk Hockemeyer, Frank Soldner 2022. Maintenance and inactivation of mouse embryonic fibroblasts (MEFs) as feeder cells for human pluripotent stem cell culture. **protocols.io**  
<https://dx.doi.org/10.17504/protocols.io.b4pbqvin>



#### FUNDERS ACKNOWLEDGEMENT

Aligning Science Across Parkinson's  
Grant ID: ASAP-000486

#### KEYWORDS

ASAPCRN

#### LICENSE

————— This is an open access collection distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

#### CREATED

Feb 03, 2022

#### LAST MODIFIED

Sep 06, 2022

#### COLLECTION INTEGER ID

57795

## MATERIALS TEXT

Item	Vendor	Catalog #
DMEM	Corning	10-013-CV
FB Essence	Avantor	10803-034
FBS	Gibco	10437028
200mM L-Glutamine	Sigma	G8540
Penicillin & Streptomycin	Gibco	15140-122
MEM Non-Essential Amino Acids	Gibco	11140-050
DMSO	Thermo Fisher	BP231-100
Gelatin powder	Sigma	G2625
0.25% Trypsin/EDTA	Gibco	25200-056
DNase	Roche	4536282001
DPBS w/o Ca & Mg	Corning	MT21031CV
50ml centrifuge tubes	Corning	1495949A
Nunc 1.8 ml cryovials	Thermo Fisher	377267
10cm petri dish	Fisher	08757100D
10ml serological pipet	Corning	7200574
15cm tissue culture dish	Corning	0877224
No5. forceps	Roboz	RS-5010
Microdissecting curved surgical scissors	Roboz	RS-5881
Styrofoam microtube freezer box	Labnet	R8000

## ABSTRACT

This collection describes the maintenance and inactivation of mouse embryonic fibroblasts (MEFs) as feeder cells for human pluripotent stem cell (hPSC) culture.

## Collection overview

Thawing of mouse embryonic fibroblasts (MEFs) for hPSC cultures  
 Expansion of mouse embryonic fibroblasts (MEFs) for hPSC cultures  
 Freezing of mouse embryonic fibroblasts (MEFs) for hPSC cultures  
 Harvesting and irradiation of mouse embryonic fibroblasts (MEFs) for hPSC cultures  
 Mitomycin C inactivation of mouse embryonic fibroblasts (MEFs) for hPSC cultures  
 Preparation of mouse embryonic fibroblast (MEF) feeder plates for hPSC cultures  
 A. Starting with frozen irradiated or Mitomycin C inactivated MEFs (optional)  
 B. Starting with fresh irradiated or Mitomycin C inactivated MEFs

## General notes

1. Throughout these protocols, the term hPSC is used to collectively refer to both hiPSCs and hESCs. All described procedures have been tested and work equally well

for hiPSCs and hESCs.


2. Either fresh or frozen stocks of irradiated or Mitomycin C inactivated MEFs can be used to prepare hPSC feeder cells.

3. The indicated MEF densities are recommended starting densities and might have to be adjusted for each hPSC line and hPSC media formulation (KSR, serum-free versus serum-containing media).

4. MEFs were obtained as described in *Manipulating the Mouse Embryo: A Laboratory Manual*, Third Edition (ISBN: 0879695919)

Andras Nagy, Marina Gertsenstein, Kristina Vintersten, & Richard Behringer. *Manipulating the Mouse Embryo: A Laboratory Manual*, 3rd ed.. Cold Spring Harbor Laboratory Press.

## FILES

-   Thawing of mouse embryonic fibroblasts (MEFs) for hPSC cultures  
**Version 1**  
by Devin E Snyder
-   Expansion of mouse embryonic fibroblasts (MEFs) for hPSC cultures  
**Version 1**  
by Devin E Snyder
-   Freezing of mouse embryonic fibroblasts (MEFs) for hPSC cultures  
**Version 1**  
by Devin E Snyder
-   Harvesting and irradiation of mouse embryonic fibroblasts (MEFs) for hPSC cultures  
**Version 1**  
by Devin E Snyder
-   Mitomycin C inactivation of mouse embryonic fibroblasts (MEFs) for hPSC cultures  
**Version 1**  
by Devin E Snyder
-   Preparation of mouse embryonic fibroblast (MEF) feeder plates for hPSC cultures  
**Version 1**  
by Devin E Snyder