




Jun 01, 2021

DNA Short-insert Library Construction Protocol for Illumina HiSeq 2500/4000/X Ten or NovaSeq

Hongfang Zhang¹¹GigaScience Press

1 Works for me

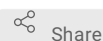

Share


dx.doi.org/10.17504/protocols.io.bvewn3fe

BGI GIGA 2 more workspaces

Hongfang Zhang
GigaScience Press

ABSTRACT

HiSeq DNA pair-end library construction is used for DNA sequencing on Illumina HiSeq 2500/HiSeq 4000/HiSeq X Ten/ and NovaSeq sequencers. The insert size of the DNA fragments can be 300bp, 350bp and 500bp.

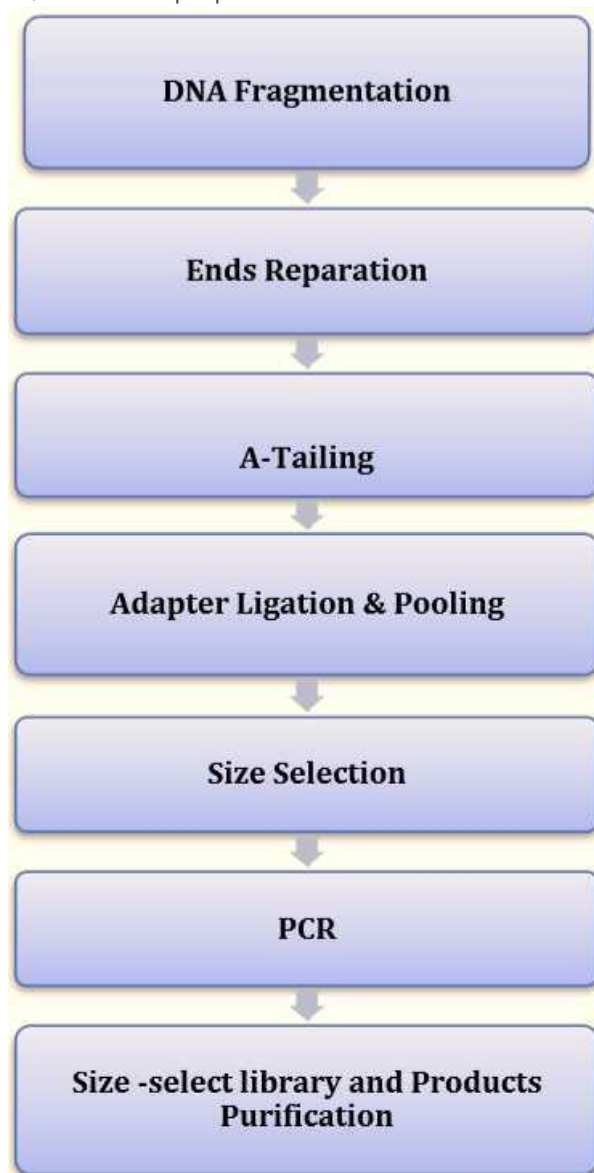

DNA pair end library construction workflow

DOI

dx.doi.org/10.17504/protocols.io.bvewn3fe

PROTOCOL CITATION

Hongfang Zhang 2021. DNA Short-insert Library Construction Protocol for Illumina HiSeq 2500/4000/X Ten or NovaSeq. **protocols.io**
<https://dx.doi.org/10.17504/protocols.io.bvewn3fe>

LICENSE

————— This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

CREATED

May 31, 2021


LAST MODIFIED

Jun 01, 2021

PROTOCOL INTEGER ID



50358

1 Fragment Genomic DNA

 **1 µl** genomic DNA was randomly fragmented using a Covaris ultrasonicator. The fragmented DNAs were tested by Gel-electrophoresis, then purified using an paramagnetic bead based AxyPrep Mag PCR clean up Kit.



2 End Repair

30m

The fragmented DNAs were combined with the End Repair Mix, incubated at  **20 °C** for  **00:30:00**. The end-repaired DNA is then purified with a AxyPrep Mag PCR clean up Kit.



3 Add A-Tails to the 3' end

30m

The repaired DNAs are combined with the A-Tailing Mix, incubated at  **37 °C** for  **00:30:00**. The DNA needs to be purified again with an AxyPrep Mag PCR clean up Kit.

4 Adapter Ligation and Pooling

16h

Illumina adapters are ligated to the Adenylated 3'Ends of the DNA, incubated at  **16 °C** for  **16:00:00**. Purify the Adapter-ligated DNAs with AxyPrep Mag PCR clean up Kit. The concentration is measured by a Caliper high-throughput Bioanalyzer.

5 Fragments Selection

The PCR products are selected by agarose gel electrophoresis with target fragments. Purify the gel with a QIAquick Gel Extraction kit (QIAGEN). Re-dissolve the appropriate amount of Elution Buffer and attach the library tag to complete the library construction.

6 PCR

Several rounds of PCR amplification with PCR Primer Cocktail and PCR Master Mix are performed to enrich the Adapter-ligated DNA fragments.

7 Products Purification

PCR products are selected by running another 2% agarose gel to recover the target fragments. Purify the gel with a QIAquick Gel Extraction kit (QIAGEN). Re-dissolve the appropriate amount of Elution Buffer and attach the library tag to complete the library construction.

8 Library QC and Sequencing

Library was qualified by the Agilent Technologies 2100 bioanalyzer and ABI StepOnePlus RealTime PCR System. Libraries that pass these QC steps can then be sequenced using pair-ends on the Illumina HiSeq/NovaSeq System.