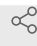


Sep 15, 2022

Plasma Exosome Isolation

Silvia Cerri¹¹Laboratory of Functional Neurochemistry, IRCCS Mondino Foundation, Pavia, Italy

1 Works for me

 Sharedx.doi.org/10.17504/protocols.io.81wgb7xnqvpk/v1

Sofia Koletsi

University College London, University of London

ABSTRACT

This protocol details methods for the isolation of small extracellular vesicles from plasma.

DOI

dx.doi.org/10.17504/protocols.io.81wgb7xnqvpk/v1

PROTOCOL CITATION

Silvia Cerri 2022. Plasma Exosome Isolation. **protocols.io**
<https://protocols.io/view/plasma-exosome-isolation-bqcdmss6>



KEYWORDS

exosome, isolation, plasma, EV, extracellular vesicles, ASAPCRN

LICENSE

————— This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

CREATED

Dec 03, 2020

LAST MODIFIED

Sep 15, 2022

OWNERSHIP HISTORY

Dec 03, 2020



Liz Brydon

Protocols.io

Jan 29, 2021



Sofia Koletsi

University College London, University of London

PROTOCOL INTEGER ID

45157

MATERIALS TEXT

Materials

- 13mm sterile syringe filters with a 0.8 µm Supor (PES) membrane
- 13mm sterile syringe filters with 0.2 µm Supor (PES) membrane
- 1X PBS

Equipment

- Ultracentrifuge

SAFETY WARNINGS

Please refer to the Safety Data Sheets (SDS) for health and environmental hazards.

BEFORE STARTING

Keep samples **⚡ On ice** during the entire procedure!

Plasma Exosome Isolation

2h 30m

1



Incubate plasma samples (starting volume: 1-1.5ml) **⚡ On ice** for complete thaw.

2

Filter samples through a 13mm sterile syringe filters with a **0.8 µm** Supor (PES) membrane.

3



Centrifuge at **⚡ 20000 x g, 4°C, 01:00:00** to remove large extracellular vesicles.

4

Collect the supernatant and filter it through 13mm sterile syringe filters with **0.2 µm** Supor (PES) membrane.

5



Centrifuge at **⚡ 100000 x g, 4°C, 01:00:00** to pellet small extracellular vesicles.

6



Remove the supernatant, resuspend the pellet in **1 mL 1X PBS** and repeat the previous step (repeated again in the next step for convenience).

This step is optional, according the final application of the pellet.

7



Centrifuge at **100000 x g, 4°C, 01:00:00** to pellet small extracellular vesicles.

8

Process the pellet according to the type of analysis.