



APR 12, 2023

🌐 Cell culture and Western blot

Andrea Dickey¹, rabrisch¹

¹UCSD Department of Cellular and Molecular Medicine



Yu Xuan Lin

ABSTRACT

Protocol for detecting Rab7a phosphorylation at S72 by LRRK1 and LRRK1 mutants in HEK 293T cells.

OPEN ACCESS

Protocol Citation: Andrea Dickey, rabrisch 2023. Cell culture and Western blot. **protocols.io** <https://protocols.io/view/cell-culture-and-western-blot-csqcwdsu>

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working
We use this protocol and it's working

Created: Apr 11, 2023

Last Modified: Apr 12, 2023

PROTOCOL integer ID:
80356

Cell culture and transfection

1d 0h 15m




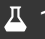
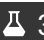


1 Day - 1

1d

Split 293T cells in 6-well dishes 24:00:00 before transfection at 200K/well in 2 mL DMEM (with FBS and with P/S)


2 Day - 2








Transfection

- 2.1** Warm PEI and Opti-MEM at RT for  00:15:00 15m
Label a tube for each well.
- 2.2** Add DNA ( 500 ng of **GFP-Rab7** +  1 µg of **LRRK1 construct of interest**) to  150 µL of Opti-MEM and gently vortex for 5 sec
- 2.3** Add  3 µL room-temperature PEI to each tube containing Opti-MEM/DNA mix. Mix gently by pipetting and incubate DNA/DMEM/PEI mixture at room temperature in the hood for  00:15:00 . 15m
- 2.4** In the meantime...remove media from 6-well dish containing cells and replace with fresh 1 mL DMEM (with FBS and withoutP/S)...put back in 37 deg incubator until ready to add transfection mixture
- 2.5** After 15 min incubation, add DNA/Optimem/PEI mixture dropwise  150 µL onto each well
Give a bit of a light swirl before putting back in incubator

Cell lysis

1d 12h

- 3**  36:00:00 after transfection, begin cell lysis 1d 12h
- 3.1** Wash plate on ice with cold PBS 1x



- 3.2 Add  300 μL RIPA buffer (0.5% Triton, 50 mM Tris pH7.5, 150 mM NaCl, 0.1%SDS) with protease and phosphatase inhibitors (cOmplete mini EDTA free + PhosSTOP tablets)
- 3.3 Lift with cell lifters on ice
- 3.4 Pipet up lysate, put in eppendorf tube, and shake 15 mins in the cold room
- 3.5 Spin at MAX at 4 deg for 15 mins
- 3.6 Remove supernatant and make sample (boil at  95 $^{\circ}\text{C}$ for  00:10:00), store in  -80 $^{\circ}\text{C}$. *I like to store the lysate and take some out to make a sample –I use the 4x NuPage LDS sample buffer:*  65 μL lysate +  10 μL 10x Reducing Agent +  25 μL 4x LDS buffer

10m

Western blot

1h

4 SDS-PAGE with Bis-Tris gel and MOPS running buffer

- 4.1 Load a 4-12% Bis-Tris gel with  25 μL of prepared lysate in sample buffer and run at 180V for ~  00:50:00 or until dye front has reached the bottom of the gel.

50m

- 4.2 Assemble gel with Immobilon-FL PVDF membrane for transfer according to instructions from your western blot transfer apparatus. When using fluorescence detection for Western blot it is important to use low fluorescence background membrane (Immobilon-FL or equivalent).

We activate membrane with MeOH and rinse with water.

Transfer in Western transfer buffer with Tris/Glycine and 20% MeOH at 200 mA for 4 hr at

🌡️ 4 °C

4.3 After transfer is complete, rinse membrane with water and allow to dry between sheets of Whatman paper.

4.4 Block in 5% milk in TBS (no Tween 20)

4.5 Dilute primary antibody in 5% milk with TBST (with Tween 20)
-rabbit anti Rab7 phospho-S72 (MJF-38) at 1:1000
-mouse anti-GFP at 1:2500 (Santa Cruz) (for total Rab quantification)
-rabbit anti-LRRK1 (ab228666) at 1:500
-rabbit anti-GAPDH at 1:3000 (Cell signaling technologies)

4.6 Rock 🕒 Overnight at 🌡️ 4 °C

4.7 Rinse 3x with TBST for 🕒 00:05:00

5m

4.8 Rinse 1x with 5% milk in TBST

4.9 Add secondary antibodies in 5% milk with TBST and incubate at 🌡️ Room temperature for 1 hr
-LiCor mouse and rabbit secondary IRdye antibodies at 1:5000

4.10 Rinse 3x with TBST for 🕒 00:05:00

5m

4.11 Image on LiCor Odyssey CLx

