



Mar 24,  
2020

# Protocol for Differentiation of Blood-Brain Barrier Endothelial Cells from Human Pluripotent Stem Cells v.1

Ethan Lippmann<sup>1</sup>, Hannah Wilson<sup>2</sup>, Emma Neal<sup>1</sup>


<sup>1</sup>Department of Chemical Engineering, Vanderbilt University, Nashville, TN, USA, <sup>2</sup>Department of Biomedical Engineering, Georgia Institute of Technology, Atlanta, GA, USA

**1** Works for me dx.doi.org/10.17504/protocols.io.8g3htyn

Neurodegeneration Method Development Community

## ABSTRACT

Human induced pluripotent stem cell (iPSC)-derived developmental lineages are key tools for in vitro mechanistic interrogations, drug discovery, and disease modeling. iPSCs have previously been differentiated to endothelial cells with blood-brain barrier (BBB) properties, as defined by high transendothelial electrical resistance (TEER), low passive permeability, and active transporter functions. Typical protocols use undefined components, which impart unacceptable variability on the differentiation process. We demonstrate that replacement of serum with fully defined components, from common medium supplements to a simple mixture of insulin, transferrin, and selenium, yields BBB endothelium with TEER in the range of 2,000-8,000  $\Omega \times \text{cm}^2$  across multiple iPSC lines, with appropriate marker expression and active transporters. The use of a fully defined medium vastly improves the consistency of differentiation, and co-culture of BBB endothelium with iPSC-derived astrocytes produces a robust in vitro neurovascular model. This defined differentiation scheme should broadly enable the use of human BBB endothelium for diverse applications.

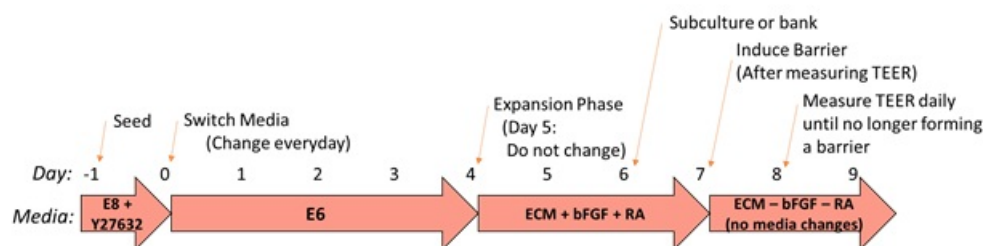

Schematic of E6 method for BBB differentiation

## THIS PROTOCOL ACCOMPANIES THE FOLLOWING PUBLICATION

A Simplified, Fully Defined Differentiation Scheme for Producing Blood-Brain Barrier Endothelial Cells from Human iPSCs. Neal EH, Marinelli NA, Shi Y, McClatchey PM, Balotin KM, Gullett DR, Hagerla KA, Bowman AB, Ess KC, Wikswo JP, Lippmann S. Stem Cell Reports. 2019 Jun 11;12(6):1380-1388. doi: [10.1016/j.stemcr.2019.05.008](https://doi.org/10.1016/j.stemcr.2019.05.008)

## ATTACHMENTS

[BBB\\_differentiation\\_from\\_hPSCs\\_-\\_B27\\_update.pdf](#)

## GUIDELINES

We recommend the following antibodies to monitor BMEC differentiation:

Target antigen	Antibody species	Vendor	Clone or product number	Dilution
PECAM-1	Rabbit	Thermo Scientific	RB-10333P	1:10 (FC) 1:25 (ICC)
GLUT-1	Mouse	Thermo Scientific	SPM498	1:50 (ICC&FC)
Occludin	Mouse	Life Technologies	OC-3F10	1:100 (ICC) 1:50 (FC)
Claudin-5	Mouse	Life Technologies	4C3C2	1:100 (ICC) 1:50 (FC)
VE-Cadherin	Mouse	Santa Cruz Biotechnologies	F8	1:25 (ICC) 1:500 (FC)
E-cadherin	Goat	R&D Systems	AF648	1:100 (ICC&FC)
P-glycoprotein	Mouse	Life Technologies	F4	1:25 (ICC) 1:50 (FC)
Breast cancer resistance protein (BCRP)	Mouse	Millipore	5D3	1:25 (ICC) 1:50 (FC)
Multidrug resistance protein 1 (MRP1)	Mouse	Millipore	QCRL-1	1:100 (ICC) 1:50 (FC)
Glial fibrillary acidic protein (GFAP)	Rabbit	Dako	Z0334	1:500 (ICC)
$\beta$ III tubulin	Rabbit	Sigma	T2200	1:1000 (ICC)
Nestin	Mouse	Millipore	10C2	1:500 (ICC)
$\alpha$ smooth muscle actin (SMA)	Mouse	American Research Products	1A4	1:100 (ICC)
Platelet-derived growth factor $\beta$ (PDGFR $\beta$ )	Rabbit	Cell Signaling	28E1	1:100 (ICC)

Lippmann, E. S.; Al-Ahmad, A.; Azarin, S. M.; Palecek, S. P.; Shusta, E. V. A retinoic acid-enhanced, multicellular human blood- brain barrier model derived from stem cell sources. *Sci. Rep.* 2014, 4, 4160.



The VE-Cadherin antibody listed above is no longer appropriate.

#### MATERIALS

NAME ▾	CATALOG # ▾	VENDOR ▾
B-27 Supplement	17504044	Gibco - Thermo Fischer
Essential 8™ Medium	A1517001	Gibco, ThermoFisher
Insulin solution human	I9278	Sigma Aldrich
Recombinant Human FGF-basic (154 a.a.)	100-18B	peprotech
Human Endothelial-SFM	11111044	Thermo Fisher
Versene Solution	15040066	Thermo Fisher
Essential 6™ Medium	A1516401	Thermo Fisher Scientific
Human Holo-Transferrin Protein CF	2914-HT-001G	R&D Systems
Retinoic acid	R2625-50MG	Sigma Aldrich

#### MATERIALS TEXT

***If desired, E8 and E6 may be purchased commercially rather than prepared in-house.***

***If purchasing E8 and E6 commercially, human holo-transferrin and human insulin solution are not needed.***

Plasticware:



FISHER

Corning Tissue Culture Plates (6- or 12-well, 3513 or 3516)  
500 ml filter-top bottles (S2GPT05RE)

#### SAFETY WARNINGS

Please see SDS (Safety Data Sheet) for hazards and safety warnings.

#### BEFORE STARTING

#### REAGENT/MEDIUM PREPARATION:

##### **E4 (prepared according to E4 large scale basal media production protocol)**

Large batch previously prepared and stored at  $-80^{\circ}\text{C}$ . If the preparation is not urgent, remove a bottle of E4 from the  $-80^{\circ}\text{C}$  and place it in the fridge overnight, which will allow the bottle slowly thaw. However, if preparation is desired for the same day, remove a bottle of E4 from the  $-80^{\circ}\text{C}$ , place it on the countertop for  $\sim 00:20:00$ , then place it in the  $37^{\circ}\text{C}$  water bath until it thaws completely ( $\sim 01:00:00 - 02:00:00$ ). Make sure to follow these steps precisely, as premature addition of a frozen bottle of E4 to a water bath can result in rupture of the plastic bottle.

##### **Insulin**

Pre-made solution provided by Sigma (catalog #I9278) that needs no additional preparation. Bottles are stored at  $4^{\circ}\text{C}$ .

##### **Transferrin**

Comes as a powder from R&D Systems (Human Holo-Transferrin, CF; catalog #2914-HT-001G). Add  $50\text{ mg}$  of transferrin to  $5\text{ ml}$  of phosphate-buffered saline (PBS), aliquot at  $500\text{ }\mu\text{l}$ /vial, and store at  $-80^{\circ}\text{C}$ . This mixture does not need to be sterile-filtered.

##### **E6 media (prepared according to E6 and E8 media preparation protocol)**

Dispense the thawed bottle of E4 into a bottletop filter attached to a 500 ml glass bottle. Add  $100\text{ }\mu\text{l}$  of insulin solution and  $500\text{ }\mu\text{l}$  of transferrin solution. Vacuum filter and store at  $4^{\circ}\text{C}$ . E6 media is stable indefinitely.

##### **B27 Supplement**

Thaw 10 ml bottle and mix thoroughly. Aliquot into sterile microcentrifuge tubes at  $280\text{ }\mu\text{l}$ /tube and store at  $-20^{\circ}\text{C}$ . Upon thawing, unused portions of an aliquot may be stored at  $4^{\circ}\text{C}$  for up to 1 week for further media preparation.

##### **bFGF, 100 $\mu\text{g}/\text{ml}$ (prepared according to E8 media protocol)**

Thaw a  $500\text{ }\mu\text{l}$  aliquot of bFGF and dilute 1:5000 in EC medium for a final concentration of  $20\text{ ng}/\text{ml}$  as described below. Divide the remaining bFGF in  $100\text{ }\mu\text{l}$  aliquots and re-freeze at  $-80^{\circ}\text{C}$ . These remaining aliquots can be thawed and used for EC medium but cannot be refrozen a second time.

##### **Retinoic acid (RA)**

Dilute  $50\text{ mg}$  RA in  $16.6\text{ ml}$  DMSO to create a stock solution of  $10\text{ mM}$  and store  $1\text{ mL}$  aliquots at  $-80^{\circ}\text{C}$ . To prepare working stocks, divide a  $1\text{ mL}$  stock tube into  $50\text{ }\mu\text{l}$  aliquots and store at  $-20^{\circ}\text{C}$ . Dilute working stocks 1:1000 in EC medium for a final concentration of  $10\text{ Micromolar } (\mu\text{M})$ .

##### **EC medium w/ 200X B27 + 20 ng/ml bFGF**

For 50 ml: add  $250\text{ }\mu\text{l}$  of B27 and  $10\text{ }\mu\text{l}$  bFGF to  $50\text{ ml}$  of hESFM.

Good for up to two weeks at  $4^{\circ}\text{C}$ .

## EC medium w/ 200X B27

For 50 ml: add **250 µl** of B27 to **50 ml** of hESFM.

### BBB differentiation (Day 0–4)

1 



*Note: Cells are seeded for differentiation in E8 medium according to the standardized single cell seeding protocol*

On **day 0**, aspirate E8 medium and add **2 ml** of E6 per well.

2 

Change medium every day using **2 ml** of E6 per well.

### BBB expansion (Day 4–6)

3 

At **day 4** of E6 treatment, aspirate and add **2 ml** of EC medium with bFGF (basic fibroblast growth factor) and **10 Micromolar (µM)** RA to each well.



Medium is NOT changed during expansion phase.

4 BBB subculturing:

On **day 6**, subculture BBB onto plates and Transwell filters according to the following protocol:



Protocol for Subculture of Differentiated Blood-Brain Barrier Endothelial Cells onto Plates and Filters

by Emma Neal,

Department of Chemical Engineering, Vanderbilt University, Nashville, TN, USA

PREVIEW

RUN

4.1 Please select between subculturing onto **plates** or **filters**.

step case

#### Plates


Subculturing onto Plates using Accutase.

#### 4.2

Coat plates with ECM plate solution for at least  **01:00:00** at  **37 °C** . Volume depends on plate type (see Table):

Plate type for subculture phase	Volume of ECM solution for coating	Working volume of EC media for cell culture
6-well	800 µl	2 ml
12-well	250 µl	1 ml
24-well	200 µl	500 µl
48-well	100 µl	400 µl
96-well	50 µl	200 µl



If desired, plates may be coated  **Overnight** . If coating overnight, add necessary volume of ECM and an equal volume of ddH<sub>2</sub>O to each well to prevent excessive evaporation. If using glass plates, overnight incubation is needed to achieve adequate protein adsorption.

#### 4.3 Aspirate plates and allow to dry in sterile hood (place the plate in the back of the hood and leave the lid slightly ajar).



Plates only need to dry for  **00:05:00** (can be aspirated during accutase incubation). Do not over dry!


#### 4.4

Retrieve cells from incubator and transfer equal volume of spent media to 15 ml conical corresponding to the number of wells being accutased.



For example, if accutasing 4 wells, save 4 ml of spent media and discard the rest.

#### 4.5

Wash each well once with  **2 ml** PBS.

#### 4.6

Add  **1 ml** accutase (warmed to  **Room temperature** ) to each well.

4.7 

Incubate at **37 °C** , length of time depends on cell treatment:

step case

#### If cells have not been treated with RA

4.8 

If cells have not been treated with RA, incubate at **37 °C** for **00:20:00** , or until cells are dissociated from plate (whichever comes first).

4.9 Using p1000, collect cells, and spray gently over surface 2–3x to dislodge any remaining cells. Triturate briefly to break up cell clumps.

4.10 

Add cells to 15 ml conical containing spent media.

4.11 

Spin down cells at **1000 rpm 00:04:00** .

4.12 

Aspirate media, and resuspend cells in appropriate volume of EC media. For 6- and 12-wells, cells are seeded based on a split ratio:

- 1 well of a 6-well plate is split to 1 well of a 6-well plate [1:1]
- 1 well of a 6-well plate is split to 3 wells of a 12-well plate [1:3]
- For smaller plates (24-, 48-, or 96-wells), seed 1 million cells/cm<sup>2</sup>.
- Multiply split ratio by the working volume found in the table to arrive at total volume of EC media in which to resuspend cells.

4.13 Thoroughly triturate 3 – 4 times to yield single cell suspension.

4.14 

Add appropriate volume of cells to each well.

4.15 

Place plate in incubator, shaking plate back and forth to distribute cells evenly (do not swirl).

4.16 

24 hours later (i.e., day 7), aspirate spent media and add appropriate volume of EC medium (**without bFGF or RA**).



This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited