

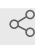


Jul 20, 2022

🌐 Mass-Spectrometry analysis of ATP13A2 samples

Sue Sim¹, eunyong_park¹¹University of California, Berkeley

1 Works for me

 Sharedx.doi.org/10.17504/protocols.io.n92ldzdyxv5b/v1

Sue Sim

ABSTRACT

Preparing mass-spectrometry samples to analyze polyamine content in purified ATP13A2 samples

DOI

dx.doi.org/10.17504/protocols.io.n92ldzdyxv5b/v1

PROTOCOL CITATION

Sue Sim, eunyong_park 2022. Mass-Spectrometry analysis of ATP13A2 samples. **protocols.io**
<https://dx.doi.org/10.17504/protocols.io.n92ldzdyxv5b/v1>



LICENSE

————— This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

CREATED

Jul 17, 2022

LAST MODIFIED

Jul 20, 2022

PROTOCOL INTEGER ID

66883

MATERIALS TEXT

MS Buffer

25 mM Tris pH 7.5

100 mM NaCl

1 mM EDTA

1 mM DTT

0.03% DDM/ 0.006% CHS

- 1 Pool purified ATP13A2 sample after SEC (sample should be in MS buffer after SEC)
- 2 Concentrate the protein to ~1.2 mg/mL using an Amicon Ultrafilter (cut-off 100kDa) at **4 °C**
- 3 Prepare standards of spermidine (Sigma-Aldrich) and spermine (Sigma-Aldrich) at a concentration of 10 μ M in MS buffer
- 4 Flash-freeze samples in liquid nitrogen and store at -80 °C until use
- 5 Dilute each sample 1:1 into acetonitrile with 1% formic acid (volume/volume)
- 6 Analyze by nanoelectrospray ionization (nanoESI) high-resolution mass spectrometry, using an LTQ-Orbitrap-XL mass spectrometer (Thermo Fisher Scientific)
 - 6.1 Performed by QB3/Chemistry Mass Spectrometry Facility (University of California, Berkeley)
 - 6.2 Mass spectrometer was equipped with a nanoESI source and operated in the positive ion mode
 - 6.3 Mass spectra were acquired at a mass resolution setting of 100,000, as measured at mass-to-charge ratio (m/z) = 400, full width at half-maximum

peak height

- 6.4 Mass spectrometry data acquisition and processing were performed using Xcalibur software (version 2.0.7, Thermo)