





Apr 21, 2022

Engineering brain assembloids to interrogate human neural circuits

Yuki Miura^{1,2}, Min-Yin Li^{1,2}, Omer Revah^{1,2}, Se-Jin Yoon^{1,2}, Genta Narazaki¹, Sergiu Pasca^{1,2}

¹Department of Psychiatry and Behavioral Sciences, Stanford University, Stanford, CA, USA;

²Stanford Brain Organogenesis, Wu Tsai Neurosciences Institute, Stanford University, Stanford, CA, USA.

1 «

dx.doi.org/10.17504/protocols.io.36wgq4xxkvk5/v1

Neurodegeneration Method Development Community



DISCLAIMER - FOR INFORMATIONAL PURPOSES ONLY; USE AT YOUR OWN RISK

The protocol content here is for informational purposes only and does not constitute legal, medical, clinical, or safety advice, or otherwise; content added to protocols.io is not peer reviewed and may not have undergone a formal approval of any kind. Information presented in this protocol should not substitute for independent professional judgment, advice, diagnosis, or treatment. Any action you take or refrain from taking using or relying upon the information presented here is strictly at your own risk. You agree that neither the Company nor any of the authors, contributors, administrators, or anyone else associated with protocols.io, can be held responsible for your use of the information contained in or linked to this protocol or any of our Sites/Apps and Services.

The development of neural circuits involves wiring of neurons locally following their generation and migration, as well as establishing long-distance connections between brain regions. Studying these developmental processes in the human nervous system remains difficult because of limited access to the tissue that can be maintained as functional over time in vitro. We have previously developed a method to convert human pluripotent stem cells into regionalized neural organoids that can be fused and integrated to form assembloids and study neuronal migration. In this protocol, we describe approaches to model long-range neuronal connectivity in human brain assembloids. We present how to generate 3D spheroids resembling specific domains of the nervous system and then how to integrate them physically to allow axonal projections and synaptic assembly. In addition, we describe a series of assays including viral labeling and retrograde tracing, 3D live imaging of axon projection and optogenetics combined with calcium imaging and electrophysiological recordings to probe and manipulate the circuits in assembloids. We anticipate that these approaches will be useful in deciphering human-specific aspects of neural circuit assembly and in modeling neurodevelopmental disorders with patient-derived cells.

DOI

dx.doi.org/10.17504/protocols.io.36wgq4xxkvk5/v1

Yuki Miura, Min-Yin Li, Omer Revah, Se-Jin Yoon, Genta Narazaki, Sergiu Pasca 2022. Engineering brain assembloids to interrogate human neural circuits. **protocols.io** https://dx.doi.org/10.17504/protocols.io.36wgq4xxkvk5/v1

protocol	
https://www.nature.com/articles/s41596-021-00632-z	

human neural circuits, brain assembloids, viral labeling, retrograde tracing, 3D live imaging of axon projection, optogenetics

_____ protocol ,

Nov 01, 2021

Apr 21, 2022

Nov 01, 2021 Renuka Suresh

Sergiu Pasca

Nov 11, 2021



1

Citation: Yuki Miura, Min-Yin Li, Omer Revah, Se-Jin Yoon, Genta Narazaki, Sergiu Pasca Engineering brain assembloids to interrogate human neural circuits https://dx.doi.org/10.17504/protocols.io.36wqq4xxkvk5/v1

Troubleshooting

Α	В	С	D
Step	Problem	Possible reason	Solution
Throughout	Cultures were contaminated	Inadequate sterile technique	Improve sterile technique. Wear laboratory coat. Use 70% ethanol to spray your hands and any equipment or reagents that go into the safety cell culture cabinet. Autoclave all non-sterile equipment. Change HEPA filters in ncubators regularly. Autoclean incubators regularly. Identify sources of contamination in the incubator or cell culture room. Clean the water bath in the incubator and in the cell culture room on a weekly basis.
29	Individual spheroids do not form	Overfilled individual wells of the AggreWell plate	Count cell number carefully and add the required amount of cells to each well. Centrifuge the plate with optimal speed
31	Spheroids do not form	hPS cells are	Always check hPS cell quality
		hPS cell colonies are too small	Allow hPS cell colonies to grow larger. Do not use until the optimal day for splitting (80–85% confluency)
		Over-crowding of hPS cell Colonies	Start differentiation by using hPS cell colonies with the optimal confluency (80–85%)
		Mycoplasma or other contamination	Perform mycoplasma testing regularly (ideally, once per month) and discard the cells and use new hPS cell stock if mycoplasma positive
38	Difficulty in inducing forebrain fate (low FOXG1 expression)	hPS cells are challenging to differentiate into forebrain cells	Consider adding XAV-939 (1.25 μ M) from days 1–6 of differentiation
	Spheroids stick to one another	Too many spheroids cultured per plate	Transfer some spheroids to a new or larger low-attachment plate. If already attached, separate them by gently pipetting up and down by using a 10-or 25-ml pipette (for spheroids up to day 25) or two P20 pipette tips at later stages
	Spheroids are too small	hPS cells are differentiating hPS cell colonies are	See solutions for Step 31
		too small Mycoplasma or other	See solutions for Step 31 See solutions for Step 31
	Spheroids are disintegrated	contamination Excessive pipetting	Be gentle when handling spheroids at all times

38 and 39	Growth of	Abnormal genome	Check the genome integrity of
	spheroids is fast	integrity in the	the hPS cell line by using SNP
		hPS cell line	array, and confirm that there is no
			de novo
			large duplication or deletion
39 and 40	Medium turns	Too many spheroids	In the initial stages of
	acidic overnight	cultured per	differentiation, 12 ml of medium
		plate	per 10-cm
			plate is sufficient. As spheroids
			grow, use 12–15 ml of medium
			per plate.
			Maintaining <30 spheroids per
40.4	0	The state of the s	10-cm plate is recommended
40A	Sparse or no viral	The viral titer was	Make sure that viral titers are
	labeling	too low	sufficiently high
		Insufficient time	Wait ≥1-2 weeks for viral gene
		between	expression
		infection and	
		fluorescent protein	
		expression	
		Incorrect serotype	Use optimal AAV serotypes
40B	Unable to assemble		Try to fuse smaller or
	hCS and	stages of	early-stage spheroids (days 47-
	hStrS	differentiation	80). After day 80, the success rat
			is lower
40F	Unable to find	The signal of	Check expression of GCaMP6 by
	GCaMP6+ cells	GCaMP6 is weak	using an epifluorescence
	inassembloids		microscope (e.g.,
	under a confocal		EVOS or REVOLVE). Increase
	microscope		laser intensity
		Imaging only a few	Increase thickness of slices by
		z-slices	changing the pinhole size and/or
			image more slices
40G	Cells are hard to	Cell bodies are not	Identify cells by using
	find under	always easy	fluorescence imaging
	bright field	to identify by using	
		bright field	
	Cells are blown	The pressure is too	Apply minimal positive pressure
	away by pressure	high	
		Pipette tips are too large	Use a pipette with smaller tips
	It is dif ficult to	Pinette tins are too	Use a ninette with smaller tins
	It is dif ficult to	Pipette tips are too	Use a pipette with smaller tips
	form a tight	Pipette tips are too large	Use a pipette with smaller tips
	form a tight seal. Cells cannot		Use a pipette with smaller tips
	form a tight seal. Cells cannot form a fully tight		Use a pipette with smaller tips
	form a tight seal. Cells cannot form a fully tight seal even though		Use a pipette with smaller tips
	form a tight seal. Cells cannot form a fully tight seal even though cells are partially		Use a pipette with smaller tips
	form a tight seal. Cells cannot form a fully tight seal even though cells are partially sucked into the		Use a pipette with smaller tips
	form a tight seal. Cells cannot form a fully tight seal even though cells are partially sucked into the pipette	large	
	form a tight seal. Cells cannot form a fully tight seal even though cells are partially sucked into the pipette Photobleaching of	large The LED power is too	Lower the LED power and
	form a tight seal. Cells cannot form a fully tight seal even though cells are partially sucked into the pipette Photobleaching of fluorescence when	large	Lower the LED power and increase
	form a tight seal. Cells cannot form a fully tight seal even though cells are partially sucked into the pipette Photobleaching of	large The LED power is too	Lower the LED power and increase gain and/or exposure time in the
	form a tight seal. Cells cannot form a fully tight seal even though cells are partially sucked into the pipette Photobleaching of fluorescence when checking cells	The LED power is too high	Lower the LED power and increase gain and/or exposure time in the camera software setting
	form a tight seal. Cells cannot form a fully tight seal even though cells are partially sucked into the pipette Photobleaching of fluorescence when checking cells Cells are hard to	The LED power is too high For some cells, it is	Lower the LED power and increase gain and/or exposure time in the camera software setting Try 'zap' in the 700B panel for
	form a tight seal. Cells cannot form a fully tight seal even though cells are partially sucked into the pipette Photobleaching of fluorescence when checking cells	The LED power is too high For some cells, it is hard to	Lower the LED power and increase gain and/or exposure time in the camera software setting Try 'zap' in the 700B panel for multiple trials; start from short
	form a tight seal. Cells cannot form a fully tight seal even though cells are partially sucked into the pipette Photobleaching of fluorescence when checking cells Cells are hard to	The LED power is too high For some cells, it is hard to break in with just	Lower the LED power and increase gain and/or exposure time in the camera software setting Try 'zap' in the 700B panel for multiple trials; start from short durations such as 50 µs and there
	form a tight seal. Cells cannot form a fully tight seal even though cells are partially sucked into the pipette Photobleaching of fluorescence when checking cells Cells are hard to	The LED power is too high For some cells, it is hard to	Lower the LED power and increase gain and/or exposure time in the camera software setting Try 'zap' in the 700B panel for multiple trials; start from short durations such as 50 µs and the increase
	form a tight seal. Cells cannot form a fully tight seal even though cells are partially sucked into the pipette Photobleaching of fluorescence when checking cells Cells are hard to	The LED power is too high For some cells, it is hard to break in with just	Lower the LED power and increase gain and/or exposure time in the camera software setting Try 'zap' in the 700B panel for multiple trials; start from short durations such as 50 µs and thei increase durations step by step; try strong
	form a tight seal. Cells cannot form a fully tight seal even though cells are partially sucked into the pipette Photobleaching of fluorescence when checking cells Cells are hard to	The LED power is too high For some cells, it is hard to break in with just	Lower the LED power and increase gain and/or exposure time in the camera software setting Try 'zap' in the 700B panel for multiple trials; start from short durations such as 50 µs and thei increase durations step by step; try strong suctions immediately after zaps.
	form a tight seal. Cells cannot form a fully tight seal even though cells are partially sucked into the pipette Photobleaching of fluorescence when checking cells Cells are hard to	The LED power is too high For some cells, it is hard to break in with just	Lower the LED power and increase gain and/or exposure time in the camera software setting Try 'zap' in the 700B panel for multiple trials; start from short durations such as 50 µs and ther increase durations step by step; try strong suctions immediately after zaps. Use
	form a tight seal. Cells cannot form a fully tight seal even though cells are partially sucked into the pipette Photobleaching of fluorescence when checking cells Cells are hard to break in	The LED power is too high For some cells, it is hard to break in with just suction	Lower the LED power and increase gain and/or exposure time in the camera software setting Try 'zap' in the 700B panel for multiple trials; start from short durations such as 50 µs and ther increase durations step by step; try strong suctions immediately after zaps. Use pipettes with slightly bigger tips
	form a tight seal. Cells cannot form a fully tight seal even though cells are partially sucked into the pipette Photobleaching of fluorescence when checking cells Cells are hard to break in	The LED power is too high For some cells, it is hard to break in with just suction Pipette tips are too	Lower the LED power and increase gain and/or exposure time in the camera software setting Try 'zap' in the 700B panel for multiple trials; start from short durations such as 50 µs and ther increase durations step by step; try strong suctions immediately after zaps. Use pipettes with slightly bigger tips Use pipettes with slightly bigger
	form a tight seal. Cells cannot form a fully tight seal even though cells are partially sucked into the pipette Photobleaching of fluorescence when checking cells Cells are hard to break in	The LED power is too high For some cells, it is hard to break in with just suction	Lower the LED power and increase gain and/or exposure time in the camera software setting Try 'zap' in the 700B panel for multiple trials; start from short durations such as 50 µs and ther increase durations step by step; try strong suctions immediately after zaps. Use pipettes with slightly bigger tips

Osmolarity of the	Check osmolarity of both the
external	external solution and the internal
solution and/or the	solution
internal solution is	
not appropriate	
The relative position	Slightly adjust the position of
of the	pipette tips
pipette tip and the	
cell changes in	
whole-cell mode	

Media

CRITICAL: When ready for use, distribute aliquots of the required volume of medium into **50 mL** sterile tubes and pre-warm in a water bath at **37°C** or at **Room temperature**. Avoid multiple cycles of refrigeration and warming. Add growth factors and small molecules before use and only after pre-warming aliquots of media.

```
Biological materials
hPS cells
REAGENT
⊠ DMEM/F-12, HEPES Thermo Fisher
Scientific Catalog #11330032

    Essential 8™ Medium Thermo

Fisher Catalog #A1517001

    Essential 6™ Medium Thermo Fisher

Scientific Catalog #A1516401
⊠ Neurobasal™-A Medium Thermo
Fisher Catalog #10888022
⊠B-27<sup>™</sup> Supplement (50X), minus vitamin A Thermo
Fisher Catalog #12587010
CRITICAL The B-27 supplement must be without vitamin A.
⊠B-27<sup>™</sup> Plus Supplement (50X) Thermo
Fisher Catalog #A3582801
 ⊠ GlutaMAX™ Supplement Thermo Fisher
Scientific Catalog #35050061

⊠ Penicillin-Streptomycin (10,000 U/mL) Thermo Fisher

Scientific Catalog #15140122
CRITICAL: Differentiations can also be performed without penicillin-streptomycin, but the risk of contamination, especially
in long-term cultures, will be higher.
              ⊠ Y-
●ROCK inhibitor 27632 Selleckchem Catalog #S1049
               ⊠ Dorsomorphin Sigma
●SMAD inhibitor Aldrich Catalog #P5499-5MG
⊗SB
431542 Tocris Catalog #1614
⊠XAV
939 Tocris Catalog #3748


    □ Recombinant Human EGF Protein CF R&D

Systems Catalog #236-EG
Recombinant Human FGF basic/FGF2/bFGF (146 aa) Protein R&D
Systems Catalog #233-FB
```



4

Citation: Yuki Miura, Min-Yin Li, Omer Revah, Se-Jin Yoon, Genta Narazaki, Sergiu Pasca Engineering brain assembloids to interrogate human neural circuits https://dx.doi.org/10.17504/protocols.io.36wgg4xxkvk5/v1

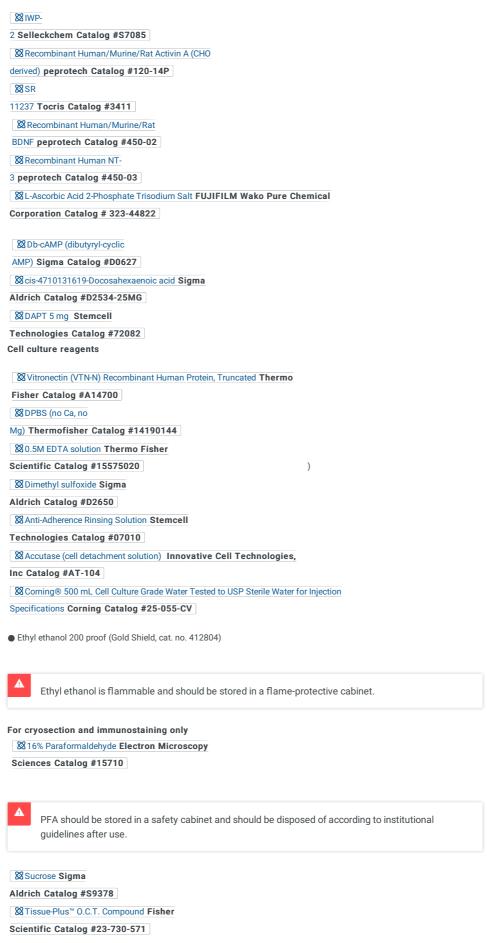

👸 protocols.io

5

Citation: Yuki Miura, Min-Yin Li, Omer Revah, Se-Jin Yoon, Genta Narazaki, Sergiu Pasca Engineering brain assembloids to interrogate human neural circuits https://dx.doi.org/10.17504/protocols.io.36wqq4xxkvk5/v1

```
XTriton X-100 Sigma
Aldrich Catalog #T8787-50ML
 ⊠ Donkey Serum Emd
Millipore Catalog #S30-100ML

    Recombinant Anti-DARPP32 antibody [EP720Y]

(ab40801) Abcam Catalog #ab40801
X Anti-GAD67 Antibody clone 1G10.2 Merck
Millipore Catalog #MAB5406
X Anti-Ctip2 antibody [25B6]
(ab18465) Abcam Catalog #ab18465
 X Anti-NeuN antibody [1B7] - Neuronal
Marker Abcam Catalog #ab104224
antibody Genetex Catalog #GTX13970
⊠ mCherry
antibody Genetex Catalog #GTX128508
MAP 2 antibody - 188 004 Synaptic
Systems Catalog #188 004
 XAnti-SATB2 antibody [SATBA4B10] - C-terminal
(ab51502) Abcam Catalog #ab51502
 ⊠ Donkey anti-Rabbit IgG (H L) Highly Cross-Adsorbed Secondary Antibody Alexa Fluor 488 Thermo Fisher
Scientific Catalog #A-21206
■ Donkey anti-Mouse IgG (H L) Highly Cross-Adsorbed Secondary Antibody Alexa Fluor 568 Thermo Fisher
Scientific Catalog #A10037

    ⊠ Alexa Fluor® 647 AffiniPure Donkey Anti-Rat IgG (H L) Jackson

Immunoresearch Catalog #712-605-153
 Immunoresearch Catalog #703-545-155
 ⊗Donkey anti-Rabbit IgG (H L) Highly Cross-Adsorbed Secondary Antibody Alexa Fluor 568 Thermo Fisher
Scientific Catalog #A10042
🛮 🛠 Alexa Fluor® 647 AffiniPure Donkey Anti-Guinea Pig IgG (H L) (min X Bov Ck Gt Sy Hms Hrs Hu Ms Jackson
Immunoresearch Catalog #706-605-148
 🛭 AF647 Donkey anti-mouse
 lgG Thermofisher Catalog #A-31571

    ⊠ Hoechst 33258 Thermo Fisher

Scientific Catalog #H3569
⊗ Aqua-
Poly/Mount Polysciences Catalog #18606-100
For clearing and 3D staining of assembloids only
Chemicals Catalog #T3740

    ▼Tissue-Clearing Reagent CUBIC-R (M) [for Animals] TCI

Chemicals Catalog #T3741
For viral labeling only (see Table 1 for full details)

    AAV-DJ-hSvn1::eYFP (Stanford University Neuroscience Gene Vector and Virus Core, cat. no. GVVC-AAV-16)

 AAV-DJ-hSyn1::mCherry (Stanford University Neuroscience Gene Vector and Virus Core, cat. no. GVVC-AAV-17)

    Rabies-ΔG-Cre-eGFP (Salk Institute Viral Vector Core)

    ⊗pAAV-EF1a-CVS-G-WPRE-
 ■ pGHpA addgene Catalog #67528
 AAV-DJ-EF1-DIO-mCherry (Stanford University Neuroscience Gene Vector and Virus Core, cat. no. GVVC-AAV-14)
```

protocols.io

6

Citation: Yuki Miura, Min-Yin Li, Omer Revah, Se-Jin Yoon, Genta Narazaki, Sergiu Pasca Engineering brain assembloids to interrogate human neural circuits https://dx.doi.org/10.17504/protocols.io.36wgg4xxkvk5/v1

- AAV1-Syn::ChrimsonR-tdT (Addgene, cat. no. 59171-AAV1)
- AAV-DJ-Ple94 (GPR88)-iCre (Addgene, cat. no. 49125)
- AAV-DJ-EF1a-DIO-GCaMP6s (Stanford University Neuroscience Gene Vector and Virus Core, cat. no. GVVC-AAV-91)


Table 1:

Α	В	С	D	E	F
	Virus and	Gene of	Target	Stock	Working
	serotype	interest	cell type	concentration	concentration
AAV-DJ-	AAV-DJ	eYFP	Neuron	4.56 × 10^12	1.14 × 10^10
Syn1::eYFP				Vg/ml	Vg/ml
AAV-DJ-	AAV-DJ	mCherry	Neuron	1.50 × 10^13	3.75 × 10^10
Syn1::mCherry				Vg/ml	Vg/ml
AAV-DJ-EF1-DIO-	AAV-DJ	mCherry	Any	1.60 × 10^13	4.00 × 10^10
mCherry				Vg/ml	Vg/ml
AAV-DJ-EF1a::G	AAV-DJ	Rabies	Any	3.00 × 10^12-	7.50 × 10^9-
		glycoprotein		2.08 × 10^13	5.20 × 10^10
				Vg/ml	Vg/ml
AAV1-hSyn-	AAV1	ChrimsonR-	Neuron	2.10 × 10^13	5.25 × 10^10
ChrimsonR-tdT		tdTomato		Vg/ml	Vg/ml
AAV-DJ-Ple94-	AAV-DJ	iCre	Neuron	1.09 × 10^12	2.73 × 10^9
iCre				Vg/ml	Vg/ml
AAV-DJ-EF1a-	AAV-DJ	GCaMP6s	Any	1.60 × 10^13	4.00 × 10^10
DIO-GCaMP6s				Vg/ml	Vg/ml
ΔG-Rabies-Cre-	G-deletedrabies	Cre-eGFP	Any	6.52 × 10^8	1.63 × 10^5
eGFP	virus			TU/ml	TU/ml

For whole-cell patch-clamp recording only **⊠LOW MELT AGAROSE IBI** Scientic Catalog #IB70051 Sodium chloride Sigma Aldrich Catalog #S7653 **⊠KCI Millipore** Sigma Catalog #P9333 Sodium phosphate monobasic monohydrate Sigma Aldrich Catalog #S9638

 Magnesium sulfate heptahydrate (MgSO4) Millipore Sigma Catalog #63138 **⊠**CaCl2 solution **Sigma** Aldrich Catalog #21115 Sodium Bicarbonate Sigma Aldrich Catalog #S6297 **⊠**D-()-glucose **Sigma** Aldrich Catalog #G7528 Aldrich Catalog #P1847 Aldrich Catalog #A9187 Sodium GTP Sigma Aldrich Catalog #G8877 (MilliporeSigma, cat. no. G8877) **XHEPES Sigma** Aldrich Catalog #H3375 ⊠ Ethylene glycol-bis(2-aminoethylether)-NNN'N'-tetraacetic acid Sigma Aldrich Catalog #03777 Scientific Catalog #A13100



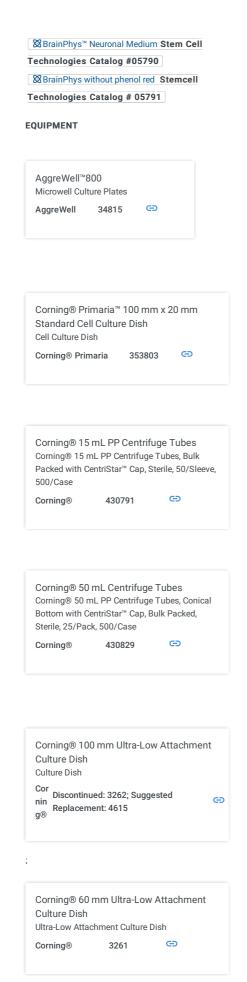



8

Costar® 24-well Clear Flat Bottom Ultra-Low Attachment Multiple Well Plates 24-well Clear Flat Bottom Ultra-Low Attachment Multiple Well Plates

Corning® 3473 G

Fisherbrand™ Premium Microcentrifuge Tubes: 1.5mL Microcentrifuge Tubes

Fisherbrand™ 05408129

Autoclave tubes before use.

5 mL Serological Pipet, Polystyrene, 0.1 Increments, Individually Packed, Sterile, 50/Bag, 200/Case Serological Pipet, Polystyrene

Falcon® 356543 👄

10 mL Serological Pipet, Polystyrene, 0.1 Increments
Serological Pipet, Polystyrene
Falcon® 356551



Falcon® 25 mL Serological Pipet, Polystyrene, Space Saver, 0.25 Increments, Individually Packed, Sterile, 50/Bag, 200/Case Serological Pipet, Polystyrene G Falcon® 356525

Falcon® 50 mL Serological Pipet, Polystyrene, 1.0 Increments, Individually Wrapped, Sterile, 25/Pack, 100/Case 50 mL Serological Pipet, Polystyrene, 1.0 Increments ⊕

356550 Falcon®

Low Retention Aerosol Barrier Pipette Tip MBP® ART® SoftFit-L™ 1,000 μL Sterile TIP, PIPETTE ART SOFTFIT-L HINGED RACK

MBP® ART® SoftFit-L™ 2779-HR 👄

Aerosol Barrier Pipette Tip MBP® ART® SoftFit-L™ 200 µL Sterile TIP, PIPETTE W/FILTER HINGED COVER 200UL-200L (5PK/CS)

MBP® ART® SoftFit-L™ 2769-HR 👄

Aerosol Barrier Pipette Tip MBP® ART® SoftFit-L™ 20 µL Sterile TIP, PIPETTE W/FILTER HINGED COVER 200UL-20L (5PK/CS)

MBP® ART® SoftFit-L™ 2749-HR 👄

Sterilized scissors

VWR® Air Jacketed CO2 Incubators Air Jacketed CO2 Incubator ⊕ **VWR**® 10810-902

- Biological safety cabinet (Labconco, cat no. 97000-862)
- Water bath



10

Citation: Yuki Miura, Min-Yin Li, Omer Revah, Se-Jin Yoon, Genta Narazaki, Sergiu Pasca Engineering brain assembloids to interrogate human neural circuits https://dx.doi.org/10.17504/protocols.io.36wgg4xxkvk5/v1

- Culture microscope (Olympus, cat. no. CKX41)
- EVOS FL cell imaging system (Life Technologies, model AMF4300)
- REVOLVE 4 upright, inverted, bright field, fluorescent microscope (Echo Laboratories, cat. no. FJSD2001)
- Confocal microscope (Leica, model TCS SP8)
- Stainless steel surgical blade, sterile, no. 10 (Medicon, cat. no. 01.22.10)

Falcon® 40 µm Cell Strainer, Blue, Sterile, Individually Packaged, 50/Case Cell Strainer
Falcon® 352340 GD

Hemocytometer

Nylon Syringe Filter, 0.22µm, 30mm, Bulk Packed, Non-sterile Syringe Filter Nylon 229775 🖘

• Centrifuge 5810 R (Eppendorf, cat. no. 022628102)

For viral labeling and assembloid formation only

- Sterile 1.5-ml microcentrifuge
- Scissors and/or razor blades


For cryosection and immunostaining only

VWR® Superfrost® Plus Micro Slide,
Premium
Superfrost slides

VWR® Superfrost® 48311-703

Humidified chamber



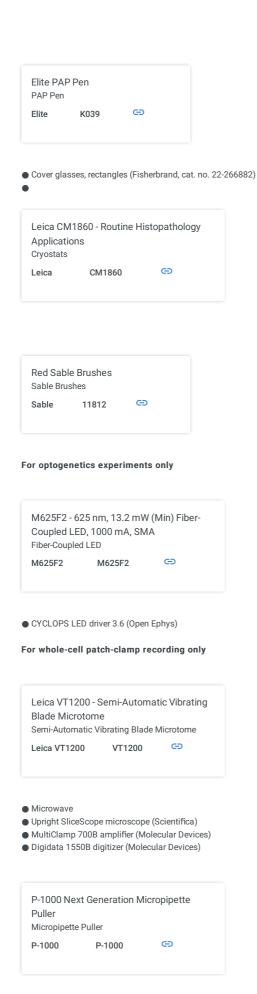



• Thick/standard wall borosilicate glass with filament (Sutter Instrument, cat. no. BF150-86-10)



• Infinity2 CCD camera (Teledyne Lumenera)

SOFTWARE

- Fiji (ImageJ, version 2.1.0, National Institutes of Health)
- MATLAB (version R2018a, 9.4.0, MathWorks)
- pCLAMP software (version 10.6, Molecular Devices)
- Mini Analysis software (version 6.0.3, Synaptosoft)
- Infinity Capture software (Teledyne Lumenera) compatible with the CCD camera

REAGENT SETUP

Growth factors and chemical stock solutions

BiologyResuspend **□25 mg** of dorsomorphin in **□12.52 mL** of sterile DMSO to obtain a [M]**5 millimolar (mM)** stock solution and use at a final concentration of [M]2.5 micromolar (µM) (1:2,000). Resuspend 110 mg of SB-431542 in **2.60 mL** of 100% ethanol to obtain a [M]10 millimolar (mM) stock solution and use at a final concentration of 10 µM (1:1,000). Resuspend EGF, FGF2, BDNF and NT-3 growth factors in sterile 0.1% (wt/vol) BSA/PBS to 20 µg/ml stock concentration and use at **20 ng/ml** (1:1,000). Resuspend **20 ng** of IWP-2 in **28.58 mL** of sterile DMSO to obtain a [M]2.5 millimolar (mM) stock solution; use at a final concentration of [M]2.5 micromolar (µM) (1:1,000). Resuspend $\blacksquare 100~\mu g$ of activin A in $\blacksquare 2~mL$ of cell culture-grade water to obtain a $\blacksquare 50~\mu g/ml$ stock solution; use at a final concentration of 50 ng/ml (1:1,000). Resuspend 10 mg of SR11237 in 5.26 mL of sterile DMSO to obtain a [m]5 millimolar (mM) stock solution and further dilute the stock solution to obtain a [M]5 micromolar (μM) solution and use at a final concentration of [M]5 nanomolar (nM) (1:1,000). Resuspend □322 mg of AA in □5 mL of cell culture-grade water to obtain a [M]200 millimolar (mM) stock solution; use at a final concentration of [M]200 micromolar (μM) (1:1,000). Resuspend ■250 mg of dibutyryl cAMP sodium salt in ■5.088 mL of cell culture-grade water to obtain a [M]100 millimolar (mM) stock solution and use at a final concentration of [M]50 micromolar (μM) (1:2,000). Resuspend \Box 5 mg of DAPT in \Box 2.31 mL of sterile DMSO to obtain a [M]5 millimolar (mM) stock solution and use at a final concentration of [M]2.5 Millimolar (mM) (1:2,000). Aliquots of stock solutions are stored at either § -20 °C or § -80 °C for 2 months. Once thawed, stocks can be kept at § 4 °C for \leq \bigcirc 336:02:00 .

Essential 8 medium (for use in Steps 1 –30). Prepare the medium as shown in the table below. 314 This medium can be stored for 2 weeks at 4 °C.

Α	В	С
Composition	Volume (500 ml	Final
	total) (ml)	concentration
Essential	490	-
8 medium		
Essential	10	2% (vol/vol)
8 supplement		

Essential 6 medium (for use in Steps 31–38). Prepare the medium as shown in the table below. This medium can be stored for 2 weeks at § 4 °C.

Α	В	С	D
Composition	Volume (500 ml total) (ml)	Final concentration	Comments
Essential 6 medium	495	_	
Penicillin- streptomycin	5	1% (vol/vol)	Optional
Dorsomorphin	-	2.5 µM (1:2,000)	Add just before use
SB-431542	-	10 μM (1:1,000)	Add just before use
XAV-939 (optional)	-	1.25 µM (1:2,000)	Add just before use

Neural differentiation medium (NM) for hCS (for use in Steps 39 and 40, on days 6 and onward of differentiation). Prepare the medium as detailed in the table below. Store stock for $\leq \bigcirc 168:02:00$ at $\S 4 ^{\circ}C$ and add growth factors just before use if required for the culture.

Α	В	С	D
Composition	Volume (500 ml	Final	Comments
	total) (ml)	concentration	
Neurobasal A	480		-
B-27	10	2% (vol/vol)	Needed from day 6 to day
supplement			46
B-27 Plus	10	2% (vol/vol)	Needed from day 46 onward
supplement			
GlutaMax	5	1% (vol/vol)	-
Penicillin-	5	1% (vol/vol)	-
streptomycin			
EGF (20-µg/ml	-	20 ng/ml (1:1,000)	Needed from day 6 to day
411 Q25			22. Add
stock)			just before use
FGF2 (20-	_	20 ng/ml (1:1,000)	Needed from day 6 to day
μg/ml stock)			22. Add
			just before use
BDNF (20-	_	20 ng/ml (1:1,000)	Needed from day 22 to day
μg/ml stock)			46. Add
			just before use
NT-3 (20-	_	20 ng/ml (1:1,000)	Needed from day 22 to day
μg/ml stock)			46. Add
			just before use
AA (200 mM	_	200 μM (1:1,000)	Needed from day 22 to day
stock)			46. Add
			just before use
Dibutyryl	_	50 μM (1:2,000)	Needed from day 22 to day
cAMP sodium			46. Add
salt (100 mM stock)			just before use
DHA (10 mM	_	10 μM (1:1,000)	Nooded from day 22 to day
stock)	_	10 μΝΙ (1.1,000)	Needed from day 22 to day 46. Add
SIUCK)			just before use
			Just before use

NM for hStrSs (for use in Steps 39 and 40, on days 6 and onward of differentiation). Prepare the medium as detailed in the table below. Store the stock for \leq 1 week at 4 °C and add growth factors just before use if required for the culture.



Α	В	С	D
Composition	Volume (~500 ml	Final	Comments
	total) (ml)	concentration	
Neurobasal A	480	-	-
B-27	10	2% (vol/vol)	Needed from day 6 to day
supplement			46
B-27 Plus	10	2% (vol/vol)	Needed from day 46 onward
supplement			
GlutaMax	5	1% (vol/vol)	-
Penicillin-	5	1% (vol/vol)	Needed from day 6 to day
streptomycin			22. Add just before use
IWP-2 (2.5 mM	-	2.5 µM (1:1,000)	Needed from day 6 to day
stock)			22. Add just before use
Activin A (50-	-	50 ng/ml (1:1,000)	Needed from day 6 to day
μg/ml stock)			22. Add just before use
SR11237 (100	-	100 nM (1:1,000)	Needed from day 22 to day
μM stock)			46. Add just before use
BDNF (20-	-	20 ng/ml (1:1,000)	Needed from day 22 to day
μg/ml stock)			46. Add just before use
NT-3 (20-µg/ml	-	20 ng/ml (1:1,000)	Needed from day 22 to day
stock)			46. Add just before use
AA (200 mM	-	200 μM (1:1,000)	Needed from day 22 to day
stock)			46. Add just before use
Dibutyryl cAMP	_	50 μM (1:2,000)	Needed from day 22 to day
sodium salt			46. Add just before use
(100 mM			
stock)			
DHA (10 mM	-	10 μM (1:1,000)	Needed from day 22 to day
stock)			46. Add just before use
DAPT (5 mM	-	2.5 µM (1:2,000)	Needed from days 42 to 46.
stock)			Add just before use

PFA

Prepare 4% (vol/vol) PFA by mixing □10 mL of 16% PFA with □30 mL of DPBS and store at § 4 °C for ≤1 week.

Embedding solution

Blocking solution for immunostaining

Prepare the blocking solution for immunostaining by mixing 10% (vol/vol) normal donkey serum with 0.1% (wt/vol) BSA and 0.3% (vol/vol) Triton-X100 in PBS. Keep this solution at 8.4° C for \leq 1 week.

Antibody solution for immunostaining

Prepare the antibody solution for immunostaining by mixing 2% (vol/vol) normal donkey serum with 0.1% (vol/vol) Triton-X100 in PBS. Keep this solution at δ 4 °C for \leq 1 week.

Artificial cerebrospinal fluid (aCSF)

Prepare concentrated (10×) aCSF stock solution containing [M]1.26 Molarity (M) NaCl, [M]25 millimolar (mM) KCl, [M]12.5 millimolar (mM) NaH2PO4, [M]260 millimolar (mM) NaHCO3 and [M]100 millimolar (mM) D-(+)-glucose. Keep this solution at § 4 °C for ≤1 month. Dilute this solution and add MgSO4 and CaCl2 to prepare aCSF working solution containing [M]126 millimolar (mM) NaCl, [M]2.5 millimolar (mM) KCl, [M]1.25 millimolar (mM) NaH2PO4, [M]1 millimolar (mM) MgSO4, [M]2 millimolar (mM) CaCl2, [M]26 millimolar (mM) NaHCO3 and [M]10 millimolar (mM) D-(+)-glucose before use.

Internal solution for whole-cell recording

Prepare internal solution containing [M]127 millimolar (mM) K-gluconate, [M]8 millimolar (mM) NaCl,



[M]4 millimolar (mM) MgATP, [M]0.3 millimolar (mM) Na2GTP, [M]10 millimolar (mM) HEPES and [M]0.6 millimolar (mM) EGTA, pH 7.2, adjusted with KOH (290 mOsm). Keep this solution at & -20 °C for ≤1 year.

DISCLAIMER - FOR INFORMATIONAL PURPOSES ONLY; USE AT YOUR OWN RISK


The protocol content here is for informational purposes only and does not constitute legal, medical, clinical, or safety advice, or otherwise; content added to protocols.io is not peer reviewed and may not have undergone a formal approval of any kind. Information presented in this protocol should not substitute for independent professional judgment, advice, diagnosis, or treatment. Any action you take or refrain from taking using or relying upon the information presented here is strictly at your own risk. You agree that neither the Company nor any of the authors, contributors, administrators, or anyone else associated with protocols.io, can be held responsible for your use of the information contained in or linked to this protocol or any of our Sites/Apps and Services.

CRITICAL: Before starting hPS cells and the long-term organoid/assembloid cultures, carefully cleaning the incubator is important to avoid potential contamination. Because most culture plates will be in the incubator for several hundred days, the surface of the incubator and culture dishes should be regularly cleaned. The amount of water in the incubator should be monitored frequently to maintain the required humidity. All cell culture-related procedures should be performed in a sterile environment with biological safety cabinets and sterilized materials. The biological safety cabinets should be regularly certified.

Maintenance of hPS cells in feeder-free conditions Timing 4-6 d


1 Culture hPS cells on human recombinant vitronectin with Essential 8 medium. When hPS cell cultures reach ~80% confluency, passage hPS cultures by using [MIO.5 millimolar (mM)] EDTA in PBS.

CRITICAL: To ensure consistency and obtain high-quality differentiation, the hPS cell culture should not grow denser than 90%. Ideally, passages can be performed at 80% confluency

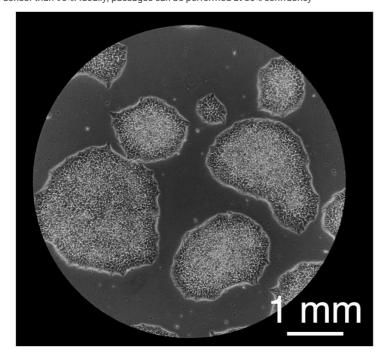


Representative image of hiPS cell culture.



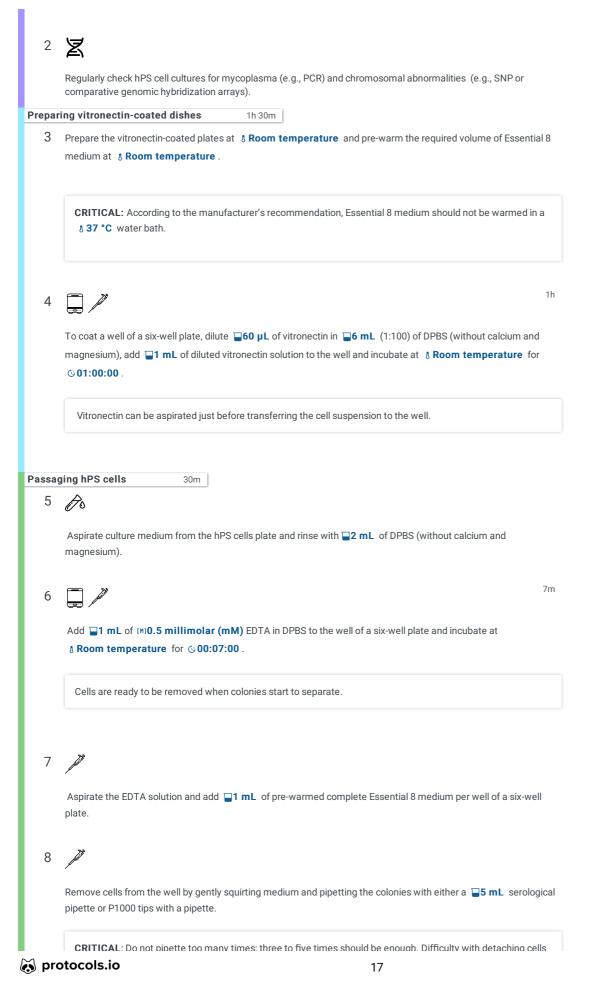

Citation: Yuki Miura, Min-Yin Li, Omer Revah, Se-Jin Yoon, Genta Narazaki, Sergiu Pasca Engineering brain assembloids to interrogate human neural circuits https://dx.doi.org/10.17504/protocols.io.36wqq4xxkvk5/v1

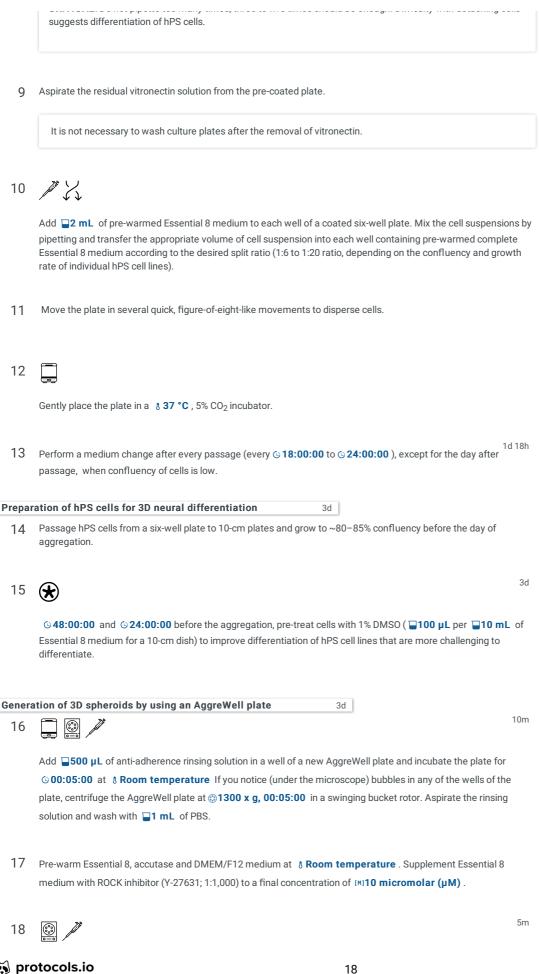

m protocols.io

Add **_0.5** mL of Essential 8 medium supplemented with Y-27632 to each well of an Aggrewell plate. Before starting aggregation, an Aggrewell plate can be centrifuged at **@2000** x g, 00:05:00 in a swinging bucket rotor that is fitted with a plate holder to remove any air bubbles in the microwells.

CRITICAL: Remaining air bubbles in the Aggrewell plate interrupt the aggregation of hPS cells.

19 Set the plate aside in an incubator while preparing a single-cell suspension of hPS cells for aggregation.

20

22

23

Aspirate the maintenance medium from the hPS cell culture plates and rinse cells twice with DPBS (without calcium and magnesium).

Add **4 mL** of accutase per 10-cm plate and incubate at **37 °C** for **© 00:07:00** in a 5% CO₂ incubator until cells detach by gentle shaking from the plate.

Gently pipette the cell suspension two to three times with a serological pipette to ensure that any remaining clumps

are fully dissociated.

Transfer the cells to a **50 mL** conical tube.

CRITICAL: The use of **□50 mL** conical tubes for cell culture (instead of **□15 mL** conical tubes) is highly recommended to avoid the risk of contamination by accidentally attaching the side of the pipette to a wall of the conical tube when pipetting.

24

25

26

Add $\Box 6$ mL of DMEM/F12 to rinse any remaining cells off the plate, transfer to a conical tube, and mix the well with the cell suspension in a total volume of $\Box 10$ mL.

4m

Count cells and centrifuge the conical tube containing the cell suspension at **200 x g, 00:04:00**.

Resuspend the cells in Essential 8 medium supplemented with Y-27632 by adjusting the concentration of the single-cell suspension to 2.5–3 million in **1 mL**, and transfer the **1 mL** of cell suspension into a well of an AggreWell plate to achieve a final volume of **1.5 mL** per well

protocols.io

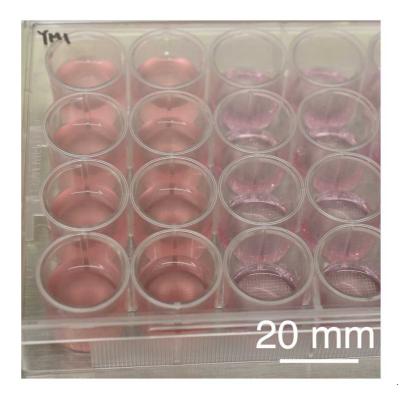


Image of detached cells in an Aggrewell plate during the aggregation process.


27 Pipette gently several times to re-establish an even distribution of cells throughout the well.

28

Centrifuge the AggreWell plate at **100** x g, **00:03:00** to distribute the cells in the micro-wells.



3m

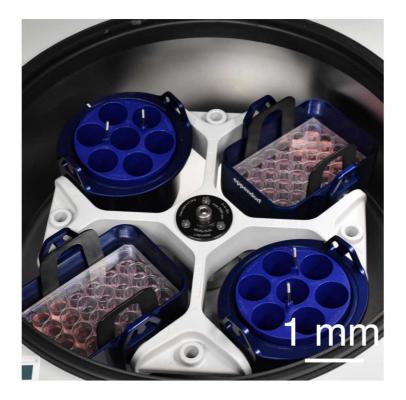


Centrifugation of an Aggrewell plate.

29



Check under a microscope to verify that cells are evenly distributed among the microwells.

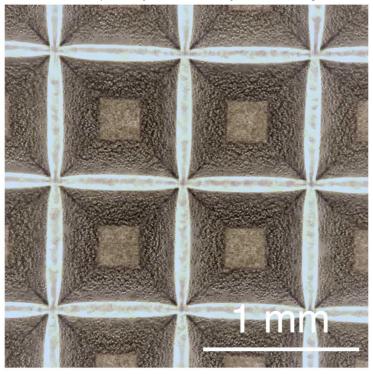

hiPS cells in the Aggrewell plate after centrifugation.

CRITICAL: Overfilled individual wells of the AggreWell plate after centrifugation will facilitate aggregation across wells and prevent the formation of individual spheroids.

1d

30

Incubate the cells at $\ 8\ 37\ ^{\circ}\text{C}\ ,\,5\%\ \text{CO}_2\ \text{for}\ \ \textcircled{0}\ 24:00:00\ .$


TROUBLESHOOTING

Dislodging 3D spheroids from the AggreWell plate ● Timing 1 h 1w 0d 1h 30m

31 Pre-warm DMEM/F12 and Essential 6 medium at & Room temperature for © 00:30:00 to © 01:00:00 before dislodging.

32

Collect the spheroids from the microwells by gently pipetting the medium in the well up and down two to three times with a micropipetter with a $\Box 1000~\mu L$ tip that was cut with scissors to allow a wider opening.

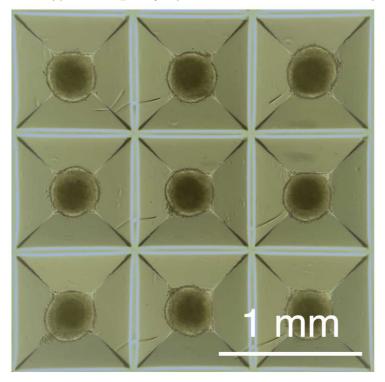

Aggregated hiPS cells days after aggregation in the Aggrewell plate.

33 Place a □40 μm cell strainer on top of a □50 mL conical tube and pass the suspension of spheroids through the strainer. The aggregates will remain on the filter of the strainer.

34

protocols.io

Pipette 1 mL of DMEM/F-12 medium across the entire surface of the well to dislodge any remaining spheroids. Collect washes and pass over the strainer used in the previous step.


35 / 60 6

Invert the strainer and place it over a new **30** mL conical tube. Collect the spheroids by adding **12** mL of Essential 6 medium on top of the inverted strainer. Observe the AggreWell plate under the microscope to ensure that all aggregates have been removed from the wells. Repeat the wash if necessary.

CRITICAL: Difficulty with detaching cells from the Aggrewell plate suggests differentiation of hPS cells.

Place spheroids in ultra-low-attachment 10-cm plates in Essential 6 medium supplemented with dorsomorphine (
[M]5 micromolar (μM)), SB-431542 ([M]10 micromolar (μM)).

If forebrain induction (e.g., low FOXG1 expression) is challenging with spheroids derived from a specific hPS cell line, we recommend testing whether addition of XAV-939 (IM1.25 micromolar (μ M)) from days \odot 24:00:00 to \odot 144:00:00 of differentiation may facilitate this process.

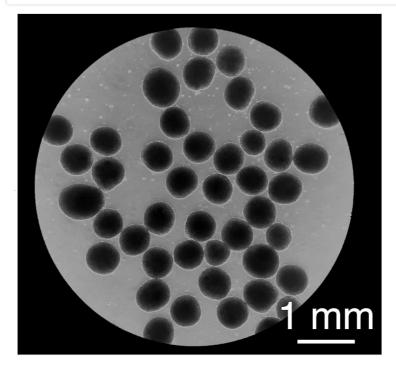

Representative images of spheroids

Neural induction 1w 3d

1w 3d

Add dorsomorphine and SB-431542 to Essential 6 medium for the first \odot **120:00:00** . XAV-939 (\blacksquare **1.25 micromolar (\muM)**) can also be added for the first \odot **120:00:00** of differentiation if forebrain differentiation is not highly efficient.

38

23

Citation: Yuki Miura, Min-Yin Li, Omer Revah, Se-Jin Yoon, Genta Narazaki, Sergiu Pasca Engineering brain assembloids to interrogate human neural circuits https://dx.doi.org/10.17504/protocols.io.36wqq4xxkvk5/v1

Perform a medium change daily (but not on day 1) at approximately the same time of the day by gently transferring the neural spheroids to a **50 mL** Falcon tube (if they are too small) and aspirating the medium once they sit at the bottom of the tube. When spheroids are large, perform the medium change by gently tilting the plate (until spheroids settle at the bottom) and removing the remaining medium. This can minimize the transferring time of the spheroids in and out of the plate at this critical differentiation period.

CRITICAL: The time that spheroids are kept outside of the incubator should be minimized at all stages of differentiation. After they reach 4 mm in diameter, keep no more than 30 spheroids per 10-cm dish. **TROUBLESHOOTING**

Patterning and differentiation

- 39 From day 6, replace the medium daily with NM with supplements listed in option A for hCSs and option B for hStrSs.
 TROUBLESHOOTING
- 40 (A) Generation of hCSs
 - From day 6, replace the medium daily with NM supplemented with [M]20 ng/ml FGF2 and [M]20 ng/ml EGF for generating hCSs. Change the medium daily from day 6 to day 15.
 - 40.2 From day 16 onward, change the medium every other day until day 22. Continue to use NM supplemented with [M]20 ng/ml FGF2 and [M]20 ng/ml Q29 EGF.
 - 40.3 From day 22 to day 46, replace FGF-2 and EGF with [M]20 ng/ml BDNF, [M]20 ng/ml NT-3, [M]200 micromolar (μM) AA, [M]50 micromolar (μM) dibutyryl cAMP sodium salt and [M]10 micromolar (μM) DHA. Perform medium changes with □12 mL to □15 mL of supplemented NM per 10-cm plate every 2 d.
- 41 (B) Generation of hStrSs
 - 41.1 From day 6, replace the medium daily with NM supplemented with M 50 ng/ml activin A and 704 [M]2.5 micromolar (μM) IWP-2. Change the medium daily from day 6 to day 11.
 - 41.2 From day 12 to day 15, supplement NM with №150 ng/ml activin A, ■2.5 μm IWP-2 and №1100 nanomolar (nM) SR11237. Perform medium changes with ■12 mL to ■15 mL of NM per 10-cm plate every day.
 - 41.3 From day 16 onward, change the medium every other day until day 22. Continue to use NM supplemented with [M]50 ng/ml activin A, 2.5 μm IWP-2 and [M]100 nanomolar (nM)

 SR11237. Perform medium changes with 12 mL to 15 mL of NM per 10-cm plate every 2 d.
 - 41.4 From day day 22 to day 41, supplement NM with [M]20 ng/ml BDNF, [M]20 ng/ml NT-3, [M]200 micromolar (μM) AA, [M]50 micromolar (μM) dibutyryl cAMP sodium salt and [M]10 micromolar (μM) DHA. Perform medium changes with □12 mL to □15 mL of NM per 10-cm plate every 2 d.
 - 41.5 From day 42 to day 46, supplement NM with [M]20 ng/ml BDNF, [M]20 ng/ml NT-3,



[M]200 micromolar (µM) AA, [M]50 micromolar (µM) dibutyryl cAMP sodium salt, [M]10 micromolar (µM) DHA and [M]2.5 micromolar (µM) DAPT. Perform medium changes with 12 mL to 15 mL of NM per 10-cm plate every 2 d.

Long-term culture and phenotyping Timing indefinite

18h 40m

After day 46, perform medium changes with 12-15 ml of NM without growth factors every 4-5 d.

Examples of additional experiments that can be carried out by using these spheroids include viral labeling of spheroids (option A), generation of cortico-striatal assembloids (option B), 3D clearing with CUBIC (option C), axon projection imaging (option D), retrograde labeling with ΔG rabies virus (option E), optogenetics coupled with calcium imaging (option F) and electrophysiological recording (option G)

- 43 (A) Viral labeling of neural spheroids Timing 3 d
 - 43.1 For labeling neural spheroids, dilute $\Box 0.5 \, \mu L$ of virus solution (recommended virus particle volumes are listed in Table 1) in ■200 µL of NM and add mixed NM to a ■1.5 mL Eppendorf tube. Then, transfer one to three spheroids in the tube containing NM and viruses. The tube can be placed overnight in an incubator at § 37 °C, 5% CO2.

We note that transduction efficiency depends on the target cell types, the promoter, stage of differentiation and brain region. CRITICAL The optimal serotype and titer of each virus depends on the gene construct of interest, as well as the target cell type. AAV-1, AAV-8, and AAV-DJ serotypes have been more extensively tested in hCSs and hStrSs. Virus at a titer higher than $1.0 \times 10^{12}-1.0 \times 10^{13}$ genome copies/ml was obtained from the Stanford Neurosciences Institute Gene Vector and Virus Core or from Addgene and stored at 8 -80 °C. ΔG -rabies virus at a titer higher than 1.0 × 100.8 transduction units/ml was obtained from the Salk Institute, diluted to 1:10 by using NM and stored at 8-80 °C.

43.2

1d

The next day, add an extra $\blacksquare 800~\mu L$ of fresh NM to the tube and incubate for @24:00:00.

43.3

On day 3, wash the spheroids with 📮 1 mL of fresh NM (two to three times) and then transfer the spheroids back to their culture plate.

Expression of the gene of interest is usually observed 7-14 d after AAV infection and 3-5 d after ΔG -rabies virus infection. ? TROUBLESHOOTING


(B) Generating cortico-striatal assembloids ● Timing 3-4 d

44.1

To generate cortico-striatal assembloids, transfer one hCS and one hStrS into a 1.5 mL

Eppendorf tube containing **1 mL** of NM.



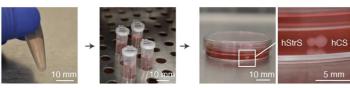

Scheme illustrating viral labeling and assembly of an hCS and an hStrS, and images showing the process of assembly.

CRITICAL The size of the spheroids for assembly should be <3 mm. Once in the tube, each spheroid should be placed horizontally next to each other at the bottom of the tube (but not stacked vertically) (Fig. 2m).

44.2

Incubate the spheroids at 8 37 °C for 3-4 d, completely replacing the medium on day 2.

CRITICAL The medium change should be performed gently to avoid breaking the fusion.

44.3

After spheroids are assembled, transfer them into an ultra-low-attachment plate by using a P1000 pipette with the tip cut for a larger bore opening. Assembly of hCS and hStrS between days 46 and 80 is more efficient.

Assembly is usually complete after 3-4 d. ? TROUBLESHOOTING

45 (C) Cryosection and immunostaining ● Timing 6-7 d

45.1

When ready to perform immunocytochemistry, transfer spheroids to a \Box 1.5 mL microcentrifuge tube by using a cut P1000 tip. Gently remove the medium and add \Box 500 μ L to \Box 1 mL of cold 4% PFA/PBS (or warm 4% PFA/4% sucrose/PBS). Leave spheroids at &4 °C \odot 0 Overnight.

45.2

After fixation, gently remove PFA and wash with PBS at 8 Room temperature.

protocols.io

1d

45.3

Add 11 mL of 30% (wt/vol) sucrose/PBS for cryopreservation. Leave spheroids or assembloids at § 4 °C until they sink to the bottom of tube ($\geq 24-48$ h).

Prepare square disposable molds that fit with the sample size and number and fill the molds with 45.4 embedding solution (1:1, OCT/30% sucrose; see Reagent setup).

> CRITICAL To avoid creating bubbles while filling the mold, use a P1000 pipette in which the tip has been slightly cut to widen the opening.

45.5

45.7

molds at 8-80 °C.

Transfer spheroids from tubes into the mold by using a P1000 cut pipette tip. Place 4-10 spheroids in the mold. If spheroids move while embedding, use a P20 tip to gently place them in the embedding solution.

CRITICAL Avoid carrying sucrose when transferring spheroids.

- $45.6 \quad \text{Place the spheroids in a cryomold} \quad \textbf{\& On ice} \quad \text{for } \odot \textbf{00:20:00} \quad \text{and allow the spheroids to sink to the}$ bottom of the mold.
- (11) Snap-freeze the spheroids by placing the mold directly on dry ice. Once completely frozen, store

Frozen molds can be stored indefinitely at 8-80 °C until ready to cryo-section.

- 45.8 When ready for sectioning, transfer the mold from & -80 °C to the cryostat chamber at & -20 °C for ~ © 00:30:00 before sectioning.
- Remove the block from the mold, paste the specimen on the stage by using OCT and section by using standard techniques. Use a brush to prevent crumbling of the sections. Spheroids and assembloids can be cryo-sectioned at → 10 μm - → 30 μm thickness, depending on the purpose of the experiments.
- 45.10 Collect sections on Superfrost Plus slides. Slides can be stored at & -20 °C to & -80 °C until ready to immunostain.
- For immunostaining, remove the slide from the freezer and leave at 8 **Room temperature** to thaw for © 00:05:00 - © 00:10:00


45.12

Wash the slide three times with PBS to remove excess OCT/sucrose on the slides.

45.13 Remove PBS and block for **© 01:00:00** at **§ Room temperature** with 10% normal donkey serum (NDS), 0.3% Triton- X and 0.1% BSA in PBS.

Parafilm can be placed on top of the slide glass to avoid drying of the antibody solution. Alternatively, a hydrophobic PAP pen can be used to draw circles around the sections, and staining can be performed inside these circles.

45.14 Perform the immunostaining steps in a humidified chamber to avoid drying of sections.

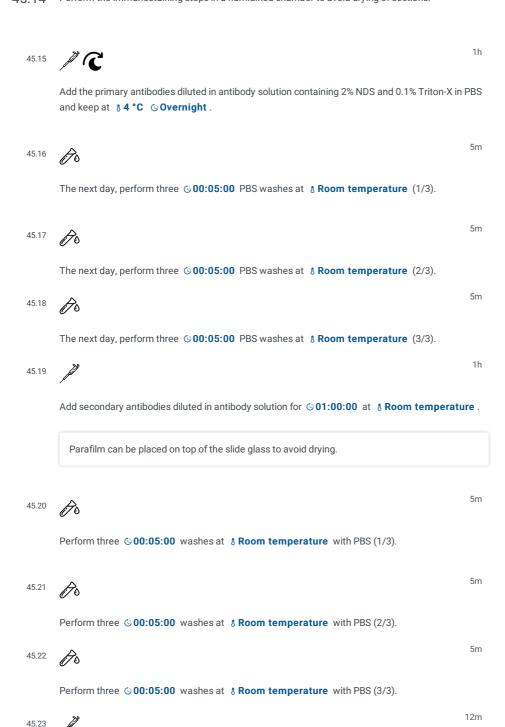

protocols.io

Add Hoechst in PBS (1:10,000 dilution) for © 00:05:00 - © 00:07:00 at § Room temperature.

45.24



Aspirate the Hoechst solution and wash twice with PBS.


45.25 Aspirate PBS and mount the sections by adding one drop of Aquamount solution on top of sections and by placing a coverglass on top.

45.26




Use a fluorescent microscope to image.


DARPP32 GAD67 CTIP2

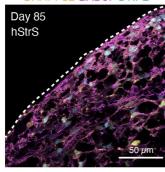


DARPP32 NeuN CTIP2

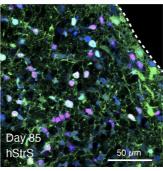

Left: Image showing immunostaining for DARPP32 (yellow), GAD67 (magenta) and CTIP2 (cyan) in an hStrS at day 85. Right: Image showing immunostaining for DARPP32 (green), NeuN (cyan) and CTIP2 (magenta) in an hStrS at day 85.

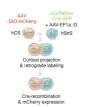

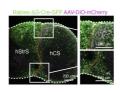







Scheme illustrating retrograde viral labeling of cortical projecting neurons in cortico-striatal assembloids, and representative retrograde viral images of Daf21 cortico-striatal assembloids.

 ${\tt CRITICAL\ Ideally,\ cryosections\ should\ be\ sampled\ from\ several\ layers\ of\ the\ spheroids\ or\ assembloids}$

46 (D) CUBIC clearing and imaging ● Timing 4-7 d

46.1


40m

Transfer spheroids to a \square 1.5 mL microcentrifuge tube by using a cut P1000 tip. Gently remove the medium and add \square 500 μ L to \square 1 mL of cold 4% PFA/PBS (or add to \square 1 mL of warm 4% PFA/4% sucrose/PBS and leave at 8 Room temperature for \bigcirc 00:20:00). Leave spheroids at 8 4 °C \bigcirc 0 Overnight .



29

Citation: Yuki Miura, Min-Yin Li, Omer Revah, Se-Jin Yoon, Genta Narazaki, Sergiu Pasca Engineering brain assembloids to interrogate human neural circuits https://dx.doi.org/10.17504/protocols.io.36wqq4xxkvk5/v1

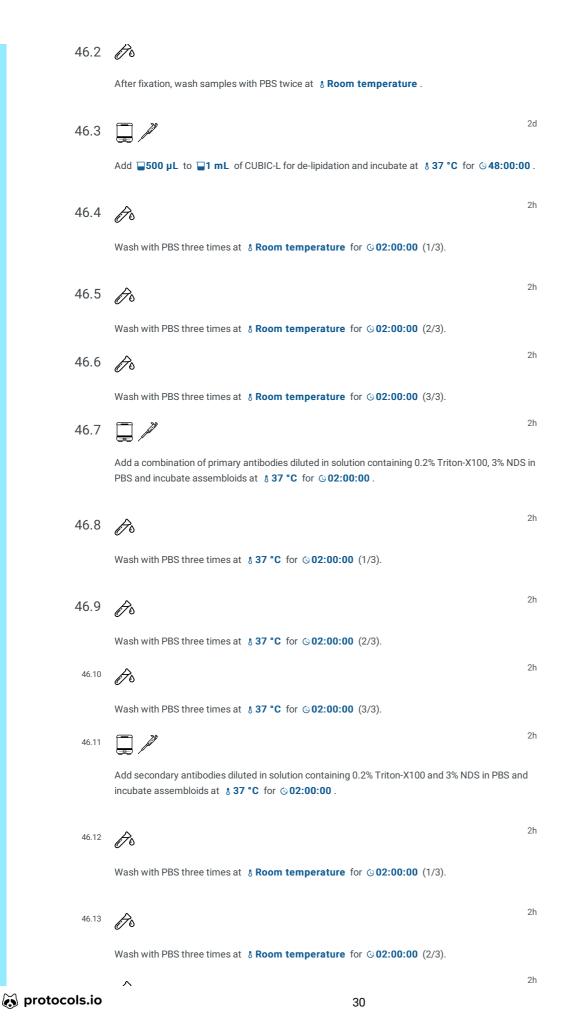

Citation: Yuki Miura, Min-Yin Li, Omer Revah, Se-Jin Yoon, Genta Narazaki, Sergiu Pasca Engineering brain assembloids to interrogate human neural circuits https://dx.doi.org/10.17504/protocols.io.36wgq4xxkvk5/v1

46.14 10

Wash with PBS three times at & Room temperature for © 02:00:00 (3/3).

46.15

2h

For refractive index-matching, remove PBS from the $\square 1.5$ mL tubes, add $\square 1$ mL of CUBIC-R in the $\square 1.5$ mL tube and leave for $\bigcirc 02:00:00$ at & Room temperature.


46.16




Once assembloids become transparent, transfer them to a well of a 96-well imaging- compatible plate filled with $\square 150~\mu L$ of CUBIR-R+ solution and image them with a confocal microscope

3D clearing and imaging



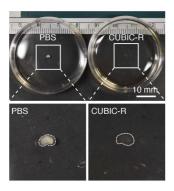

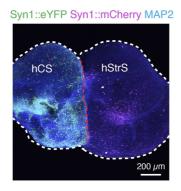



Scheme illustrating 3D clearing and imaging, and images of cleared assembloids and 3D stained cortico-striatal assembloids (eYFP: green; mCherry: magenta; MAP2: blue). RI, refractive index.


47 (E) Live imaging of axon growth ● Timing 2-8 h

CRITICAL For live imaging of axon growth in assembloids, viral labeling and fusion are performed at days 60–65, and assembloids are imaged at 8, 14 and 21 d after fusion.

47.1 Using a cut P1000 tip, transfer assembloids to a well of a 96-well imaging-compatible plate in $\Box 150~\mu L - \Box 200~\mu L$ of NM.

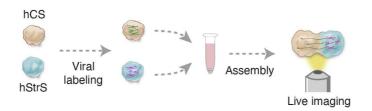

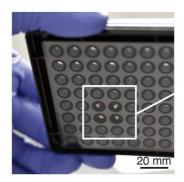


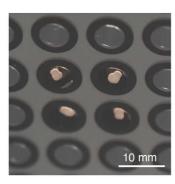

Scheme and images illustrating viral labeling, assembly and live imaging of cortico-striatal assembloids after their transfer into 96-well plates.

CRITICAL Avoid adding bubbles at the bottom of the wells. The distance from the lens to the sample is critical for imaging.

- 47.2 Place the 96-well plate into an inverted microscope with a motorized stage under environmentally controlled conditions (§ 37 °C, 5% CO2), and keep the assembloids for © 00:15:00 © 00:30:00 before experiments to let assembloids settle at the bottom of the plate.
- 47.3 Set up imaging positions with z-stacks and take images by using an $\times 10$ objective lens at a depth of $0-500~\mu m$.
- 47.4 Process fluorescent images with Fiji and quantify the percentage of eYFP coverage in hStrSs and the percentage of mCherry coverage in hCSs.

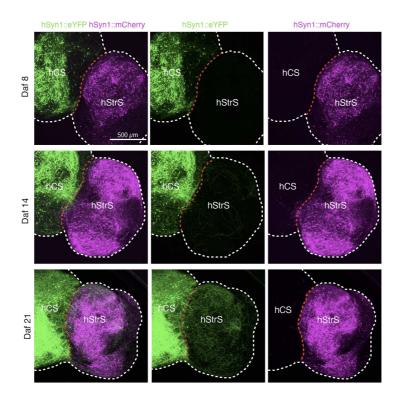

Representative images of Daf8-21 cortico-striatal assembloids labeled with eYFP and mCherry.


48 **(F) Optogenetic stimulation and calcium imaging ●** Timing 2-6 h

CRITICAL: For optogenetic stimulation and calcium imaging in assembloids, viral labeling and fusion are performed at days 60–65, and assembloids are imaged at 30–90 d after fusion.

48.1


Transfer an assembloid onto a +200 mm glass coverslip in a +35 mm plate in NM and image by using an ×10 objective on a Leica TCS SP8 confocal microscope.

48.2 For optogenetic stimulation, activate ChrimsonR-tdTomato+ cells in hCSs with +-625 nm light by using an optical fiber-coupled LED. Our stimulation experiments included 1,500 frames, and one +-625 nm pulse of LED light (68 ms) was applied every 150 frames as generated by a Cyclops LED connected to the LEICA TCS SP8 microscope.

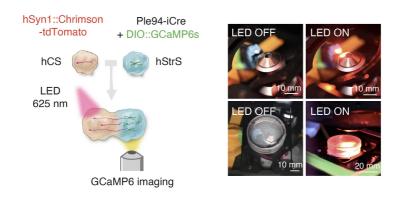

Scheme illustrating viral labeling and assembly for optogenetics coupled with calcium imaging, representative images of imaging setup during optogenetic stimulation using a 625-nm LED.

48.3 Image GCaMP6 at a frame rate of 14.7 frames/s.

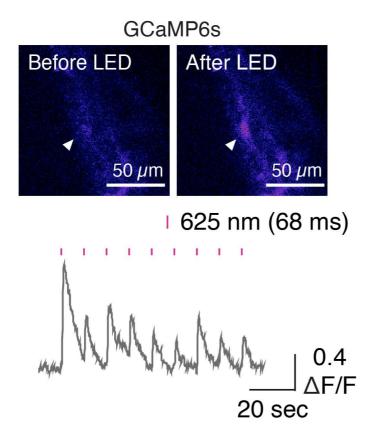


Representative image and trace of GCaMP6s signal in response to LED stimulation.

- 48.4 Process raw GCaMP6 signals with Fiji.
- 48.5 After registration of ROIs, transform raw time series to relative changes in fluorescence by using the following formula: $\Delta F/F = (F_t F_0)/F_0$.
- 48.6 To verify whether responses are time-locked to LED stimulation, compare $\Delta F/F$ responses to $\Delta F/F$ values obtained at randomly selected time points.



For quantification of time-locked $\Delta F/F$ values, calculate the amplitude of $\Delta F/F$ values from each cell as the maximum $\Delta F/F$ values within 20–30 frames (1,360–2,040-ms window) after LED stimulation (minus the mean of the baseline 1 s before stimulation).

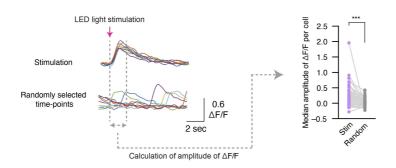

Quantification of stimulation-dependent response.n=180 cells from 10 assembloids with three hiPS cell lines; two-tailed Wilcoxon test. ***P=0.0002. Stim, stimulation.

TROUBLESHOOTING

49 (G) Whole-cell electrophysiological recording • Timing 2-6 h

CRITICAL For whole-cell electrophysiological recording in assembloids, viral labeling and fusion are performed at days 60-65, and recordings are performed at 30-60 d after fusion.

49.1

Add **Q0.2** g of low-melt agarose into a **35 mm** plate containing **Q5 mL** of aCSF and mix well (4% agarose). Dissolve low-melt agarose by using a microwave (avoid spillover due to overheating).

CRITICAL:

- Usually 10 s is sufficient.
- STEP BrainPhys can be used in place of aCSF.

49.2

Wait for a few minutes for the agarose to cool down to § 37 °C and transfer spheroids or assembloids from the medium into agarose with a cut P1000 pipette tip.









Images showing the process of agarose embedding and slicing of cortico-striatal assembloids.


CRITICAL: Avoid transferring the medium to prevent changes in the local concentration of agarose. Cool down agarose at § 4 °C for ⊙ 00:10:00 until it becomes solid.


- 49.3 Cut a block of agarose out of the culture plate by using a blade and glue the block on the plate of the vibratome. Slice spheroids at a + 200 μm to + 250 μm thickness with cold aCSF in a vibratome.

 Transfer slices into aCSF (at δ Room temperature bubbled with 95% O₂/5% CO₂).
- 49.4 Patch-clamping on spheroid slices is usually performed on fluorescently labeled neurons. Record labeled neurons with glass pipettes with small tips (8–10 M Ω with K-gluconate internal solution, 7–9 M Ω with Cs internal solution)


CRITICAL: Use glass pipettes with small tips instead of large tips.


- 49.5 Once cells are identified through fluorescence, patch under differential interference contrast (DIC) fluorescence. For DIC, the microscope is switched to a bright field. For patching under fluorescence, there are at least two approaches:
 - Increase the bright field exposure under epifluorescence so that both fluorescence (to visualize the cell) and bright field (to visualize the pipette tip) can be seen simultaneously. Once a 'dimple' is seen under fluorescence or increased pipette resistance is achieved, release positive pressure and form a tight seal.
 - Use fluorophore-coated pipettes, patching under epifluorescence. Back-fill glass pipettes by dipping the tips in the BSA-Alexa Fluor (0.02% in PBS) in a 1.5-ml Eppendorf tube for **© 00:00:05** to **© 00:00:10** and, after drying for **© 00:10:00**, fill pipettes with internal solution and patch.

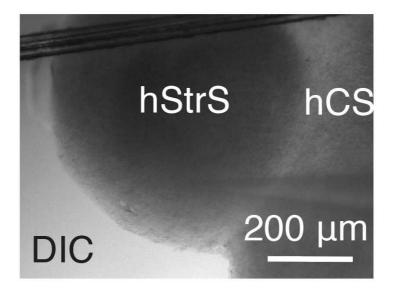

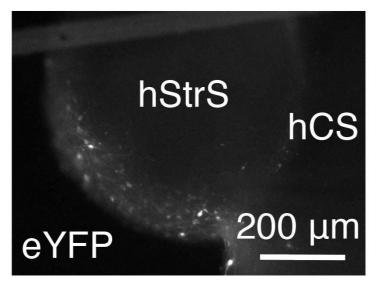

DIC (top) and fluorescent (bottom) images of sliced cortico-striatal assembloids expressing eYFP in an hStrS.

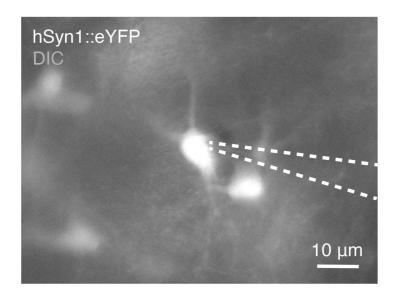

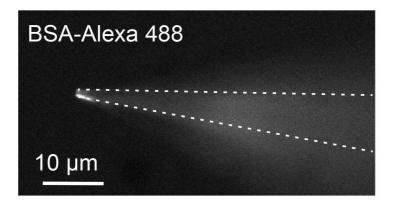


Representative image of patching from hSyn::eYFP+neurons with a non-fluorescently labeled pipette, and image showing afluorophore-coated patch pipette.

49.6 Once the pipette is approaching the target cell, provide minimal positive pressure. Once the tight seal is formed, strong suction is required to break cells.

CRITICAL: It is not always necessary to form a 'dimple' to get a Giga-seal; in most cases, an increase in the pipette resistance is suggestive of contact of the pipette tip with the cell membrane.

49.7 For recording optically evoked responses in assembloids, apply light at the maximal power through the 40× objective by using a CoolLED (++480 nm for ChR2 and ++550 nm for Chrimson). CoolLED is controlled by a protocol in pClamp.

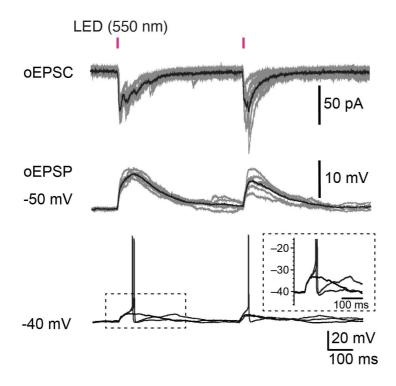

Representative traces of oEPSC, oEPSP and spikes from striatal neurons in cortico-striatal assembloids. oEPSC, optically evoked excitatory postsynaptic current; oEPSP, optically evoked excitatory postsynaptic potential.

11h

CRITICAL: Always check opsin expression before recording.

49.8 Make recordings at $\, \&\, 30\, ^{\circ}C$ to $\, \&\, 37\, ^{\circ}C$ or $\, \&\, Room\ temperature\$ within $\odot\, 05:00:00$ to

TROUBLESHOOTING

© 06:00:00 after sectioning.