



Mar 01, 2022

Human Placenta Tissue Collection and Preservation Methods - UCSD Female Reproductive TMC

Valentina Stanley¹, Scott Lindsay-Hewett¹, Louise Laurent¹, Mana Parast¹¹University of California, San Diego

1

dx.doi.org/10.17504/protocols.io.b5j9q4r6

Human BioMolecular Atlas Program (HuBMAP) Method Development Community
Tech. support email: Jeff.spraggins@vanderbilt.edu

Valentina Stanley
University of California, San Diego

Human Placenta Tissue collection and storage protocol for HuBMAP's UCSD Female Reproductive TMC.

DOI

dx.doi.org/10.17504/protocols.io.b5j9q4r6

Valentina Stanley, Scott Lindsay-Hewett, Louise Laurent, Mana Parast 2022.
Human Placenta Tissue Collection and Preservation Methods - UCSD Female Reproductive TMC. **protocols.io**
<https://dx.doi.org/10.17504/protocols.io.b5j9q4r6>



HuBMAP, UCSD Female Reproductive TMC, Placenta

protocol ,

Feb 24, 2022


Mar 01, 2022

58721


1. Sterile surgical dissecting instruments (scissors/forceps)
2. Gauze pad
3. Petri dish
4. Sterile 1X PBS (cold, from fridge)
5. Liquid Nitrogen
6. Empty microfuge tubes x9
7. Microfuge tubes pre-filled with RNA later x9
8. 10% Formalin filled tubes x2
9. MACS tissue buffer filled tubes x4 or OCT compound and cryomolds x4

Preparation


- 1 Collect placenta from delivery room within 1 hour of delivery.
- 2 Place placenta into bucket and photograph prior to sampling.
- 3 Select 2 sites (A and B) from which samples will be collected. Aim for sites to be equidistant from the cord, avoid the edges, avoid areas with fibrin deposits/infarcts/necrosis and/or large blood vessels.
- 4 From site A, cut a circle all the way down from the fetal to the maternal surface, about one inch in diameter.

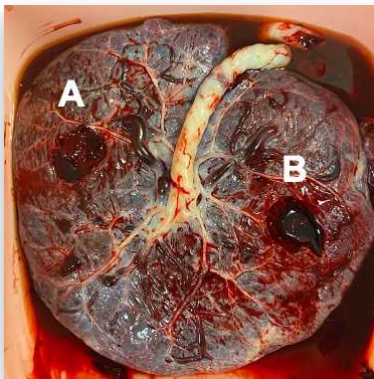

- 5 Wash the full-thickness core in cold PBS in a petri dish. Keep the fetal surface pointing up.

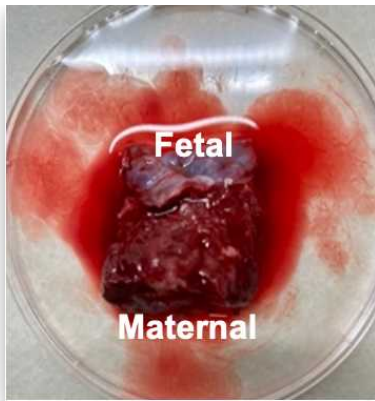


- 6 Divide the full core in half vertically, from fetal to maternal surface.

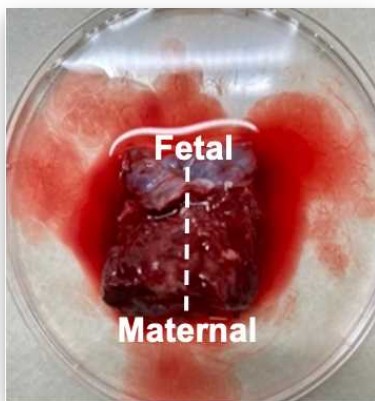


Flash freezing and storing in RNAlater

- 7 With the left half of site A, divide the core horizontally into fetal (F), central (C) and maternal (M) regions.

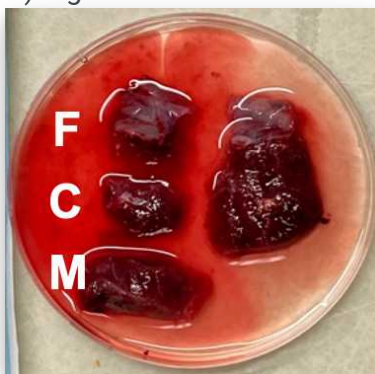

- 7.1 Cut each of the three regions from the previous step into 6 small chunks for a total of 18 chunks. Chunks should be small enough to just cover the bottom of a microfuge tube (~5mm x 5mm x 5mm). Collect tissue directly adjacent to


the fetal and maternal membranes but *do not* cut from the membranes themselves.

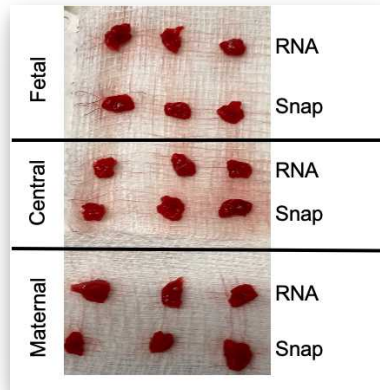

- 7.2 Dab the small chunks on a sterile gauze pad and place each one into a labeled microfuge tube.
- 7.3 Place 9 small chunks (3 fetal, 3 central, 3 maternal) into the empty tubes to be snap frozen. Drop tubes into liquid Nitrogen. Leave for ~2-10 minutes. Remove and store in -80C freezer.
- 7.4 Place the other 9 chunks (3 fetal, 3 central, 3 maternal) into the RNAlater filled tubes. Place RNAlater tubes into a 4C fridge. Allow to sit for 24-48 hours, then remove the RNAlater with a sterile transfer pipette and store in -80C freezer.


Storing in formalin and MACS

- 8 With the right half of site A, divide the tissue into 4 quadrants so that you have 2 fetal chunks and 2 maternal chunks, keeping all the membranes intact. Some of the central region tissue may need to be trimmed off if the core is very thick.

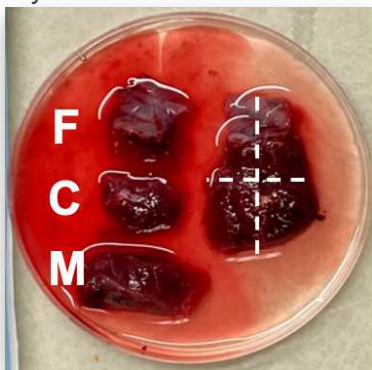

- 9 Place one fetal piece and one maternal piece into the 10% formalin tubes (to be processed into FFPE blocks). Place the other fetal and maternal piece into cold MACS tissue buffer. Store in

4C fridge.

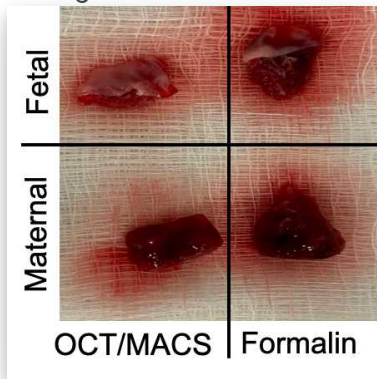

Embedding in OCT

- 10 Alternatively, if not collecting into MACS tissue buffer, freeze the other fetal and maternal piece in OCT.
 - 10.1 Fill plastic cryomold halfway with OCT compound being careful not to create bubbles.
 - 10.2 Place fetal tissue (membrane side pointing up) in one mold and maternal tissue (membrane side pointing down) in another mold.
 - 10.3 Cover tissue with a thin layer of OCT compound being careful not to create bubbles.
 - 10.4 Slowly lower the mold into liquid nitrogen until the whole block freezes, then store in -80C freezer.
- 11 Repeat Steps 4-10 with placenta Site B.
- 12 Take photos of the placenta after both sites (A and B) are removed. Measure the distance from the cord insertion to each sampling site.