



OCT 18, 2023


## Polysomnography

Núria  
Peñuelas<sup>1</sup>

<sup>1</sup>Vall d'Hebron Research Institute

Vilalab Public

Vilalab Public

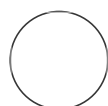

Núria Peñuelas  
Vall d'Hebron Research Institute

OPEN ACCESS



### ABSTRACT

Polysomnography in mice

**DOI:**  
[dx.doi.org/10.17504/protocols.io.rm7vzx77xgx1/v1](https://dx.doi.org/10.17504/protocols.io.rm7vzx77xgx1/v1)

**Protocol Citation:** Núria Peñuelas 2023. Polysomnography. **protocols.io** <https://dx.doi.org/10.17504/protocols.io.rm7vzx77xgx1/v1>

**License:** This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

**Protocol status:** Working  
We use this protocol and it's working

**Created:** Oct 18, 2023

**Last Modified:** Oct 18, 2023

- 1 Anesthetize the animals (ketamine and xylazine, 100/10 mg/kg. i.p.) for polysomnography.
- 2 Screw 3 electrodes to skull above parietal and frontal cortex and cerebellum as the reference for electroencephalogram (EEG).
- 3 Slip 2 wire electrodes between neck muscles for electromyogram (EMG).
- 4 Braze the leads to a miniature plug (Plastics One, Billerica, MA), cemented to skull using acrylic cement (Superbond, C&B Sun Medical) and coated using dental cement (Paladur, Heraeus Kulzer).
- 5 For pain caring, administer a subcutaneous injection of carprofen (5 mg/kg).
- 6 Once recovered, house mice in individual barrels under 12h light/dark cycle (light-on 7 a.m.) and constant temperature, tethered to the acquisition system and record 1 week out of 3 over 12 consecutive months until natural death.
- 7 Based on amplified/digitized unipolar EEG and bipolar EMG signals, visually score vigilance states using a 5s sliding window frame and assign them as Waking (W), Slow Wave Sleep (SWS) and Paradoxical Sleep (PS) to produce hypnograms and state quantifications (Sleepscore, Viewpoint).