




SEP 30, 2023

Protein Digestion and Mass Spectrometry Analysis Protocol

Anthony Iavarone¹, Annan SI Cook^{1,2}

¹University of California, Berkeley; ²ASAP CRN

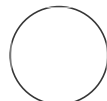

Annan SI Cook

ABSTRACT

This protocol details a protocol for the trypsin digestion of PI3KC3-C1 and subsequent analysis of digested peptides using LC-MS/MS.

ATTACHMENTS

[852-2203.pdf](#)

MATERIALS

Materials

- PI3KC3-C1 (mCherry-ATG14|VPS15-TSF) protein sample
- Urea
- tris(hydroxymethyl)aminomethane (TRIS) pH 8.0
- tris(2-carboxyethyl)phosphine (TCEP)
- Iodoacetamide (IAA)
- Trypsin (Gold Trypsin, 0.5 mg/mL Promega V5280)
- CaCl₂
- NH₄CH₃CO₂
- Acetonitrile
- Formic acid (Optima LC-MS grade, 99.9% minimum)
- Water (purified using the Milli-Q Gradient system)
- Zorbax 300SB-C8 Micro Bore Rapid Resolution column (150 mm length, 1.0 mm inner diameter, 3.5 µm particle size, Agilent)
- LC-MS system (e.g., Agilent 1200 series LC and Thermo Fisher Scientific LTQ-Orbitrap-XL mass spectrometer)
- Proteome Discoverer software (version 1.3, Thermo Fisher Scientific)

OPEN ACCESS



Protocol Citation: Anthony Iavarone, Annan SI Cook 2023. Protein Digestion and Mass Spectrometry Analysis Protocol. [protocols.io](https://protocols.io/view/protein-digestion-and-mass-spectrometry-analysis-p-c2pfydjn) <https://protocols.io/view/protein-digestion-and-mass-spectrometry-analysis-p-c2pfydjn>

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working
We use this protocol and it's working


Created: Sep 25, 2023

Last Modified: Sep 30, 2023

PROTOCOL integer ID:
88519

Protein Denaturation and Reduction

20m

- 1 Mix PI3KC3-C1 (mCherry-ATG14|VPS15-TSF) protein with [M] 8 Molarity (M) urea, [M] 50 millimolar (mM) TRIS  8.0, and [M] 10 millimolar (mM) TCEP.



- 2 Add Iodoacetamide (IAA) to achieve a final concentration of [M] 15 millimolar (mM) .



- 3 Incubate the mixture at  Room temperature for  00:20:00 to alkylate cysteine residues.



20m

Trypsin Digestion Setup

- 4 Prepare a trypsin digestion solution containing Gold Trypsin (0.5 mg/mL), [M] 50 millimolar (mM) TRIS  8, and [M] 100 millimolar (mM) CaCl₂.

Trypsin digestion solution

A	B
Gold Trypsin (0.5 mg/mL)	8 uL
50 mM TRIS pH 8	160 uL
100 mM CaCl ₂	2 uL

- 5 Mix the trypsin digestion solution with the alkylated protein sample.



Digestion Incubation

20m

- 6 Incubate the digestion mixture at  37 °C  Overnight with shaking at  200 rpm .



20m

Confirmation by SDS-PAGE

- 7 Analyze the digested proteins by SDS-PAGE to confirm successful digestion. The absence of bands corresponding to the PI3KC3-C1 components indicates successful digestion.



Liquid Chromatography (LC) Setup


- 8 Prepare LC mobile phase solvents: Solvent A and Solvent B.

Solvent A

A
99.9% water
0.1% formic acid

Solvent B

A
99.9% acetonitrile
0.1% formic acid

- 9 Set up the LC system with a Zorbax 300SB-C8 Micro Bore Rapid Resolution column (150 mm length, 1.0 mm inner diameter, 3.5 μ m particle size).
- 10 Maintain the column compartment at  50 °C.







LC-MS Analysis


2h

- 11 Load a  10 μ L sample onto the column.

- 12 Implement the following gradient elution program:

2h

- Isocratic flow at 1% (volume/volume) B for  00:02:00 .
- Linear gradient to 35% B over  01:30:00 .
- Linear gradient to 95% B over  00:01:00 .
- Isocratic flow at 95% B for  00:06:00 .
- Linear gradient to 1% B over  00:01:00 .
- Isocratic flow at 1% B for  00:20:00 .

13 Set the flow rate to  120 μL /min.

Mass Spectrometry Analysis

14 Acquire full-scan, high-resolution mass spectra in positive ion mode over the m/z range of 340 to 1800 using the Orbitrap mass analyzer with a mass resolution of 60,000 (at $m/z = 400$, FWHM).

Data-Dependent MS/MS Analysis


15 Select the ten most intense ions exceeding an intensity threshold of 10,000 raw ion counts from each full-scan mass spectrum for tandem mass spectrometry (MS/MS) analysis using collision-induced dissociation (CID).

16 Acquire MS/MS spectra using the linear ion trap.

Data Analysis

17 Use Proteome Discoverer software (version 1.3, SEQUEST algorithm) to search raw data files against the amino acid sequences of mCherry-labeled Class III phosphatidylinositol 3-kinase complex I (mCherry-PI3KC3-C1) proteins.

18 Consider tryptic peptides with up to three missed cleavages and specify post-translational modifications.



19 Validate peptide assignments by manually inspecting MS/MS spectra.