




Sep 06, 2022

Lentivirus production

Itika Saha¹, F. Ulrich Hartl¹, Mark S. Hipp^{2,3}¹Department of Cellular Biochemistry, Max Planck Institute of Biochemistry, Am Klopferspitz 18, 82152 Martinsried, Germany;²Department of Biomedical Sciences of Cells and Systems, University Medical Center Groningen, University of Groningen, Antonius Deusinglaan, 1, 9713 AV Groningen, The Netherlands;³School of Medicine and Health Sciences, Carl von Ossietzky University Oldenburg, 26129 Oldenburg, Germany

1 Works for me

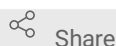

Share

dx.doi.org/10.17504/protocols.io.6qpvr4xn3gmk/v1

Felix Kraus

ABSTRACT

This protocol describes the production of lentiviruses to transduce HEK293T cells and has to be performed in a biosafety level 2 laboratory

DOI

dx.doi.org/10.17504/protocols.io.6qpvr4xn3gmk/v1

EXTERNAL LINK

<https://www.biorxiv.org/content/10.1101/2022.02.18.481043v1.full>

PROTOCOL CITATION

Itika Saha, F. Ulrich Hartl, Mark S. Hipp 2022. Lentivirus production.

protocols.io<https://protocols.io/view/lentivirus-production-cf7ntrme>

MANUSCRIPT CITATION please remember to cite the following publication along with this protocol

The AAA+ chaperone VCP disaggregates Tau fibrils and generates aggregate seeds Itika Saha, Patricia Yuste-Checa, Miguel Da Silva Padilha, Qiang Guo, Roman Körner, Hauke Holthusen, Victoria A. Trinkaus, Irina Dudanova, Rubén Fernández-Busnadiego, Wolfgang Baumeister, David W. Sanders, Saurabh Gautam, Marc I. Diamond, F. Ulrich Hartl, Mark S. Hipp *bioRxiv* 2022.02.18.481043; doi: <https://doi.org/10.1101/2022.02.18.481043>

LICENSE

————— This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

CREATED

Sep 05, 2022

LAST MODIFIED

Sep 06, 2022

PROTOCOL INTEGER ID

69582

SAFETY WARNINGS

Has to be performed in a biosafety level 2 laboratory.

Lentivirus production

1d 1h 10m

- 1 Plate $\sim 3.6 \times 10^6$ Lenti-X HEK293T cells (Takara) in 10 cm dish in 10 mL standard DMEM. Cells should be $\sim 80\%$ confluent at the time of transfection.

NOTE: Only low passage cells should be used.

- 2 Next day, remove 5 mL medium and replenish with fresh medium.
- 3 Warm up reduced serum medium e.g. Opti-MEM (Gibco) and transfection reagent to room temperature (RT). This protocol was performed with Lipofectamine 3000 transfection reagent (Thermo).
- 4 Add 24 μL Lipofectamine 3000 to 600 μL Opti-MEM, mix by vortexing and incubate 5 min at RT.^{5m}
- 5 In another tube, mix 6 μg plasmid containing gene of interest, 5 μg packaging plasmid psPAX2 (RRID:Addgene_12260), 1 μg envelope plasmid pMD2.G (RRID:Addgene_12259) and 24 μL P3000 reagent (provided by manufacturer along with Lipofectamine 3000 reagent) in 600 μL Opti-MEM, mix by vortexing and incubate 5 min at RT.
- 6 Mix contents of both tubes and incubate for 15 min at RT. 15m
- 7 Add DNA-lipid complex to cells dropwise.

8 2 days later, collect virus-containing medium and centrifuge for 5 min at 1,000 x g. 2d

9 Collect supernatant in a fresh tube and proceed with concentration.

Concentration 1d 1h 10m

10 Add Lenti-X concentrator (Takara) to clarified virus-containing medium at 1:4 dilution and mix well by gently inverting tube.

11 Incubate overnight or 2 h at 4 °C. 1d

12 Next day, centrifuge for 45 min at 1,500 x g at 4 °C followed by gently aspirating supernatant.^{45m}

13 Resuspend viral pellet in 100 - 1000 µL PBS, aliquot and store at -80 °C until use.