

A





Yingyu Hu¹, Zoe V Finkel¹

¹Dalhousie University

1



protocol

Marine Microbial Macroecology Lab Tech. support email: ruby.hu@dal.ca



Here we describe a protocol for extracting and quantifying bulk RNA and DNA from microalgae, which is adapted from Berdalet E. et al. (2005).

RNA and DNA are extracted from microalgae samples and then quantified by fluorochrome SYBR Green II.

The level of sensitivity of this method was set at ca. 40 ~300 ng RNA or 10 ~ 100 ng DNA (ml assay)-1.

Berdalet E, Roldán C, Olivar MP, Lysnes K. Quantifying RNA and DNA in planktonic organisms with SYBR Green II and nucleases. Part A. Optimisation of the assay. Scientia Marina. https://doi.org/10.3989/scimar.2005.69n11

Berdalet E, Roldán C, Olivar MP. Quantifying RNA and DNA in planktonic organisms with SYBR Green II and nucleases. Part B. Quantification in natural samples. Scientia Marina.

https://doi.org/10.3989/scimar.2005.69n117

Yingyu Hu, Zoe V Finkel 2022. Total RNA and DNA from Microalgae (24 samples per day) . **protocols.io** https://protocols.io/view/total-rna-and-dna-from-microalgae-24-samples-per-d-b6nnrdde Yingyu Hu

Explain the factors that can affect the fluorescence of stained RNA and DNA

RNA, DNA, SYBR Green II, DNase, RNase, microalgae, fluorescence

_____ protocol,

Mar 23, 2022

Mar 24, 2022

59822



1

Estimation of RNA/DNA in the collected microalgae samples:

Under replete condition, RNA and DNA is about 5.7% and 1% in total dry mass, while Chl-a is bout 1.1% in total dry mass. Therefore, RNA_ug/L = Chl-a_ug/L \times (5.7/1.1), DNA_ug/L = Chl-a_ug/L \times (1/1.1).

Common dilution from sample collected on the filter to assay is 1/40.

```
STEP MATERIALS

    ⋈ Nuclease decontamination

solution IDT Catalog #11-05-03-01 Step 2
Aldrich Catalog #R6513-50MG Step 6.1

    ■ DEOXYRIBONUCLEASE1 RNase and Protease

Free Bioshop Catalog # DRB002.10 Step 7.1

    Magnesium chloride solution Sigma

Aldrich Catalog #63069-100ML Step 40.1

    □ Calcium chloride solution Sigma

Aldrich Catalog #21115-100ML Step 42.1
SYBR™ Green II RNA Gel Stain, 10,000X concentrate in DMSO Thermo
Fisher Catalog #S7564 Step 49
Scientific Catalog #AAJ60080AK Step 3.1

    □ Deoxyribonucleic acid from calf thymus Sigma

Aldrich Catalog #D4522-1MG Step 5.1
N-Lauroylsarosine sodium salt solution (20% RNase/DNase free) Sigma
Aldrich Catalog #L744-50mL Step 10.1
⊠ EDTA buffer (0.5M DNase/RNase
free) Bioshop Catalog #EDT333.100 Step 10.2
⊠ UltraPure<sup>™</sup> DNase/RNase-Free Distilled
Water Thermofisher Catalog #10977023 Step 3.2
Scientific Catalog #AM7940 Step 4.1
```

No data is available addressing the mutagenicity or toxicity of SYBR® Green II Nucleic Acid Gel Stain. Because this reagent binds to nucleic acids, it should be treated as a potential mutagen and used with appropriate care. The DMSO stock solution should be handled with particular caution as DMSO is known to facilitate the entry of organic molecules into tissues. We strongly recommend using double gloves when handling the DMSO stock solution. As with all nucleic acid stains, solutions of SYBR Green II Nucleic Acid Gel Stain should be poured through activated charcoal before disposal or collected in waste container to be treated later. The charcoal must then be incinerated to destroy the dye.

Day 1: Freeze-dry samples

- 1 Freeze dry samples and blank filters. Freeze at 8 -80 °C until processed.
 - 1. Freeze-drying should be as short as possible to reduce sample degradation.
 - 2. The exact duration of freeze-drying depends on size of filter, quantity of sample and the size of container.

FreeZone® 2.5 L Benchtop Freeze Dryers Labconco® 700202000

Day 1: Prepare primary solutions

2 Turn on UV light in biosafety cabinet for © 00:15:00 and clean working surface with decontamination solution.

⊗ Nuclease decontamination

solution IDT Catalog #11-05-03-01

- 3 Prepare Tris buffer [M]5 mM pH8.0
 - 3.1 Pour MI M PH.0 Tris into an RNase free 15 mL Falcon tube

 Scientific Catalog #AAJ60080AK

Falcon® Centrifuge Tubes Polypropylene, Sterile, 15 mL Corning® 352096

3.2 Directly add 22.5 mL [M] M pH8.0 Tris into 500 mL RNase free water in its original package.

⊠UltraPure™ DNase/RNase-Free Distilled

Water Thermofisher Catalog #10977023

BT Barrier Pipet Tips
Pre-Sterile
Neptune® BT1250, BT100, BT10

- 4 RNA primary standard solution ([M]200 ug/ml)

https://assets.thermofisher.com/TFS-Assets/LSG/manuals/sp. 7940.ndf



4.2 Uncap the original package of E. Coli Total RNA and directly add $\blacksquare 800~\mu L$ Tris buffer (IMJ5 mM , pF8.0).

Cap the package and vortex for a thorough mix.

4.3 Aliquot 30 uL by stepper with sterile tip to 600 uL RNase free microtubes. Keep frozen at & -80 °C.

Finnpipette Stepper Pipette
Thermo Scientific™ 4540000

Finntip stepper pipette tips

500 ul (sterile)

Thermo Scientific™ Scientific 9404173

Microcentrifuge Tubes 1.7 mL/0.6 mL

Axygen Scientific MCT-175-C/MCT-060-L-C

- 5 DNA primary standard solution (≈ [M]500 ug/ml)
 - 5.1 Uncap the original package of Deoxyribonucleic acid from calf thymus and add ☐2 mL Tris buffer ([M]5 mM , [P+8.0]).

⊠Deoxyribonucleic acid from calf thymus **Sigma**

Aldrich Catalog #D4522-1MG

https://www.sigmaaldrich.com/content/dam/sigmaaldrich/docs/Sigma/Product_Information_Sheet/d4522pis.pdf

- 5.2 Cap the package. Do not vortex or sonicate.
- 5.3 Keep the solution at & 0 °C ~ & 4 °C overnight to completely solubilize the DNA. Gentle reversion is recommended.

5.4 Aliquot 10 uL by stepper with sterile tip to 600 uL RNase free microtubes. Keep frozen at 8-80 °C.

Finntip stepper pipette tips

500 ul (sterile)

Thermo Scientific™

Scientific 9404173

5.5 Dilute 5 ul primary DNA standard solution with 95 ul Tris buffer (M15 mM, PH8.0) in a microtube (600 ul).

Measure DNA concentration by using µdrop plate (sample volume: 4 ul)

Use Tris buffer ([M]5 mM , p⊦8.0) as blank.

µDrop™ Plates

Thermo Scientific N12391

Varioskan LUX Multimode Microplate

Reader

Thermo Fisher VL0L00D0

The dilution is to avoid saturated observation at 260 nm.

5.6 DNA concentration (μ g/ml) = (Abs_{260} - Abs_{260} (blank))x 50 μ g/ml x (10mm/0.5 mm) X DF Where, DF=20.

If the measured DNA concentration is not close to [M]500 ug/ml , check reverse pipetting technique.

- 6 RNase primary stock solution (M10 mg/ml)
 - 6.1 Uncap the original package of Ribonuclease A from bovin pancreas and add ☐5 mL Tris buffer ([M]5 mM , [p+8.0]).

Cap the package and vortex for a thorough mix.

Aldrich Catalog #R6513-50MG

6.2 Aliquot 30 uL by stepper with sterile tip to 600 uL RNase free microtubes. Keep frozen at 8-20 °C.

Finntip stepper pipette tips

500 ul (sterile)

Thermo Scientific™

Scientific 9404173

Finntip™ Stepper Pipette Tips

500 ul (Sterile)

Thermo Scientific 21-377-149

- 7 DNase primary stock solution ([M]5 mg/ml = 10,000 U/mL)
 - 7.1 Uncap the original package of Deoxyribonuclease1 and add 11 mL Tris buffer (1815 mM , pH8.0)

Cap the package and vortex for a thorough mix.

⊠DEOXYRIBONUCLEASE1 RNase and Protease

Free Bioshop Catalog # DRB002.10

7.2 Aliquot 100 uL by stepper with sterile tip to 600 uL RNase free microtubes. Keep frozen at 8 -20 °C

Finntip stepper pipette tips

1250 ul (sterile)

Thermo Scientific™

Scientific 9404183

Day 2: Exact RNA and DNA

- 8 Turn on UV light in biosafety cabinet for ©00:15:00 and clean working surface with decontamination solution.
- 9 Prepare falcon tubes and tube rack in biosafety cabinet

Α	В
Volume of tube (mL)	Contents in the tube
5	0.5 M EDTA
5	20% sarcosine
15 or 50	5 mM Tris
15 or 50	1% STEB

Falcon® Centrifuge Tubes
Polypropylene, Sterile, 15 mL
Corning® 352096

Falcon® Centrifuge Tubes
Polypropylene, Sterile, 50 mL
Corning® 352070

10 Prepare STEB ([M]1 %)

Use the following formula to determine the total volume of 1% STEB required: $(\# \text{ samples} + \# \text{ blank filters}) \times (500 \text{ ul}) + (500 \text{ ul}) = \text{total volume of } 1\% \text{ STEB required}$

10.1 Pour sarcosine ($\mbox{\scriptsize [M]} \mbox{\bf 20}$ %) into an RNase free 5 mL falcon tube.

⊗N-Lauroylsarosine sodium salt solution (20% RNase/DNase free) Sigma
 Aldrich Catalog #L744-50mL

- 10.2~ Pour EDTA ($\mbox{\scriptsize IM}\mbox{\scriptsize 0.5~M}$) into an RNase free 5 mL falcon tube.


⊠EDTA buffer (0.5M DNase/RNase

free) Bioshop Catalog #EDT333.100

10.3 Pour Tris buffer ([M] $\mathbf{5}$ mM , [p+ $\mathbf{8.0}$]) into an RNase free 15 or 50 mL falcon tube.

For large volume of Tris buffer usage, use sterile serological pipet

10.4 Mix $\square 500 \, \mu L$ sarcosine (M120 %), $\square 10 \, \mu L$ EDTA (M10.5 M) and $\square 9 \, mL + \square 490 \, \mu L$ Tris buffer (M15 mM , P+8.0) to obtain STEB (M11 %).

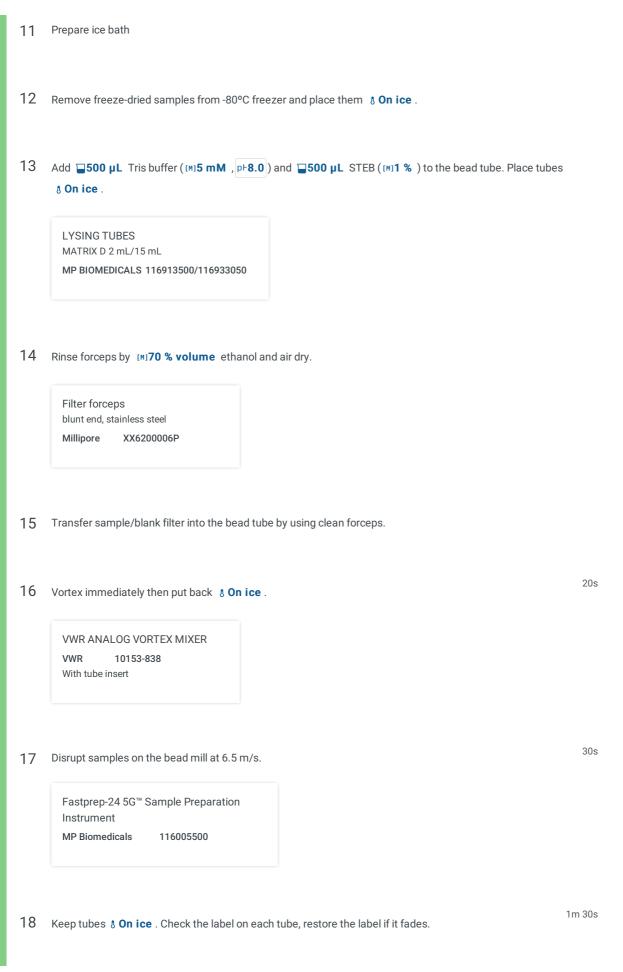

19	Disrupt samples on the bead mill at 6.5 m/s.	30s
20	Keep tubes & On ice . Check the label on each tube, restore the label if it fades.	1m 30s
21	Disrupt samples on the bead mill at 6.5 m/s	30s
22	Keep tubes & On ice . Check the label on each tube, restore the label if it fades.	1m 30s
23	Disrupt samples on the bead mill at 6.5 m/s.	30s
24	Continuously shake homogenate in a multi-head vortex at the highest speed for © 01:00:00 & Room temperature	1h
	Votex mixer should be able to remain stable on the bench under this vortex speed.	
25	In the biosafety cabinet, transfer 150 uL of homogenate into RNase free 600 uL tube.	
26	Freeze both aliquote (150 uL) and the rest of homogenate (in bead mill tube) at 8-80 °C until analyzed.	
Day 3:	Run the assay	
27	Prepare ice bath.	
28	Turn on UV light in biosafety cabinet for $©$ 00:15:00 and clean working surface with decontamination sol	lution.
29	Prepare falcon tubes, microtubes and tube racks in biosafety cabinet * Maximum number of samples (including blanks) per assay is 20.	

Α	В	С
Number of tubes	Type of tubes	Contents
3	5 mL falcon tubes	1 M MgCl2, 1 M CaCl2, Sybr Green II working solution (SG-II WS)
1	50 mL falcon tube	5 mM Tris buffer
4	15 mL falcon	0.05% STEB, Working solution A (WS-A), Working solution B (WS-B), Working
	tubes	solution (WS-C)
5	1.7 mL RNase free	RNase working solution, Secondary RNA standard solution, Secondary DNA
	tubes	standard solution, 900mM MgCl2, 900 mM CaCl2
33	1.7 mL RNase free	RNA standard solutions for RNA standard curves, DNA standard soutions for
	tubes	DNA standard curves
N= total number of	1.7 mL RNase free	Samples and blanks
samples and blanks	tubes	
3XN	1.7 mL RNase free	Diluted samples and blanks
	tubes	
5	Microtube racks	Tubes of 1.7 mL
1	Tube racks	Falcon tubes

Screw-Cap Centrifuge Tube

5 mL

VWR 10002-738

Day 3: Run the assay

30 Organize and label the tubes as shown below

Set 1:

In microtube rack, label 1.7 mL tubes for samples and blanks to be further diluted.

1	2	3	4	5	6	7	8	9		Blank
---	---	---	---	---	---	---	---	---	--	-------

Set 2:

In microtube rack, label 1.7 mL tubes for RNA (marked in pink) and DNA (marked in blue) standard solutions to be used as standard curves.

Tubes A is for standard solutions treated with working solution A (WS-A)

Tubes B is for standard solutions treated with working solution B (WS-B)

Tubes C is for standard solutions treated with working solution C (WS-C)

Tubes A	R1A	R2A	R3A	R4A	R5A	R6A	D1A	D2A	D3A	D4A	D5A
Tubes B	R1B	R2B	R3B	R4B	R5B	R6B	D1B	D2B	D3B	D4B	D5B
Tubes C	R1C	R2C	R3C	R4C	R5C	R6C	D1C	D2C	D3C	D4C	D5C

Set 3:

In microtube rack, label 1.7 mL tubes for diluted samples and blanks.

Tubes A is for diluted samples and blanks treated with working solution A (WS-A)

Tubes B is for diluted samples and blanks treated with working solution B (WS-B)

Tubes C is for diluted samples and blanks treated with working solution C (WS-C)

Tubes A	1A	2A	3A	4A	5A	6A	7A	8A	9A	 BlankA
Tubes B	1B	2B	3B	4B	5B	6B	7B	8B	9B	 BlankB
Tubes C	1C	2C	3C	4C	5C	6C	7C	8C	9C	 BlankC



31 Label tubes for reagents as following.
Follow the sheet, add Tris buffer ([M]5 mM , pF8.0) to the reagent tubes:

Α	В
Reagent type	5 mM Tris (uL)
0.05% STEB	9X1000 + 500
RNA	990+495
DNA	998
RNase	380
900 mM MgCl2	40
900 mM CaCl2	40
WS-A	5X1000 + 640
WS-B	5X1000 +640
WS-C	5X1000 + 880
SG-II WS	2X1000 + 500

32 Add $\square 900 \ \mu L$ Tris buffer ([M]5 mM , p+8.0) to each tube in Set 1

Depending on the dilution of extracted sample

Follow the sheet, add Tris buffer ([M]5 mM , p+8.0) to each tube in Set 2. The unit of volume is uL.

Tubes A										112027	AR ASSESSED
Tubes B	650	640	625	600	550	500	640	630	610	580	550
Tubes C	600	590	575	550	500	450	590	580	560	530	500

Follow the sheet, add Tris buffer ([M]5 mM, pF8.0) to each tube in Set 3. The unit of volume is ul.

Tubes A	650	650	650	650	650	650	650	650	650	 650
Tubes B	650	650	650	650	650	650	650	650	650	 650
Tubes C	600	600	600	600	600	600	600	600	600	 600

35 Prepare STEB ([M]**0.05** %)

Add $\blacksquare 500~\mu L$ STEB ([M]1 %) to 0.05% STEB tube, and vortex.

- 36 Add **250 μL** STEB ([M]**0.05** %) to each tube in **Set 2** by reverse pipetting.
- 37 Place RNase and DNase primary stock solutions, RNA and DNA primary standard solutions and samples & On ice.
- Turn on refrigerated centrifuge and set the temperature to $~8\,4\,^{\circ}\text{C}$.

CENTRIFUGE 5430 R
Eppendorf MP2231000510

Turn on shaker/incubator and set temperature to § 37 °C.

SHAKING INCUBATOR
71L
Corning® LSE™ 6753

- 40 Prepare [M]900 mM MgCl₂

 - 40.2 Transfer ■360 µL [M]1 M MgCl₂ solution into 900 mM MgCl₂ tube
- 41 Add 120 μL [M]900 mM MgCl₂ to WS-A and WS-B
- 42 Prepare [M]900 mM CaCl₂
 - 42.1 Pour M1 M CaCl₂ solution into 5 mL RNase free Falcon tube

 8 Calcium chloride solution Sigma

 Aldrich Catalog #21115-100ML
 - 42.2 Transfer \blacksquare 360 μ L [M]1 M CaCl₂ solution into 900 mM CaCl₂ tube
- 43 Add 120 μL [M]900 mM CaCl₂ to WS-A and WS-B

44

Lunch break!



45 Prepare RNase working solution [M]0.5 mg/ml Add 20 µL RNase primary stock solution ([M]10 mg/ml) to RNase tube 46 Add **120 μL** [M]**0.5 mg/ml** RNase to WS-B and WS-C. Keep WS-B and WS-C & On ice . 47 Add ■120 µL DNase primary stock solution ([M]5 mg/ml) to WS-A. Keep WS-A & On ice. Two microtubes of aliquot (70 ul/tube) 48 Centrifuge extracted samples @10000 x g, 4°C, 00:04:00 Take one tube of 49 SYBR™ Green II RNA Gel Stain, 10,000X concentrate in DMSO Thermo Fisher Catalog #S7564 out of the freezer, keep it at 8 Room temperature . If open a new package, wrap 1.7 mL microtube with foil and aliquot 1000 ul to each tube, store at § -20 °C. Prepare RNA secondary standard solution [M]2 ug/ml Add 15 µL RNA primary standard solution to RNA standard tube and vortex. Keep & On ice . 51 Prepare DNA secondary standard solution [M]1 ug/ml Add 22 µL DNA primary standard solution to DNA standard tube and vortex. Keep & On ice. 52 Load ⊒50 µL WS-A to Tubes A in Set 2 and Set 3. From Step 51 to 54: Reverse pipetting Decontaminate pipet between different WS. 53 Load $\Box 50 \mu L$ WS-A to Tubes C in **Set 2** and **Set 3**. 54 Load ⊒50 µL WS-B to Tubes B in Set 2 and Set 3. 55 Load **□50 µL** WS-C to Tubes C in **Set 2** and **Set 3**.

- 56 Add □100 μL centrifuged samples to its corresponding tubes in **Set 1**. Vortex each tube.
- 57 From Set 1, transfer **250 μL** of diluted samples to each corresponding tubes in **Set 3**.

In order to avoid cross contamination from RNase or DNase, use one tip for each dispensing. Pipette solution in the tube up and down for mixing.

58 Follow the sheet:

Add RNA secondary standard to tubes (marked in pink) in Set 2. Add DNA secondary standard to tubes (marked in blue) in Set 2. The unit of volume is uL.

Tubes A	0	10	25	50	100	150	10	20	40	70	100
Tubes B	0	10	25	50	100	150	10	20	40	70	100
Tubes C	0	10	25	50	100	150	10	20	40	70	100

Vortex each tube for @00:00:02 and place all tubes into the shaker/incubator at &37 °C, continuously shaking at 200 RPM for @00:20:00.

Incubation time is critical. Temperature might be disturbed by door open/close. Don't start the timer until temperature returns to 37° C.

Read fluorescence

- 60 Prepare SYBR Green II working solution (SG-II WS)
 - 60.1 Each 96-well microplate requires 1 mL of SG-II WS.

For 24 samples: 2.5 mL Tris 17.5 uL SG-II WS

- 60.2 Centrifuge one tube of SG-II concentrate at & Room temperature @13000 rpm, 00:05:00 to deposit DMSO.
- 60.3 Wrap SG-II WS tube with foil, transfer [M]17.5 ul supernatant of SYBR Green II 10,000X concentrate to SG-II WS tube in biosafety cabinet.

Any step involving SYBR Green II should be operated in dark room or at least dim light.

60.4 Adhere black film on the top of a microplate lid.

Black Vinyl Films for Fluorescence and Photoprotection

VWR

89087-692

Microplate Lids Polystyrene

Greiner Bio-One 07000288

60.5 Load $\Box 10~\mu L$ SG-II WS to each well in the microplate with 0.5 mL tip of stepper, and cover the plate with the black-film lid.

Finntip™ Stepper Pipette Tips

500 uL

Thermo Scientific™ 9404170

96-Well Black Microplates

Polystyrene

Greiner Bio-One 655076

Wipe or dab the liquid drop on the outside of the tip, avoid wiping the tip open before dispensing the liquid.

- After incubation, vortex each tube for © 00:00:02 and then place into the fridge to stop the reaction.
- 62 Allow samples to reach § Room temperature for © 00:02:00 before loading the microplate.

Since fluorescence decreases with increasing temperature, with percentage changes depending on the fluorophore (Bashford, 1987), the SG-II WS must be kept dark at RT (22°C) and the samples must be equilibrated at RT (c. 2 min).

63 Label the microplate.

64 Organize tubes in 96-well microtube rack in the same order as how microplates are loaded.

45m

65 Load **190 μL** working sample to the microplate by reverse pipetting.

Blank must be included in each plate.

Pink area: RNA standard solutions for RNA standard curves Blue area: DNA standard solutions for DNA standard curves

Yellow area: Samples and blanks

	1	2	3	4	5	6	7	8	9	10	11	12
Α	R1A	R2A	R3A	R4A	R5A	R6A	1A	1A	1B	1B	1C	1C
В	R1B	R2B	R3B	R4B	R5B	R6B	2A	2A	2B	2B	2C	2C
C	R1C	R2C	R3C	R4C	R5C	R6C	3A	3A	3B	3B	3C	3C
D	D1A	D2A	D3A	D4A	D5A	1	4A	4A	4B	4B	4C	4C
E	D1B	D2B	D3B	D4B	D5B		5A	5A	5B	5B	5C	5C
F	D1C	D2C	D3C	D4C	D5C		6A	6A	6B	6B	6C	6C
G	9A	9A	9B	9B	9C	9C	7A	7A	7B	7B	7C	7C
Н	Blank A	Blank A	Blank B	Blank B	Blank C	Blank C	8A	8A	8B	8B	8C	8C

	1	2	3	4	5	6	7	8	9	10	11	12
Α	Blank A	Blank A	Blank B	Blank B	Blank C	Blank C	17A	17A	17B	17B	17C	17C
В	10A	10A	10B	10B	10C	10C	18A	18A	18B	18B	18C	18C
С	11A	11A	11B	11B	11C	11C	19A	19A	19B	19B	19C	19C
D	12A	12A	12B	12B	12C	12C	20A	20A	20B	20B	20C	20C
Ε	13A	13A	13B	13B	13C	13C	21A	21A	21B	21B	21C	21C
F	14A	14A	14B	14B	14C	1C	22A	22A	22B	22B	22C	22C
G	15A	15A	15B	15B	15C	15C	23A	23A	23B	23B	23C	23C
Н	16A	16A	16B	16B	16C	16C	24A	24A	24B	24B	24C	24C

Wipe or dab the liquid drop on the outside of the tip, avoid wiping the tip open before dispensing the liquid.

66 Shake black film covered microplate at 8 Room temperature for © 00:10:00

10m

- (1) The fluorescence of the RNA-bound SG-II is highly stable during the 0 180 min period
- (2) The fluorescence of the DNA-bound SG-II varies during the first 10 min and remains stable during the 10-60 min period, but tended to decrease afterwards.
- (3) In summary, readings are performed within the 10-60 min period following the SG-II addition; meanwhile the samples are kept dark at RT (22°C).
- 67 Setup microplate reader:

Plate: Greiner F bottom chimney well PP 96 well;

Shake: Continuous 5s at 600 rpm

Endpoint reading: Ex 490 nm/Em 520 nm;

mprotocols.io

Varioskan LUX Multimode Microplate

Reader

Thermo Fisher VL0L00D0

68 Read fluorescence and export data to excel sheet.

Calculate

- 69 RNA standard curve
 - 69.1 Concentrations of RNA standards in the microplate

Standard	2 ug/mL (uL)	Tris+WS (uL)	0.05% STEB (uL)	SG II WS (uL)	Final (ng/mL)
R1	0	700	250	50	0
R2	10	690	250	50	20
R3	25	675	250	50	50
R4	50	650	250	50	100
R5	100	600	250	50	200
R6	150	550	250	50	300

- 69.2 Slope of fluorescence in Tube A vs concentration of RNA standard gives $m_{\text{RNA+DNase}}$ (\approx 0.03) Slope of fluorescence in Tube B vs concentration of RNA standard gives $m_{\text{RNA+RNase}}$
- 69.3 Calculate ρ

$$\rho = \frac{m_{RNA+RNase}}{m_{RNA+DNase}}$$

70 Total RNA of the samples

ug Total RNA (ml assay)-1

$$= 0.001 \times \frac{\left(RFU_A - RFU_{A_{Blank}}\right) - \left(RFU_C - RFU_{C_{Blank}}\right)}{(1 - \rho)} / m_{RNA + DNase}$$

Where,

 RFU_A and RFU_C are the fluorescence in Tube A and Tube C of the same sample.

 $\mathsf{RFU}_\mathsf{ABlank}$ and $\mathsf{RFU}_\mathsf{CBlank}$ are the fluorescence in Tube A and Tube $_\mathsf{C}$ of the blank.

- 71 DNA standard curve
 - 71.1 Concentrations of DNA standards in the microplate

Standard	1 ug/mL (uL)	Tris+WS (uL)	0.05% STEB (uL)	SG II WS (uL)	Final (ng/mL)
D1	10	690	250	50	10
D2	20	680	250	50	20
D3	40	660	250	50	40
D4	70	630	250	50	70
D5	100	580	250	50	100

- 71.2 Slope of fluorescence in Tube A vs concentration of DNA standard gives $m_{\text{DNA+DNase}}$ Slope of fluorescence in Tube B vs concentration of DNA standard gives $m_{\text{DNA+RNase}}$ (\approx 0.12)
- 71.3 Calculate δ

$$\delta = \frac{m_{DNA+DNase}}{m_{DNA+RNase}}$$

72 Total DNA of the samples

ug Total DNA (ml assay)-1

$$=0.001\times\frac{\left(RFU_B-RFU_{B_{Blank}}\right)-\left(RFU_C-RFU_{C_{Blank}}\right)}{(1-\delta)}/m_{DNA+RNase}$$

Where

 RFU_B and RFU_C are the fluorescence in Tube B and Tube C of the same sample

RFU_{BBlank} and RFU_{CBlank} are the fluorescence in Tube B and Tube _C of the blank.

73 Dilution factor=40

If,

- Sample is extracted by 1 mL extraction reagent
- In Set 1, sample is diluted to 1/10
- In Set 3, diluted by Tris and all working solutions to 250/950
- In microplate, diluted by SG-II WS to 190/200