



VERSION 1

FEB 27, 2023

OPEN  ACCESS

**Protocol Citation:** Annabel J Curle, Sarah Howlett, Joanne Jones 2023. MoDC immunogenicity assay. **protocols.io** <https://protocols.io/view/modc-immunogenicity-assay-cp3vvqn6>

**License:** This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

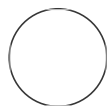
**Protocol status:** Working  
We use this protocol and it's working

**Created:** Feb 27, 2023**Last Modified:** Feb 27, 2023

**PROTOCOL integer ID:**  
77653

**Keywords:** Immunogenicity, T cell, In vitro, monocyte-derived DC

## MoDC immunogenicity assay V.1

Annabel J Curle<sup>1</sup>, Sarah Howlett<sup>1</sup>, Joanne Jones<sup>1</sup><sup>1</sup>University of Cambridge

Annabel J Curle  
University of Cambridge

### ABSTRACT

The monocyte-derived dendritic cell (moDC) immunogenicity assay can be utilised to test immunogenicity of different cell types. Similar to the mixed lymphocyte reaction (MLR) assay, cells of interest are co-cultured with T cells. Professional antigen presenting cells (APCs; moDCs) are added here to ensure optimal T cell activation in the case of immunogenic cells and improve limited sensitivity of the MLR-like assays.

Autologous CD14+ monocytes and T cells are isolated from PBMCs. Monocytes are differentiated to mature dendritic cells over 7 days, then co-cultured with T cells and cells of interest for 5 days. Recommended readouts of the assay are ELISA (for IFN $\gamma$ /TNF $\alpha$  release) and flow cytometry (for T cell proliferation and CD25 expression).

Positive controls wells of allogeneic moDC/T cells and bead-activated T cells are utilised. Biological replicates can be run in parallel. Well-well variation is minimal, but samples can be run in technical duplicates if desired.

## MATERIALS

### Cell preparations:

PBMC isolation materials

MACS equipment and LS tubes

CD14+ microbeads for MACS (we use Miltenyi CD14+ microbeads)

PanT cell isolation kit for MACS (we use Miltenyi PanT cell isolation kit)

96 well tissue culture plates

Cell proliferation dye (we use eBioscience™ Cell Proliferation Dye eFluor™ 450)

Cell freezing media (FBS + 10% DMSO)

### moDC differentiation

RPMI

Penicillin/streptomycin

Human AB serum

rhIL-4

rhGM-CSF

LPS

PBS

### Co-culture

Cell of interest media

T cell activation beads (we use Gibco Dynabeads™ Human T-Activator CD3/CD28 for T Cell Expansion and Activation)

### Readouts

ELISA materials

Flow cytometry materials

## BEFORE START INSTRUCTIONS

Fresh blood must be used each time as monocytes do not survive the freeze/thaw process

Optimisation experiments must be performed to ensure all cell types survive without compromise to survival, proliferation and differentiation state in the 1:1 mix of medias (T cell/moDC media: cell of interest media)

## Cell preparations and moDC differentiation - Day 0-7

1w 0d 3h

- 1 Make PBMCs using preferred method from whole blood
- 2 Take whole PBMCs and enrich for CD14+ cells using magnetic activated cell sorting (MACS) with CD14+ microbeads

- 3 Count monocytes, then seed in a 96 well plate at 150,000 cells/well in 200 ul RPMI + 10% P/S + 5% human AB serum (HAB) + 50 ng/ml rhIL-4 and 50 ng/ml rhGM-CSF (day 0)

Incubate in 37°C incubator for 7 days refreshing half media on day 3 and 6

- 3.1 Each media refresh should be done with double concentration IL-4 and GM-CSF to account for half media change: remove 100 ul media, add 100 ul RPMI + 10% P/S + 5% HAB + 100 ng/ml rhIL-4 and 100 ng/ml rhGM-CSF

Day 6 media change, plus 100 ng/ml LPS for activation : remove 100 ul media, add 100 ul RPMI + 10% P/S + 5% HAB + 100 ng/ml rhIL-4 and 100 ng/ml rhGM-CSF + 200 ng/ml LPS

- 3.2 **Optional:** collect monocytes (day 0), immature moDCs (day 6) and mature moDCs (day 7) to check purity of MACS enrichment by flow cytometry

Monocytes: CD14+

Immature DCs: CD14-, CD80+, CD86+, HLA-ABC+, HLA-DR+

Mature DCs: CD14-, CD80+, CD86+, HLA-ABC+, HLA-DR+, CD83+

- 4 Using the CD14- population, enrich for T cells using MACS with a PanT negative selection kit

- 4.1 If specific T cell population is required, perform MACS or FACS (fluorescence activated cell sorting) for required population

- 4.2 **Optional:** collect PanT to check purity of MACS enrichment by flow cytometry

- 5 Stain T cells with proliferation dye as per preferred protocol

- 6 Count then freeze the stained-T cells in FBS + 10% DMSO for the duration of the monocyte-DC differentiation

- 7 Prepare your "cells of interest" (COI) ready to have enough for 150,000 cells per well on day 7

## Seeding the immunogenicity assay - Day 7 (assay Day 0) 2h

- 8 Thaw proliferation dyed-T cells as per preferred protocol
- 9 Remove all media from moDCs and wash once with PBS (moDCs are adherent but aspirate gently to avoid disturbing cells)
- 10 **Experimental wells:**  
Seed 150,000 autologous T cells per well in 150 ul RPMI + 10% P/S + 5% HAB  
Seed 150,000 COI per well with 150 ul media specific to cells of interest (+ ROCKi if required)  
  
Final: 1:1:1 ratio of moDC:T cell:COI in 300 ul of 50:50 media RPMI media: COI media
- 10.1 **Note:** optimisation will be required to ensure both cell types survive and proliferate as expected in the 1:1 media mix
- 11 **Positive control wells:**  
  
(1) *Allogeneic Tcell-moDC*  
Seed 150,000 allogeneic T cells into moDC wells in 150 ul RPMI + 10% P/S + 5% HAB  
Add 150 ul media specific to cells of interest  
  
(2) *Polyclonally activated T cells*  
Seed 150,000 T cells into an empty well (no moDC) in 150 ul RPMI + 10% P/S + 5% HAB  
Add activation beads at recommended ratio in 150 ul media specific to cells of interest

## Immunogenicity assay - Day 7-12 (assay Day 0-5) 5d

- 12 Incubate co-culture for 5 days, refreshing half media every 48 hours (day 2 and day 4)  
  
Each change: remove 150 ul media and add 150 ul prepared 1:1 media

### 13 For cytokine assays:

Collect 150 ul supernatant per well into a v-bottom 96 well plate

**13.1** Centrifuge supernatants plate 1500 g x 10 minutes then re-collect supernatants avoiding any pellet/debris that have collected

**13.2** Freeze supernatants at -80°C until use.

Perform ELISA or luminex cytokine assays including IFN $\gamma$  and TNF $\alpha$  as readouts of immunogenicity.

### 14 For flow cytometry:

Collect remaining supernatants containing non-adherent T cells by agitating gently, washing the well 1-2 times with the supernatant.

**Note:** moDCs are adherent so will not be collected, consider that COI may be adherent or non-adherent (it is ok if some are collected as they will be gated out during flow cytometry).

**14.1** Perform flow cytometry immediately, including a minimum panel of: CD3, CD4 and CD8 (T cell markers), CD25 (activated T cell marker) and live/dead stain.

**Note:** cells are already pre-stained with proliferation dye

Activated T cells (CD25+, proliferation dye low/negative) are used as readout of immunogenicity of COI