



OCT 20, 2023

OPEN ACCESS



**DOI:**  
[dx.doi.org/10.17504/protocols.io.bp2l6x7zdlqe/v1](https://dx.doi.org/10.17504/protocols.io.bp2l6x7zdlqe/v1)

**Protocol Citation:** J Bons, J P Rose, M A Watson, B Schilling 2023. Protein Digestion with S-trap Spin Columns. **protocols.io** <https://dx.doi.org/10.17504/protocols.io.bp2l6x7zdlqe/v1>

**License:** This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

**Protocol status:** Working  
 We use this protocol and it's working


**Created:** Oct 10, 2023

## Protein Digestion with S-trap Spin Columns

M A

J Bons<sup>1</sup>, J P Rose<sup>1</sup>, Watson<sup>1</sup>, B Schilling<sup>1</sup>

<sup>1</sup>Buck Institute for Research on Aging

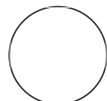

M A Watson

### ABSTRACT

Trypsin digestion of isolated proteins using S-trap Spin columns in preparation for downstream proteomic profiling.

### MATERIALS

- S-trap spin columns (Protifi)
- Centrifuge
- Centrifugal vacuum concentrator
- Small hot box/incubator
- HPLC-grade water
- 10% SDS solution
- 1 M triethylammonium bicarbonate (TEAB) solution, pH 8
- 100 mM triethylammonium bicarbonate (TEAB) solution, pH 8 in water
- Dithiothreitol (DTT; 250 mM in 100 mM TEAB, pH 8)

**Last Modified:** Oct 20, 2023

**PROTOCOL integer ID:**  
89096

**Keywords:** Trypsin  
Digestion, S-trap spin  
columns, Proteomics, Mass  
Spectrometry

- Iodoacetamide (IAA; 250 mM in 100 mM TEAB, pH 8)
- S-trap buffer (90% methanol in 100 mM TEAB)
- Sequencing-grade trypsin
- 12% phosphoric acid in water
- Digestion buffer/Elution buffer 1: 50 mM triethylammonium bicarbonate (TEAB) solution, pH 8 in water
- Elution buffer 2: 0.5% formic acid (FA) in water
- Elution buffer 3: 50% acetonitrile (ACN), 0.5% formic acid (FA) in water
- 0.2% formic acid (FA) in water

- 1 In a 2-mL microcentrifuge tube add up to 300 µg of the protein lysate, 10% SDS for a final concentration of 4% SDS, 1M TEAB pH 8 solution for a final concentration of 50 mM TEAB, and HPLC-grade water if necessary to bring the volume up to a minimum of 50 µL.
- 2 Add 250 mM DTT for a final concentration of 20 mM and incubate for 10min at 50°C to reduce the proteins. Then immediately leave for 10 min on the bench at room temperature (RT).
- 3 Add 250 mM IAA for a final concentration of 40 mM and incubate for 30min at RT in the dark to alkylate the proteins.

- 4 Acidify the sample with 12% phosphoric acid for a final concentration of 1.2%.
- 5 Add 7 volumes of S-Trap buffer to the acidified lysate and mix immediately by inversion. Formation of protein colloid may be observed.
- 6 Use S-trap micro spin columns for  $\leq 100$   $\mu\text{g}$  of protein or S-trap mini spin columns for 100-300  $\mu\text{g}$  of protein. Ensure that the S-Trap spin column is in a 2.0-mL flow-through catch tube.
- 7 Add 100  $\mu\text{L}$  (S-trap micro spin column) or 200  $\mu\text{L}$  (S-trap mini spin column) of the acidified lysate/S-Trap buffer mix into the S-Trap spin column. Centrifuge at  $4,000 \times g$  for 10 seconds or until all the solution has passed through the column. Discard the flow-through.
- 8 Repeat step 7 until the entire acidified lysate and S-Trap buffer mix has passed through the column.
- 9 Add 200  $\mu\text{L}$  (S-trap micro spin column) or 400  $\mu\text{L}$  (S-trap mini spin column) of S-Trap buffer to wash the column. Centrifuge at  $4,000 \times g$  for 10 seconds or until all solution has passed through the column. Discard the wash solution.
- 10 Add 200  $\mu\text{L}$  (S-trap micro spin column) or 400  $\mu\text{L}$  (S-trap mini spin column) of S-Trap buffer and set aside.
- 11 Prepare a 1:25 (wt/wt) trypsin/protein ratio solution of sequencing-grade trypsin.
- 12 Centrifuge the S-trap spin column at  $4,000 \times g$  for 10 seconds or until the column is dry.

- 13** Place the S-Trap spin column into a clean 2.0-mL elution tube.
- 14** Add 125  $\mu\text{L}$  of the sequencing-grade trypsin solution to the column at a 1:25 (wt/wt) trypsin/protein ratio. Be careful not to pierce the trap material.
- 15** Loosely cap the S-trap micro spin column or close the S-trap mini spin column and incubate for 1 hour at 47°C with no agitation.
- 16** After 1 hour, add another 125  $\mu\text{L}$  of the sequencing-grade trypsin solution to the column and incubate overnight at 37°C with no agitation.
- 17** After overnight incubation take the samples out of the incubator and elute sequentially in the same 2.0-mL elution tube as follows:
  - 17.1** a. Add 80  $\mu\text{L}$  of elution buffer 1 (50 mM TEAB, pH 8) and centrifuge at 1,000 x *g* for 1 minute.
  - 17.2** b. Add 80  $\mu\text{L}$  of elution buffer 2 (0.5% FA) and centrifuge at 1,000 x *g* for 1 minute.
  - 17.3** c. Add 80  $\mu\text{L}$  of elution buffer 3 (50% ACN, 0.5% FA) and centrifuge at 4,000 x *g* for 1 minute.

- 18 Vacuum dry the eluted peptide solution in a centrifugal vacuum concentrator.
- 19 Reconstitute the dried peptides in 0.2% FA and thoroughly mix the solution before desalting.