



HEK293T Landing Pad Recombination Protocol with Fugene (Based on a 24-well plate) V.1


Sarah Roelle¹, Kenneth Matreyek¹

¹Case Western Reserve University

Version 1

Jun 28, 2021

1 Works for me

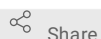

Share

dx.doi.org/10.17504/protocols.io.bcdeis3e

MatreyekLab

Sarah Roelle

ABSTRACT

For using Fugene 6 to recombine Bxb1 attB plasmids into Bxb1 attP Landing Pads already integrated into HEK 293T cells.

DOI

dx.doi.org/10.17504/protocols.io.bcdeis3e

PROTOCOL CITATION

Sarah Roelle, Kenneth Matreyek 2021. HEK293T Landing Pad Recombination Protocol with Fugene (Based on a 24-well plate). **protocols.io**
<https://dx.doi.org/10.17504/protocols.io.bcdeis3e>

LICENSE

This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

CREATED

Feb 10, 2020

LAST MODIFIED

Jun 28, 2021

PROTOCOL INTEGER ID

32902

SAFETY WARNINGS

Engineering Controls: use of biosafety cabinet
PPE: lab coat and gloves

Day 0: Day of transfection 5m

- 1 Trypsinize and count the cells. 5m
- 2 Make 2 transfection mixtures per sample, 5m
Opti1 and **Opti2** (DNA volumes should be kept between 1.0 and 5.0 uL, if possible):

2.1 Opti1: 23.8 µL Opti-MEM + 1 µL Fugene6

2.2 For the Bxb1-containing sample - Opti2: 23 µL Opti-MEM + µL corresponding to 16 ng Bxb1 expression plasmid + µL corresponding to 238 ng attB plasmid ---
OR----

2.3 For the No-Bxb1 sample - Opti2-noB: 23 µL Opti-MEM + µL corresponding to 254 ng attB plasmid

3 Add 300 µL media to each well.

(Be sure media contains doxycycline to a final concentration of 2 µg / mL .)

4 Plate 120,000 cells per well within a 24-well plate.

5 Once at least 15 minutes have past since the mixtures were completed, add them dropwise throughout the well. Rock the plate back and forth to mix it around the well. 15m

Day 1: Adding more media to cells 5m

6 Add at least 500 µL media to each well.

5m

Day 3+: Negative and/or positive selection 5m

7 If using a landing pad that encodes iCasp9, then add AP1903 to a final concentration of 10 nM to induce apoptosis of un-recombined cells. Dying cells can be observed as soon as 2 hours after AP1903 addition, and can be removed by changing the media. 5m

If you want to estimate how many cells were recombined (this is recommended, especially for library -based experiments), make sure you leave some cells WITHOUT selection, so you can see what percent of the cells are reporter (eg. mCherry) positive. Then, you can multiply that percentage with the total number of cells you recombined (eg. 120,000 cells if only recombining a single well) to get the total number of recombined cells.

8 If using an attB vector encoding a positive selection marker, add the small molecule. For puromycin resistant recombination plasmids (this is most often the case), add puromycin to a final concentration of 1 µg / mL . 5m

If performing both negative AND positive selection, we typically perform the negative selection first, and only perform the positive selection after the well has reached greater than 20% confluency.

Day 7+: Read out the recombination

1h

- 9 Run the samples through the flow cytometer. Watch the cells on the days in between. They may need to be split at ^{1h} some point. It's probably easiest to simply plate all the cells into a 6-well plate when they need to be split.