



Aug 13, 2020

Target Guide Sequence Cloning Protocol Version 2

Forked from [Target Guide Sequence Cloning Protocol](#)Skye Waterland¹, Yang Li¹¹Baylor College of Medicine1 Works for me dx.doi.org/10.17504/protocols.io.bjp3kmaq

Skye Waterland

ABSTRACT

Create single gRNA vectors for targeted cloning utilizing CRISPR or CRISPR-based systems.

ATTACHMENTS

[Lentivirus_Protocol.pdf](#) [Addgene_Protocol - Bacterial Transformation.pdf](#)

DOI

dx.doi.org/10.17504/protocols.io.bjp3kmaq

PROTOCOL CITATION

Skye Waterland, Yang Li 2020. Target Guide Sequence Cloning Protocol Version 2. **protocols.io**
<https://dx.doi.org/10.17504/protocols.io.bjp3kmaq>

FORK FROM

Forked from [Target Guide Sequence Cloning Protocol](#), Skye Waterland

KEYWORDS

Lentivirus vector, cloning, vector digestion, oligo annealing, CRISPR

LICENSE

————— This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

IMAGE ATTRIBUTION

<https://www.addgene.org/52963/>

CREATED

Aug 13, 2020

LAST MODIFIED

Aug 13, 2020

PROTOCOL INTEGER ID

40411

MATERIALS

NAME	CATALOG #	VENDOR
NEBuffer 3.1 - 5.0 ml	B7203S	New England Biolabs
T4 DNA Ligase - 20,000 units	M0202S	New England Biolabs
Agar		
lentiGuide-Puro	52963	addgene

NAME	CATALOG #	VENDOR
double distilled water (ddH2O)		
SOC Media		
1X TAE Buffer		
10X NEB T4 DNA ligase buffer		New England Biolabs
10X T4 PNK Reaction Buffer		New England Biolabs
ethanol		
10X PCR Buffer	10966-034	Life Technologies
LB-Broth Miller (= LB mix)	LMM0104	Formedium
One Shot™ TOP10 Chemically Competent <i>E. coli</i>	C404010	Thermo Fisher
BsmBI-v2	R0739L	New England Biolabs
HotStarTaq Plus DNA Polymerase (1000)	203605	Qiagen
dNTP Set (100mM each A C G T)	95038-256	Ge Healthcare

STEPS MATERIALS

NAME	CATALOG #	VENDOR
ddH2O		
NEBuffer 3.1 - 5.0 ml	B7203S	New England Biolabs
lentiGuide-Puro	52963	addgene
double distilled water (ddH2O)		
10X NEB T4 DNA ligase buffer		New England Biolabs
10X T4 PNK Reaction Buffer		New England Biolabs
BsmBI-v2	R0739L	New England Biolabs
Agar		
1X TAE Buffer		
T4 DNA Ligase - 20,000 units	M0202S	New England Biolabs
One Shot™ TOP10 Chemically Competent <i>E. coli</i>	C404010	Thermo Fisher
ethanol		
10X PCR Buffer	10966-034	Life Technologies
dNTP Set (100mM each A C G T)	95038-256	Ge Healthcare
HotStarTaq Plus DNA Polymerase (1000)	203605	Qiagen
LB-Broth Miller (= LB mix)	LMM0104	Formedium
SOC Media		

MATERIALS TEXT

lentiGuide-Puro: RRID:Addgene_52963
Sigma-Aldrich: RRID:SCR_008988

EQUIPMENT

NAME	CATALOG #	VENDOR
Snap Cap Microcentrifuge Tube or equivalent	07200210	
Mini-centrifuge	S67601B	
Oven	15-103-0510	
SimpliAmp Thermal Cycler	A24811	

DISCLAIMER:

This protocol is a modified version of the Zhang Lab's *GeCKOv2* Target Guide Sequence Cloning Protocol attached below based off of Joung, J., Konermann, S., Gootenberg, J. *et al.* Genome-scale CRISPR-Cas9 knockout and

More information about the specific lentiGuide-puro plasmid can be found here: <https://www.addgene.org/52963/>.

Design and order gRNA oligos from [Sigma-Aldrich](#) (RRID:SCR_008988).



ddH₂O

to a  1.5 mL



Snap Cap Microcentrifuge Tube or
equivalent
Polypropylene Microcentrifuge Tube
Corning Costar Snap Cap Microcentrifuge
Tube
2 mL snap cap polypropylene micro tube



07200210 [↗](#)

1.2 Add  5 µl of



NEBuffer 3.1 - 5.0 ml
by New England Biolabs
Catalog #: [B7203S](#)

to solution.

1.3 Add  3 µl of



BsmBI-v2
by New England Biolabs
Catalog #: [R0739L](#)

to solution.



Note: we used BsmBI Cat #R0580, NEB, but this one has been discontinued. #R0739L is considered as an effective replacement.

1.4 Add  2 µl of



lentiGuide-Puro

by addgene

Catalog #: 52963

to solution.

1.5 Close cap on microcentrifuge tube and place in



Mini-centrifuge Centrifuge

Fisher S67601B [↗](#)

Any standard mini centrifuge with adapters for different tube sizes will suffice



for [🕒 00:00:10](#) on until all of the solution is at the bottom of the tube.

1.6 Place microcentrifuge tube with vector digestion mixture in a



Oven

Oven forced-air convection

Fisher Isotemp 15-103-0510 [↗](#)



on [🌡 55 °C](#)

1.7 Close lid and set a timer for [🕒 01:00:00](#)

Gel purify the digested plasmid from Step 1

2 Prepare gel

2.1 Create gel concentration of 1.2-1.5%



Agar

for  **100 mL** of



1X TAE Buffer

solution.

3 Run gel

3.1 Isolate 2kb and 8kb band. Collect 8kb (8318 bp) band for gel purification.

Preparing the gRNAs

4 Design gRNA sequence for CRISPR strategy using [CRISPR direct](#).


5 Order oligos from [Sigma-Aldrich](#) (RRID:SCR_008988).

6 Create a dilution from stock oligos at a 1:10 ratio in




double distilled water (ddH₂O)



Diluted oligo concentration should be  **10 Micromolar (μM)** .

Phosphorylate and anneal each pair of oligos

7 Prepare phosphorylation/annealing reaction

7.1 Add  **6.5 μl**



double distilled water (ddH₂O)

to a microcentrifuge tube.

7.2 Add  **1 µl** each of **10 Micromolar (µM)** Oligo 1 (F), **10 Micromolar (µM)** Oligo 2 (R), and



10X NEB T4 DNA ligase buffer
by New England Biolabs

7.3

Add  **0.5 µl**



10X T4 PNK Reaction Buffer
by New England Biolabs

7.4 Vortex and microcentrifuge

8 Place the phosphorylation/annealing reaction in a



SimpliAmp Thermal Cycler
PCR

Applied Biosystems A24811 [↗](#)

Any standard PCR thermocycler will suffice



Settings should be as follows:  **37 °C** for  **00:30:00** ,  **95 °C** for  **00:05:00** , and then ramp down to  **25 °C** at  **5 °C** /  **00:01:00** .

Setting up and incubating the ligation reaction

9 Making the ligation reaction

9.1 Place  **4.8 µl** of



double distilled water (ddH₂O)

in a microcentrifuge tube.

9.2 Add  **2.2 µl**



BsmBI-v2

by New England Biolabs

Catalog #: R0739L

9.3 Add  **1 µl** each of



10X NEB T4 DNA ligase buffer

by New England Biolabs

, diluted oligo duplex from [go to step #8](#), and





T4 DNA Ligase - 20,000 units

by New England Biolabs

Catalog #: M0202S


9.4 Lightly vortex and microcentrifuge

10 Incubate the ligation reaction at room temperature for  **02:00:00** -  **03:00:00**

Transformation into E. coli bacteria

11 Prepare LB Agar plates

11.1 Wipe down your bench with at least 70%

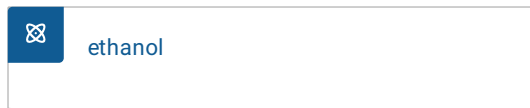


and light a bunsen burner.

11.2 Remove Agar Ampicillin Plates 250 µl from 4 °C and let warm up to room temperature.


12 E. coli competent cell transfection

12.1 Take competent cells

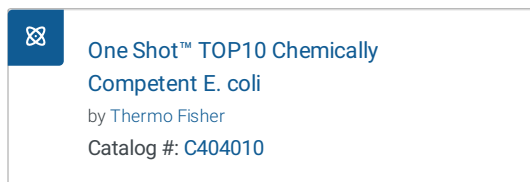

out of  -80 °C and thaw on ice ( 00:20:00 -  00:30:00).

12.2 Add 100 µl of E.coli cells to 10 µl of DNA in a microcentrifuge tube next to the bunsen burner.

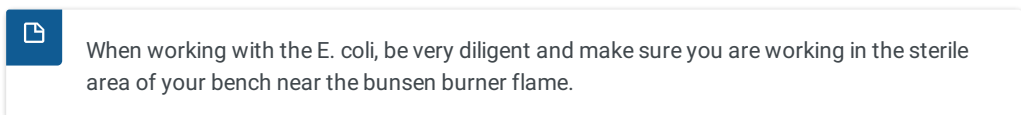

12.3 Gently flick tube a few times with your finger to mix.

12.4 Incubate the competent cell/DNA mixture on ice for 00:30:00 .

12.5 Place transformation tube(s) into water bath at 42 °C for 00:00:30 - 00:01:00 to heat shock E. coli cells.

12.6 Place the transformation tube(s) back on ice for 00:02:00 .

12.7 Add 250 µl



LB-Broth Miller (= LB mix)
by Formedium
Catalog #: [LMM0104](#)

(without antibiotic) or



SOC Media
[View](#)


to the tube(s).

12.8 Place tube(s) in 37 °C shaking incubator for 00:45:00 - 01:00:00 .


12.9 Plate all of transformation onto LB agar plate(s) with ampicillin. Incubate plates at 37 °C overnight.

Picking Colony for Suspension Growth

13 Before you get started




Start this whole section (picking colony and suspension growth) sometime in the late afternoon (~4-5pm), so that you can run PCR for the bacterial plasmids you collect the next morning (~16 hours later). More than 24 hours of incubation can cause other non-ampicillin resistant bacteria to grow in the suspension tubes.




Set up an aseptic area for handling the bacterial colonies.

14 Picking colony


14.1 Pick 10-15 colonies from each agar plate to suspend in an LB solution.



We originally selected only 4 colonies from each plate, but did not have a successful PCR. To increase chances of successfully amplifying the plasmid vector, we suggest picking 10-15 colonies.

14.2 Place each colony in a tube with  **3 mL** of LB.

15 Culture the bacteria

15.1 Place tubes into shaking incubator at  **37 °C** overnight.

Run PCR to Identify Positive Clones

16 Prepare Master PCR Mix by adding each of the below reagents to a microcentrifuge tube.



Before you start: multiply each volume below by the same amount (x 10, 15, etc.) according to how much master mix you need to run your sets.

16.1

 **15.875 µl**



ddH₂O

16.2

 **2.5 µl**



10X PCR Buffer

by Life Technologies

Catalog #: 10966-034

16.3

 **0.5 µl** ^[M] **10 Milimolar (mM)**



dNTP Set (100mM each A C G T)

by Ge Healthcare

Catalog #: 95038-256

16.4  2 µl  10 Micromolar (µM) Forward Primer hU6-02



Forward Primer hU6-02 sequence: TAATTAGAATTAATTTGACT
Ordered from [Sigma-Aldrich](#) (RRID:SCR_008988).

16.5  2 µl  10 Micromolar (µM) Reverse Primer (gRNA reverse primer oligo)


16.6  0.125 µl



HotStarTaq Plus DNA Polymerase
(1000)
by Qiagen
Catalog #: 203605

17 Loading and Prepping PCR wells

17.1 Aliquot equal proportions of master mix to each PCR well (we used  23 µl).

17.2 Add  2 µl of each bacterial suspension sample from [go to step #15](#) to each PCR well.



Total volume of Master Mix and Sample should be  25 µl in each well.

17.3 Seal off PCR wells tightly with a clear plastic cover or tube tops.

17.4 Centrifuge  1200 rpm 00:01:00 .

18 Place PCR reaction in thermocycler and run PCR.



Settings should be as follows:

- 1) 95C for 5mins (ramp up).
- 2) 95C for 30sec.
- 3) 55C for 45 sec.
- 4) 72C for 45 sec.
- x35 cycles steps 2-4.
- 5) 72 for 10mins.

Running Gel for PCR Product Verification

- 19** Prepare a gel of 2% concentration. Run gel and examine bands. Desired band length is about 200bp with the gRNA insertion.

Congrats!

- 20** You have successfully transformed a lentiviral vector with your gRNA sequence of interest! For confirmation, feel free to sequence your vector.