

MAR 27, 2024

## 🌐 CMT-93 Cell Culture Protocol

Laura Gómez<sup>1</sup>

<sup>1</sup>CBM Severo Ochoa



Laura Gómez  
CBM Severo Ochoa

### ABSTRACT

CMT-93 is a cell line exhibiting epithelial morphology that was isolated from the rectum of a mouse with polyploid carcinoma.

### PROTOCOL REFERENCES

<https://www.atcc.org/products/ccl-223#detailed-product-information>

### IMAGE ATTRIBUTION

ATCC - <https://www.atcc.org/products/ccl-223>

### GUIDELINES

Working with cell cultures requires a laminar flow cabinet. It has to be radiated with UV light, cleaned with any highly effective terminal disinfectant (such as Tego® 2000 or Suredis®) and 70% ethanol. All material introduced into the cabinet must also be sprayed with ethanol.

Once the work is finished, we must clean the cabinet with the detergent, then with 70% etOH and turn on the UV light for 30 min.

OPEN  ACCESS



#### DOI:

[dx.doi.org/10.17504/protocols.io.bp2l62k1dgqe/v1](https://dx.doi.org/10.17504/protocols.io.bp2l62k1dgqe/v1)

#### External link:

<https://www.atcc.org/products/ccl-223>

**Protocol Citation:** Laura Gómez 2024. CMT-93 Cell Culture Protocol. **protocols.io**  
<https://dx.doi.org/10.17504/protocols.io.bp2l62k1dgqe/v1>

#### MANUSCRIPT CITATION:

**License:** This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

**Protocol status:** In development  
We are still developing and optimizing this protocol

**Created:** Mar 18, 2024

**Last Modified:** Mar 27, 2024

**PROTOCOL integer ID:** 96872

**Keywords:** CMT-93, cell line, epithelial cell line

## MATERIALS

### Plasticware:

p60 cell culture plates

p100 cell culture plates

Cell culture flasks, 75 cm<sup>2</sup>, treated for cell attachment.

15 and 50 mL centrifuge tubes

Cryovials

### To prepare the complete medium:

DMEM 1X

Fetal bovine serum heat inactivated (FBS)


Glutamine 200 mM

### To subculture the cells:

Trypsin-EDTA

PBS 1X





## SAFETY WARNINGS

 Every reagent must be sterile in order to avoid contaminations.

## BEFORE START INSTRUCTIONS

Clean and prepare the laminar flow cabinet, turn on the water bath and warm up the culture media.


## Preparation of complete growth medium (DMEM+)

- 1 Add  445 mL 1X DMEM,  50 mL FBS and  5 mL glutamine  200 millimolar (mM) to a sterile 500 mL bottle and homogenize
- 2 Label the bottle with name, group, phone number, date and additions.


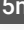

3 Close with parafilm and store at  4 °C .

## Cell thawing procedure



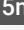
4 Remove one vial of cell stock from the liquid nitrogen tank with gloves and forceps. Transfer them to the cell culture laboratory in an appropriate container or a box with ice.

5 Thaw the vial by gently shaking it in a  37 °C water bath. Thawing should be rapid (approximately 2 min).

6 Remove the vial from the water bath as soon as the contents are thawed, and decontaminate by dipping in or spraying with 70% ethanol




7 Transfer the contents of the vial to a centrifuge tube containing  9 mL of complete culture medium and  5m  
 1200 rpm, Room temperature, 00:05:00

8 Resuspend with  10 mL DMEM+ and distribute on 2 P60 plates.

9 Incubate cultures at  37 °C , 5% CO<sub>2</sub>  Overnight  5m

## Subculturing procedure

- 10 Remove and discard culture medium.
- 11 Rinse with PBS 1X solution and discard
- 12 Briefly rinse the cell layer with 0.25% (w/v) Trypsin-0.53mM EDTA solution to remove all traces of serum which contains trypsin inhibitor. Discard.



- 13 Add  2 mL of Trypsin-EDTA solution to flask and incubate  00:10:00 at  37 °C to facilitate detachment from the plate. 10m


#### Note

To avoid clumping, **do not agitate** the cells by hitting or shaking the flask while waiting for the cells to detach.

- 14 Observe cells under an inverted microscope until cell layer is dispersed.

#### Note

If the cells are not detached already, incubate  00:05:00 more at  37 °C .

- 15 Add  5 mL of DMEM+ and aspirate cells by **gently** pipetting. Pour the existing volume down the walls of the flask in order to drag and collect as many cells as possible.

- 16 Collect the cell suspension in a centrifuge tube and  1200 rpm, Room temperature, 00:05:00 . 5m

17 Discard the supernatant into a beaker with 70% EtOH or 10% bleach.

18 Resuspend in medium according to the dilution to be made.

#### Note

A subcultivation ratio of 1:4 to 1:10 is recommended

19 Add  1 mL of the cell suspension to new culture vessels containing  14 mL DMEM+.


20 Incubate cultures at  37 °C , 5% CO<sub>2</sub>  Overnight


5m




## Cryopreservation and storage procedure

1d

21 Repeat steps of *Subculturing procedure* until the "Resuspend in medium according to the dilution to be made" step.

22 Resuspend in medium taking into account that for every p100 we can storage up to 2 cryovials of cells, containing  1 mL .

23 Prepare the cryovials with  50 µL DMSO.

- 24 Add  950  $\mu\text{L}$  of the cell suspension to every cryovial.
- 25 Label the cryovials with cell line, passage, date and lab number or phone number.
- 26 Store the cryovials in a slow freezing container at   $-80\text{ }^{\circ}\text{C}$  for  24:00:00 .
- 27 Transfer the cryovials to the liquid nitrogen tank.

1d