



Version 2 ▾

Oct 15, 2020

🌐 BGISEQ-500 WGS library construction V.2

Jie Huang¹, Xinming Liang², Yuankai Xuan³, Chunyu Geng², Yuxiang Li², Haorong Lu², Shoufang Qu¹, Xianglin Mei³, Hongbo Chen¹, Ting Yu¹, Nan Sun¹, Junhua Rao², Jiahao Wang⁴, Wenwei Zhang², Ying Chen², Sha Liao², Hui Jiang², Xin Liu², Zhaopeng Yang¹, Feng Mu², Shangxian Gao¹

¹National Institutes for food and drug control (NIFDC); ²BGI-Shenzhen;

³State Food and Drug Administration Hubei Center for Medical Equipment Quality Supervision and Testing; ⁴BGI-Qingdao

1

Works for me

This protocol is published without a DOI.

BGI



Hongling Zhou

ABSTRACT

BGISEQ-500 is a desktop sequencer developed by BGI. Using DNA nanoball and combinational probe anchor synthesis developed from Complete Genomics™ sequencing technologies, it generates short reads at a large scale. Library construction on the platform includes fragmentation, size selection, end repair and A-tailing, adaptor ligation, PCR amplification, and splint circularization.

EXTERNAL LINK

<https://doi.org/10.1093/gigascience/giy107>

EXTERNAL LINK

<https://doi.org/10.1093/gigascience/giy107>

PROTOCOL CITATION

Jie Huang, Xinming Liang, Yuankai Xuan, Chunyu Geng, Yuxiang Li, Haorong Lu, Shoufang Qu, Xianglin Mei, Hongbo Chen, Ting Yu, Nan Sun, Junhua Rao, Jiahao Wang, Wenwei Zhang, Ying Chen, Sha Liao, Hui Jiang, Xin Liu, Zhaopeng Yang, Feng Mu, Shangxian Gao 2020. BGISEQ-500 WGS library construction.

protocols.io

<https://protocols.io/view/bgiseq-500-wgs-library-construction-bngwmbxe>

Version created by Hongling Zhou



EXTERNAL LINK

<https://doi.org/10.1093/gigascience/giy107>

LICENSE

———— This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

CREATED

Oct 15, 2020

LAST MODIFIED

Oct 15, 2020

PROTOCOL INTEGER ID

43254

MATERIALS

NAME	CATALOG #	VENDOR
Fresh 80% ethanol		XILONG SCIENTIFIC
ERAT Buffer		

NAME	CATALOG #	VENDOR
ERAT Enzyme		
Adapter Mix		
Ligation Buffer		
Ligation Enzyme		
TE buffer		Ambion
PCR Enzyme Mix		
Splint Buffer		
Digestion Buffer		
Digestion Enzyme		

STEPS MATERIALS


NAME	CATALOG #	VENDOR
Fresh 80% ethanol		XILONG SCIENTIFIC
ERAT Buffer		
ERAT Enzyme		
Adapter Mix		
Ligation Buffer		
Ligation Enzyme		
TE buffer		Ambion
PCR Enzyme Mix		
Splint Buffer		
Digestion Buffer		
Digestion Enzyme		

ABSTRACT

BGISEQ-500 is a desktop sequencer developed by BGI. Using DNA nanoball and combinational probe anchor synthesis developed from Complete Genomics™ sequencing technologies, it generates short reads at a large scale. Library construction on the platform includes fragmentation, size selection, end repair and A-tailing, adaptor ligation, PCR amplification, and splint circularization.

Overview

1

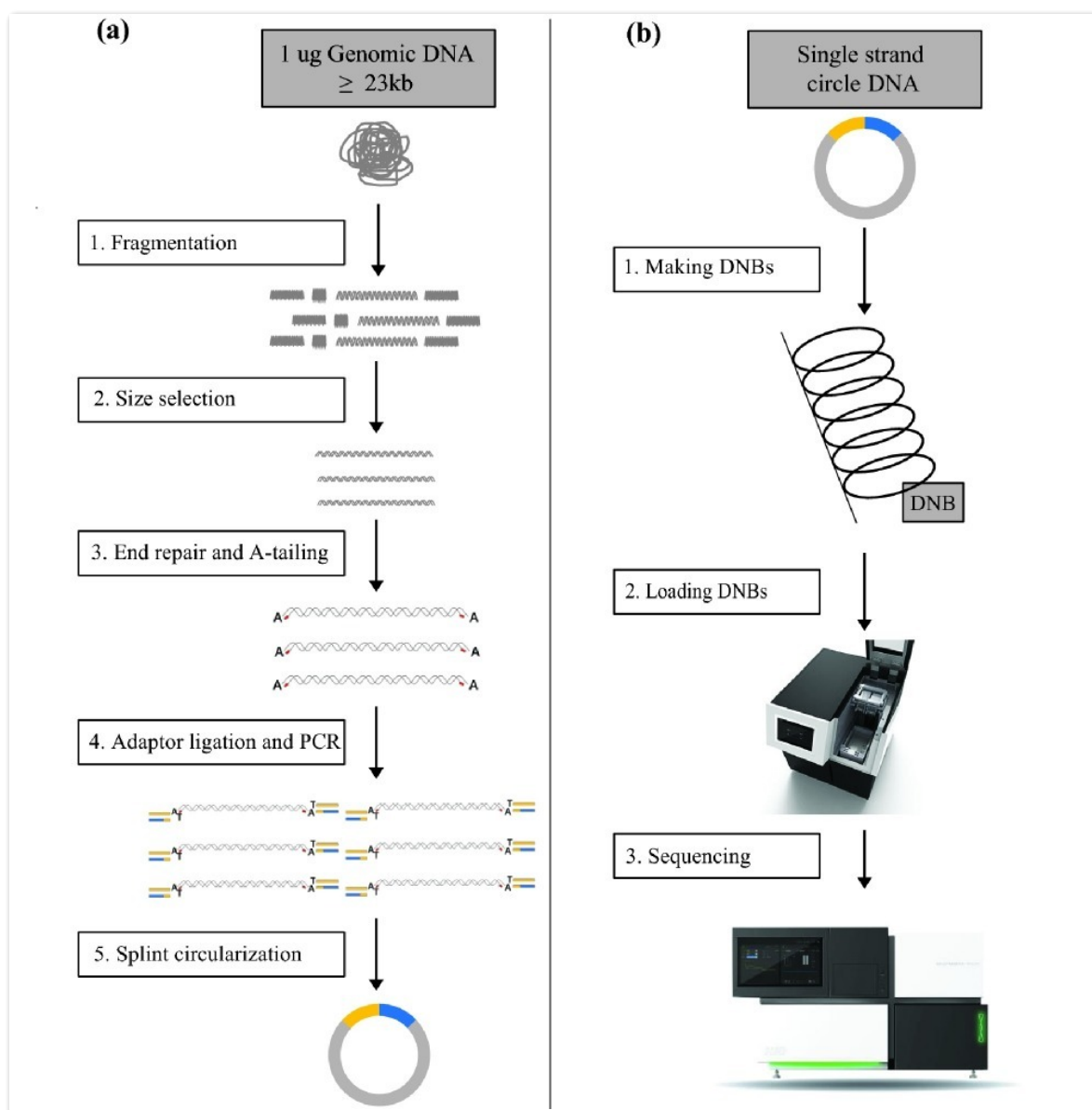

DNA fragmentation

2 1) Input genomic DNA sample

Genomic DNA Sample Recommendation

Nucleic Acid	*High-quality genomic DNA
Molecular Weight	>23k bp
Amount	1 µg
Concentration	$\geq 12.5\text{ng}/\mu\text{L}$
Purity	OD260/OD280=1.8~2.0

*High-quality genomic DNA should be free of salt or organics. It could run as an intact band with DNA length >23kb during 1% agarose gel electrophoresis

2) Fragmentation

Use the Covaris Focused-ultrasonicator for genomic DNA fragmentation following the instructions of the instrument. Optimisation should be performed on DNA prior to experiment and analyzed with agarose electrophoresis or an Agilent 2100 BioAnalyzer

Sequencing with	Input Amount	Reaction Volume	Derived Fragments
PE 100	1 µg	80 µL	100-700 bp (main band ≈ 200-300 bp)

PE 50	1μg	80μL	100-500 bp (main band≈200 bp)
PE 150	1μg	80μL	100-700 bp (main band≈400 bp)

3) Bead-based Cleanup

1. Place **AMPure XP magnetic beads** at RT for 30 min, fully thaw before use.
2. **PE100**: Pipette 48 μL AMPure XP magnetic beads to 80 μL shearing product, and mix well by gently pipetting 10 times, incubate for 5min at room temperature. (**PE50**: Pipette 80 μL AMPure XP magnetic beads to 80 μL shearing product, and mix well by gently pipetting 10 times, incubate for 5min at room temperature. **PE150**: Pipette 44 μL AMPure XP magnetic beads to 80 μL shearing product, and mix well by gently pipetting 10 times, incubate for 5min at room temperature.)
3. After brief centrifugation, place the non-stick tube on the magnet for 2min until the liquid clears, carefully transfer the supernatant to a new non-stick tube with a pipette.
4. **PE100**: Pipette 16 μL AMPure XP magnetic beads to 128 μL supernatant, mix well by gently pipetting 10 times, and incubate at room-temperature for 5 min. (**PE50**: Pipette 40 μL AMPure XP magnetic beads to 160 μL supernatant, mix well by gently pipetting 10 times, and incubate at room-temperature for 5 min. **PE150**: Pipette 12 μL AMPure XP magnetic beads to 124 μL supernatant, mix well by gently pipetting 10 times, and incubate at room-temperature for 5 min.)
5. After brief centrifugation, place the non-stick tube on the magnet for 2min until the liquid clears, remove and discard the supernatant with a pipette.
6. Add 500 μL of fresh 80% ethanol, while the tube remains on the magnet, then, rotate the tubes in the rack by half turns 4 times to wash the beads then carefully remove and discard the supernatant after 1 min.
7. Repeat step 6 once, and remove all liquid from tube without disrupting the beads.

4) Homogenization

1. Use double-strand DNA quantification kit such as Qubit® dsDNA HS Assay Kit or Quant-iT™ PicoGreen® dsDNA Assay Kit, and quantify the sample as per the instructions of the quantification kit.
2. Remove 50 ng of sample (calculated based on its concentration) to a new 0.2 mL PCR tube, and add NF water to final volume of 40 μL.

🕒 00:30:00 311

🕒 00:05:00 312

🕒 00:02:00 313

🕒 00:05:00 314

🕒 00:02:00 315

🕒 00:01:00 316

🕒 00:01:00 317



Fresh 80% ethanol
by XILONG SCIENTIFIC

End Repair and Tailing

- 3 1) Prepare the mixture as follows in PCR tube (do not vortex enzymes):

Components	Volume
DNA	40μL
ERAT Buffer	7.1μL
ERAT Enzyme	2.9μL
Total	50μL


- 2) Mix well by gently pipetting (Do not mix by vortexing), concentrate the reaction liquid to tube bottom by brief centrifugation.
- 3) Place the PCR tube containing the reaction mixture of above step in a Thermal Cycler, and initiate the reaction as per the following conditions:

Temperature	Time
-------------	------

Heated lid	On
37°C	30 min
65°C	15 min
4°C	Hold

🕒 00:45:00 3


 ERAT Buffer


 ERAT Enzyme

Ligate Adapters

- 4
- 1) Add **5 µL of Adapter Mix** to above PCR tube, and mix well by pipetting. Now 16 Adapter Mix are available, 8 libraries in one lane strategy, every sample with 4 different barcodes.
 - 2) Prepare the following reaction mixture (Note: Ligation Buffer II is viscous, pipette slowly):

Components	Volume
Ligation Buffer	23.4 µL
Ligation Enzyme	1.6 µL
Total	25 µL


- 3) Add 25 µL of above reaction mixture to the reaction solution containing adapters from above step.
- 4) Place the tube in a Thermal Cycler, then initiate reaction as per following condition:


Temperature	Time
Heated lid	On
23°C	30 min
4°C	Hold

- 5) After ligation, add **20 µL TE to final volume of 100 µL**, then transfer entire volume to a non-stick tube containing 50 µL of room temperature AMPure beads and mix by slow pipetting 10 times to avoid bubble formation.

🕒 00:30:00 4


 Adapter Mix


 Ligation Buffer


 Ligation Enzyme


 TE buffer
by Ambion

Purify Ligated DNA

- 5
 - 1) Incubate at room temperature for 5 min.
 - 2) After brief centrifugation, place the non-stick tube on the magnet for 2 min until the liquid clears, remove and discard the supernatant with a pipette:
 - 3) Add 500 µL of fresh 80% ethanol, while the tube remains on the magnet, then, rotate the tubes in the rack by half turns 4 times to wash the beads. Carefully remove and discard the supernatant after 1 min.
 - 4) Repeat step 4) once, remove all liquid from tube without disrupting the beads.
 - 5) Open the cap of non-stick tube, while the tube remains on the magnet, and dry at room temperature for 3 min.
 - 6) Remove the non-stick tube from the magnet, add 46 µL of TE for DNA elution, mix well by pipetting and incubate at room temperature for 5 min.
 - 7) After brief centrifugation, place the non-stick tube on the magnet for 2 min until the liquid clears, transfer all 44 µL of supernatant to a new 0.2 mL PCR tube ready for PCR in next step, or store at -20°C.

🕒 00:05:00 1

🕒 00:02:00 2

🕒 00:01:00 3

🕒 00:03:00 5

🕒 00:05:00 6

🕒 00:02:00 7

PCR


- 6
 - 1)

Components	Volume
DNA	44 µL
PCR Enzyme Mix	50 µL
PCR Primer Mix	6 µL
Total	100 µL

- 2) Place above PCR tube in a Thermal Cycler, and the initiate the reaction as per following conditions:

Temperature	Time	Cycles
Heated lid	On	
95°C	3 min	
98°C	20 sec	
60°C	15 sec	8
72°C	30 sec	
72°C	10 min	
4°C	Hold	

🕒 00:15:50 2

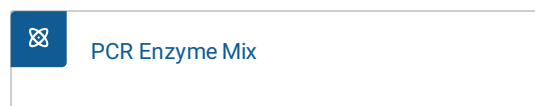

Purify PCR Product

- 7
 - 1) Place AMPure XP magnetic beads at room temperature 30 min in advance, mix well by vortexing before use.
 - 2) Add 100 µL of AMPure XP magnetic beads to 100 µL of PCR product, mix well by gently pipetting 10 times, and incubate at room temperature for 5 min.
 - 3) After brief centrifugation, place the non-stick tube on the magnet for 2 min until the liquid clears, remove and discard the supernatant with a pipette.

- 4) Add 500 μ L of fresh 80% ethanol, while the tube remains on the magnet, then, rotate the tubes in the rack by half turns 4 times to wash the beads. Carefully remove and discard the supernatant after 1 min.
- 5) Repeat step 4) once, and try to suck up all liquid from tube bottom.
- 6) Open the cap of non-stick tube, while the tube remains on the magnet, and dry at room temperature for 3 min.
- 7) Remove the non-stick tube from the magnet, add 32 μ L of TE water for DNA elution, mix well by pipetting and incubate at room temperature for 5 min.
- 8) After brief centrifugation, place the non-stick tube on the magnet for 2 min until the liquid clears, transfer the supernatant to a new non-stick tube. Proceed next step reaction or store at -20°C.

🕒 00:30:00 1

🕒 00:05:00 2

🕒 00:02:00 3

🕒 00:01:00 4

🕒 00:01:00 5

🕒 00:03:00 6

🕒 00:05:00 7


🕒 00:02:00 8

Homogenization

- 8 1 Use double-strand DNA quantification kit such as Qubit® dsDNA HS Assay Kit or Quant-iT™ PicoGreen® dsDNA Assay Kit, and quantify the sample as per the instructions of the quantification kit.
- 2 It is recommended to mix samples of different Barcodes here.
- 3 Add mixed sample (calculated based on its concentration) to a PCR tube, and add NF water to final volume of 48 μ L.

Circularization

9

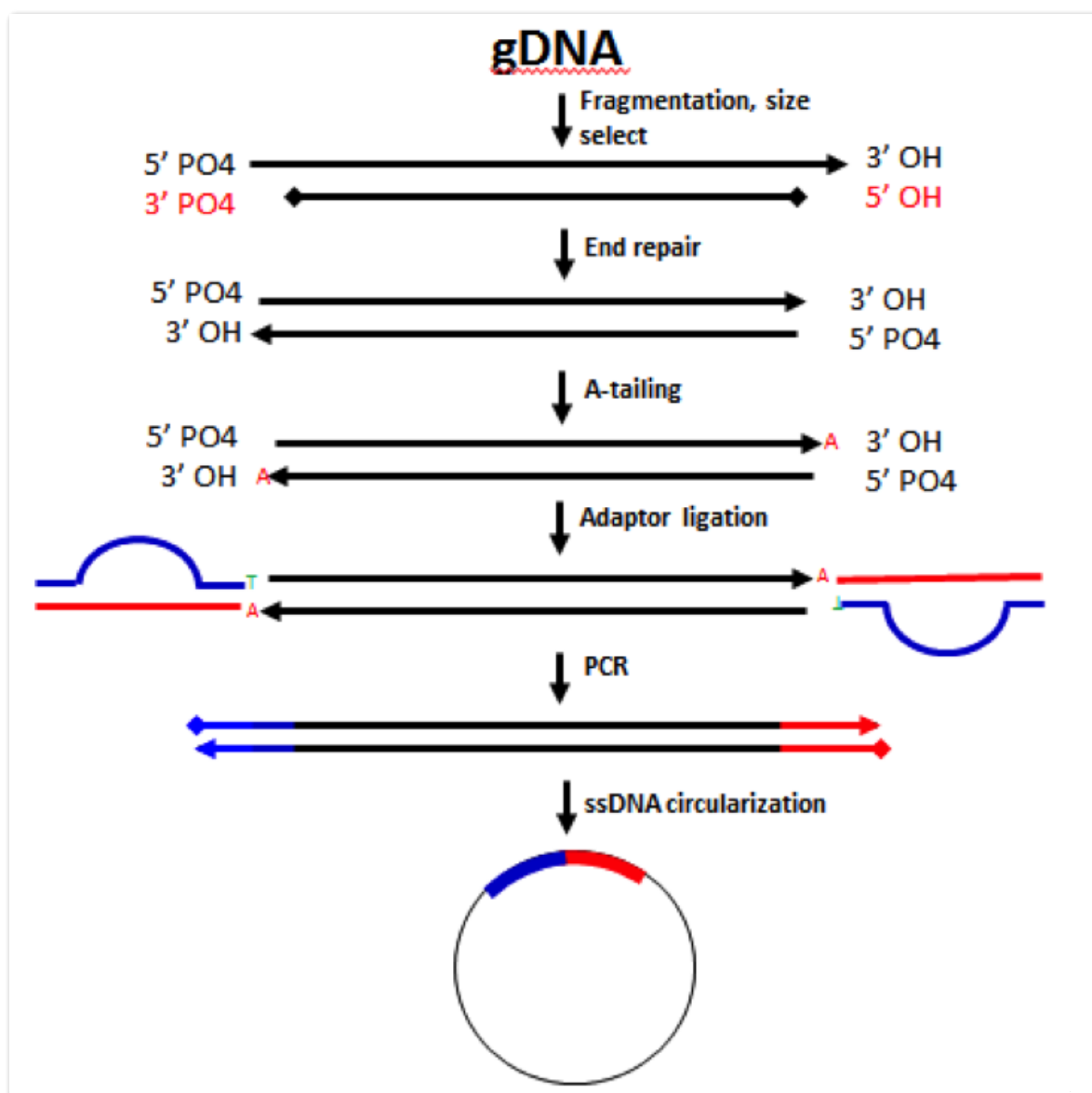

1) Denature the homogenized PCR product on a Thermal Cycler at 95°C for 3 min, then immediately transfer to ice batch.

2) Prepare reaction mixture on ice as per following system:

Components	Volume
Splint Buffer	11.6 µL
Ligation Enzyme	0.2 µL
Total	11.8 µL

3) Add 11.8 µL of above reaction mixture to 48 µL of denatured DNA.

4) Place above PCR tube in a Thermal Cycler, and initiate the reaction as per following conditions:

Temperature	Time
Heated lid	on
37°C	30 min
4°C	Hold

🕒 00:03:00 1

🕒 00:30:00 4



Splint Buffer

Digestion

- 10 1) Prepare digestion reaction solution on ice as per following system:

Components	Volume
Digestion Buffer	1.4 μ L
Digestion Enzyme	2.6 μ L
Total	4 μ L

2) After the circularization reaction is finished, directly add 4 μ L of digestion reaction solution into circularized DNA solution, mix well and briefly centrifuge, then place the PCR tube in a Thermal Cycler, and initiate the reaction as per following conditions:

Temperature	Time
Heated lid	on
37°C	30 min
4°C	Hold

- 3) Add 7.5 μ L of Digestion Stop Buffer to each reaction, mix well to terminate the reaction.
 4) Transfer all the reaction solution to a new non-stick tube, ready for purification.

🕒 00:30:00 2



Digestion Buffer



Digestion Enzyme

Purify Digestion Product

- 11 1) Place AMPure XP magnetic beads and place at room temperature for 30 min in advance. Mix well by vortexing before use.
 2) Pipette 168 μ L AMPure XP magnetic beads to digestion product, mix well by pipetting 10 times, and incubate at room temperature for 10 min.
 3) After transient centrifugation, place the non-stick tube on the magnet for 2 min until the liquid clears, remove and discard the supernatant with a pipette:
 4) Add 500 μ L of fresh 80% ethanol, while the tube remains on the magnet, then, rotate the tubes in the rack by half turns 4 times to wash the beads. Carefully remove and discard the supernatant after 1 min.
 5) Repeat step 4) once, and try to suck up all liquid from tube bottom.
 6) Open the cap of non-stick tube, while the tube remains on the magnet, and dry at room temperature for 3 min.
 7) Remove the non-stick tube from the magnet, add 32 μ L of TE for DNA elution, mix well by pipetting and incubate at room temperature for 10 min.
 8) After brief centrifugation, place the non-stick tube on the magnet for 2 min until the liquid clears, transfer the supernatant to a new non-stick tube. Store at -20°C, ready for preparation of DNB.

🕒 00:30:00 1

🕒 00:10:00 2

🕒 00:02:00 3

🕒 00:01:00 4

🕒 00:01:00 5



🕒 00:03:00 6%

🕒 00:10:00 7%

🕒 00:02:00 8%