




DEC 01, 2023

🌐 FlowFISH with PrimeFlow

Ronghao Zhou¹, Jesse Engreitz¹

¹Stanford University

Ronghao

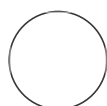

Ronghao Zhou

OPEN  ACCESS



ABSTRACT

Detect RNA expression in single cell with PrimeFlow

MATERIALS

ThermoFisher PrimeFlow RNA Assay Kit (88-18005-210)
positive control probe: RPL13A (Thermo VA4-13187-PF)
probe against gene of interest: MESDC1 (Thermo VA1-3010837-PF)

DOI:
dx.doi.org/10.17504/protocols.io.36wgq3q4olk5/v1

Protocol Citation: Ronghao Zhou, Jesse Engreitz 2023. FlowFISH with PrimeFlow. **protocols.io** <https://dx.doi.org/10.17504/protocols.io.36wgq3q4olk5/v1>

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working
We use this protocol and it's working

Created: Nov 30, 2023

Last Modified: Dec 01, 2023

Before start:

- 1 Make sure the heat block temp is set to 40°C using the digital monitor. The heat block should be ON at least a day in advance to stabilize temperature
- 2 Turn on Hyb oven (use to pre-warm target probe diluent)
- 3 Make FACS staining buffer: 0.5% BSA in PBS, filter sterilized
 - 10% BSA stock: 3g BSA (light sensitive powder in 4°C MISC) + 30mL PBS
 - Store at 4°C
- 4 Set vortex to ~1500RPM

A. Fixation and permeabilization

- 5
 - leave ~100µL when removing buffer, vortex to resuspend pellet
 - invert to mix when adding 1mL and tap the bottom if some cells stuck
 - pipette up and down to mix for fixation, permeabilization, and hybridization
 - may be in bulk, use volumes that cells don't exceed 10M cells/mL
 - pre-warm PrimeFlow RNA **Wash Buffer** to room temperature
- 6 Aliquot 5-10M cells in PBS per sample
- 7 Add 1mL PBS per sample, pipette to mix, and spin down at 500g at 4°C for 5min, discard supernatant, resuspend in ~100µL by vortexing gently

- 8 Prepare **Fixation Buffer 1** by mixing equal parts of PrimeFlow RNA **Fixation Buffer 1A** and PrimeFlow RNA **Fixation Buffer 1B**
- need 1mL per sample: 500µL Buffer 1A + 500µL Buffer 1B
 - mix gently by inverting, don't vortex or vigorously shake, prepare fresh
 - keep at 4°C
- 9 Add 1mL **Fixation Buffer 1**, **pipette** to mix, incubate in dark (put cardboard box on it) rotating on rotator for 30min at 2–8°C in cold room
- fixation time is critical, do NOT overshoot since the cells will lose integrity
- 10 Prepare 1X RNA **Permeabilization Buffer** with RNase Inhibitors by diluting PrimeFlow RNA **Permeabilization Buffer** (10X) and **RNase Inhibitors** (100X) with RNase-free water
- need 2mL per sample: 200µL Perm Buffer + 1.78mL H₂O + 20µL RNase Inhibitor
 - mix gently by inverting, don't vortex or vigorously shake, prepare fresh
 - keep at 4°C
- 11 Spin down cells at 800g at 4°C for 5min, discard supernatant
- speed is critical here: slower spin will lose a lot of cells
- 12 Add 1mL RNA **Permeabilization Buffer** with RNase Inhibitors, **pipette** to mix, and spin down at 800g at 4°C for 5min, discard supernatant
- 13 Repeat wash with 1mL RNA **Permeabilization Buffer**
- set centrifuge to room temp
- 14 Prepare 1X RNA **Fixation Buffer 2** by diluting PrimeFlow RNA **Fixation Buffer 2** (8X) with PrimeFlow RNA **Wash Buffer**
- need 1mL per sample: 125µL Fix Buffer 2 + 875µL RNA Wash Buffer
 - mix gently by inverting, don't vortex or vigorously shake, prepare fresh
 - keep at room temp
- 15 Add 1mL RNA **Fixation Buffer 2**, **pipette** to mix, and incubate for 60min in the dark at room temperature while rotating

- 16 Spin down cells at 800g at room temp for 5min, discard supernatant, resuspend (~100µL) by vortexing gently
- 17 Add 1mL PrimeFlow RNA **Wash Buffer**, invert to mix, and spin down at 800g at room temp for 5min, discard supernatant, resuspend (~100µL) by vortexing
- if in bulk, the cells should be transferred to the 1.5mL tubes from kit
- 18 Repeat wash with 1mL RNA **Wash Buffer**
- * can store samples overnight at 4°C, then the last wash with RNase Inhibitors

B. Target Probe hybridization

- 19 Thaw Probe Sets (20X), including positive control (RPL13A, Type 4) on ice; pre-warm PrimeFlow **RNA Target Probe Diluent** to 40°C
- 20 take 2µL cells from “unstained” sample into 18µL PBS -> measure on countess with trypan blue and record cell amount and % live
- **critical that the residual volume after all washes be as close to 100µL as possible, use the markings on the 1.5mL tubes to assist**
- 21 Dilute **Probe Sets** (20X) 1:20 in PrimeFlow RNA **Target Probe Diluent**, mix thoroughly by **pipetting** up and down
- need 100µL per sample: 5µL Target Probe + 5µL RPL13A Probe + 90µL Diluent
 - for unstained sample: 100µL Diluent
 - keep at 40°C
- 22 Add 100µL diluted **Target Probe(s)** into the cell suspension (~100µL), pipette to mix, briefly vortex, then incubate for 2h at 40°C
- do not pipette solutions onto the walls of the tubes, and samples should be mixed well before incubating
 - vortex samples to mix every 30min
 - temp is critical for hybridization
- 23 Add 1mL PrimeFlow RNA **Wash Buffer**, invert to mix, and spin down at 800g for 5min, discard supernatant, resuspend (~100µL) by vortexing gently

24

Prepare PrimeFlow RNA **Wash Buffer** with RNase Inhibitors by diluting **RNase Inhibitors** (100X) with RNA **Wash Buffer**

- need 1 mL per sample: 10µL RNase Inhibitor + 990µL Wash Buffer
- mix gently by inverting, prepare fresh
- keep at room temp

25

Add 1 mL PrimeFlow RNA **Wash Buffer** with RNase Inhibitors, **invert** to mix, and spin down at 800g for 5min, discard supernatant, resuspend (~100µL) by vortexing gently

- Store samples overnight in the dark at 4°C

C. Signal amplification

26

- Pre-warm samples and PrimeFlow RNA **Wash Buffer** to room temperature
- Pre-warm PrimeFlow RNA **PreAmp Mix**, PrimeFlow RNA **Amp Mix**, and PrimeFlow RNA **Label Probe Diluent** to 40°C

27

Add 100µL PrimeFlow RNA **PreAmp Mix** into the cell suspension (~100µL), **pipette** to mix, briefly vortex, then incubate for 1.5h at 40°C

- do not pipette solutions onto the walls of the tubes, and samples should be mixed well before incubating

28

Add 1 mL PrimeFlow RNA **Wash Buffer**, **invert** to mix, and spin down at 800g for 5min, discard supernatant, resuspend (~100µL) by vortexing gently

29

Repeat wash two times with 1 mL RNA **Wash Buffer**, for a total of three washes

30

Add 100µL PrimeFlow RNA **Amp Mix** into the cell suspension (~100µL), **pipette** to mix, briefly vortex, then incubate for 1.5h at 40°C

- do not pipette solutions onto the walls of the tubes, and samples should be mixed well before incubating

31

Add 1 mL PrimeFlow RNA **Wash Buffer**, **invert** to mix, and spin down at 800g for 5min, discard supernatant, resuspend (~100µL) by vortexing gently

- 32 Repeat wash with 1mL RNA **Wash Buffer**
- 33 Dilute PrimeFlow RNA **Label Probes** (100X) 1:100 in PrimeFlow RNA **Label Probe Diluent**
- need 100µL per sample: 1µL Label Probes + 99µL Diluent
 - keep at 40°C
- 34 Add 100µL diluted **Label Probes** into the cell suspension (~100µL), **pipette** to mix, briefly vortex, then incubate for 1h at 40°C
- do not pipette solutions onto the walls of the tubes, and samples should be mixed well before incubating
- 35 Add 1mL PrimeFlow RNA **Wash Buffer** at room temp, **invert** to mix, and spin down at 800g for 5min, discard supernatant, resuspend (~100µL) by vortexing gently
- 36 Repeat washes with warm (35°C) 1mL RNA **Wash Buffer** 5 times
- 37 Add 1mL **FACS staining buffer** (or PrimeFlow RNA **storage buffer**), **invert** to mix, and spin down at 800g for 5min, discard supernatant, resuspend (~100µL) by vortexing gently
- Samples can be stored in the dark at 4°C for up to three days before analysis
 - QC: take stained vs unstained samples and check fluorescence with microscope


D. Flow cytometric analysis

- 38 Add 100µL FACS buffer to cells, transfer all to filter, quick spin
- 39 Add 100µL FACS buffer to cap to wash remaining cells from filter, quick spin

40


FACS:

- 400-500µL of an optimally concentrated sample (20M cells/mL) take ~30min
- when low volume left, add ~200µL staining buffer can recover more cells
- setup:
 - o sort population should be <10% of unstained sample
 - o exclusion should be <5% of stained sample
 - § if more than 5%, increase voltage or decrease compensation
 - o compensate for each gene: AF647 over AF488
 - § RPL13A (AF488) expression = total RNA, compensate gene of interest expression (AF647) relative to total RNA
 - § average AF488 for Top & Bottom 25% AF647 should be within 10% of average AF488 for all sort population

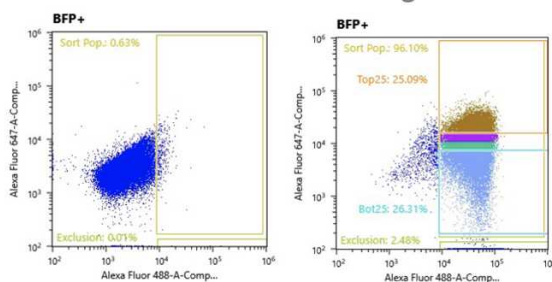

Sort Pop.	64,557	96.10%	64.56%	
Alexa Fluor 488-A-Comp...				44,080
Alexa Fluor 647-A-Comp...				12,052
Top25	16,858	25.09%	16.86%	
Alexa Fluor 488-A-Comp...				45,768
Alexa Fluor 647-A-Comp...				21,142
Bot25	17,678	26.31%	17.68%	
Alexa Fluor 488-A-Comp...				43,948
Alexa Fluor 647-A-Comp...				4,742

- sort:

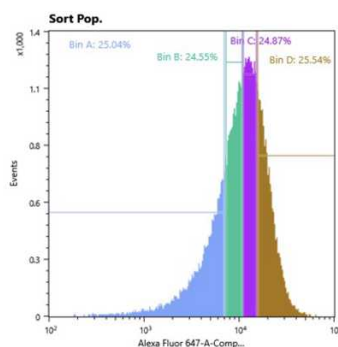

E. gDNA Extraction after FlowFISH

41

- Homemade **ChIP Lysis Buffer**:
 - o 1% SDS, 10mM EDTA, 50mM Tris-HCl, pH7.5
 - o store at 4°C
 - o for 50mL
- Make Lysis Buffer fresh or warm at room temp 30-60min to solubilize precipitate
- Always include a no cell control until PCR to check contamination
- Set thermocycler:

o 65°C Hold -> 65°C 10min -> 37°C Hold -> 37°C 30min -> 65°C Hold -> 65°C 2h -> 95°C 20min -> 4°C Hold

- Protocol for **1M** cells

42 Spin cells down for 10min at 800*g* at 4°C, remove supernatant

43 Resuspend cells in 70µL **ChIP Lysis Buffer**, transfer to 96-well plate

44 Incubate at 65°C for 10min

45 When samples cool to 37°C, add 2µL **RNase Cocktail** (Invitrogen AM2286) mix well by pipetting

46 Incubate at 37°C for 30min

47 Add 10µL **Proteinase K** (Thermo 25530049), mix well by pipetting

48 Incubate at 65°C for 2h, then 95°C for 20min

49 Store at 4°C (can store overnight)

- 50** Beads clean with 0.7X AMPure beads: 70µL AMPure beads for ~100µL samples
- a. warm AMPure beads to room temp
 - b. add beads, pipette to mix, and bind DNA for 2min
 - c. wash 3 times with 150µL 70% EtOH
 - d. let beads dry for 5-10min after 3rd wash
 - e. elute with ~40-50µL H₂O, elute for 2min, transfer to a new plate/tube