



JUL 10, 2023

OPEN ACCESS

**DOI:**  
[dx.doi.org/10.17504/protocols.io.dm6gp3ww1vzp/v1](https://dx.doi.org/10.17504/protocols.io.dm6gp3ww1vzp/v1)

**Protocol Citation:** Ana C Reis, Daniela Pinto, Mónica V. Cunha 2023. RT-PCR amplification for the discrimination between wild type and mutation alleles in SARS-CoV-2 genomes from wastewater. **protocols.io** <https://dx.doi.org/10.17504/protocols.io.dm6gp3ww1vzp/v1>

**License:** This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

**Protocol status:** Working  
 We use this protocol and it's working

**Created:** Jul 10, 2023


**Last Modified:** Jul 10, 2023

**PROTOCOL integer ID:**  
 84761

# RT-PCR amplification for the discrimination between wild type and mutation alleles in SARS-CoV-2 genomes from wastewater

Daniela Ana C Reis<sup>1</sup>, Pinto<sup>1</sup>, Mónica V. Cunha<sup>1</sup>

<sup>1</sup>Centre for Ecology, Evolution and Environmental Changes (cE3c) & CHANGE - Global Change and Sustainability Institute, Faculdade de Ciências, Universidade de Lisboa, Lisboa, Portugal

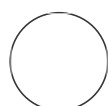

Daniela Pinto

Centre for Ecology, Evolution and Environmental Changes (cE3c)

## ABSTRACT

This protocol describes the procedure of RT-PCR for the allelic discrimination between wild type and mutation alleles in specific positions of the SARS-CoV-2 genome. Starting with isolated RNA from SARS-CoV-2 positive wastewater samples, this workflow combines Applied Biosystems TaqMan SNP Genotyping Assays with a one-step RT-PCR reaction to detect whether there are known SARS-CoV-2 mutations present in the wastewater samples.

Sequence-specific forward and reverse primers amplify the target sequence region in each mutation assay. The reverse primer also primes reverse transcription of the SARS-CoV-2 genomic RNA sequences. Each mutation assay includes two TaqMan minor groove binder probes with nonfluorescent quenchers: one VIC dye labelled probe to detect the reference sequence and one FAM dye labelled probe to detect the mutation sequence.

## MATERIALS

- TaqPath 1-Step RT-qPCR Master Mix CG, ThermoFisher Scientific, Catalog A15299
- AcroMterix SARS-CoV-2 Control, ThermoFisher Scientific, Catalog 954519
- Mutation Assay: S.G339D.GGT.GAT, ThermoFisher Scientific, Catalog CV47VRX
- Mutation Assay: S.Q493R.CAA.CGA, ThermoFisher Scientific, Catalog CVH49P2
- Mutation Assay: S.L452R.CTG.CGG, ThermoFisher Scientific, Catalog CVAAAAD
- Mutation Assay: S.T547K.ACA.AAA, ThermoFisher Scientific, Catalog CVYMJGA
- Mutation Assay: M.D3N.GAT.AAT, ThermoFisher Scientific, Catalog CVAAAAK
- Mutation Assay: ORF7b.L11F.TTG.TTT, ThermoFisher Scientific, Catalog CVCE3VH

## BEFORE START INSTRUCTIONS

Wipe bench surfaces with RNase Away and keep working environment clean.  
Maintain the RNA samples on ice, and prepare the RT-PCR mixture on ice.

### RT-PCR amplification

#### 1 Prepare the amplification mixture as follows:

Reagent	Concentration in the reaction
TaqPath 1-Step RT-qPCR Master Mix CG (4X)	1X
Mutation assay (40X)	2X
DNase/RNase free water	up to 15 µL

To the 15 µl of reaction mix, add 5 µl of wastewater RNA sample, or no-template control (DNase/RNase free water), or a wild-type AcroMetrix Coronavirus 2019 (COVID-19) RNA control, to complete 20 µl of reaction volume.

#### 2 Thermocycling conditions

Step	Temperature (°C)	Time	No. Cycles
Reverse transcription	45	15'	
Taq polymerase activation	95	2'	
Denaturation	95	15"	45
Annealing	58	45"	
Post-read	60	30'	

Set the detection module for FAM and VIC labelled probes.

### Data analysis

#### 3 The allelic discrimination results can be graphed on a scatter plot contrasting reporter dye fluorescence (*i.e.*, allele X versus allele Y), and analyzed with the QuantStudio Design and Analysis Software, using the genotyping analysis module with real-time data.