



APR 06, 2023

OPEN  ACCESS

**DOI:**  
[dx.doi.org/10.17504/protocols.io.4r3l27j4qg1y/v1](https://dx.doi.org/10.17504/protocols.io.4r3l27j4qg1y/v1)

**Protocol Citation:** Kazumasa Wakamatsu 2023. Synthesis and preparation of DOPA pheomelanin and eumelanin.  
**protocols.io**  
<https://dx.doi.org/10.17504/protocols.io.4r3l27j4qg1y/v1>

**License:** This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

**Protocol status:** Working  
We use this protocol and it's working

**Created:** Apr 06, 2023

**Last Modified:** Apr 06, 2023


**PROTOCOL integer ID:**  
80143

**Keywords:** ASAPCRN

## Synthesis and preparation of DOPA pheomelanin and eumelanin

Kazumasa Wakamatsu<sup>1</sup>

<sup>1</sup>Institute for Melanin Chemistry, Fujita Health University, Toyoake, Japan

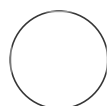

Daniel Choo

### ABSTRACT

This is the protocol for synthesizing and preparing DOPA pheomelanin and DOPA eumelanin using tyrosinase.

### ATTACHMENTS

[pcmr12121-sup-0001-appendixs1.docx](#)

- 1 As precursor(s) of melanin, 1) l-dopa, 2) DHI, 3) a mixture of DHI and DHICA, 4) DHICA, 5) a mixture of l-dopa and l-cysteine, and 6) 5SCD in the presence of l-dopa is used.

- 2 To a solution of the precursor(s) (1 mmol) in 0.1 M sodium phosphate buffer, pH 6.8 (98 ml), add a solution of mushroom tyrosinase in the same buffer (2 ml) at 25°C.
- 2.1 The amount of mushroom tyrosinase (Sigma-Aldrich, 4270 U/mg) is either 50,000 U or 100,000 U.
- 3 Carry out oxidation with vigorous (maximal, but avoiding splashing) stirring in a 500 ml Erlenmeyer flask.
- 4 Stop oxidation (1-4) at 4 h by adding 6 M HCl (2 ml) to make pH 1 and keep the mixture at 4°C overnight to ensure complete precipitation of the pigment.
- 5 Collect the precipitate by centrifugation at 3,000 rpm for 10 min and wash with 0.1 M HCl (40 ml x 3).
- 6 For the preparation of pheomelanins (5 and 6), acidify the oxidation mixture with acetic acid (3 ml) to make pH 3, and wash the precipitate with 1% acetic acid (40 ml x 3).
- 7 Dry the melanin powder by lyophilization and equilibrate with moisture in a desiccator with a saturated CaCl<sub>2</sub> solution.
- 8 Prior to treatment, suspend the synthetic DOPA pheomelanin (DOPA:Cys = 1:1) and DOPA eumelanin (DOPA:Cys = 1:0) in sterile PBS and sonicate overnight.

