




MAR 22, 2023

Phosphorus Extraction - Bray Method

In 2 collections

maggie.bowman¹

¹Environmental Molecular Sciences Laboratory

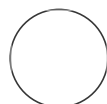

nicholas.sconzo

ABSTRACT

The method estimates the relative bioavailability of inorganic ortho-phosphate (PO₄-P) in soils with acid to neutral pH.

MATERIALS

- 50 mL polypropylene centrifuge tube (3 per site and depth)
- 15 mL polypropylene centrifuge tube (3 per site and depth)
- ~ 12 g of fresh soil
- Bray extract solution (See OSU Soil health lab protocol for making extraction solution)
- 25 mL pipet
- Electronic pipettor
- Centrifuge
- Shaker

OPEN ACCESS

DOI:
dx.doi.org/10.17504/protocols.io.3byl4jkbolo5/v1

Protocol Citation: maggie.bowman 2023. Phosphorus Extraction - Bray Method. protocols.io
<https://dx.doi.org/10.17504/protocols.io.3byl4jkbolo5/v1>

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working

Created: Oct 06, 2022



Last Modified: Mar 22, 2023

PROTOCOL integer ID:
70965


- 1 Print the required labels needed for Bray extraction for each site, depth, and replicate.

1.1 1000S_SITE_P_BRAY_DEPTH_REPLICATE

1.2 1000S_SITE_Pext_BRAY_DEPTH_REPLICATE


2 Weigh out  2 g of soil into labeled  50 mL centrifuge tubes.

2.1 For example: 1000S_PRS1_P_Bray_TOP_1


3 Add  14 mL of Bray Extract solution (in cabinet under hood in 1521) using the pipet and electronic pipettor.

4 Cap the samples and shake to mix. Add samples to foam block on shaker.


5 Shake samples at  1000 rpm, 00:05:00 .

6 Centrifuge samples (using centrifuge on west side of 1521) for  4000 rpm, 00:08:00 .

8m

7 Pour off liquid into a clean and labelled  15 mL centrifuge tube.

7.1 For example: 1000S_PRS1_Pext_Bray_TOP_1

8 Freeze samples on their side to prevent tube breaking in the  -20 °C freezer (Staging Freezer 2) bottom shelf if unable to filter immediately.

9 Follow sample filtering protocol below when ready to filter extracts.