




AUG 25, 2023

# Protocol to isolate and fix nuclei from flash frozen mouse adrenal glands for IGVF

Elisabeth Rebboah<sup>1</sup>

<sup>1</sup>University of California, Irvine

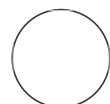

Elisabeth Rebboah  
University of California, Irvine

## ABSTRACT

This protocol describes isolation of nuclei from 10 week old mouse adrenal glands from 8 founder strains (B6J, AJ, 129S1J, NZOJ, WSBJ, NODJ, PWKJ, and CASTJ), preparation of a single nucleus suspension, and fixation for single nucleus RNA-seq using Parse Biosciences. We process 1 rep from each strain per day; e.g. female rep 1 across all 8 strains. The main products we use are Parse Biosciences Nuclei Fixation Kit (v2) and Miltenyi Biotec's gentleMACS Octo Dissociator with accessories. This protocol takes about 3.5 hours from start to finish.

The results are 2 aliquots of fixed single-nucleus suspensions for Parse per each of the 8 samples.

## GUIDELINES

1. We recommend using a 5 ml pipette for aspirations and resuspensions > 1 ml.
2. Record everything in the [IGVF spreadsheet](#), "Samples into experiment" tab.

## MATERIALS

Name	Manufacturer	Cat. #
Nuclei Fixation Kit v2	Parse Biosciences	ECF2003
Nuclei Extraction Buffer	Miltenyi Biotec	130-128-024
RNase Inhibitor, murine	New England Biolabs	M0314L
PBS	HyClone	SH30256.02
7.5% BSA	Life Technologies	15260037
gentleMACS C Tube	Miltenyi Biotec	130-093-237
gentleMACS Octo Dissociator	Miltenyi Biotec	130-095-937
MACS SmartStrainers (30 um)	Miltenyi Biotec	130-098-458

OPEN ACCESS



**Protocol Citation:** Elisabeth Rebboah 2023. Protocol to isolate and fix nuclei from flash frozen mouse adrenal glands for IGVF. [protocols.io](https://protocols.io/view/protocol-to-isolate-and-fix-nuclei-from-flash-froz-cy7dxzi6) <https://protocols.io/view/protocol-to-isolate-and-fix-nuclei-from-flash-froz-cy7dxzi6>

**License:** This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

**Protocol status:** Working  
We use this protocol and it's working

**Created:** Aug 25, 2023

**Last Modified:** Aug 25, 2023

**PROTOCOL integer ID:**  
86981

Name	Manufacturer	Cat. #
NucBlue Fixed Cell ReadyProbes	Thermo Fisher	R37606
Hemacytometer	Fisher Scientific	02-671-51B
Mr. Frosty	Sigma-Aldrich	635639

### Reagents/equipment, manufacturer and catalog number

Name	Reagent	Volume (for 8 samples)	Final Concentration
Lysis buffer	Nuclei Extraction Buffer	35 ml	NA
	40 U/ul RNase inhibitor	175 ul	0.2 U/ul
NB + BSA + RNase inhibitor	Nuclei Buffer (Parse Biosciences)	3.15 ml	NA
	7.5% BSA	350 ul	0.75%
	RNase inhibitor (Parse Biosciences)	44.1 ul	

### Buffers

## Setup

- 1 Label tubes.
- 2 Pre-chill centrifuge to 4C.
- 3 Prepare ice buckets.

- 4 Prepare lysis buffer in a 50 ml conical tube on ice. Distribute 2 ml into 8 gentleMACS C Tubes on ice. Add RNase inhibitor to the lysis buffer aliquot the day of the experiment.
- 5 Prepare NB + BSA + RNase inhibitor. Add RNase inhibitor included in Parse Biosciences fixation kit the day of the experiment.
- 6 Prepare 1.5 ml Nuclei Buffer + RNase inhibitor for final resuspension. Add 18.6 ul RNase inhibitor to 1.5 ml Nuclei Buffer.
- 7 Thaw components of 1 Parse Biosciences Nuclei Fixation v2 kit at room temperature, then place on ice.
- 8 Distribute 20 ul NucBlue Fixed Cell ReadyProbes into 16 PCR strip tubes for cell counting. Need 8 tubes for counting after nuclei extraction, and another 8 tubes for final fixed nuclei.

## **Tissue lysis and nuclei extraction**

- 9 Keep flash frozen tissue samples on dry ice until lysis.
- 10 Drop whole frozen tissue into a chilled gentleMACS C Tube with 2 ml lysis buffer. Close tubes firmly and invert immediately, ensuring tissue is not stuck to the bottom or side. Keep tubes on ice and proceed immediately to dissociation. There should be 2 adrenal glands.
- 11 Run the gentleMACS Program 4C\_nuclei\_1 on the Octo Dissociator (~5 minutes).


- 12 Remove tubes, ensuring tissue did not get stuck on the sides, and spin down in a 4C centrifuge for ~10 seconds to bring liquid to the bottom, then place tubes back on ice.
- 13 Filter nuclei suspension through 30 um MACS SmartStrainer into a 5 ml tube. Fit a tube rack in ice for extra stability while filtering.
- 14 Wash 30 um MACS SmartStrainer with 2 ml additional lysis buffer. Add 2 ml to C tubes, cap, and swish to recover any nuclei stuck to the sides and cap of the C tubes, then wash the strainer.
- 15 Discard strainer and centrifuge the 4 ml nuclei suspension at 4C, 350g for 5 minutes.
- 16 Discard supernatant and resuspend nuclei pellet in 375 ul NB-BSA + RNase inhibitor.

## Parse Nuclei Fixation

- 17 Filter the nuclei suspension (375 ul) NB-BSA + RNase inhibitor through a 40 um strainer (provided in Parse Biosciences kit) into a new 5 ml tube. Count nuclei. Use 1:6 dilution factor, e.g. 4 ul + 20 ul dye.
- 18 Add 125 uL Nuclei Fixation Solution and mix 3 times. Do not over-mix.
- 19 Incubate nuclei for 10 minutes on ice. Set 1 P200 pipette to 40 ul and keep the other P200 at 125 ul.
- 20 Add 40 uL Nuclei Permeabilization Solution and mix by pipetting 3 times with the P200 still set to

125 uL. Do not over-mix.

- 21 Incubate 3 minutes with nuclei on ice.
- 22 Add 2 ml Nuclei Neutralization Solution and invert the tube once to mix.
- 23 Centrifuge at 4C, 750g for 10 minutes.
- 24 Aspirate and discard supernatant.
- 25 Resuspend the samples in 150 ul Nuclei Buffer with RNase inhibitor without BSA and filter through a 40 um filter into a labeled 1.5 ml tube.
- 26 Count nuclei. Use a 1:11 dilution factor, e.g. 2 ul + 20 ul dye.
- 27 Add Nuclei DMSO: For 150 ul samples: add 2.5 ul and gently flick tubes to mix. One minute later, add another 2.5 ul and flick to mix, then after another minute add a final 2.5 ul for a total volume of 7.5 ul. Mix by gently pipetting 5x with a P200 set to 50 ul.
- 28 Split nuclei suspension into 2 aliquots of equal volume in labeled 1.5 ml tubes.

- 
- 29** Place tubes in a Mr. Frosty at -80C. The next day, move tubes to boxes in -80C racks.
  - 30** Move leftover nuclei suspension to labeled 2 ml tubes and spin at 4C, 750g for 5 minutes. Remove supernatant and flash-freeze nuclei in liquid nitrogen as dry pellets. Store at -80C.