



Version 3

Apr 22, 2021

# Concentration of viruses from sewage using HA filters

## V.3

Ad  la  de Roguet<sup>1</sup>, Shuchen Feng<sup>1</sup><sup>1</sup>UWM

In Development

[dx.doi.org/10.17504/protocols.io.bue9nth6](https://dx.doi.org/10.17504/protocols.io.bue9nth6)

mclellan lab



McLellan Lab

### ABSTRACT

#### Concentration of viruses from sewage using HA filters

### DOI

[dx.doi.org/10.17504/protocols.io.bue9nth6](https://dx.doi.org/10.17504/protocols.io.bue9nth6)

### PROTOCOL CITATION

Ad  la  de Roguet, Shuchen Feng 2021. Concentration of viruses from sewage using HA filters. **protocols.io**  
<https://dx.doi.org/10.17504/protocols.io.bue9nth6>  
Version created by [McLellan Lab](#)



### LICENSE

————— This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

### CREATED

Apr 22, 2021

### LAST MODIFIED

Apr 22, 2021

### PROTOCOL INTEGER ID

49345

### GUIDELINES

Each sewage samples is processed in triplicate: two technical duplicates and one archive.

## MATERIALS TEXT

Bovine coronavirus solution prepared according to

[Bovine](#)

[coronavirus](#) **Zoetis Catalog #CALF-GUARD**

[Cellulose esters filters 0.8 µm 47 mm diameter](#) **Emd**

**Millipore Catalog #AAWP04700** Step 6.2

[RNeasy PowerMicrobiome](#)

**Kit Qiagen Catalog #26000-50** Step 6.5

[ZR BashingBead Lysis Tubes \(2mm\)](#) **Zymo**

**Research Catalog #S6003-50** Step 6.5

Gast DOA-P704-AA  
High-Capacity Vacuum Pump, with Gauge,  
Regulator, and Relief

Cole-Parmer EW-07061-40 [↗](#)

Magnetic Filter Funnel  
300 mL, 47 mm

Pall 4242 [↗](#)

Pipette controller  
25 ml serological pipette

## SAFETY WARNINGS

The filtration has to be performed in a biosafety cabinet to avoid any contact with raw sewage.

## DISCLAIMER:

DISCLAIMER – FOR INFORMATIONAL PURPOSES ONLY; USE AT YOUR OWN RISK

The protocol content here is for informational purposes only and does not constitute legal, medical, clinical, or safety advice, or otherwise; content added to [protocols.io](https://protocols.io) is not peer reviewed and may not have undergone a formal approval of any kind. Information presented in this protocol should not substitute for independent professional judgment, advice, diagnosis, or treatment. Any action you take or refrain from taking using or relying upon the information presented here is strictly at your own risk. You agree that neither the Company nor any of the authors, contributors, administrators, or anyone else associated with [protocols.io](https://protocols.io), can be held responsible for your use of the information contained in or linked to this protocol or any of our Sites/Apps and Services.

## BEFORE STARTING

Before starting the filtration:

1. Clean the working area and all equipments: wipe down with 10% bleach and let stand for 10 min, then rinse with water and dry. Next, wipe down using 70% ethanol and let dry.
2. Verify that the carboy to receive filtration waste is not full or near to be full.
3. Keep forceps sterile by placing them in a small beaker containing 70% ethanol.
4. Prepare an ice bucket.
5. Thaw the bovine coronavirus (BCoV) recovery control on ice. Do not thaw the same aliquot more than twice.
6. Add  $\text{MgCl}_2$  stock solution in 50-mL tubes for a final concentration of 25 mM.
7. Warm Solution PM1 (provided in RNeasy PowerMicrobiome kit) at 55°C for at least 10 min prior to dissolve precipitates. Shake to mix before transferring 650  $\mu\text{L}$  into each 2-mL bead beating tubes (3 tubes needed per sample).

When work is completed, remove equipments and supplies from the cabinet. Wipe the work area with 10% bleach, let stand for 10 min, rinse with water, then with 70% ethanol. Finally, turn on the UV lamp for 15 minutes. Empty the carboy.

#### Filtration of the sewage samples (in the biosafety cabinet)

- 1 Take the sample out from the refrigerator.
- 2 Homogenize the sample thoroughly, avoiding foaming.
- 3 Use a pipette controller, transfer 25 mL of the sample into a 50-mL tube (where  $\text{MgCl}_2$  has been previously added **25 Milimolar (mM) final concentration** ). Repeat for each replicate. Save the rest of sample at 4°C for archive or until decontamination and disposal.
- 4 Spike in each replicate 5  $\mu\text{L}$  of Bovine Coronavirus recovery control (~100,000 cp/uL), prepared as Bovine coronavirus solution prepared according to <http://dx.doi.org/10.17504/protocols.io.bpg8mjzw>
- 5 Put the 50-mL tubes on ice until further processing.
- 6 Filter the samples in the 50-mL tubes:

##### 6.1 Rinse the clean filtration funnel with sterile deionized water.

Magnetic Filter Funnel

300 mL, 47 mm

Pall

4242





Make sure some water is still on the porous plate.

## 6.2

Place the HA filter in the filtration funnel.

Cellulose esters filters 0.8  $\mu\text{m}$  47 mm diameter Emd

Millipore Catalog #AAWP04700

filter



Make sure the filter gets wet due to the residual water on the porous plate.

6.3 Homogenize the sample in the 50-mL tube. Avoid as much as possible the liquid to touch the lid.

6.4 Filter the content in the 50-mL tube through a 0.80  $\mu\text{m}$  47mm cellulose ester filter at low vacuum pressure (~50 mmHg VAC).



Gast DOA-P704-AA  
High-Capacity Vacuum Pump, with Gauge,  
Regulator, and Relief  
Cole-Parmer EW-07061-40 [Link](#)

## 6.5

As soon as the filtration is complete, transfer the filter to a 2-mL ZR BashingBead Lysis Tubes containing 650  $\mu$ L of PM1 buffer solution (previously added).

[RNeasy PowerMicrobiome](#)

Kit [Qiagen Catalog #26000-50](#)



It is important not to “over-dry” the filter on the filtration system.



Fold the filter as shown in Figure 1. Do not crease the filters with the forceps.

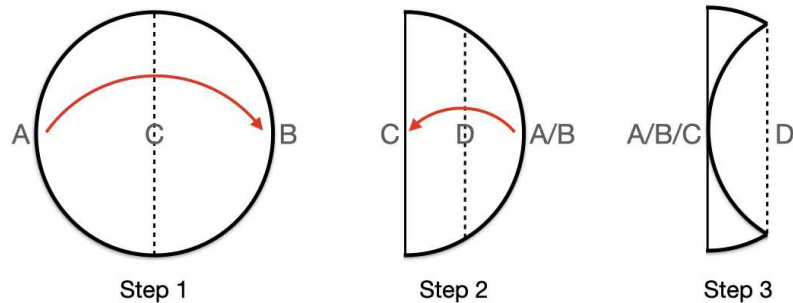


Figure 1 Filter folding procedure. Step 1: fold the filter in half. Step 2: repeat that step. Step 3: transfer to a 2-mL tube.

[ZR BashingBead Lysis Tubes \(2mm\)](#) [Zymo](#)

Research Catalog #S6003-50

## 6.6 Put the 2-mL tubes on ice until further processing.

6.7 Repeat Steps 6 for the other replicates. ↻


7 Store filters at -20°C until all samples are processed

8 ↻

Store filters at -80°C for a minimum of 2h or overnight prior to the nucleic acid extraction.  
Archive filters at -80°C.

#### Nucleic acid extraction

9

 Nucleic acid extraction - RNeasy PowerMicrobiome Kit (Qiagen)  
by McLellan Lab

PREVIEW

RUN



9.1 For **HA filter** extraction, add 6.5 µL of beta-Mercaptoethanol to each of the 2-mL ZR BashingBead Lysis tubes. Go to **step 2**.



Add the beta-Mercaptoethanol while samples are still frozen. Let thaw at room temperature.

For **BCoV/BRSV** extraction, add 5 µL of BCoV/BRSV solution to the 2-mL tube containing the warmed PM1 + beta-Mercaptoethanol solution. Vortex for 15 seconds (speed 7 out of 10) and let sit for 10 min. Skip step 2.

For **Direct extraction**, add 150 µL of wastewater to the 2-mL tube containing the warmed PM1 + beta-Mercaptoethanol solution. Vortex for 15 seconds (speed 7 out of 10) and let sit for 10 min. Skip step 2.

9.2

For **HA filter** extraction, place the 2-mL tubes in the bead beater.

Mini-Beadbeater-16  
high-energy cell disrupter

BioSpec

607



1 speed

9.2.1 Bead beat for ⌚ 00:02:30

2m



Start the bead beating when the beads start to be loose in the tubes.

9.2.2 Cool down the samples on ice for ⌚ 00:05:00 .

5m

9.2.3 Repeat Steps 9.1 and 9.2 once ↺ .

9.3 Centrifuge at maximum speed for 1 min at room temperature. ⚙️ 150000 rpm, Room temperature , 00:01:00

9.4 For **HA filter** extraction, transfer 450 µL of supernatant to a Collection Tube (provided in the RNeasy PowerMicrobiome kit).  
For **BCoV/BRSV/Direct extraction**, transfer all supernatant.

9.5 For **HA filter/BCoV/BRSV** extraction, add 150 µL of Solution IRS.  
For **Direct extraction**, add 100 µL of Solution IRS.

9.6 Vortex briefly to mix (speed 7 out of 10). Incubate at 2–8°C for 5 min. ⌚ 00:05:00

5m

Place the tubes in a cold rack (stays in the refrigerator).

9.7 Centrifuge at maximum speed for 1 min ⚙️ 150000 rpm, Room temperature , 00:01:00 . Avoiding the pellet and transfer the supernatant to a new Collection Tube.

9.8 Add 650 µL each of Solution PM3 and Solution PM4. Vortex by inverting the tubes 15 times.

Visually, check that the 3 solutions are well mixed.

9.9 Load the mixture into an MB Spin Column.

9.9.1 **Using centrifuge:**

1m

- Load 650 µL of the mixture into an MB Spin Column


- Centrifuge at max speed for 1 min 🌀 **150000 rpm, Room temperature , 00:01:00** . Discard the flow-through and repeat until all the mixture has been loaded onto the MB Spin Column.

### 9.9.2 Using manifold:

1m

- Load the all mixture into an MB Spin Column (see picture)

- Centrifuge at max speed for 1 min 🌀 **150000 rpm, Room temperature , 00:01:00** . Discard the flow-through.

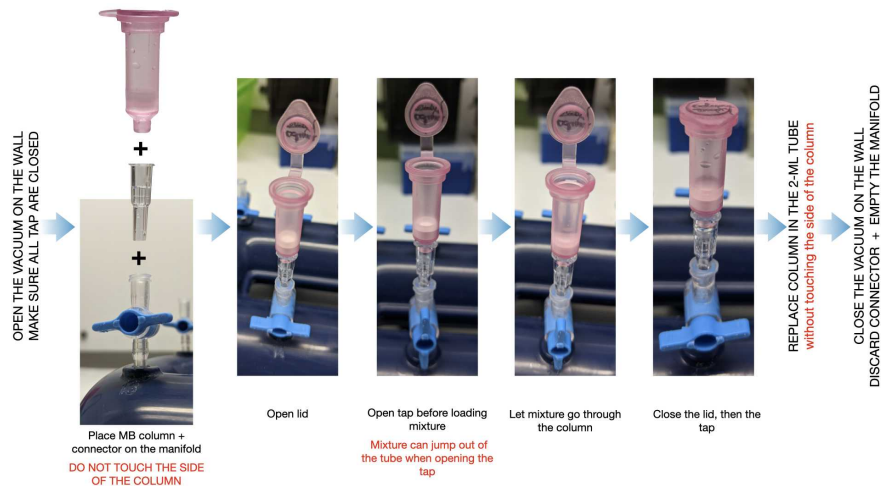

9.10 Shake well to mix Solution PM5 and add 650 µL to the MB Spin Column.

9.11 Centrifuge at max speed for 1 min 🌀 **150000 rpm, Room temperature , 00:01:00** .

9.12 Discard flow-through. Add 600 µL of Solution PM4.

9.13 Centrifuge at max speed for 1 min 🌀 **150000 rpm, Room temperature , 00:01:00** .

9.14 Discard flow-through and centrifuge filter at max speed for an additional 2 min  
🌀 **150000 rpm, Room temperature , 00:02:00** .

9.15 Place the MB Spin Column in a clean 2-ml Collection Tube (provided in the RNeasy PowerMicrobiome kit).

9.16 Add 60 µL of RNase-Free Water (warmed to 55°C) to the center of the MB Spin Column membrane. Incubate at room <sup>5m</sup>



temperature for 5 min ⌚00:05:00 .

9.17 Centrifuge at max speed for 1 min ⚙️150000 rpm, Room temperature , 00:01:00 . Discard the MB Spin Column.



The DNA/RNA is now ready for downstream applications. RNA extract may be stored in RNase-free water at -80°C for 1 year.

## Quantification

10



Quantification by Droplet Digital PCR (ddPCR)  
by McLellan Lab

PREVIEW

RUN



10.1 When all reagents are thawed on ice, vortex Supermix, Reverse transcriptase and DTT thoroughly for 30 seconds. Vortex to mix primers and probes stocks.

10.2



Plate Set up: see sub-steps below.



This step requires to keep the mixture cold/on ice.

10.2.1 Prepare the reaction matrix (for one well, beside sample RNA) according to the table below. Prepare Use a low-binding tube of appropriate volume to mix all the components according to the reaction numbers. Always include extra wells when setting up reaction to avoid potential volume shortage caused by pipetting.

Component	Volume per reaction, $\mu\text{L}$	Final concentration
Supernatant	5.5	1x
Reverse transcriptase	2.2	20 U/ $\mu\text{L}$
300 mM DTT	1.1	15 mM
Primer mix (forward + reverse)	1.1	900 nM
Probe	1.1	250 nM
RNase-free water*	5.5	/
Total	16.5	/

\* Note: Water volume can be replaced accordingly by other ingredients, such as another assay (e.g., duplex assay) or another RNA template (e.g., inhibition test).



1. Keep all reagents on ice during the process as well as the matrix.
2. Always prepare a total of  $8 \times (\text{N columns})$  wells for droplet generation, or use ddPCR Buffer Control for Probes (#1863052) to fill empty wells on the last cartridge.
3. Make sure sticky reagents are added the correct volumes and not left in the tips, i.e., Supernatant, RT.

- 10.2.2** Place a 96-well PCR plate onto a low temp PCR rack, or on ice. Disperse 16.5  $\mu\text{L}$  of reaction matrix into each well of a 96-well PCR plate. For runs with multiple columns, calculate (the matrix total volume/8) and evenly distribute the matrix into an 8-well PCR strip. Then use an appropriate multichannel pipette (e.g., 2-20  $\mu\text{L}$  range) to add 16.5  $\mu\text{L}$  of matrix into each well, column by column.



1. Keep the PCR plate cold/on ice during the process.
2. Lower the pipetting speed to avoid liquid leftover to the inner side of the tips.

- 10.2.3** Gently vortex at half speed to mix the RNA sample. Make sure no liquid is attached on the lid. Add 5.5  $\mu\text{L}$  of sample RNA into each well containing 16.5  $\mu\text{L}$  of reaction matrix, making the total volume of each reaction 22  $\mu\text{L}$ .



1. Keep the PCR plate and RNA samples cold/on ice throughout the process.
2. Do not over vortex the RNA samples.
3. Pipetting robot may be used to add RNA samples to the matrix.

- 10.2.4** Seal the PCR plate. Centrifuge down gently at 1000 rpm for 30s. Take out the plate and vortex on a 96-well plate mixer at 1600 rpm for 30s. Centrifuge again at 1000 rpm for 30s to settle down the plate.

## 10.3



Droplet generation handling.



It is NOT required to keep the plate cold/on ice during this process.

- 10.3.1** Prepare materials/reagents on the working bench top, e.g., cartridges, gaskets, ddPCR 96-well plate, droplet generation oil, foil cover. Label the cartridges with corresponding column numbers (e.g., 1 to 12 for a full plate).
- 10.3.2** Remove the sealing on the PCR plate. Place the cartridge in the cartridge holder. Align well the PCR plate and the cartridge on the bench.  
Use an appropriate multichannel pipette (e.g., 2- 20 µL) for liquid transfer. Adjust the pipette at 20 µL and make sure the tips are well positioned. Gently mix the liquid by aspirating up to 2/3 of the tip height and then releasing to a lower level of the tip height, i.e., not to the end of the tip, to avoid creating bubbles. Repeat this mixing step 10 times. In the last movement, slowly aspirate to the full volume of 20 µL.
- 10.3.3** Transfer the 20 µL reaction matrix to the middle column of the cartridge. Position the tip end to the ridge in well (where the well wall connects to the bottom) at 15° angle. Avoid creating bubbles when releasing liquid from tips into the cartridge; this can be realized by only pressing the plunger to position 2 and not position 3 before pull the tips out from the cartridge. Make sure the cartridge wells are in the same direction as on the PCR plate.
- 10.3.4** Fill in 70 µl of Droplet Generation Oil into bottom wells of the cartridge and cover the cartridge with a red gasket. Loop outer holes of red gasket around hooks on left and right sides of cartridge holder. Place the gasket equipped cartridge into the droplet generator.



You can start preparing the next column of droplet generation while waiting for the previous column to be done.

- 10.3.5** When droplet generation is done, take out the cartridge from the droplet generator and remove the red gasket. Using a Rainin multichannel with recommended Rainin tips, in a leaning position, count to 5 to aspirate all the liquid (i.e., 40 µl) from the droplets column, and press against side of wells of the corresponding column in the ddPCR 96-well plate (i.e., not the previous PCR plate), count to 5 to expel the droplets into the wells.



1. Avoid multiple times of liquid transfer.
2. Eye ball the ddPCR plate when droplet generation and transferring are done.

## 10.4

Turn on the PX1 PCR plate sealer and let heat to 180°C. Correctly place the plate support block, the ddPCR plate, foil cover (i.e., red line up) and the metal holder. Seal the plate at 180°C for 5s and remove the plate immediately from the sealer.



Sealed plate should have indentations around walls. Always check the sealing before loading the plate on the

Sealed plate should have indentations around wells. Always check the sealing before loading the plate on the PCR thermal cycler.

- 10.5 Load the plate on to a PCR thermocycler. Our lab's assay conditions are shown as below. We use Eppendorf Mastercycler Pro and the ramp speed is set to 50% for RT-ddPCR.

Assay	Step 1	Step 2	Step 3 (40 cycles)	Step 4 (40 cycles)	Step 5	Step 6
N1/N2, BCoV, HepG	50°C 60 min	95°C 10 min	94°C 30s	55°C 1 min	98°C 10 min	4°C 30 mins and hold
PMMoV	50°C 60 min	95°C 10 min	94°C 30s	60°C 1 min	98°C 10 min	4°C 30 mins and hold

- 10.6 List of assay primers and probes.

A	B	C	D
CDC N1	GACCCCAAAATCAGCGAAAT	FAM- ACCCCGCATTACGTTTGGTGGACC- BHQ1	TCTGGTTACTGCCAGTTGAATCTG
CDC N2	TTACAAACATTGGCCGCAAA	HEX-ACAATTTGCCCGAGCGCTTCAG- IowaBHQ	GCGCGACATTCCGAAGAA
BCoV	CTGGAAGTTGGTGGAGTT	FAM - CCTTCATATCTATACACATCAAGTTGTT- BHQ1	ATTATCGGCCTAACATACATC
HepG	CGGCCAAAAGGTGGTGGATG	HEX- AGGTCCCTCTGGCGCTTGTGGCGAG- BHQ1	CGACGAGCCTGACGTCGGG
PMMoV	GAGTGGTTTGACCTTAACGTTGA	FAM-CCTACCGAAGCAAATG- BHQ1	TTGTCGGTTGCAATG CAA GT