

Jul 25, 2024 Version 2

T cell purification and activation V.2

DOI

dx.doi.org/10.17504/protocols.io.81wgbz431gpk/v2

Moustafa Nouh Elemeery^{1,2}, Salix Boulet³, louis-eric.trudeau⁴, Nathalie Labrecque⁵

¹Department of neurosciences, Faculty of Medicine, Université de Montréal, Canada;

²Medical Biotechnology Department, National Research Centre, Dokki, Cairo, Egypt;

³Centre de recherche de l'hôpital Maisonneuve-Rosemont (CRHMR), Department of Medicine, Faculty of Medicine, Université de Montréal, Canada;

⁴Department of Pharmacology and Physiology, Faculty of Medicine, Université de Montréal, Canada;

⁵Institut de recherches cliniques de Montréal (IRCM), Department of Microbiology, Infectious Diseases and Immunology, Faculty of Medicine, Université de Montréal, Canada

ASAP Collaborative Rese...



Lilia Rodriguez

Université de Montréal

OPEN  ACCESS



DOI: dx.doi.org/10.17504/protocols.io.81wgbz431gpk/v2

Protocol Citation: Moustafa Nouh Elemeery, Salix Boulet, louis-eric.trudeau, Nathalie Labrecque 2024. T cell purification and activation. **protocols.io** <https://dx.doi.org/10.17504/protocols.io.81wgbz431gpk/v2> Version created by **Lilia Rodriguez**

Manuscript citation:

Elemeery, M. N., Chung, A., Boulet, S., Mukherjee, S., Giguere, N., Daudelin, J. F., ... & Trudeau, L. E. (2024). Adoptive transfer of mitochondrial antigen-specific CD8⁺ T-cells in mice causes parkinsonism and compromises the dopamine system. *bioRxiv*, 2024-02.

License: This is an open access protocol distributed under the terms of the **Creative Commons Attribution License**, which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working

We use this protocol and it's working

Created: July 02, 2024



Last Modified: July 25, 2024

Protocol Integer ID: 104114

Keywords: ASAPCRN

Funders Acknowledgement:

Aligning Science Across

Parkinson's (ASAP)

Grant ID: ASAP 000525

Abstract

This protocol details the purification and activation of Mouse Naïve CD8+ T Cell using an isolation kit from STEMCELL. This kit is designed to isolate naïve CD62L+CD44-CD8+ T cells from single-cell suspensions of splenocytes by negative selection. Unwanted cells are targeted for removal with biotinylated antibodies that are directed against non-naïve CD8+ T cells (CD4, CD11b, CD11c, CD19, CD44, CD45R/B220, CD49b, TCR γ/δ , TER119) and streptavidin-coated magnetic particles. Labeled cells are separated using an EasySep™ magnet without the use of columns.

Materials

Reagents and solution:

A	B	C
Reagent or solution	Supplier	Catalogue #
Anti-mouse CD3 (Clone: 145-2C11)	BioXcell	BE0001
Anti-mouse CD28 (clone: 37.51)	Leinco	C379-5.0 mg
RPMI media	VWR (Corning)	CA45000-396
Phosphate Buffer saline (PBS)	Gibco	14190144
EasySep buffer	STEMCELL	20144
EasySep™ Mouse CD8+ T Cell Isolation Kit	STEMCELL	19858
ACK (Ammonium-Chloride-Potassium) Lysing Buffer	ThermoFisher	A1049201
L-glutamine 200 mM	VWR	CA45000-676
HEPES 1M	Fisher	MT25060CI
Sodium Pyruvate	VWR	CA45000-710
2-Mercaptoethanol	ThermoFisher	21985023
Non-Essential Amino acids	VWR (Corning)	CA45000-700
Fetal Bovin Serum	Gibco	12483020
24 wells plate flat bottom suspension plates	Sarstedt	83.1836.500
96 wells plate flat bottom suspension plates	Sarstedt	82.1581.001
Cell strainers (70um)	Fisher	08-771-2

⊗ InVivoMAb anti-mouse CD3 **BioXcell Catalog #BE0001**

⊗ Anti-Mouse CD28 (Clone 37.51) **Leinco Catalog #C379**

⊗ RPMI **VWR International Catalog #45000-396**

⊗ Phosphate buffered saline (PBS) without Ca/Mg **Thermo Fisher Scientific Catalog #14190144**

⊗ EasySep™ Buffer **STEMCELL Technologies Inc. Catalog #20144**

⊗ EasySep™ Mouse Naïve CD8+ T Cell Isolation Kit **STEMCELL Technologies Inc. Catalog #19858**

⊗ ACK Lysing Buffer **Thermo Fisher Scientific Catalog #A1049201**

⊗ L(+)-Glutamine solution 200 mM **VWR International Catalog #45000-676**

⊗ Sodium pyruvate solution 100 mM **VWR International Catalog #45000-710**

⊗ 2-mercaptoethanol **Gibco - Thermo Fisher Catalog #21985023**

⊗ Non essential amino acids **VWR International Catalog #45000-700**

⊗ Fetal Bovine Serum, qualified, Canada **Thermo Fisher Catalog #12483020**

⊗ Falcon™ Cell Strainers **Fisher Scientific Catalog #08-771-2**

**RPMI complete (RPMIc):**

A	B
RPMI	500 mL
FBS 10%	50 mL (decomplemented)
L-Glutamine	5 mL
Sodium pyruvate	5 mL
Antibiotic (Pen-Strep)	5 mL
Non-essential amino acids	5 mL
2-Mercaptoethanol	50 mmol/L (final)


















Note

Note! Very important, 2-Mercaptoethanol is an essential growth factor for mouse T-lymphocytes.







Purification and activation

4d 0h 30m

- 1 Dilute CD3 antibody (145-2C11) (clone KT3 can also be used) to  1 undetermined in PBS.
- 2 Coat plates with anti-CD3 antibody.
- 3 Use 24 or 96 wells flat bottom suspension plates. If these plates are not used, there is a risk of partial stimulation due to low absorbance of the antibody on the plate.
- 3.1 Add  100 μ L of antibody/well if 96 wells or  1 mL per well if 24 wells. 
- 3.2 Incubate the plate for  24:00:00 at  4 °C or  01:00:00 at  37 °C . 
1d
- 4 Remove the antibody (aspirate) and wash 2 times with PBS. 
- 5 Add  100 μ L of RPMIc if 96 wells or  1 mL if 24 wells and incubate the plate at  37 °C for  00:15:00 . 
15m
- 6 Collect spleen from mice 6-8 weeks in complete RPMI media (RPMIc), sex matched with recipient mice.
- 7 Purify CD8+ T cells using EasySep mouse naïve purification kit (STEMCELL, catalo # 19858) as follow:
- 7.1 Use a frosted microscope slide to homogenize spleens in PBS or Hanks' Balanced Salt Solution (HBSS) containing 2% fetal bovine serum (FBS).
- 7.2 Remove aggregates and debris by passing cell suspension through a 70 μ m mesh nylon strainer. Collect cells in a 15 mL tube.
- 7.3 Centrifuge at  1300 rpm, 00:05:00 and discard the supernatant. 
5m




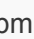
7.4 Red blood cells lysis is done by adding  5 mL /spleen of 0.83% ammonium chloride and incubate for  00:05:00 (or  00:02:00 for ACK lysing buffer) at  Room temperature while continuously shaking tubes.



5m



7.5 Quench by filling up tube with RPMIc.

7.6 Centrifuge at  1300 rpm,  00:05:00 and discard supernatant.


5m





7.7 Resuspend in EasySep buffer at 1×10^8 cells/ml (Easysep mouse naïve CD8+ T cell isolation kit (#19858A) and follow the protocol provided by STEMCELL.


8 Wash cells with RPMIc medium.





9 Resuspend cells in  5 mL of RPMIc, count cells and adjust the concentration to 2×10^6 cells/mL

10 Add  100 μ L of cells (2×10^5 cells) per well if 96 wells or  1 mL of cells (2×10^6 cells) per well if 24 wells).



11 Add purified anti-CD28 antibodies to reach a concentration of  5 undetermined .




12 Incubate the plate at  37 °C and 5% CO₂ for  24:00:00 .

1d



13 Add 20 Unit/ml of IL-2.



14 Incubate the plate at  37 °C and 5% CO₂ for another  48:00:00 .

2d



**Note******* Don't forget the non-stimulated controls.**

- Resuspend the cells well throughout the experiment (the cells quickly settle to the bottom of the tube)
- Check purity of CD8+ T cells after STEMCELL isolation by staining with anti-CD8 antibody followed by flow cytometry.
- Check activation after stimulation using anti-CD8 and anti-CD44 staining followed by flow cytometry.