

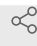



Oct 05, 2022

# Single-cell RNA-seq for mDA neurons

gurvir.virdi<sup>1</sup><sup>1</sup>UCL Institute of Neurology

1 Works for me

 Share[dx.doi.org/10.17504/protocols.io.6qpvr4dxpgmk/v1](https://dx.doi.org/10.17504/protocols.io.6qpvr4dxpgmk/v1) gurvir.virdi

## ABSTRACT

Harvesting and performing single-cell RNA-seq on hiPSC-derived mDA neurons.

## DOI

[dx.doi.org/10.17504/protocols.io.6qpvr4dxpgmk/v1](https://dx.doi.org/10.17504/protocols.io.6qpvr4dxpgmk/v1)

## PROTOCOL CITATION

gurvir.virdi 2022. Single-cell RNA-seq for mDA neurons. **protocols.io**  
<https://dx.doi.org/10.17504/protocols.io.6qpvr4dxpgmk/v1>



## LICENSE

————— This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

## CREATED

Oct 05, 2022

## LAST MODIFIED

Oct 05, 2022

## PROTOCOL INTEGER ID

70858

### Harvesting cells for RNA-seq 5m

- 1 At the desired age of mDA neurons, they are harvested for single-cell RNA-seq:

#### 1.1 mDA neurons are washed 1x in PBS.

1.2 They are incubated with Accutase (Thermo Fisher Scientific) for 5m 00:05:00 at 37 °C .

1.3 mDA neurons are collected as a single cell suspension and diluted 1:3.

#### Quality and concentration of cells

- 2 The quality and concentration of each single-cell suspension was measured using Trypan blue and the Eve automatic cell counter.
- 3 Each sample presented a concentration between a 1200-1700 cell/μl and viability ranged between 55-68%, samples with a viability above 57% were used for sequencing.

#### single-cell RNA-seq

- 4 Approximately 10000 cells were loaded for each sample into a separate channel of a Chromium Chip G for use in the 10X Chromium Controller (cat: PN-1000120).
- 5 The cells were partitioned into nanoliter scale Gel Beads in emulsions (GEMs) and lysed using the 10x Genomics Single Cell 3' Chip V3.1 GEM, Library and Gel Bead Kit (cat: PN-1000121).
- 6 cDNA synthesis and library construction were performed as per the manufacturer's instructions.
- 7 The RNA was reversed transcribed and amplified using 12 cycles of PCR.
- 8 Libraries were prepared from 10μl of the cDNA and 13 cycles of amplification. Each library was prepared using Single Index Kit T Set A (cat: PN-1000213) and sequenced on the HiSeq4000 system (Illumina) using 100 bp paired-end run at a depth of 65-100 million reads. Libraries were generated in independent runs for the different samples.