



JUL 04, 2023

OPEN ACCESS

Protocol Citation: Koen Theunis, Suresh Poovathingal, Sarah Geurs 2023. Nuclei isolation and permeabilisation of fresh frozen human brain samples for 10X Genomics Multiome. **protocols.io** <https://protocols.io/view/nuclei-isolation-and-permeabilisation-of-fresh-frozen-human-brain-samples-for-10x-genomics-multiome-cwn7xdhn>

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working
We use this protocol and it's working

Created: Jul 03, 2023

Last Modified: Jul 04, 2023

PROTOCOL integer ID:
84415

Keywords: ASAPCRN, single nuclei preparation, multiome, snATAC-seq, snRNA-seq, 10X Genomics

Nuclei isolation and permeabilisation of fresh frozen human brain samples for 10X Genomics Multiome

Koen


Theunis¹, Suresh Poovathingal², Sarah Geurs³

¹Laboratory of Computational Biology, VIB-KU Leuven, Center for Brain & Disease Research, Leuven, Belgium;

²VIB-KU Leuven/VIB Center for Brain & Disease Research, Leuven, Belgium;

³VIB Tech Watch, VIB Technologies, Ghent, Belgium

Aertslab

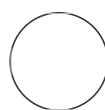

Koen Theunis

ABSTRACT

This protocol details the procedure of nuclei isolation and debris removal from fresh frozen brain tissue in preparation for the 10X Genomics multiome ATAC/GEX assay.

MATERIALS

HOW MANY SAMPLES TO PROCESS? 2

A	B	C	D	E	F
2X salt-Tris solution	Stock conc		Final conc		V (µl)
NaCl	5000	mM	292	mM	1226,4
Tris pH7.5	1000	mM	20	mM	420,0
CaCl	1000	mM	2	mM	42,0
MgCl ₂	1000	mM	42	mM	882,0
Water					1842,6
Total Volume					2100,0

TST+NP Homogenization buffer	Stock conc		Final conc		V (µl)
2X salt-Tris solution	2	X	1	X	1100,0
Tween-20	10	%	0,03	%	6,6

NP-40	10	%	0,03	%	6,6
BSA	7,5	%	1	%	293,3
KCl	3000	mM	25	mM	18,3
Sucrose	2000	mM	250	mM	275,0
DTT	1000	mM	1	mM	2,2
cOmplete protease inhibitor	50	X	1	X	44,0
Rnase In plus (promega)	40	U/u	1	U/u	55,0
Water					398,9
Total Volume					2200,0

Wash buffer 1	Stock conc		Final conc		V (μl)
2X salt-Tris solution	2	X	1	X	561,0
BSA	7,5	%	1	%	149,6
KCl	3000	mM	25	mM	9,4
Sucrose	2000	mM	250	mM	140,3
DTT	1000	mM	1	mM	1,1
cOmplete protease inhibitor	50	X	1	X	22,4
Rnase In plus (promega)	40	U/u	1	U/u	28,1
Water					210,2
Total Volume					1122,0

GM (Gradient Medium)	Stock conc		Final conc		V(μl)
CaCl	1000	mM	1	mM	1,3
Optiprep	60	%	50	%	1083,3
MgCl ₂	1000	mM	5	mM	6,5
Tris pH 7.5	1000	mM	10	mM	13,0
Sucrose	2000	mM	75	mM	48,8
DTT	1000	mM	1	mM	1,3

cOmplete protease inhibitor	50	X	0,5	X	13,0
Rnase In plus (promega)	40	U/u 	1	U/u 	32,5
Water					100,3
Total Volume					1300

ODM (Optiprep Diluent Medium)	Stock conc		Final conc		V(μl)
KCl	3000	mM	150	mM	62,0
MgCl ₂	1000	mM	30	mM	37,2
Tris pH 8	1000	mM	60	mM	74,4
Sucrose	2000	mM	250	mM	155,0
RNAse in plus (promega)	40	U/u L	0,5	U/u L	15,5
Water					895,9
Total Volume					1240,0

29% Cushion	Stock conc	Final conc	V(μl)		
Optiprep	60	%	29	%	966,7
ODM					1033,3
Total Volume					2000

1x Lysis	Stock conc		Final conc		V (μl)
Tris-HCl (pH 7.4)	1000	mM	100	mM	50,0
NaCl	5000	mM	100	mM	10,0
MgCl ₂	1000	mM	30	mM	15,0
Tween-20	10	%	0,1	%	5,0
NP-40	10	%	0,1	%	5,0
Digitonin	5	%	0,02	%	2,0
Water					413,0
Total Volume					500,0

0.1x Lysis	Stock conc		Final conc		V (μl)
1x lysis	1	x	0,1	x	100,0
BSA	7,5	%	1	%	133,3
DTT	1000	mM	1	mM	1,0
Rnase Inhibitor (Promega)	40	U/μl	1	U/μl	25,0
Water					740,7
Total Volume					1000,0

Wash Buffer 2	Stock conc		Final conc		V (μl)
Tris-HCl (pH 7.4)	1000	mM	10	mM	22,0
NaCl	5000	mM	10	mM	4,4
MgCl ₂	1000	mM	3	mM	6,6
BSA	7,5	%	1	%	293,3
Tween-20	10	%	0,01	%	2,2
DTT	1000	mM	1	mM	2,2
Rnase Inhibitor (SIGMA)	40	U/μl	1	U/μl	55,0
Water					1814,3
Total Volume					2200,0

1X Diluted Nuclei buffer	Stock conc		Final conc	V(μl)	
Nuclei buffer	20	X	1	X	7,5
DTT	100	mM	1	mM	1,5
Rnase Inhibitor (SIGMA)	40	U/μl	1	U/μl	3,75
Water					137,25
Total Volume					150

Preparation

- 1 Prepare the 2X salt-Tris solution and aliquot to 1 ml tube and freeze at -20°C
- 2 Prepare the BSA in PBS and filter through 0.2um filter.
- 3 Prepare the Tween-20 Homogenization lysis buffer fresh for each preparations

- 4 Filter sucrose 2M solution through 0.2um/0.45um filter
- 5 Prepare cOmplete protease inhibitor (50X) by dissolving a tablet in 1 mL of water.
- 6 Filter protease inhibitor solution through 0.2um filter
- 7 Place the homogenizer at -80°C at the start and leave 10 min before homogenization on ice.

Homogenization (Perform all steps on ice)

- 8 Cut brain piece on dry ice (if a large section), then transfer immediately in homogenizer containing 250 uL of TST+NP Homogenization Lysis Buffer. Immediately add 500uL TST+NP Homogenization Lysis Buffer
- 9 Allow the tissue to thaw for 2 mins in the homogenization buffer.
- 10 Homogenize tissue with pestle A (10X), wait 1min, and with pestle B (10X).
- 11 Put the homogenate through 70µm cell strainer and place the falcon on ice. Incubate on ice for 5 mins. Rinse dounce and filter with 250uL TST+NP HB.
- 12 During the incubation (in about 4 mins), use LUNA-fl counter to assess the viability of the nuclei.

- 13 Transfer the contents from the 50 mL falcon to a 2 mL protein lo-bind tube.
- 14 Centrifuge @ 500xg for 5 mins. Discard the supernatant.
- 15 Resuspend the pellet in 200 uL of Wash Buffer 1, and transfer to a 2mL DNA lobind tube.
- 16 Add additional 320 uL Wash Buffer 1 to a final volume of 520 uL.
- 17 Add 520uL of Gradient medium to sample ($V_f = 1040$ uL).

Isolation by centrifugation

- 18 Layer 770 uL of 29% Cushion in the ultracentrifuge tube.
- 19 Layer the sample on top of cushion using a P1000, without disturbing the cushion.
- 20 Optional: Check tube weight, adjust weight to counter differences.

- 21 Centrifuge at least 3,000 rcf in a swinging bucket, at 4°C for 20 minutes with brake off.
- 22 Remove supernatant, remove the lower supernatant with P200, leaving about 50-100 uL.

Resuspension and permeabilisation

- 23 Gently resuspend nuclei in ultracentrifuge tube and transfer to 1.5 mL DNA lobind tube.
- 24 Rinse ultracentrifuge tube with Wash buffer 2 and transfer also to 1.5 mL DNA lobind tube.
- 25 Centrifuge @ 350-450xg for 5 mins. Discard the supernatant.
- 26 Resuspend in 0.1x lysis buffer (200 uL) and gently pipetmix 5x
- 27 Incubate on ice for 2 min
- 28 Add Wash buffer 2 (1mL) & gently pipetmix 5x

29 Centrifuge @ 350-450xg for 5 mins. Discard the supernatant.

30 Resuspend in 1X Diluted Nuclei buffer. Count