

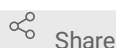


Jul 12, 2022

Detection of JEV and WNV by quantitative real-time reverse-transcription (qRT)-PCR

Wenjing Liu¹, Shihong Fu²¹Department of Basic Medicine, School of Medicine, Qingdao University, Qingdao, People's Republic of China;²Department of Arbovirus, NHC Key Laboratory of Biosafety, National Institute for Viral Disease Control and Prevention, Chinese Center for Disease Control and Prevention, Beijing, People's Republic of China

1 Works for me

dx.doi.org/10.17504/protocols.io.bdayi2fw Wenjing Liu

ABSTRACT

you can use this protocol to detect the RNA of your sample

DOI

dx.doi.org/10.17504/protocols.io.bdayi2fw

PROTOCOL CITATION

Wenjing Liu, Shihong Fu 2022. Detection of JEV and WNV by quantitative real-time reverse-transcription (qRT)-PCR. **protocols.io**
<https://dx.doi.org/10.17504/protocols.io.bdayi2fw>

LICENSE



————— This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

CREATED

Mar 06, 2020

LAST MODIFIED

Jul 12, 2022

PROTOCOL INTEGER ID

33848

MATERIALS TEXT

RNA of mosquito grinding supernatant , serum sample and CSF samples,AgPath-IDTM One-step RT-PCR Kit (AM1005, Thermo Fisher Scientific, USA) ,Stratagene real-time PCR instrument (model MX300; Thermo Fisher Scientific, Waltham, MA, USA) .

BEFORE STARTING

RNA was extracted from all collected samples (serum, CSF, and mosquito grinding supernatant) using a Tianlong Nucleic Acid Automatic Extractor with the Tianlong Nucleic Acid Extraction Kit(EX-RNA/DNA virus, Suzhou Tianlong Biotechnology Co., Ltd, Suzhou, China).

Detection of JEV and WNV by quantitative real-time reverse-transcription (qRT)-PCR

- 1 The instrument used was a Stratagene [™], real time PCR machine (manufacturer: Thermo Scientific, model: MX300) for quantitative fluorescence qRT-PCR detection of JE virus and West Nile virus. AgPath-IDTMOne-step RT-PCR Kit (REF: AM1005, manufacturer: Thermo Fisher Scientific.) was the detection kit used, the total volume of the detection system was 20 µl, containing 10 µl of 2Xbuffer4 µl of RNase Free Water4 µl, 1 µl of primers, probes, and enzyme mixture (Mix), and 2 µl of template RNA.
- 2 The reaction program was: 45 ° C for 10 min, 95 ° C for 10 min for 1 cycle, 95 ° C for 15 s and 60 ° C for 1 min for 40 cycles. Using primers and probes specific for the JE virus gene and the West Nile virus gene.
- 3 The final result is checked on the real time PCR machine,the sample with CT value over 35 thought to be negative, sample with CT value lower than 35 thought to be positive.
- 4 The positive sample will be sent to do PCR specific to JEV genes to be further tested.