




NOV 02, 2023

🌐 TEM

Alexander  
Röntgen<sup>1</sup>

<sup>1</sup>University of Cambridge

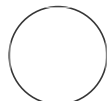

Alexander Röntgen

ABSTRACT

TEM protocol.

OPEN  ACCESS



**DOI:**

[dx.doi.org/10.17504/protocols.io.q26g7pd79gwz/v1](https://dx.doi.org/10.17504/protocols.io.q26g7pd79gwz/v1)

**Protocol Citation:** Alexander Röntgen 2023. TEM.

**protocols.io**

<https://dx.doi.org/10.17504/protocols.io.q26g7pd79gwz/v1>

**License:** This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

**Protocol status:** Working  
We use this protocol and it's working

**Created:** Nov 01, 2023

**Last Modified:** Nov 02, 2023

## Preparation

- 1 Recover samples by pipetting out of aggregation plate
- 2 Glow-discharge grids (3-mm 300-mesh (EM Resolutions Ltd))

## Sample spotting and staining

- 3 Take up grid with tweezers
- 4 Spot 2.3  $\mu$ L sample on grid
- 5 Stain sample with 2.3  $\mu$ L 2% uranyl acetate for 40 s at room temperature

## Imaging

- 6 Image grids on an electron microscope (FEI Tecnai G2 transmission electron microscope)