

Oct 28, 2021

# Purification protocol of Mouse (*Mus Musculus*) E1-like enzyme ATG7

Dorotea Fracchiolla<sup>1</sup>

<sup>1</sup>Sascha Martens lab, University of Vienna, Max Perutz Labs - Vienna (AT)

1



[dx.doi.org/10.17504/protocols.io.bsennbde](https://dx.doi.org/10.17504/protocols.io.bsennbde)



Dorotea Fracchiolla

Sascha Martens lab, University of Vienna, Max Perutz Labs - ...

This protocol outlines the procedures for expression and purification of Mouse (*Mus Musculus*) E1-like enzyme ATG7 (AuTophagy-related protein) of the ATG8 ubiquitin-like conjugation system in autophagy.

[ASAP\\_Team\\_Hurley\\_Martens\\_lab\\_Dorotea\\_Fracchiolla\\_Exp\\_purif\\_mATG7.pdf](#)

DOI

[dx.doi.org/10.17504/protocols.io.bsennbde](https://dx.doi.org/10.17504/protocols.io.bsennbde)

Dorotea Fracchiolla 2021. Purification protocol of Mouse (*Mus Musculus*) E1-like enzyme ATG7.

**protocols.io**

<https://dx.doi.org/10.17504/protocols.io.bsennbde>

protocol

<https://doi.org/10.1083/jcb.201912098>

expression, purification, E1-like enzyme ATG7, autophagy, ASAPCRN

protocol ,

© Dorotea Fracchiolla

Feb 14, 2021

Oct 28, 2021

Feb 14, 2021



Emily Hasser

University of Washington

Feb 25, 2021




Dorotea Fracchiolla

Sascha Martens lab, University of Vienna, Max Perutz Labs - Vienna (AT)



47278

**General information:** expression system: Sf9 cells, medium: SF921+Pen/Strep; plasmid origin: Gammoh Noor (University of Edinburgh); internal lab Construct Database number: SMC911; backbone: pFast-BacHT(B); resistance: Ampicillin/Gentamicin; insert: *Mus musculus* ATG7, Isoform 2\_NP\_00124647.1; tags & cleavage sites: N-term 6xHis, followed by Tobacco Etch Virus (TEV) cleavage site, 6xHis-TEVcs-mATG7; protein length without tags: 698 aa; Epsilon (all Cysteines reduced, with tags) = 88810 M<sup>-1</sup> cm<sup>-1</sup>; MW with tags= 83 kDa; ORF (only mATG7):


MGDPGLAKLQFAPFNSALDVGFWHELTQKKLNEYRLDEAPKDIKGYYYNGDSAGLPTRLTLEFSAFDMSA  
STPAHCCPAMGTLHNTNTLEAFKTADKKLLLEQSANIWEAIKSGAALNPMLLNKFLLLTFADLKKYHF  
YYWFCCPALCLPESILIRGPVSLDQRLSPKQIQALEHAYDDL CRAEGVTALPYFLFKYDDDTVLVSLLK  
HYSDDFFQGQRKITVGVYDPCNLAQYPGWPLRNFLVLAHRWSGSFQSVEVLCFRDRTMQGARDVTHSII  
FEVKLP EMAFSPDCPKAVGW EKNQKGGMGPRMVNLSGC MDPKRLAESSVDLNLKLMCWRLVPTLDLDKVV  
SVKCLLLGAGTLGCNVARTLMGWGVVRHVTFVDNAKISYSNPVRQPLYEFEDCLGGGKPKALAAAERLQKI  
FPGVNARGFNMSIPMPGHPVNFSDVTMEQARRDVEQLEQLIDNHDVIFLLMDTRESRWLPTVIAASKRKL  
VINAALGFDTFVVMRHGLKKPKQQGAGDLCPSHLVAPADLGSSLFANIPGYKLGCFYFCNDVVAPGDSTRD  
RTL DQQCTVSRPGLAVIAGALAVELMVSVLQHPEGGYAIASSSDDRMNEPPTSLGLVPHQIRGFLSRFDN  
VLPVSLAFDKCTACSPKVL DQYEREGFTFLAKVFNSSH SFLDLTGLTLLHQETQAAEIWDMSDEETV

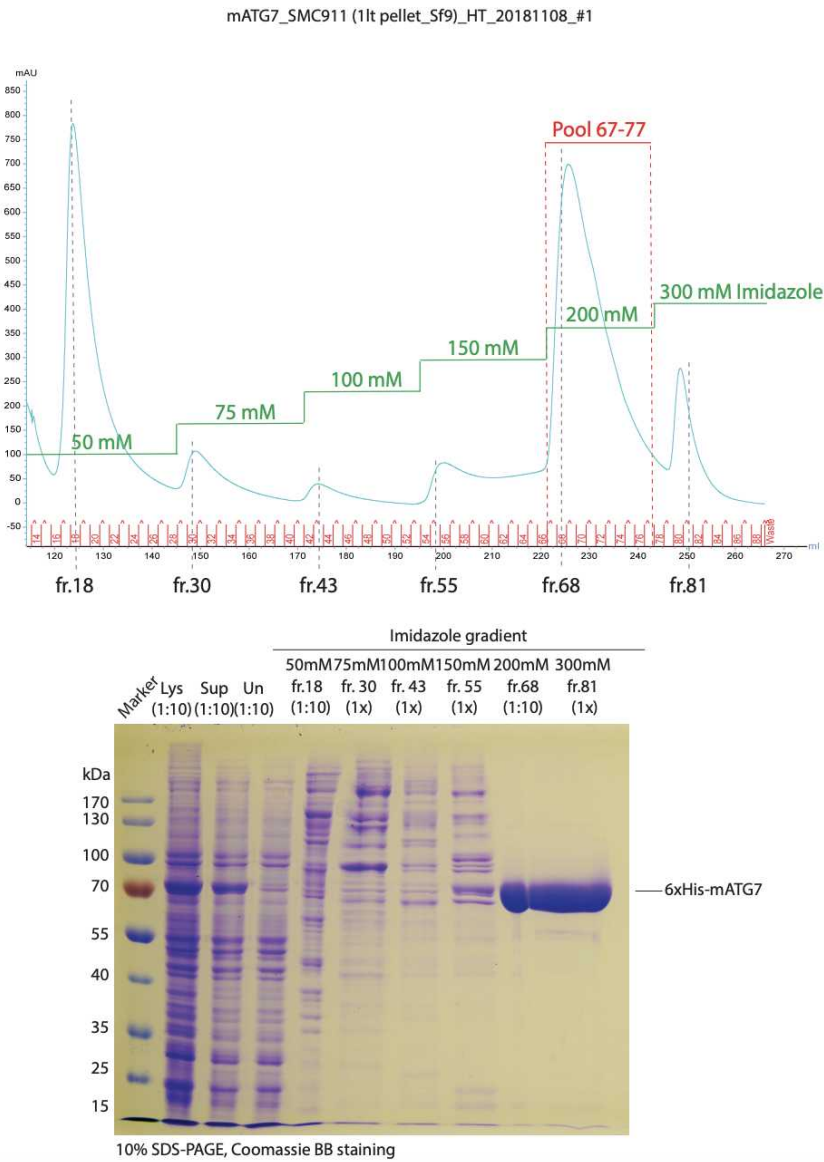

Chromatograph and Coomassie BB stained gel of His-tag affinity purification.

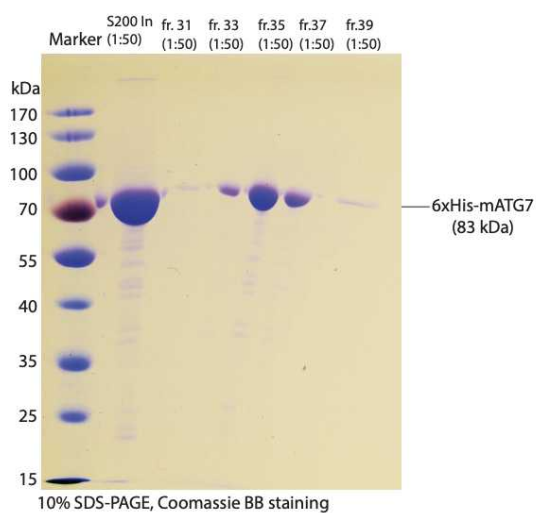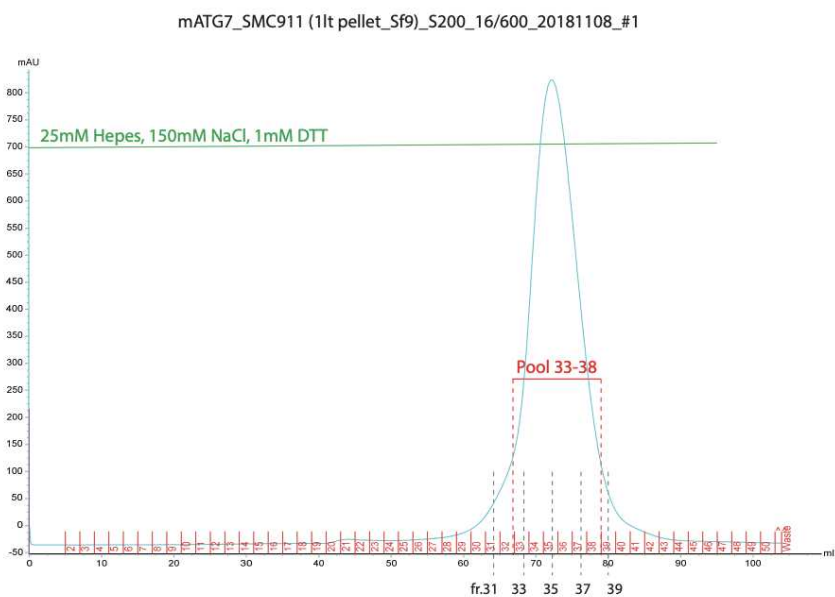

Chromatograph and Coomassie BB stained gel of Size Exclusion purification.

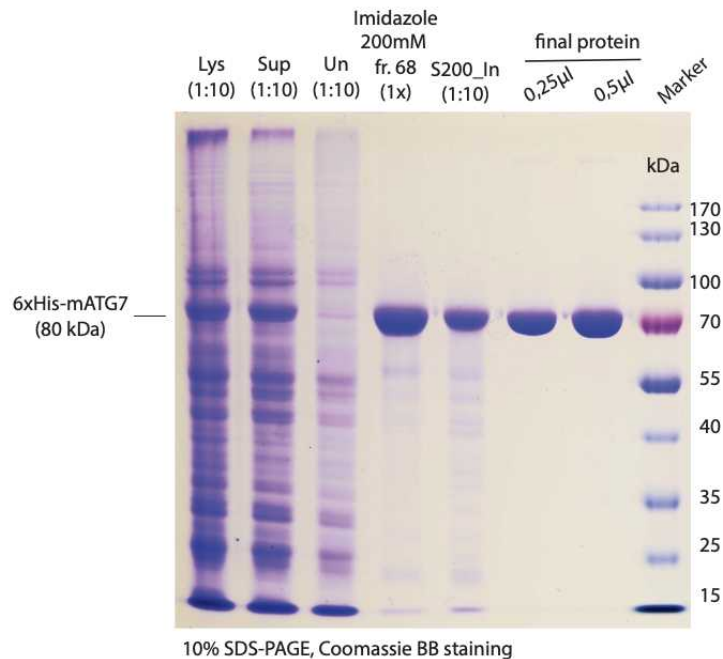

Summary gel of mATG7 purification.

#### Materials and Reagents

- Sf9 insect cells
- SF921 medium with antibiotics 100 IU/ml Penicillin and 100 µg/ml Streptomycin
- sterile cell culture hood
- 27°C shaker incubator
- sterile flasks/sterile pipettes
- Baculovirus for mATG7 (SMC911)
- douncer 40 mL

#### Materials for Protein Purification

- **Lysis Buffer:** 50mM Hepes pH=7.5; 300mM NaCl, Benzonase (1 µl/50ml lysis buffer), 10mM Imidazole, 2mM MgCl<sub>2</sub>, 1mM Dithiothreitol, 1xProtease Inhibitors/50 ml lysis buffer (EDTA-free CIP tablet, Roche).
- **Buffer A:** 50mM Hepes pH=7.5, 300mM NaCl, 10mM Imidazole (filtered and degassed) + 1mM β-mercaptoethanol
- **Buffer B:** 50mM Hepes pH=7.5, 300mM NaCl, 300mM Imidazole (filtered and degassed) + 1mM β-mercaptoethanol
- **SEC Buffer:** 25mM Hepes pH=7.5, 150mM NaCl (filtered and degassed) + 1mM Dithiothreitol.


Note: all purification buffers are filtered and degassed. Reducing agents (β-mercaptoethanol and Dithiothreitol) are added after degassing step.

**Columns:** - HT 5ml column (GE Healthcare)  
- S200\_16/60 (GE Healthcare)

**Gels:** 10% Poly-acrylamide SDS-gels

Please refer to the Safety Data Sheets (SDS) for health and environmental hazards.

## Infection/expression/harvest

- 1 Infect  **1 L culture** of Sf9 cells growing in Sf921 medium with Penicillin/Streptomycin at 1-1.5 mil/ml cells/volume at 99-100% viability in log phase with **Virus 1 (V1)** according to viral titer.  
(NOTE: protein yield from 1 lt culture is high, e.i. 32 mg/1lt culture).
- 2 Monitor infection and harvest cells when viability goes to 97-98%.




Always check under microscope: when all alive cells are brightly fluorescent and only few dead -> harvest!

## Protein Purification

45m

3

This section describes the procedure for His-Trap affinity purification followed by Size Exclusion Chromatography. All steps are to be executed at 4°C or on ice.

- 4 Resuspend cell pellet corresponding to  **1 L culture** in  **50 mL ice cold lysis buffer** ; gently stir at  **4 °C** until pellet dissolves avoiding bubbling.
- 5 Additionally, mechanically lyse the cells passing them through a pre-cooled douncer for 3x (10x pestle A followed by 10x pestle B).

6



45m

Clear the lysate by spinning it in a Beckman centrifuge, at 25K using a Ti45 Rotor for  **00:45:00** at  **4 °C** .

7

Inject the supernatant onto a 5ml HT column pre-equilibrated with Buffer A at 1ml/min flow rate to allow protein binding.

8




Wash the column for 5CV with Buffer A at 2 ml/min flow rate to remove unspecific bound proteins.

9

Elute protein through a step elution gradient in 50mM, 75mM, 100mM, 150mM, 200mM and 300mM Imidazole concentration. Perform the elution at 1ml/min flow rate.




10



Collect a sample from peak fractions of each elution step and check them on a SDS-PAGE (see gel in the "Guidelines" section). Pool and concentrate fractions corresponding to the peak containing the protein of interest (usually 200mM and/or 300mM Imid.) by spinning them down at  **4 °C** using a 30 kDa cut-off Amicon Filter to 2 ml

in a 5810R centrifuge (Eppendorf).


Keep centrifugation steps short (🕒 **00:05:00**) to avoid protein local concentration/aggregation on the filter.

- 11 Inject  **2 mL protein** onto a S200\_16/60 column at 4  **4 °C** pre-equilibrated in buffer containing 25mM Hepes pH=7.5, 150mM NaCl and 1mM DTT (see profile in the "Guidelines" section).
- 12 Check fractions on a 10% SDS-PAGE (see gel in the "Guidelines" section), pool and concentrate down those containing the protein of interest by spinning them at  **4 °C** using a 30kDa cut-off Amicon Filter in a 5810R centrifuge (Eppendorf). Protein elutes at around 72.5ml, a retention volume corresponding to that of a dimeric globular protein of ~150-160 kDa.

Keep centrifugation steps short to avoid protein precipitation.

- 13 

Measure protein absorbance at  $A_{280}$  is measured with Spectrophotometer against Size Exclusion Chromatography buffer.

- 14 Aliquot the protein, snap freeze it in liquid Nitrogen and store it at  **-80 °C**.

Usually protein activity is kept for 18 months when stored at -80°C.