



NOV 09, 2023

OPEN ACCESS



**DOI:**  
[dx.doi.org/10.17504/protocols.io.bwq4pdyw](https://dx.doi.org/10.17504/protocols.io.bwq4pdyw)

**Protocol Citation:** Steven B. Wells, Peter A. Szabo 2023. Preparation of Single Cell Suspension from Human Spleen Tissue. **protocols.io** <https://dx.doi.org/10.17504/protocols.io.bwq4pdyw>

**License:** This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

**Protocol status:** Working  
 We use this protocol and it's working

**Created:** Jul 20, 2021


## Preparation of Single Cell Suspension from Human Spleen Tissue

Peter A.  
 Steven B. Wells<sup>1</sup>, Szabo<sup>2</sup>

<sup>1</sup>Department of Systems Biology, Columbia University Irving Medical Center, New York, NY 10032, USA;

<sup>2</sup>Department of Microbiology and Immunology, Columbia University Irving Medical Center, New York, NY 10032, USA

Columbia

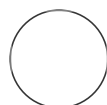

Steven B. Wells

### ABSTRACT







This protocol describes a method for the isolation of pan-lymphocytes, pan-myeloid cells, and progenitors from human spleen tissue. By providing defined media formulations, volumes at each step, and a defined dilution factor for density centrifugation, it yields consistent single-cell suspensions across samples.

### ATTACHMENTS

[dzhbkb587.pdf](#)

### MATERIALS










#### Materials:

-  Fisherbrand™ Sterile Syringes for Single Use Fisher Scientific Catalog #14955459
-  DPBS no calcium no magnesium Thermo Fisher Scientific Catalog #14190144
-  Penicillin-Streptomycin-Glutamine (100X) Thermo Fisher Catalog #10378016
-  Thermo Scientific™ Nunc™ 50mL Conical Sterile Polypropylene Centrifuge Tubes Fisher Scientific Catalog #12-565-271
-  Gibco™ IMDM (Iscoves Modified Dulbeccos Medium) Fisher Scientific Catalog #12-440-053
-  Gibco™ Fetal Bovine Serum qualified Australia Fisher Scientific Catalog #10-099-141

**Last Modified:** Nov 09, 2023

**PROTOCOL integer ID:** 51708

**Keywords:** Spleen, CD45, Lymphocytes, Myeloid, Isolation, Density gradient, Ficoll, Immune, 10x, scRNAseq, Flow cytometry, Leukocyte, Single cell suspension, T cell

-  UltraPure™ 0.5M EDTA pH 8.0 Fisher Scientific Catalog #15575020
- 100µM cell strainer (Fisher Scientific, Cat. No.: 50-146-1428)
-  Ficoll-Paque™ PLUS Media Fisher Scientific Catalog #45-001-749
-  Mr. Frosty™ Freezing Container Fisher Scientific Catalog #5100-0001
-  CryoStor CS10 100ML Fisher Scientific Catalog #NC9930384
-  Corning™ Externally Threaded Cryogenic Vials Fisher Scientific Catalog #09-761-71
-  5mL Falcon™ Round-Bottom Polypropylene Test Tubes Fisher Scientific Catalog #14-959-11A
-  Solution 13 AO – DAPI Chemometec Catalog #910-3013
-  NC-Slide A8™ box with 25 Slides Chemometec Catalog #942-0003
-  Falcon™ Plastic Disposable Transfer Pipets Fisher Scientific Catalog #1368050

#### Equipment:

- Centrifuge
- Cell Counter - NC-3000
- Surgical scissors
- Scale

## Preparing Medium and Buffer

- 1 Create the following IMDM-FBS-PSQ Media in a  500 mL bottle of IMDM by using the table below:

A	B	C	D
Component	Volume (mL)	Starting Conc.	Final Conc.*
IMDM	500	-	-
Penicillin-Streptomycin-Glutamine	5	100X	1X
FBS	50	100%	10%

Table 1.

\*Final Concentration is approximate.


- 2 Create the following DPBS-FBS-EDTA Solution in a bottle of DPBS without calcium and magnesium by using the table below:

A	B	C	D
Component	Volume (mL)	Starting Conc.	Final Conc.*
DPBS	500	-	-
FBS	25	100%	5%
EDTA	1	0.5M	1mM

Table 2.


\*Final Concentration is approximate.

## Tissue Dissociation

- 3 Add  $2 \pm 10\%$  grams of spleen tissue to a  50 mL centrifuge tube and record below.

\_\_\_\_\_g



- 4 Add  5 mL of DPBS-FBS-EDTA Solution to the spleen tissue and use a scissors to chop the tissue into a fine “mash”.



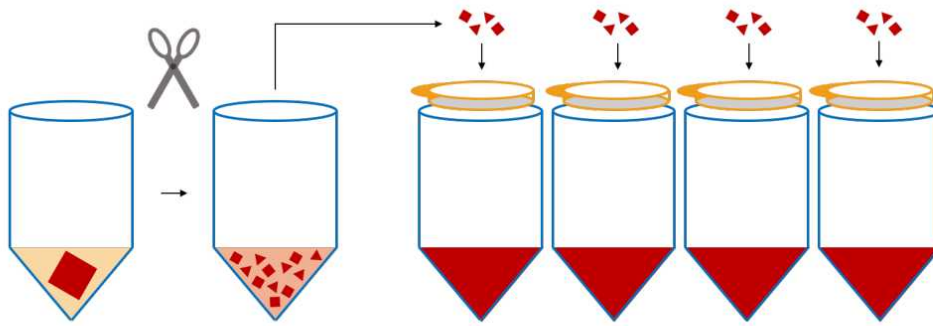






Figure 1. Steps 4.2.2 through 4.2.4.

- 5 Add  35 mL of DPBS-FBS-EDTA Solution to the mash of tissue, and distribute and filter the tissue over  100 micromolar ( $\mu\text{M}$ ) cell strainers above  50 mL tubes (about 4 filters/2 grams of tissue).

#### Note

NOTE: Cell yields and ease of pushing through the filter are increased by using multiple filters/gram of tissue, default to using more filters to decrease processing time, and increase yields.




- 6 Apply pressure with the black rubber bottom or the plastic end of a  10 mL syringe plunger to any remaining, partially digested tissue on the cell strainers, and intermittently wash through with DPBS-FBS-EDTA Solution from a transfer pipet – the aim is to push and wash through the tissue until only pink/white connective tissue remains. When finished, combine the tubes of cell suspension and proceed to the next section.

## Ficoll-Paque

50m

- 7 Centrifuge the cell suspensions for  00:10:00 at  400 x g at  20 °C .

10m



- 8 Remove the supernatants and combine the cell pellets down to a single  50 mL tube, top to  50 mL with  Room temperature DPBS-FBS-EDTA Solution.

- 9 Filter the cell suspension through a 100 micromolar (µM) cell strainer.
- 10 In two 50 mL tubes, layer 25 mL of cell suspension on top of 15 mL of Ficoll-Paque Media PLUS.
- 11 Spin for 00:20:00 at 1200 x g at 20 °C with 4 acceleration and 0 brake, evenly distribute 20m tubes across the entire rotor to prevent wobbling (use all four buckets if possible as opposed to just two).
- 12 Remove the mononuclear cell layer from both tubes with a transfer pipet and combine in one 50 mL tube. Add cold DPBS-FBS-EDTA Solution to a final volume of 50 mL and centrifuge the cell suspension for 00:10:00 at 400 x g at 4 °C.
- 13 Remove the supernatants and re-suspend the cell pellets in 50 mL cold DPBS-FBS-EDTA Solution and centrifuge the cell suspension for 00:10:00 at 120 x g at 4 °C.
- 14 Remove the supernatant and re-suspend the cell pellet in cold 10 mL IMDM-FBS-PSQ Media.

## Cell Count

- 15 Count cells, and viability by using the NC-3000 cell counter. Calculate total viable cells and record below:  
 cell number: \_\_\_\_\_ cells/mL, \_\_\_\_\_ % viable  
 final volume: \_\_\_\_\_ mL  
 $cell\ number\ (cells/mL) * viability(\%) * final\ volume(mL) = total\ viable\ cells$   
 Total Viable Cells: \_\_\_\_\_

## Freeze-down and QC

- 16 (Optional QC) Aliquot  $2 \times 10^6$  cells to a 5mL Falcon tube and place on ice for subsequent flow cytometric analysis.
- 17 Aliquot cells for analysis or experimentation, and then freeze down cells in up to  $3 \times 10^7$  aliquots using Cryostor CS10 Medium, a Mr. Frosty, and a  $-80^\circ\text{C}$  freezer (  1 mL -  1.5 mL aliquots, round down to the nearest 30 million cells and discard/freeze/use any left over cells). Record the number of vials frozen: \_\_\_\_\_.