



Version 2

Mar 02, 2021

Chlamydomonas reinhardtii nuclear transformation by electroporation. V.2

João Vitor Molino¹¹University of Zurich

1

Works for me

This protocol is published without a DOI.

Joao Vitor Molino
University of Zürich

SUBMIT TO PLOS ONE

ABSTRACT

This protocols describe the steps required for nuclear transformation of *Chlamydomonas reinhardtii* by electroporation.

[Here](#) you can find a video following the protocol.

EXTERNAL LINK

<https://doi.org/10.1371/journal.pone.0192433>

THIS PROTOCOL ACCOMPANIES THE FOLLOWING PUBLICATION

Molino JVD, Carvalho JCMd, Mayfield SP (2018) Comparison of secretory signal peptides for heterologous protein expression in microalgae: Expanding the secretion portfolio for *Chlamydomonas reinhardtii*. PLoS ONE 13(2): e0192433. doi: [10.1371/journal.pone.0192433](https://doi.org/10.1371/journal.pone.0192433)

EXTERNAL LINK

<https://doi.org/10.1371/journal.pone.0192433>

PROTOCOL CITATION

João Vitor Molino 2021. Chlamydomonas reinhardtii nuclear transformation by electroporation..
protocols.io
<https://protocols.io/view/chlamydomonas-reinhardtii-nuclear-transformation-b-7q2hmye>

MANUSCRIPT CITATION please remember to cite the following publication along with this protocol

Molino JVD, Carvalho JCMd, Mayfield SP (2018) Comparison of secretory signal peptides for heterologous protein expression in microalgae: Expanding the secretion portfolio for *Chlamydomonas reinhardtii*. PLoS ONE 13(2): e0192433. doi: [10.1371/journal.pone.0192433](https://doi.org/10.1371/journal.pone.0192433)

KEYWORDS

Microalgae, Recombinant, electroporation, plasmid

LICENSE

This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

CREATED

Sep 27, 2019

LAST MODIFIED

Mar 02, 2021

PROTOCOL INTEGER ID

28154

GUIDELINES

Cell density for harvesting is important to overall transformant yields. It should be aimed to mid-log phase cells.

*Transformation tested from $1-6 \times 10^6$ cells/mL - Worked.

MATERIALS TEXT

MATERIALS

 **MAX Efficiency™ Transformation Reagent for Algae Thermo**

Fisher Catalog #A24229

Step 2

SAFETY WARNINGS

- High voltage is used in the electroporation, use EPIs and avoid contact with electrodes on the electroporator.

DISCLAIMER:

DISCLAIMER – FOR INFORMATIONAL PURPOSES ONLY; USE AT YOUR OWN RISK

The protocol content here is for informational purposes only and does not constitute legal, medical, clinical, or safety advice, or otherwise; content added to protocols.io is not peer reviewed and may not have undergone a formal approval of any kind. Information presented in this protocol should not substitute for independent professional judgment, advice, diagnosis, or treatment. Any action you take or refrain from taking using or relying upon the information presented here is strictly at your own risk. You agree that neither the Company nor any of the authors, contributors, administrators, or anyone else associated with protocols.io, can be held responsible for your use of the information contained in or linked to this protocol or any of our Sites/Apps and Services.




BEFORE STARTING

- Prepare a ice bucket
- Separate cuvettes, keep them on ICE
- Allow linearized vectors to melt
- Keep transformation buffer on ICE/Fridge
- Prepare 50 mL centrifugal tubes with 10 mL TAP medium for recover stage

DNA Preparation

6h 30m

- 1 Digest a large enough amount of plasmid. The goal is to have a concentrated digested sample in the range of 250-700 ^{6h 30m} ng/uL.

1. Select the appropriate enzymes for linearization. Usually, restrictions sites in flanking position to the expression cassette.
2. Mix all components for digestion  **40 µg uncut vector** . Digest for  **06:00:00** at  **37 °C** .
3. Column purify digestion (Avoid gel purify, since vector backbone may helps to prevent intracellular DNAses action).
4. Quantitate by absorbance measuring (Nanodrop).

Component	Amount
10X Cutsmart NEB	6.0 uL
XbaI NEB 20 U/uL	3.0 uL
KpnI HF NEB 20 U/uL	3.0 uL
Plasmid 1219.9 ng/uL	40 uL
ddH2O, Molecular grade	8.0 uL

Typical reaction setup



Result example

[M]0.641 µg/µl Cutted vector

☐30 µl Final elution volume

☐19.230 µg total mass

Cells preparation

- 2 1. Aseptically inoculate ☐250 mL Tap media with wild type cells. Either by scraping cells of a plate with a inoculating loop or from a previous cultured cells.
2. Incubate at 25 °C , under constant shaking (~150-180 RPM) and light (60-80 µmols de photons/m²s) until a cell density from ☐3 x 10⁶ cells/mL to ☐6 x 10⁶ cells/mL is reached.
3. Pellet cells in centrifuge tubes. Separate culture in sufficient amount of sterile 50mL centrifuge tubes or larger



Pellet Cells

volume tubes, and centrifuge for 2000 x g, 25°C, 00:10:00 .

4. Gently resuspend cells at 3-6-10⁸ cells/mL in Transformation Buffer.

☒MAX Efficiency™ Transformation Reagent for Algae Thermo

Fisher Catalog #A24229



Culture at 3x10⁶ cells/mL usually yield 12-13 transformations.


Transformation

- 3 1. Add cutted vector to the bottom of the electroporation cuvette. Typically from ☐250 ng cutted vector to ☐1000 ng cutted vector
2. Add ☐250 µl resuspended cells (at approximately [M]3 x 10⁸ cells/mL) to each cuvette. Pipette up and down on DNA sample. Flick cuvette to mix DNA and cells. Shake cells to the bottom of the cuvette. Also add no DNA control (Elution buffer or water).



Cell suspension ready for electroporation

3. Incubate cells with DNA ⏱ **On ice** for ⌚ **00:10:00**
4. Wipe cuvette (to remove condensed water) and electroporate (Table Electroporation).
5. Let it recover for ⌚ **00:10:00** on the cuvette
6. Add cells to 📄 **10 mL TAP/40mM sucrose, pH 7.0** inside sterile 50mL centrifuge tubes. Gently transfer cells from cuvette to TAP/40 mM sucrose. Rinse cuvette with TAP/40 mM sucrose to transfer any remaining cells.
7. Incubate at 🌡 **Room temperature** on rocker or shaker at 50 rpm ⌚ **Overnight** ambient light.

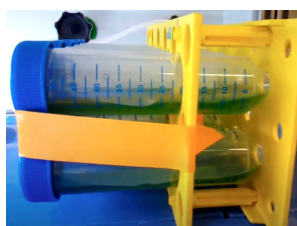

Recover step of cells in shaker, low mixing.

8. Pellet cells by centrifuging for ⌚ **2000 x g, 25°C, 00:10:00**
9. Aseptically pour off supernatant. Add 📄 **300 µl TAP/40mM sucrose** to pellet. Gently re-suspend cells and pipette onto 2 plates with appropriate antibiotics. ie. 📄 **200 µl cells resuspended** per plate, and let it dry aseptically without plate cover.
10. Spread cells evenly over the plate with an inoculation loop. Avoid spreading to the borders.
11. Use parafilm to block evaporation and place plates under constant light (60 µmol photons/m²s), 🌡 **25 °C** .
Colonies should be visible in 5-7 days.

Table Electroporation - Settings

Voltage	800 V
Time Constant	20 ms
Cuvette gap	4 mm

Gene Pulser Xcell Electroporation Systems
Electroporator

Biorad

1652660




Typical output after electroporation

Time constant (ms)	Voltage (V)	Capacitance (uF)	Resistance (Ohms)
20.1	788	50	650
20.4	789	50	600
19.8	789	50	550



Green colonies should appear in the plate as in the pictures below.

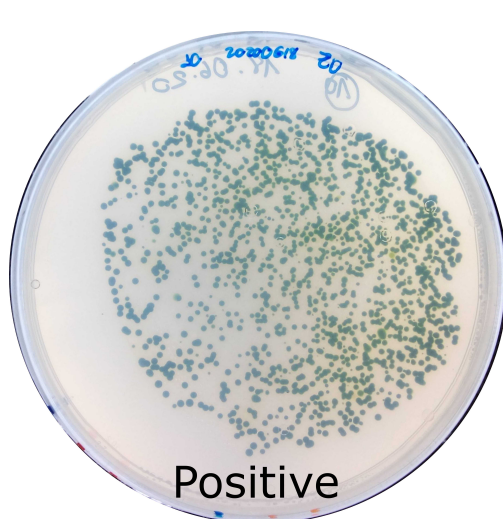

Positive

Positive result



Negative

Negative control