

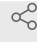


Version 1 ▾

Sep 09, 2022

# Astrocyte extraction from brain organoids

## V.1

[gustavo.parfitt](#)<sup>1</sup><sup>1</sup>Icahn School of Medicine1 *Works for me* Share[dx.doi.org/10.17504/protocols.io.261ge364wl47/v1](https://dx.doi.org/10.17504/protocols.io.261ge364wl47/v1)[Ahfeldt Lab](#) [gustavo.parfitt](#)

### ABSTRACT

Protocol for astrocyte extraction from brain organoids.

### DOI

[dx.doi.org/10.17504/protocols.io.261ge364wl47/v1](https://dx.doi.org/10.17504/protocols.io.261ge364wl47/v1)

### PROTOCOL CITATION

[gustavo.parfitt](#) 2022. Astrocyte extraction from brain organoids. **protocols.io**  
<https://protocols.io/view/astrocyte-extraction-from-brain-organoids-cgdsts6e>



### LICENSE

\_\_\_\_\_ This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

### CREATED

Sep 09, 2022

### LAST MODIFIED


Sep 09, 2022

### PROTOCOL INTEGER ID



69778

- 1 Coat a 6 well plate coat with gelatin 0.1% or  
[Geltrex LDEV Free hESC Quality 5 ml Thermo Fisher Scientific Catalog #A1413302](#)

1:100

- 2 Collect 10-20 spheres and place in a 6 well plate (day 40+ spheres)
- 3 Wash in PBS twice
- 4 Aspirate the PBS
- 5 Add  **1 mL** of 5m  

TryPLE™ Select Enzyme (1X), no phenol red Thermo  
Fisher Catalog #12563011

for  
 **00:05:00**
- 6 Triturate using glass a pipette (2-3 up and down)
- 7 Transfer the cells and tissue (even the chunks, the trituration is more effective if you have chunks of tissue) to the coated 6 well plate
- 8 Aspirate the extra TrypLe
- 9 Add  **3 mL** of astrocyte media (<https://www.sciencellonline.com/astrocyte-medium.html>)
- 10 Half media change media every 3 days or until it turns orange. 3d

- 11 Keep changing media every 3d until you observe good amount of cell attached to the plate (may vary according to the patterning)
- 12 Passage the astrocytes to a 15 cm plate in Astrocyte media once they reach 80% confluency
- 13 For the experimental protocols after passage cell will be incubated in maturation media (TBD)