


Apr 30, 2024

FUNDIS ONT Sequencing IT/Compute Pop!_OS 22.04 Setup

 Forked from [ONT Sequencing IT/Compute Pop!_OS 22.04 Setup](#)

This protocol is a draft, published without a DOI.

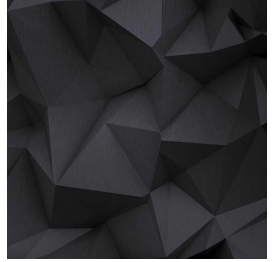

Stephen Douglas Russell^{1,2}, Harte Singer³

¹Mycota Lab; ²The Hoosier Mushroom Society; ³FUNDIS

The Hoosier Mushroom S...



Harte Singer

FUNDIS

OPEN  ACCESS



Protocol Citation: Stephen Douglas Russell, Harte Singer 2024. FUNDIS ONT Sequencing IT/Compute Pop!_OS 22.04 Setup.

protocols.io <https://protocols.io/view/fundis-ont-sequencing-it-compute-pop-os-22-04-setu-dcwy2xfw>

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working

We use this protocol and it's working

Created: January 06, 2024

Last Modified: April 30, 2024

Protocol Integer ID: 99000

Keywords: Oxford Nanopore Technologies, DNA sequencing, DNA barcoding, MinION, Flongle, system76, linux, POP! OS



Abstract

These are the steps to setting up a Linux system running Ubuntu for ONT Nanopore sequencing, basecalling, and consensus-building of fungal amplicon sequences. Adapted from dx.doi.org/10.17504/protocols.io.14egn7kzmv5d/v3

Includes an R-based run QC script



Preparing a new CPU for MinION Sequencing

- 1 Ideally you will be running a linux machine with a Ubuntu distro with an NVIDIA graphics card with at least 8GB of ram. (I am running a 4070Ti) You should have at least 32GB of system ram and a new processor.

Minimum IT requirements for MinION from ONT:



2

Install CUDA toolkit - <https://developer.nvidia.com/cuda-toolkit> :

Command

```
wget
https://developer.download.nvidia.com/compute/cuda/repos/ubuntu2204/x8
6_64/cuda-ubuntu2204.pin
sudo mv cuda-ubuntu2204.pin /etc/apt/preferences.d/cuda-repository-
pin-600
wget
https://developer.download.nvidia.com/compute/cuda/11.7.0/local_instal
lers/cuda-repo-ubuntu2204-11-7-local_11.7.0-515.43.04-1_amd64.deb
sudo dpkg -i cuda-repo-ubuntu2204-11-7-local_11.7.0-515.43.04-
1_amd64.deb
sudo cp /var/cuda-repo-ubuntu2204-11-7-local/cuda-*-keyring.gpg
/usr/share/keyrings/
sudo apt-get update
sudo apt-get -y install cuda
```

- 3 Install Boost



Command

```
sudo apt install libboost-all-dev
```

- 4 Install MinKNOW from <https://community.nanoporetech.com/downloads?from=support>

Specifically download the GPU version of the software with the following name. It is a .deb file and your computer will install it when you double-click it

MinKNOW Software for the MinION Mk1B and the PromethION 2 Solo (P2 Solo) Devices

- 5 Install Guppy from the "Archived Software" section at the bottom of this page
<https://community.nanoporetech.com/downloads?from=support>

Guppy v6.5.7

Make sure you install the Ubuntu 20 GPU version, just double click the installer and it should work.

- 6 Per this document: https://denbi-nanopore-training-course.readthedocs.io/en/latest/read_qc/MinionQC.html
Install R: <https://cran.r-project.org/>
Install MinionQC: https://github.com/roblanf/minion_qc

Install R:

**Command**

```
# update indices
sudo apt update -qq
# install two helper packages we need
sudo apt install --no-install-recommends software-properties-common
dirmngr
# add the signing key (by Michael Rutter) for these repos
# To verify key, run gpg --show-keys
/etc/apt/trusted.gpg.d/cran_ubuntu_key.asc
# Fingerprint: E298A3A825C0D65DFD57CBB651716619E084DAB9
wget -qO- https://cloud.r-project.org/bin/linux/ubuntu/marutter_pubkey.asc | sudo tee -a
/etc/apt/trusted.gpg.d/cran_ubuntu_key.asc
# add the R 4.0 repo from CRAN -- adjust 'focal' to 'groovy' or
'bionic' as needed
sudo add-apt-repository "deb https://cloud.r-project.org/bin/linux/ubuntu $(lsb_release -cs)-cran40/"
```

Command

```
sudo apt install --no-install-recommends r-base
```

Open an R terminal window for the following command

Command

```
install.packages(c("data.table",  
                  "futile.logger",  
                  "ggplot2",  
                  "optparse",  
                  "plyr",  
                  "readr",  
                  "reshape2",  
                  "scales",  
                  "viridis",  
                  "yaml"))
```

7 Install Bioconductor:

In an R command window:

Command

```
if (!require("BiocManager", quietly = TRUE))  
  install.packages("BiocManager")  
BiocManager::install(version = "3.15")
```

8 Install Anaconda:

from: <https://www.digitalocean.com/community/tutorials/how-to-install-the-anaconda-python-distribution-on-ubuntu-22-04>

**Command**

```
cd /tmp
```

Command

```
curl https://repo.anaconda.com/archive/Anaconda3-2022.05-Linux-  
x86_64.sh --output anaconda.sh
```

Command

You can now verify the data integrity of the installer with cryptographic hash verification through the SHA-256 checksum. You'll use the sha256sum command along with the filename of the script:

```
sha256sum anaconda.sh
```

Expected result

You'll receive output that looks similar to this:

```
fedf9e340039557f7b5e8a8a86affa9d299f5e9820144bd7b92ae9f7ee08ac6  
0  anaconda.sh
```

**Command**

```
bash anaconda.sh
```

Expected result

Press ENTER/yes as needed

```
Welcome to Anaconda3 2021.11
```

```
In order to continue the installation process, please review  
the license  
agreement.  
Please, press ENTER to continue  
>>>
```

Command

```
source ~/.bashrc
```

Command

```
conda list
```




Expected result

```
# packages in environment at /home/user/anaconda3:
#
# Name                      Version                      Build
Channel
_ipyw_jlab_nb_ext_conf     0.1.0                        py39h06a4308_0
_libgcc_mutex              0.1                          main
_openmp_mutex              4.5                          1_gnu
alabaster                  0.7.12                       pyhd3eb1b0_0
anaconda                   2022.05                       py39_0
```

Command

```
conda search "^python$"
```

Command

Verify Python is installed

```
python --version
```

**Command****Install bioconda**

```
conda install -c bioconda seqkit
```

- 9 Install NGSspeciesID: <https://github.com/ksahlin/NGSpeciesID>

Command

```
conda create -n NGSspeciesID python=3.6 pip  
conda activate NGSspeciesID
```

Command

```
conda install --yes -c conda-forge -c bioconda medaka==0.11.5  
openblas==0.3.3 spoa racon minimap2  
pip install NGSspeciesID
```



Command

```
conda activate NGSpeciesID
```

Command

Test the install

```
mkdir test_ngspeciesID  
cd test_ngspeciesID
```

Command

Download the test fastq file called "sample_h1.fastq" (filesize 390kb)

```
curl -LO  
https://raw.githubusercontent.com/ksahlin/NGSpeciesID/master/test/sample_h1.fastq
```



Command

Run the **NGSpecies** command on test file. Outputs will be saved in **"/test_ngspeciesID/sample_h1/"**, where the final polished consensus file ("consensus.fasta") is located in the **"/test_ngspeciesID/sample_h1/medaka_cl_id_"** directory.

```
NGSpeciesID --ont --fastq sample_h1.fastq --outfolder ./sample_h1 --  
consensus --medaka
```

- 10 You should now be ready to begin sequencing runs.

Installing Dorado

- 11 We will eventually be moving to Dorado basecalling.

Installing Dorado:



```
sudo -i

sudo apt-get update && apt-get install -y --no-install-recommends \
    curl \
    git \
    ca-certificates \
    build-essential \
    nvidia-cuda-toolkit \
    libhdf5-dev \
    libssl-dev \
    libzstd-dev \
    cmake \
    autoconf \
    automake

git clone https://github.com/nanoporetech/dorado.git dorado
cd dorado
cmake -S . -B cmake-build
cmake --build cmake-build --config Release -j
ctest --test-dir cmake-build

pip install pre-commit
pre-commit install
```