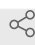




Jul 20, 2022

🌐 Fluorescence size exclusion chromatography (FSEC) from ATP13A2 microsomes

Sue Sim¹, eunyong_park¹¹University of California, Berkeley1 *Works for me* Sharedx.doi.org/10.17504/protocols.io.81wgb6bpolk/v1 Sue Sim

ABSTRACT

Using fluorescence size exclusion chromatography (FSEC) to analyze ATP13A2 expression in microsomes

DOI

dx.doi.org/10.17504/protocols.io.81wgb6bpolk/v1

PROTOCOL CITATION

Sue Sim, eunyong_park 2022. Fluorescence size exclusion chromatography (FSEC) from ATP13A2 microsomes. **protocols.io**
<https://dx.doi.org/10.17504/protocols.io.81wgb6bpolk/v1>



LICENSE

————— This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

CREATED

Jul 17, 2022

LAST MODIFIED

Jul 20, 2022

PROTOCOL INTEGER ID

66881

MATERIALS TEXT

Lysis Buffer

50 mM Tris pH 7.5

200 mM NaCl

1 mM EDTA

1 mM DTT

10% glycerol

Plus protease inhibitors (5 µg/mL aprotinin, 5 µg/mL leupeptin, 1 µg/mL pepstatin A, and 2 mM PMSF)

Running Buffer

25 mM Tris pH 7.5

100 mM NaCl

1 mM EDTA

0.03% DDM/ 0.006% CHS

- 1 Thaw 100 µg microsomes on ice
- 2 Resuspend microsomes in Lysis Buffer and final volume 1% DDM/0.2% CHS (1X pellet, 3X Lysis Buffer, 1X 5% DDM/1% CHS) at **4 °C**

2.1 DDM: n-dodecyl-β-D-maltopyranoside (Anatrace)

2.2 CHS: cholesteryl hemisuccinate (Anatrace)

- 3 Solubilize by rotating end-over-end for 2 h at 4C

- 4 Clarify lysate by spinning at **17000 x g, 4°C, 01:00:00**

1h

- 5 Equilibrate Superose 6 column, connected to an HPLC system, with Running Buffer

- 6 Inject 85 uL of lysate into HPLC
- 7 Monitor elution of GFP-tagged ATP13A2 constructs using a fluorometer (λ_{ex} = 475 nm; λ_{em} = 510 nm; with a fixed gain) connected to the system