

Jan 07, 2025

Inhibitor removal from DNA extracts

DOI

dx.doi.org/10.17504/protocols.io.rm7vzbp54vx1/v1



Dominik Buchner¹, Marie Borowski¹

¹University of Duisburg-Essen, Aquatic Ecosystem Research

Aquatic Ecosystem Resea...



Dominik Buchner

University of Duisburg-Essen, Aquatic Ecosystem Research

OPEN  ACCESS



DOI: dx.doi.org/10.17504/protocols.io.rm7vzbp54vx1/v1

Protocol Citation: Dominik Buchner, Marie Borowski 2025. Inhibitor removal from DNA extracts. **protocols.io**
<https://dx.doi.org/10.17504/protocols.io.rm7vzbp54vx1/v1>

License: This is an open access protocol distributed under the terms of the **Creative Commons Attribution License**, which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working

We use this protocol and it's working


Created: April 21, 2023

Last Modified: January 07, 2025

Protocol Integer ID: 80902



Abstract

This protocol describes how to remove inhibitory substances such as humic substances from DNA that has already been extracted. The protocol is formulated for an initial input of  100 μL of extracted DNA however it can be adjusted to any other volume as well with minor changes. A lot of the buffers can be found in the following patent <https://patents.google.com/patent/US7459548B2/en>

Guidelines

Follow general lab etiquette. Wear gloves to prevent contaminating the samples. Clean the workspace before starting with 80% EtOH.

Materials

Materials required:

Below all materials needed for the protocol are listed. Vendors and part numbers are listed but interchangeable depending on the supply situation.

Chemicals:

Sodium phosphate dibasic

☒ Sodium phosphate dibasic **Merck MilliporeSigma (Sigma-Aldrich) Catalog #S0876-100G**

Guanidinium thiocyanate ☒ Guanidinium thiocyanate **Fisher Scientific Catalog #10503345**

Sodium phosphate monobasic Sodium phosphate monobasic

☒ Sodium phosphate monobasic **Merck MilliporeSigma (Sigma-Aldrich) Catalog #S0751-100G**

SDS ultrapure ☒ Sodium dodecyl sulfate **Diagonal Catalog #A1112.0500**

Sodium chloride ☒ Sodium chloride **Fisher Scientific Catalog #10616082**

Tris ultrapure 99.9% ☒ Tris ultrapure 99.9% **Diagonal Catalog #A1086.1000**

Hydrochloric acid fuming 37%

☒ Hydrochloric acid fuming 37% **Merck MilliporeSigma (Sigma-Aldrich) Catalog #1003171011**

Ammonium acetate ☒ Ammonium acetate **Carl Roth Catalog #7869.2**

Aluminium ammonium sulfate dodecahydrate

☒ Aluminium ammonium sulfate dodecahydrate **Merck MilliporeSigma (Sigma-Aldrich) Catalog #A2140-500G**

Guanidine hydrochloride ☒ Guanidine hydrochloride **Fisher Scientific Catalog #10543325**

Acetic acid ☒ Acetic acid **Carl Roth Catalog #7332.1**

Ethanol absolute ☒ Ethanol absolute 99.8% p.a. **Carl Roth Catalog #9065.1**

Labware:

2 mL centrifuge tubes ☒ Reaction tube, 2 mL, PP **Sarstedt Catalog #72.691**

1.5 mL centrifuge tubes ☒ Reaction tube, 1.5 ml, PP **Sarstedt Catalog #72.690.001**

The EconoSpin® All-In-One DNA Only Mini Spin Column



☒ The EconoSpin® All-In-One DNA Only Mini Spin Column **Epoch Life Science Catalog #1920-250**



Stock solutions:




🧴 1 L SDS stock solution [M] 10 Mass / % volume

▪ Add 🧴 100 g SDS ultrapure to a beaker











- Adjust volume to  1 L with ddH₂O
- Sterilize by filtering and store at  Room temperature






 1 L sodium chloride stock solution  5 Molarity (M)

- Add  292.2 g sodium chloride to a beaker
- Adjust volume to  1 L with ddH₂O
- Sterilize by filtering and store at  Room temperature






 1 L Tris stock solution  1 Molarity (M)  8

- Add  121.14 g Tris ultrapure 99.9% to a beaker
- Adjust volume to  800 mL with ddH₂O
- Adjust pH to  8 with HCl
- Adjust volume to  1 L with ddH₂O
- Sterilize by filtering and store at  Room temperature



 500 mL sodium acetate stock solution  3 Molarity (M)  5

- Add  123 g sodium acetate to a beaker
- Adjust volume to  400 mL with ddH₂O
- Adjust pH to  5 with acetic acid
- Adjust volume to  500 mL with ddH₂O
- Sterilize by filtering and store at  Room temperature

 1 L Tris stock solution  1 Molarity (M)  7.5

- Add  121.14 g Tris ultrapure 99.9% to a beaker
- Adjust volume to  800 mL with ddH₂O
- Adjust pH to  7.5 with HCl
- Adjust volume to  1 L with with ddH₂O
- Sterilize by filtering and store at  Room temperature

 1 L Tris stock solution  1 Molarity (M)  8.5

- Add  121.14 g Tris ultrapure 99.9% to a beaker
- Adjust volume to  800 mL with ddH₂O

- Adjust pH to $\text{pH } 8.5$ with HCl
- Adjust volume to 1 L with ddH₂O
- Sterilize by filtering and store at Room temperature

1 L wash buffer stock solution ($50 \text{ millimolar (mM) Tris}$) $\text{pH } 7.5$

- Add $50 \text{ mL Tris stock solution}$ $\text{pH } 7.5$ to a beaker
- Adjust volume to 1 L with ddH₂O
- Sterilize by filtering and store at Room temperature

Working solutions:

$500 \text{ mL bead-beating solution}$ ($180 \text{ millimolar (mM) sodium phosphate}$, $120 \text{ millimolar (mM) guanidinium thiocyanate}$) $\text{pH } 8$

- Add $12.8 \text{ g sodium phosphate dibasic}$ to a beaker
- Add $7.1 \text{ g guanidinium thiocyanate}$
- Adjust volume to 490 mL with ddH₂O
- Adjust pH to $\text{pH } 8$ by adding sodium phosphate monobasic
- Adjust volume to 500 mL with ddH₂O
- Sterilize by filtering and store at Room temperature

$500 \text{ mL lysis solution}$ ($150 \text{ millimolar (mM) sodium chloride}$, $4 \text{ Mass / \% volume SDS}$, $500 \text{ millimolar (mM) Tris}$) $\text{pH } 8$

- Add 200 mL of $10 \text{ Mass / \% volume SDS stock solution}$ to a beaker
- Add 15 mL of $5 \text{ Molarity (M) sodium chloride stock solution}$
- Add 250 mL of $1 \text{ Molarity (M) Tris stock solution}$ $\text{pH } 8$
- Adjust volume to 500 mL with ddH₂O
- Sterilize by filtering and store at Room temperature

$500 \text{ mL ammonium acetate buffer}$ ($130 \text{ millimolar (mM) ammonium acetate}$)

- Add $5 \text{ g ammonium acetate}$ to a beaker
- Adjust volume to 500 mL with ddH₂O
- Sterilize by filtering and store at Room temperature



🧪 500 mL inhibitor removal solution ([M] 120 millimolar (mM) aluminum ammonium sulfate dodecahydrate)

- Add 🧪 27.2 g aluminium ammonium sulfate dodecahydrate to a beaker
- Adjust volume to 🧪 500 mL with ddH₂O
- Sterilize by filtering and store at 🌡 Room temperature

🧪 500 mL DNA binding buffer ([M] 2.5 Molarity (M) Guanidine hydrochloride , [M] 80 % (v/v) ethanol , [M] 0.05 % (v/v) Tween 20 , [M] 120 millimolar (mM) sodium acetate) 📏 pH 5

- Add 🧪 119.4 g guanidine hydrochloride to a beaker
- Fill up to 🧪 400 mL ethanol
- Add 🧪 20 mL [M] 3 Molarity (M) sodium acetate stock solution 📏 pH 5
- Adjust volume to 🧪 500 mL with ddH₂O
- Sterilize by filtering and store at 🌡 Room temperature

🧪 1 L wash buffer ([M] 10 millimolar (mM) Tris , [M] 80 % (v/v) Ethanol) 📏 pH 7.5

- Add 🧪 200 mL was buffer stock solution
- Adjust volume to 🧪 1 L with Ethanol absolute
- Sterilize by filtering and store at 🌡 Room temperature

🧪 1 L elution buffer ([M] 10 millimolar (mM) Tris) 📏 pH 8.5

- Add 🧪 10 mL Tris stock solution 📏 pH 8.5 to a beaker
- Adjust volume to 🧪 1 L with ddH₂O
- Sterilize by filtering and store at 🌡 Room temperature


Before start

Make sure all buffers are prepared before starting.

Inhibitor removal from DNA extracts



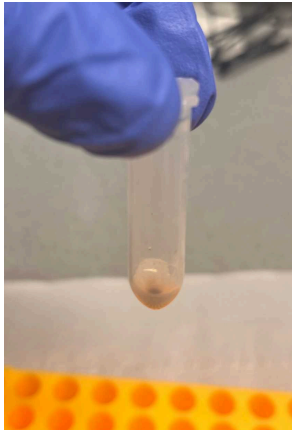

31m



- 1 Prepare  100 μL of sample in 2 mL tubes.




Sample contaminated with inhibitory substance (e.g. humic acids)

Note


This protocol is designed for a sample volume of  100 μL . Differing volumes can be used, however, all the following volumes have to be adjusted accordingly. For volumes smaller than  100 μL we recommend to simply fill up the sample to  100 μL with water and then proceed with the protocol. If a sample volume of  200 μL is to be processed, all volumes have to be doubled.

- 2 Add  345 μL bead-beating solution and  60 μL lysis solution .
Vortex shortly.




10m

- 3  10000 x g, 20°C , 00:03:00 . Transfer all of the supernatant to a new tube.

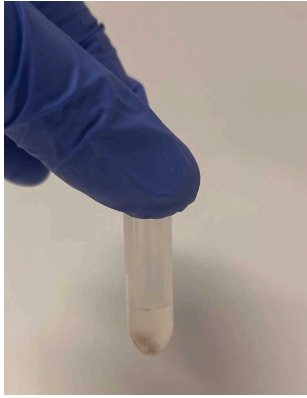
3m



Try to avoid the the pellet if any is formed.




- 4 Add  125 μ L ammonium acetat buffer , vortex shortly and incubate at  4 °C for  00:05:00 .

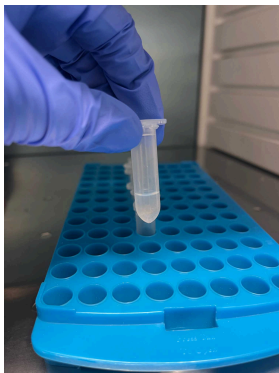5m




- 5  10000 x g, 20°C, 00:01:00 . Transfer the supernatant to a new tube.

1m

- 6 Add  100 μ L of inhibitor removal buffer . A precipitate may form. Vortex shortly, incubate at  4 °C for  00:05:00

5m

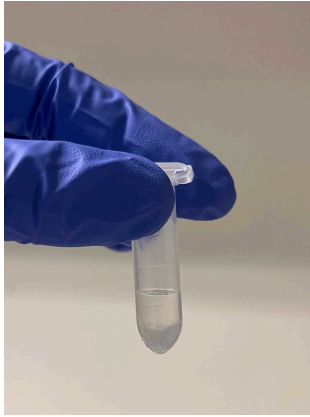

7 10000 x g, 20°C, 00:01:00 . Transfer 600 μ L of the supernatant to a new tube.

1m

8 Add 1200 μ L DNA binding buffer . Vortex to mix.

9 Load 650 μ L of the mixture to a mini spin column (e.g. Epoch Life Science).








10 10000 x g, 20°C, 00:00:30 . Discard the flow-through. Repeat two times to bind the complete sample volume.

30s

11 Add 500 μ L wash buffer . 10000 x g, 20°C, 00:00:30 to wash the column. Discard the flow-through.

30s



- 12  10000 x g, 20°C, 00:01:00 to dry the column. Transfer the spin column to a clean 1.5mL microcentrifuge tube. 1m
- 13 Add  50 µL elution buffer . Incubate for  00:03:00 at  Room temperature . 3m
- 14  10000 x g, 20°C, 00:01:00 to eluate the DNA. DNA eluate should be completely colorless and ready to go for downstream analysis. 1m

